



START YOUR DAY
WITH ENVISION



ENVISION NEXT

BASIC REPORTING

Short Manual Document

TABLE OF CONTENTS

	USER LOGIN	1
	REGISTRATION	2-3
	• NEW REGISTRATION	2
	• QUICK PATIENT REGISTRATION	3
	APPOINTMENT	4
	IMAGING	6-8
	• WORKFLOW SETTING TO SKIP IMAGING SECTION	6
	• IMAGING SECTION	7
	REPORT	9-16
	• START REPORTING	9
	• REPORTING CONFIGURATION	13
	• SEMI-STRUCTURE REPORT FINDINGS	14
	• SIMPLE REPORT FINDINGS	15
	• SAVE/ PREVIEW/ PRINT/ DOWNLOAD	16


1.USER LOGIN

1.1 SIGNING IN ENVISION NEXT

1.1.1 Log in with your username instead of your email

Login

Note : Before you can log in to Envision Next, you will need to gather your username and password from an administrator. An administrator will have to register you from within the system and provide you with your proper credentials

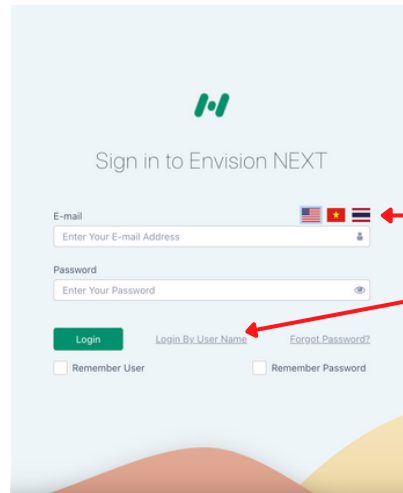

Sign in to Envision NEXT

E-mail  LANGUAGE SELECTOR

Enter Your E-mail Address


Password

Enter Your Password

[Login By User Name](#) LOGIN BY USERNAME/ EMAIL [Forgot Password?](#)

☐ Remember User ☐ Remember Password

1.1.2 After clicking the Login button, the user will get a Confirmation message for replacing another session, if there is any.

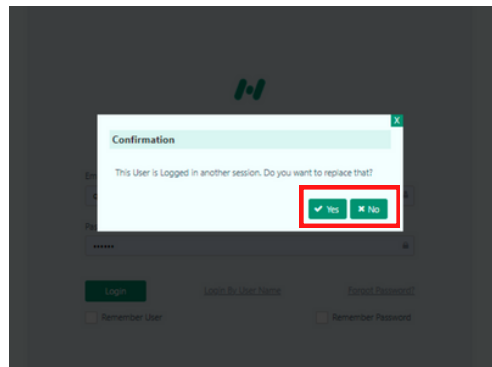


Confirmation

This User is Logged in another session. Do you want to replace that?

2. REGISTRATION

2.1 NEW REGISTRATION

2.1.1 From the Main Menu on the left, you can start a Patient section:

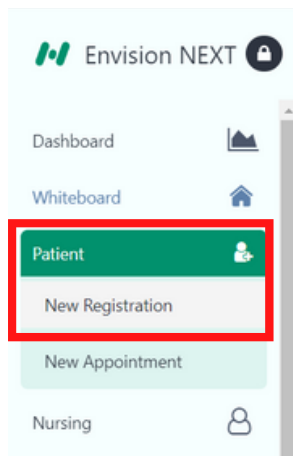


2.1.2 The New Registration form has 4 sections that you can fill up –

- Patient Information
- Additional
- Emergency Contact
- Patient Allergy


You will need to fill up the (minimum) mandatory fields in the form before you can successfully complete patient registration.

Note : Mandatory fields to be filled are marked with red asterisks (*).

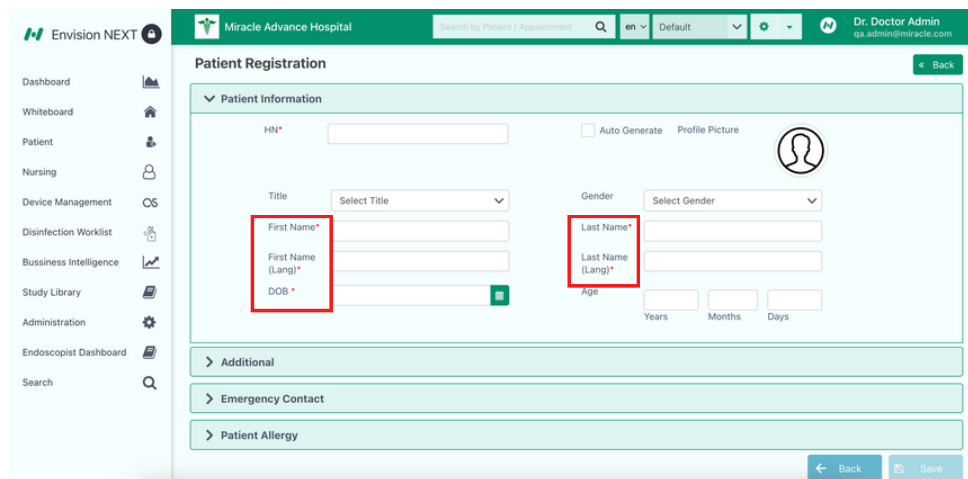

2.1.3 After that the software will show the Patient Dashboard.


The Patient Dashboard is where you can access information about the patient, edit the information and view useful details such as Invoices and History.

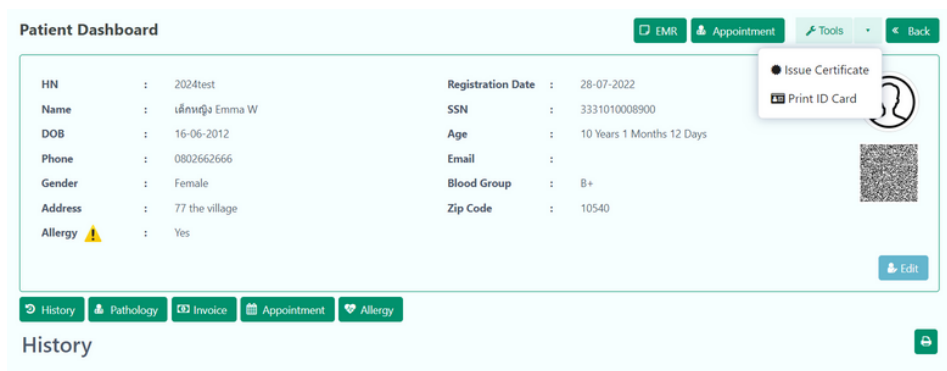


2. REGISTRATION

2.2 QUICK PATIENT REGISTRATION OPTION FROM THE APPOINTMENT SCREEN

The user can enable Quick registration by going to the Main Menu on the left, selecting Administration and selecting Update Organization, Next selecting the Org Setting tab, Then selecting the Appointment tab and clicking Enable Quick Registration.

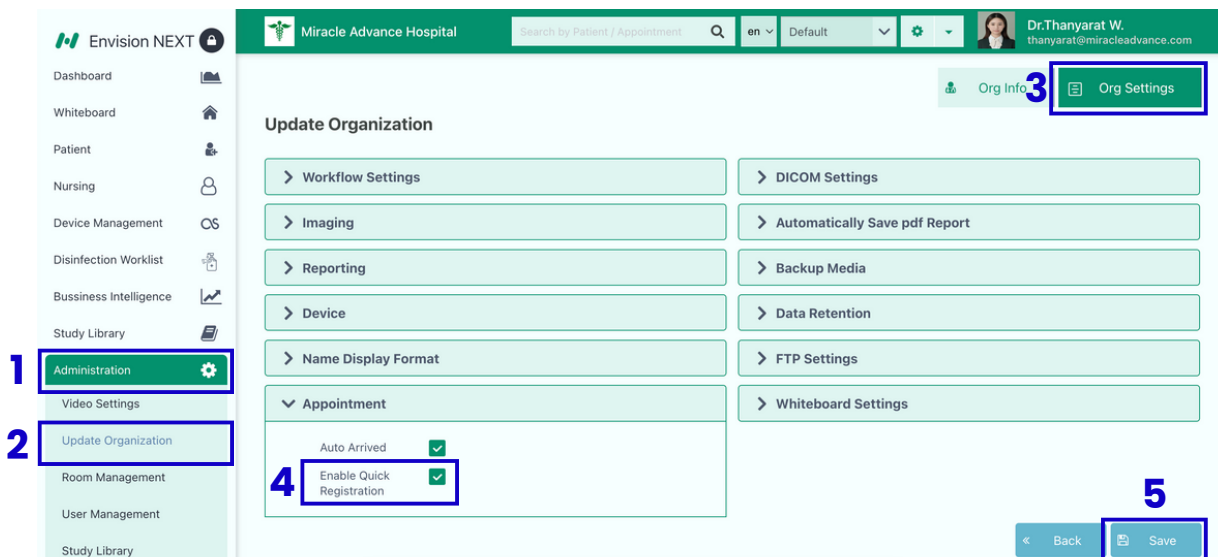


2.2.1 From the Main Menu on the left, selecting the “New Appointment”. The software will show the Appointment calendar, You can select a slot [example: red area] (by clicking on it).

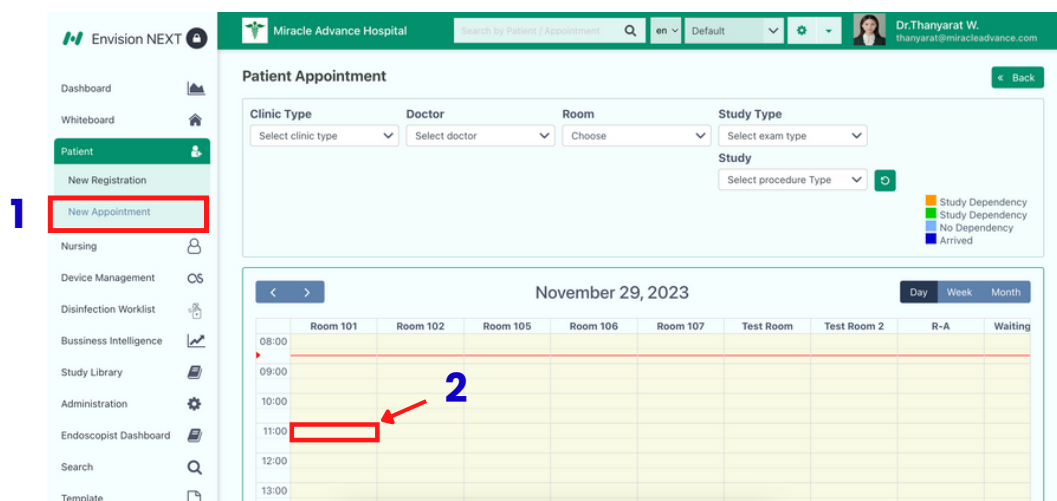

2.2.2 After select a slot, The software will show the Quick Patient Registration window. You will need to fill up the (minimum) mandatory fields in the form before you can successfully complete patient registration.


Note : Mandatory fields to be filled are marked with red asterisks (*).

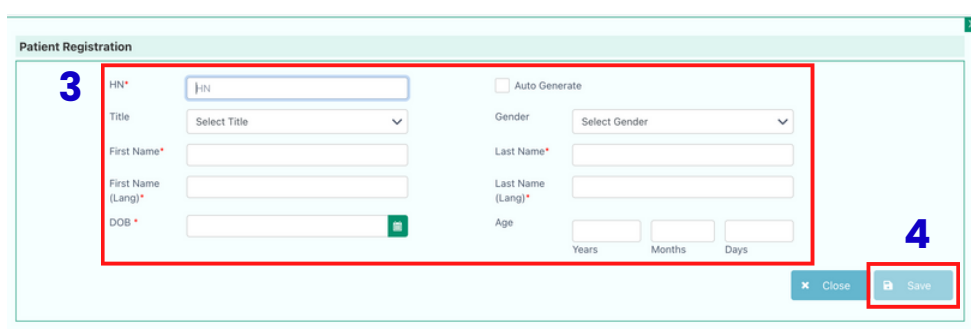


3. APPOINTMENT

3.1 NEW APPOINTMENT


3.1.1 From the Main Menu on the left, you can start a **Patient** section:
There are 2 submenus : New Registration / New Appointment. Click on **the “New Appointment”**

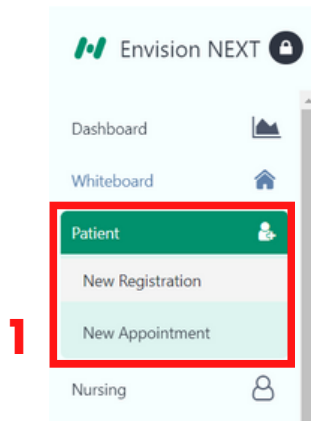

3.1.2 **The Appointment calendar** is where you can select a slot [example: red area] (by clicking on it) and set an appointment for the patient from the Patient Dashboard that you just came from.

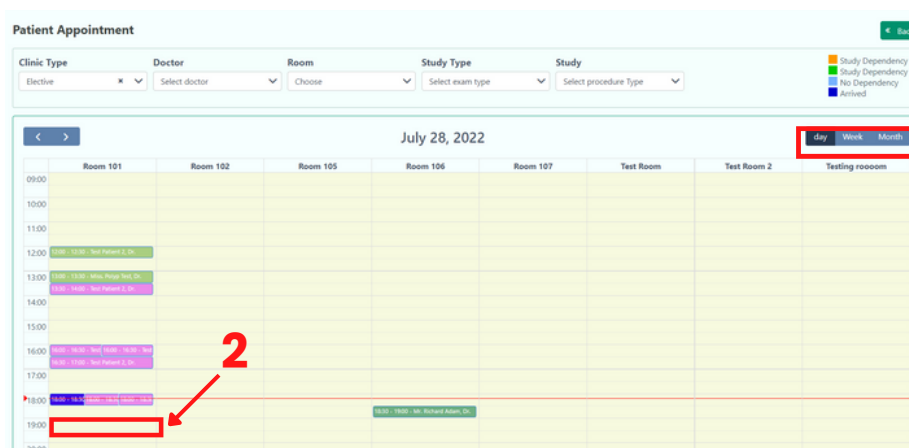

3.1.3 After select a slot, The user can **quickly register a patient**. If the patient is already registered, click Close to close the Quick Patient Registration screen and the user can also search patients in the search bar.

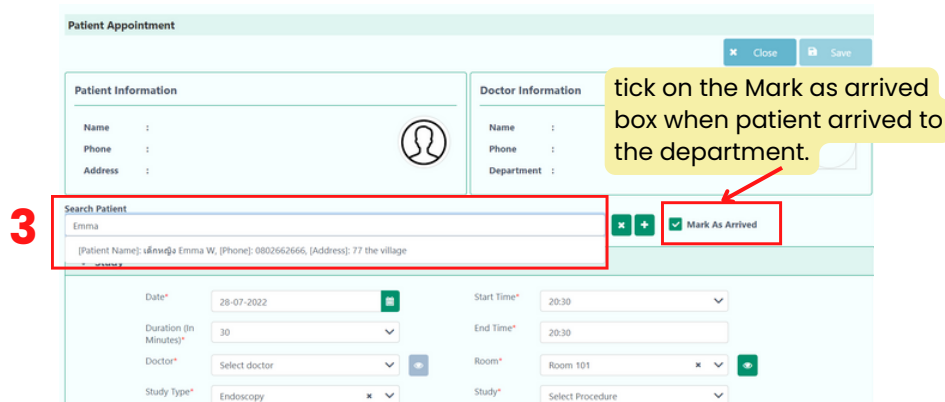

3. APPOINTMENT

3.1.4 Fill in the information in the two tabs.

Note: Please note that you will not be able to 'Save' the information until you have filled out all the required fields. The fields are highlighted in the screenshots below-

4

5

The screenshot shows the 'Patient Appointment' form with two tabs: 'Study' and 'Additional'. In the 'Study' tab, the 'Doctor' and 'Study' fields are highlighted with red boxes. In the 'Additional' tab, the 'Save' button is highlighted with a red box.


3.1.5 Once you have saved the information, you will be taken back to the calendar to show the appointment that was being made.

The screenshot shows the 'Patient Appointment' calendar view for May 9, 2023. The calendar displays a grid of rooms (Room 101, Room 102, Room 105, Room 106, Room 107, Test Room, Test Room 2, R-A) and time slots (09:00, 10:00, 11:00, 12:00, 13:00, 14:00, 15:00, 16:00, 17:00). An appointment is scheduled for 16:00 in Room 101, labeled '16:00 - 17:00 - Miss. Emma'. A red box highlights this appointment.

4. IMAGING

4.1 WORKFLOW TO SKIP THE IMAGING SECTION

4.1.1 If the user is **using DICOM or FTP function**, they do not need to use the Imaging section. the user can skip the imaging section by setting the administration option as below:

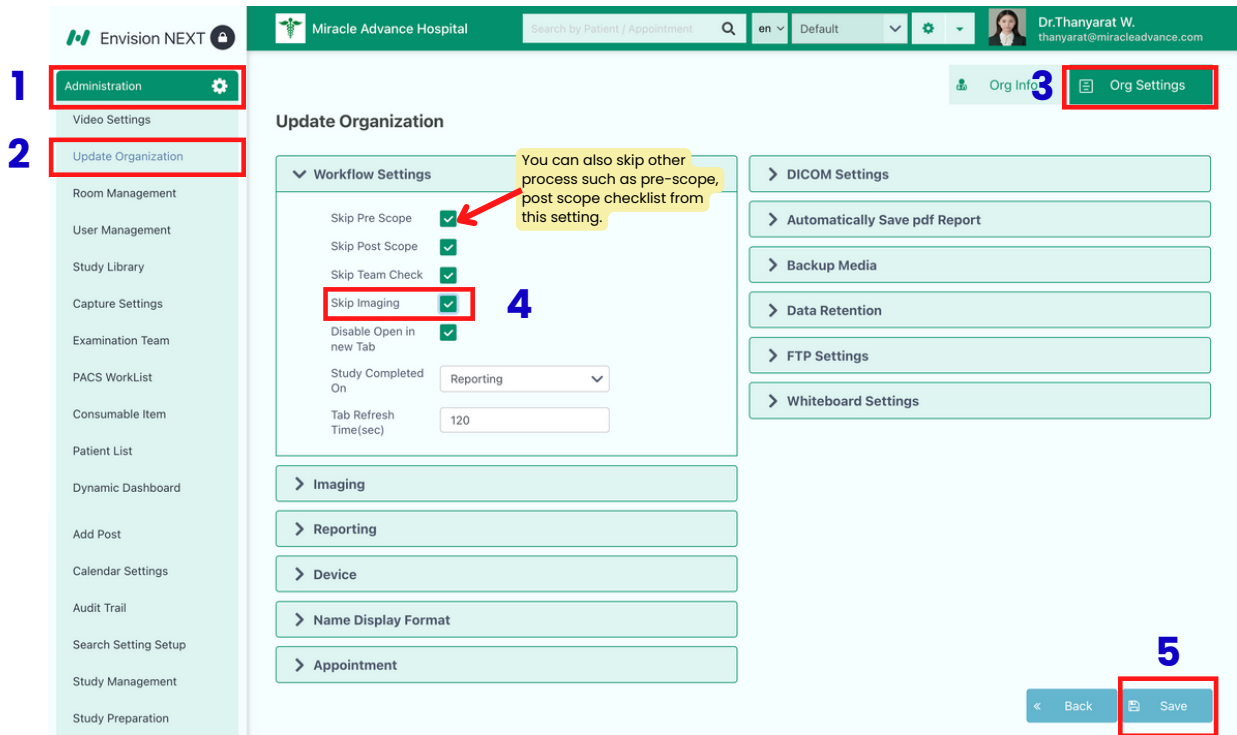

The screenshot shows the 'Update Organization' settings page in the Envision NEXT interface. The page is titled 'Update Organization' and features a sidebar on the left with various settings categories. The main content area is divided into 'Workflow Settings' and a list of other settings on the right.

- 1** Administration (gear icon)
- 2** Update Organization
- 3** Org Settings
- 4** Skip Imaging (checkbox checked)
- 5** Save

Additional details from the screenshot:


- The 'Workflow Settings' section includes checkboxes for 'Skip Pre Scope', 'Skip Post Scope', 'Skip Team Check', and 'Skip Imaging' (checked). A yellow callout box states: 'You can also skip other process such as pre-scope, post scope checklist from this setting.'
- The 'Study Completed On' dropdown is set to 'Reporting'.
- The 'Tab Refresh Time(sec)' is set to '120'.
- The right sidebar contains expandable sections for 'DICOM Settings', 'Automatically Save pdf Report', 'Backup Media', 'Data Retention', 'FTP Settings', and 'Whiteboard Settings'.
- The bottom right corner has 'Back' and 'Save' buttons, with 'Save' being the final step indicated by the number 5.

After successful the procedure in the Endoscopy room, the user can start the reporting section. (Please see page 9)


4. IMAGING

4.2 THE IMAGING SECTION

4.2.1 When the patient is ready, the user can start imaging by clicking **the start imaging button**.

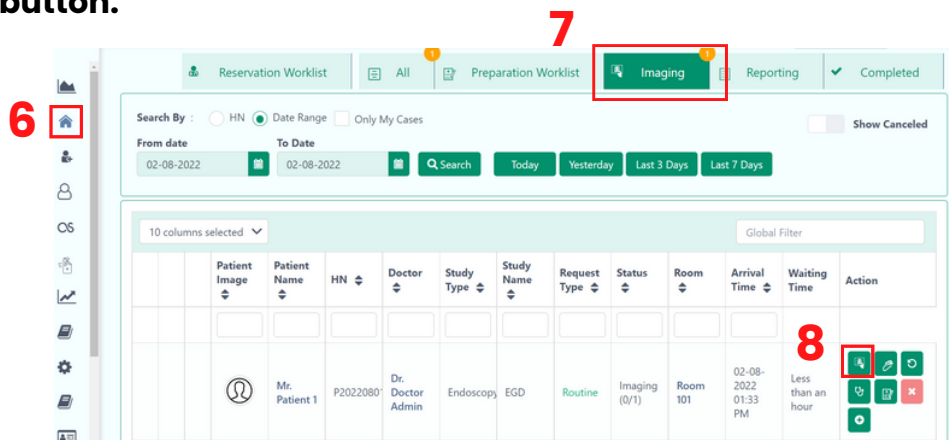


4.2.2 Then the software will open the imaging window. The user has to select the Resource Allocation for selecting the team and the scope by clicking the Resource allocation button.

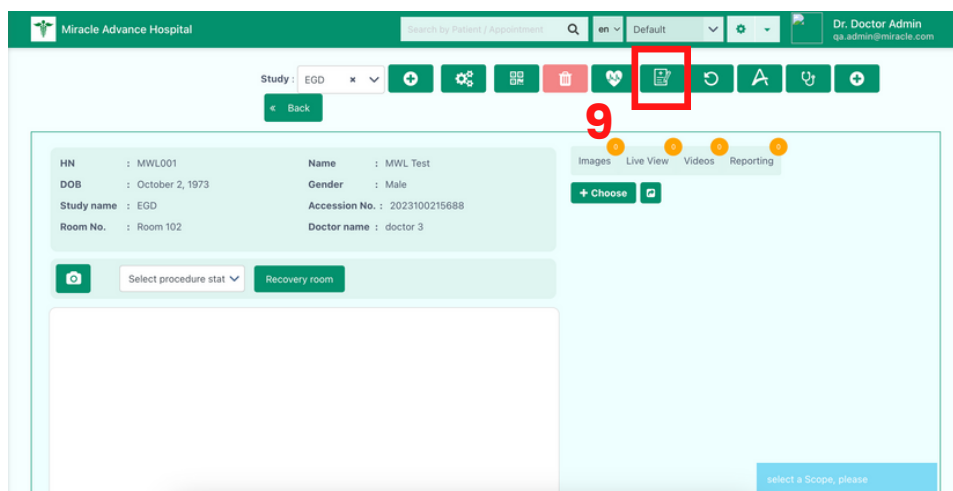


Note: If you do not select the scope in the previous step , it will continuously give notification for selecting the scope.

4.2.3 After selecting the Resource Allocation button, In the team settings, the user can set their desired team from the dropdown. They will get a confirmation notification when adding the team.

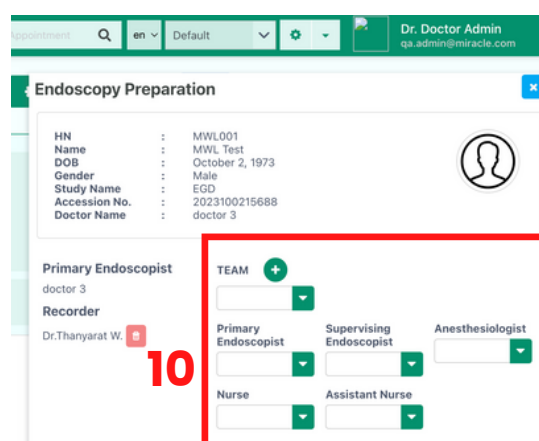


4.IMAGING


4.2.4 After that, the user has to select Include all Scopes. Then, select Scopes from the dropdown and pressing the select scope button.

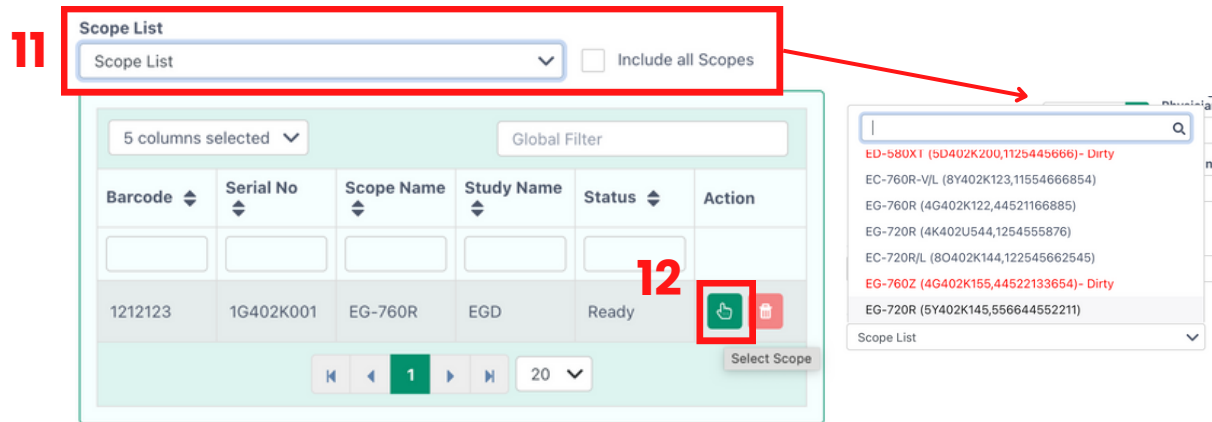

4.2.5 The user then start imaging process by trigger from Foot-switch and Endoscopy-switch. The user can notice the number of images and Live captures from the above.

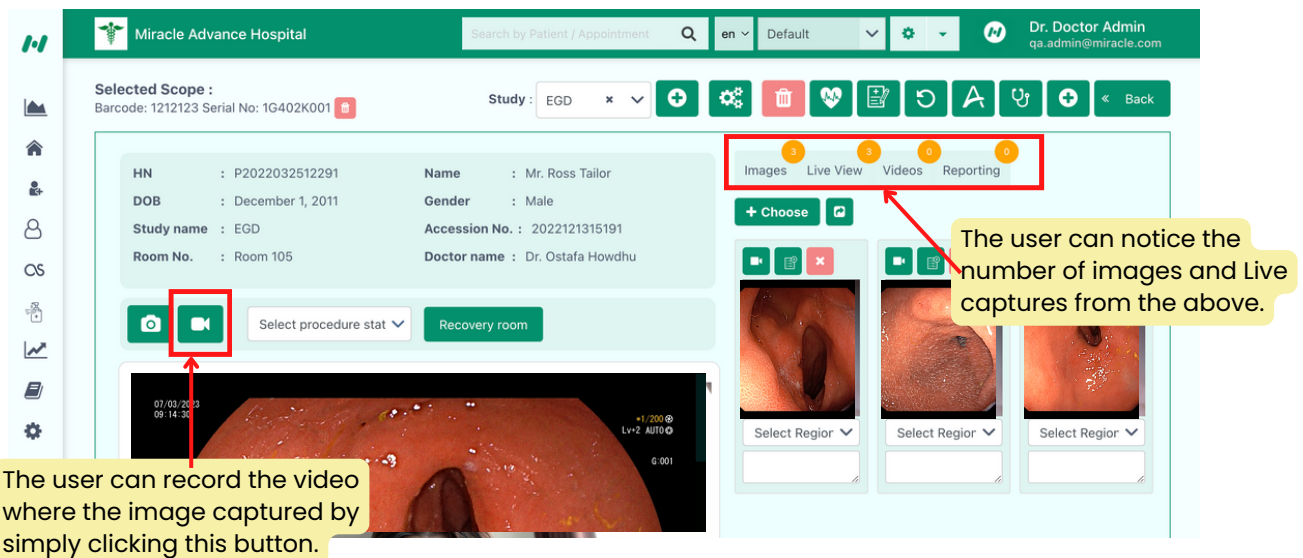

4.2.6 When the user finish the study successfully, The user can go to Reporting from Imaging by selecting the Procedure status to Success and click the Reporting button.

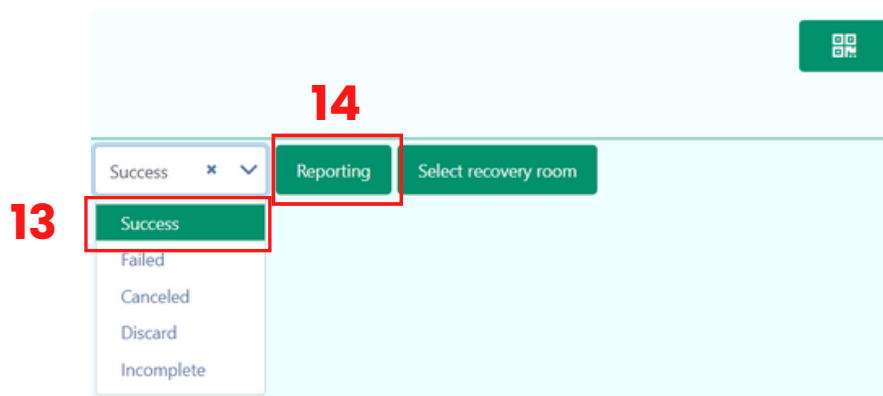


5. REPORT

5.1 START REPORTING

5.1.1 After selecting the Procedure status, The Reporting section will be automatically shown by the system.

5.1.2 Moreover the user can start writing the report by clicking the “Write report” button in the Reporting tab.

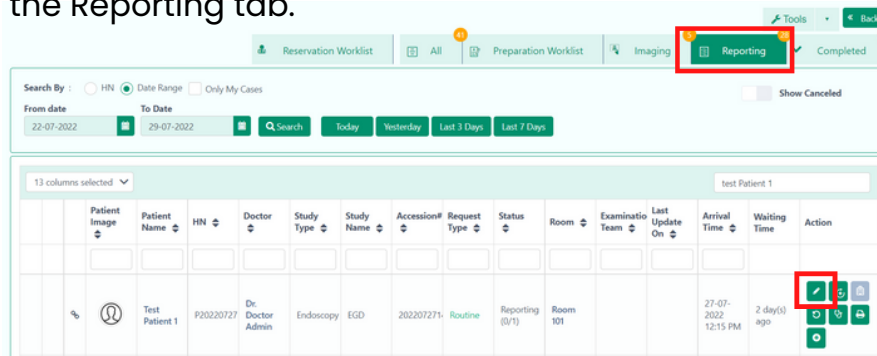


5.2 REPORTING SECTION

In this section, the user will find 3 sections.


5.2.1 The Information of the Patient.

5.2.2 Captured Files : The section contains the data from the Imaging. The user can see the images, Live captures, and Videos.

5.2.2.1 They can also choose images from their devices also. There are also some buttons where the user can control the images like- minimizing, sorting, maximizing, and selecting specifically selected images for adding to reporting.

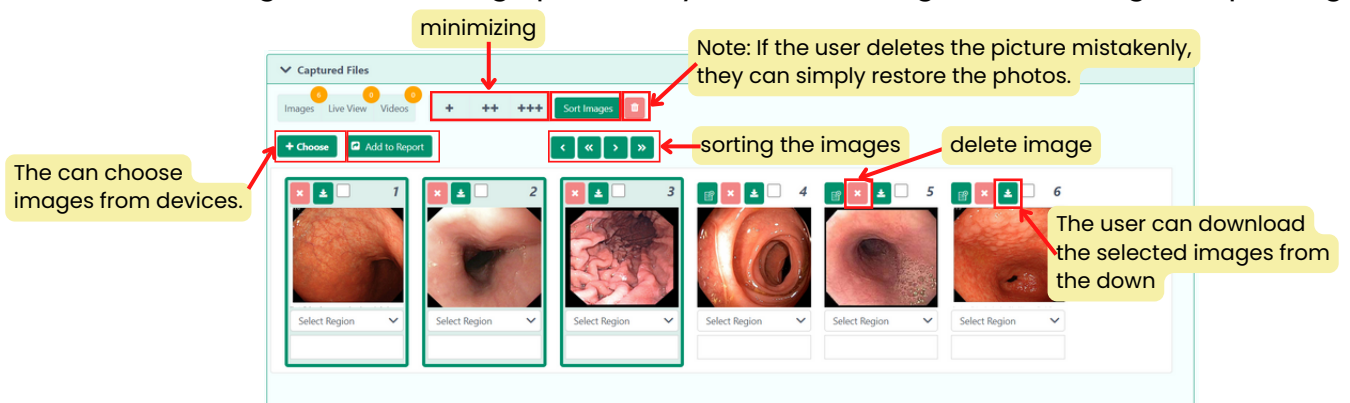


5.2.2.2 The user can edit and mark down an image by double clicking that specific image.

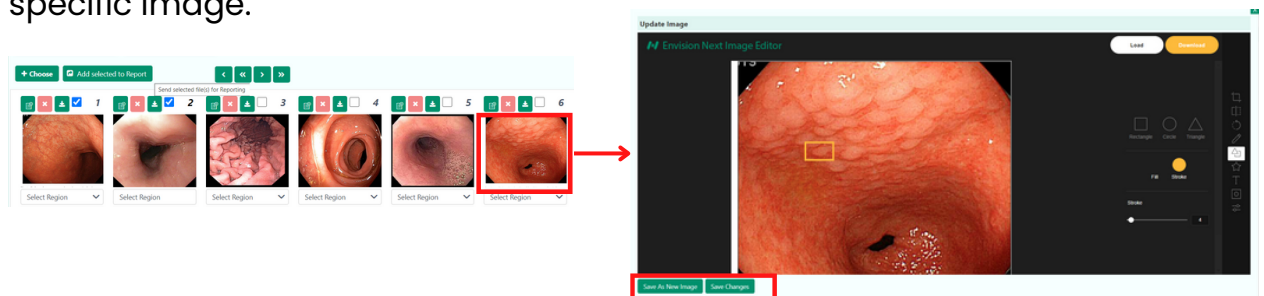

- if you want to save as new image click “save as a new image” button
 - if you want to save replace the original image click “save changes”
- The newly edited image will display on the screen.


5. REPORT

5.2.2.3 The user can inputting the name of the region in the text box under the images.

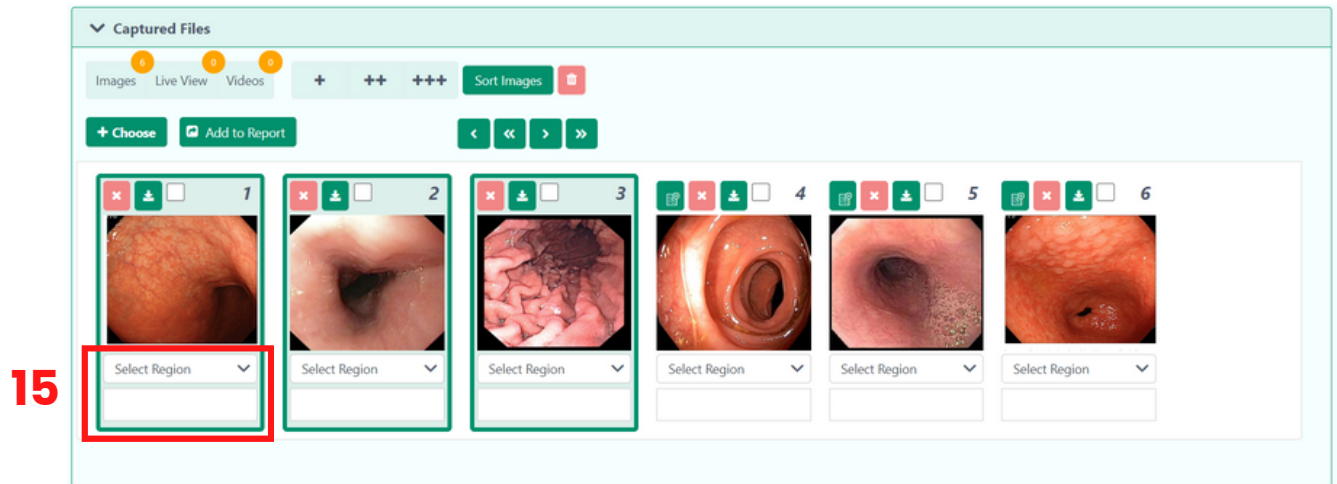


5.2.2.4 Sent the images to the reporting.

- 1) The user can **add all the images** to reporting by clicking “Add to report” button and you will get a confirmation notification.

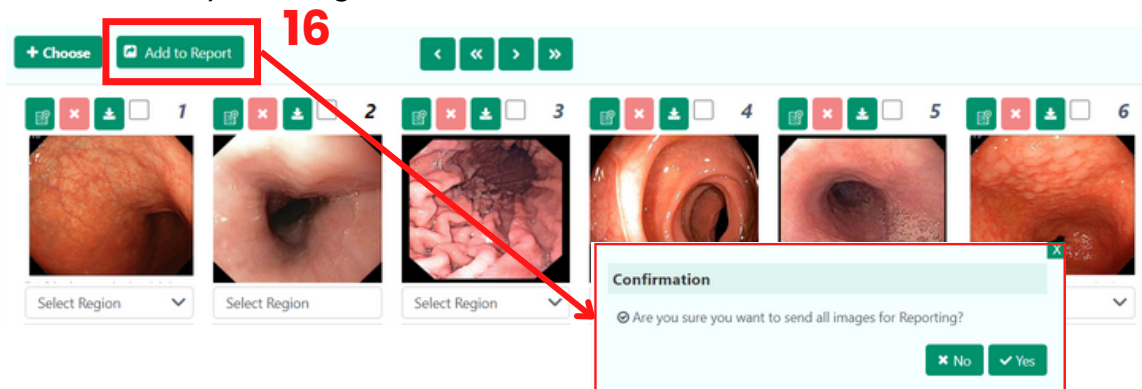

- 2) if you want to **add the selected images**, you choose by clicking on checkbox, then click “Add selected to Report”. then we will see the selected images in reporting below.

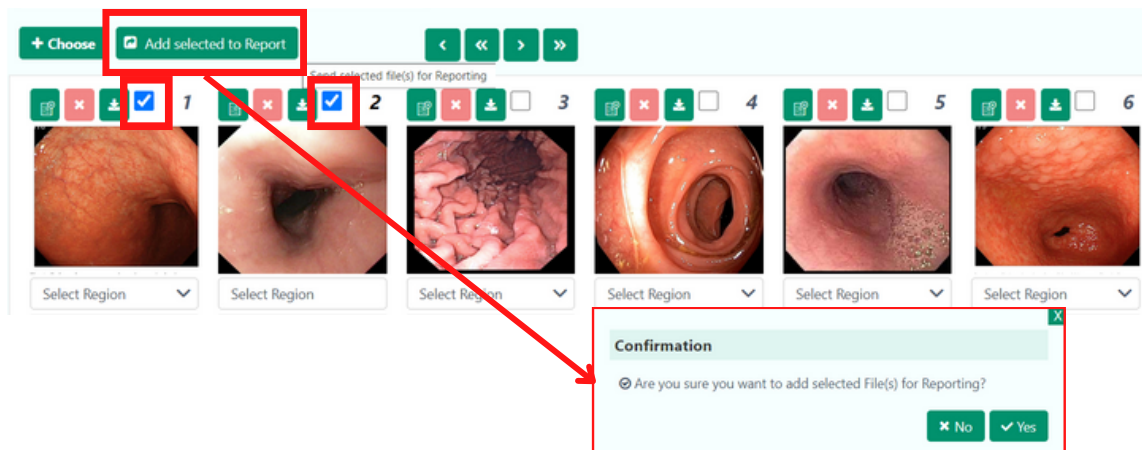


5. REPORT

5.2.3 Study metrics, There are 5 sections in study metrics :


- indication/ diagnosis
- sedation
- timing
- team
- bowel preparation section enables only use in the colonoscopy study in reporting.

5.2.3.1 indication/diagnosis section, The user will get 3 column sets. They are Chief complain, Diagnosis and Qualification.



- The user must select symptoms according to the patient's conditions and images. The user can search Chief complain, Diagnosis and Qualification from the search bar.
- After selecting from the column, then click the "Add" button to add the reporting and click the "Save" button.

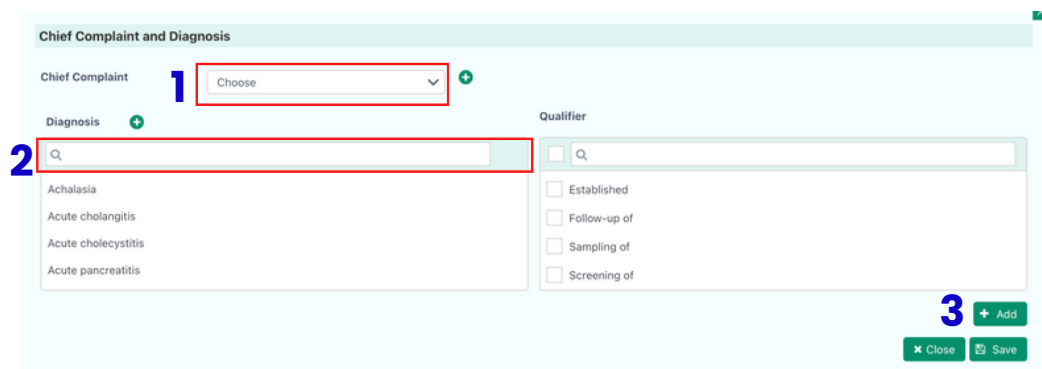

5.2.3.2 Sedation section, There is 2 types of sedation (Conscious sedation/ General sedation). The user can add sedative agent from dropdown menu and at the amount and click "Add". then click "save".

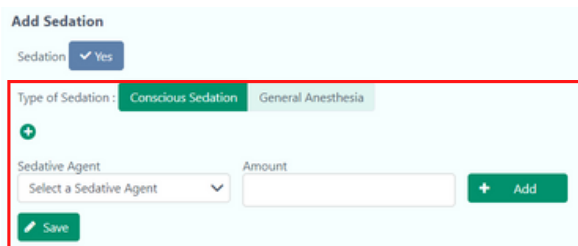


Note: "Local xylocaine throat spray" enables only use in the EGD study in reporting.

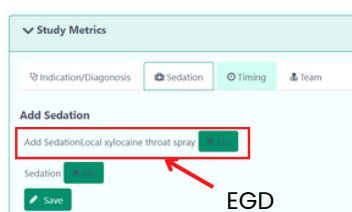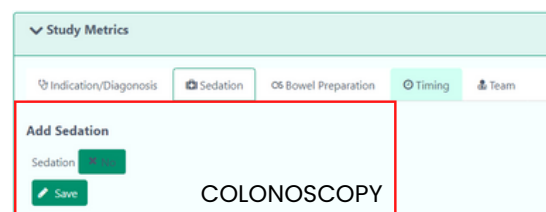
5. REPORT

5.2.3.3 Timing section, The user can view the start and end procedural time of the scope. It also provides the duration of the procedural time.

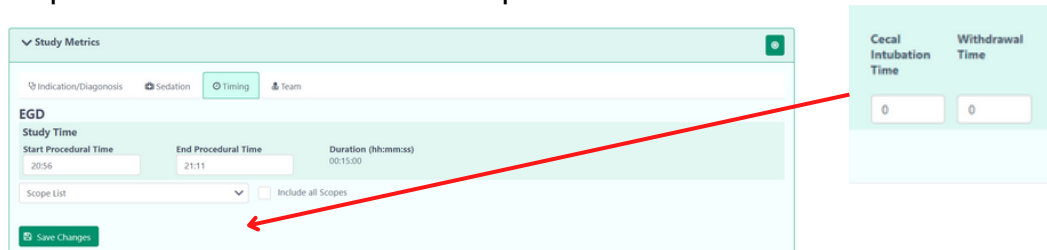


Note: the cecum intubation time and Withdrawal time enables only use in the colonoscopy study.

5.2.3.4 Team section, The user can add more teams from this section. They can select the team member from each of the dropdowns.

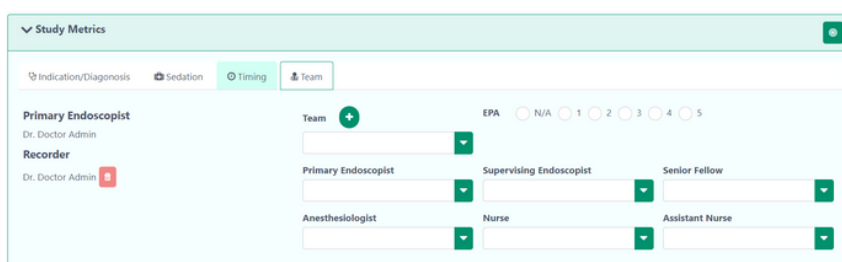

5.2.3.5 Bowel preparation section, The user will find bowel preparation for colonoscopy study.

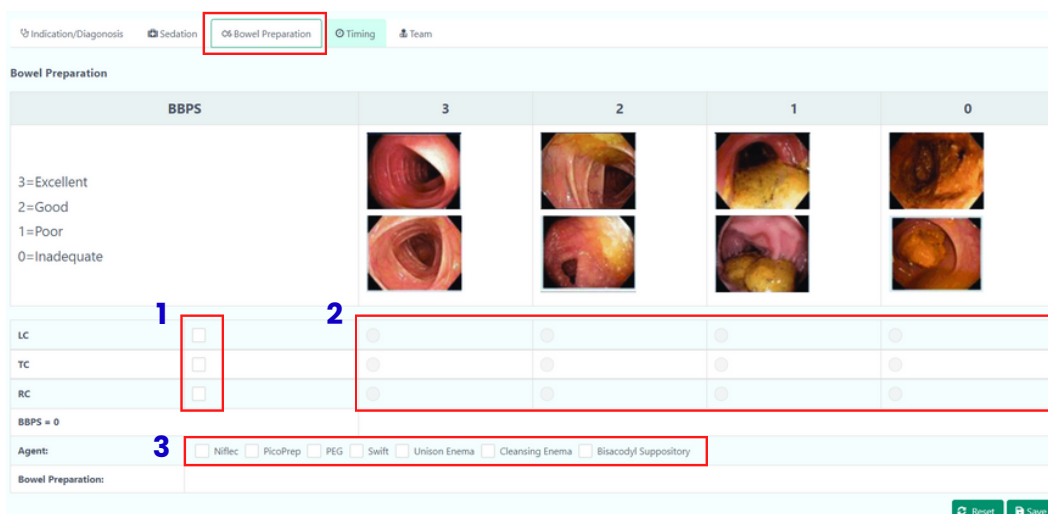

- 1) You can check the BBPS (Boston Bowel Preparation Scale) for assessment of the quality of bowel preparation each of regions of the colon.
 - LC: the right colon
 - TC: the transverse colon
 - LC: the left colon
- 2) The points are assigned as
 - 3 = Excellent
 - 2 = Good
 - 1 = Poor
 - 0 = Inadequate
- 3) Then select the laxative agent.
- 4) Finally, The user must click the "Save" button to save the data.

5. REPORT


5.3 REPORTING CONFIGURATION

We have 3 types of reporting, there are also different findings report formats.

- Simple report
- Semi-structure report
- Structure report

In this Short manual, We will show the basic report including: the Simple report and the semi-structure reporting.

5.3.1 the user can customize report settings by going to the Main Menu on the left, selecting Administration and selecting Update Organization, Next selecting the Org Setting tab, Then selecting the Reporting tab and Select the report type as below.

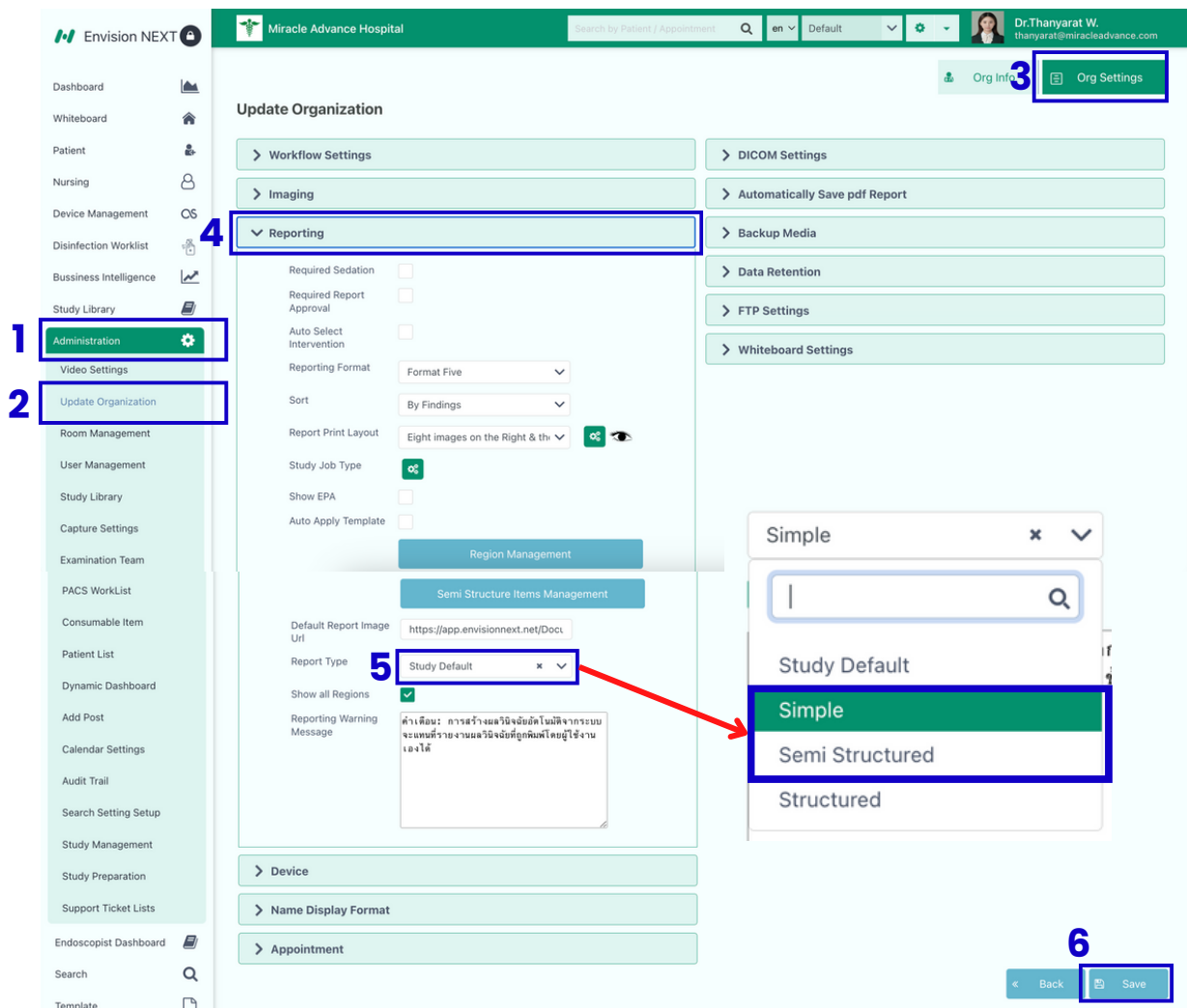

The screenshot shows the 'Update Organization' page in the Envision NEXT system. The interface includes a left sidebar with a menu, a top header with user information, and a main content area with various settings tabs. The steps are numbered 1 through 6:

1. Click on 'Administration' in the left sidebar.
2. Click on 'Update Organization' in the left sidebar.
3. Click on 'Org Settings' in the top right corner.
4. Click on 'Reporting' in the 'Update Organization' section.
5. Click on 'Study Default' in the 'Report Type' dropdown menu.
6. Click on 'Save' in the bottom right corner.

The 'Reporting' section shows various settings including 'Required Sedation', 'Required Report Approval', 'Auto Select Intervention', 'Reporting Format' (Format Five), 'Sort' (By Findings), 'Report Print Layout' (Eight images on the Right & th), 'Study Job Type', 'Show EPA', and 'Auto Apply Template'. A dropdown menu for 'Report Type' is open, showing options: 'Simple', 'Semi Structured', and 'Structured'. The 'Simple' option is selected.


5.3.2 when you click the "Save" button , You will get a notification.

5. REPORT

5.4 TYPE OF REPORTING

5.4.1 Semi-structure reporting: Findings

5.4.1.1 First, The user can click the **Start** button to add all organs or select each organ.

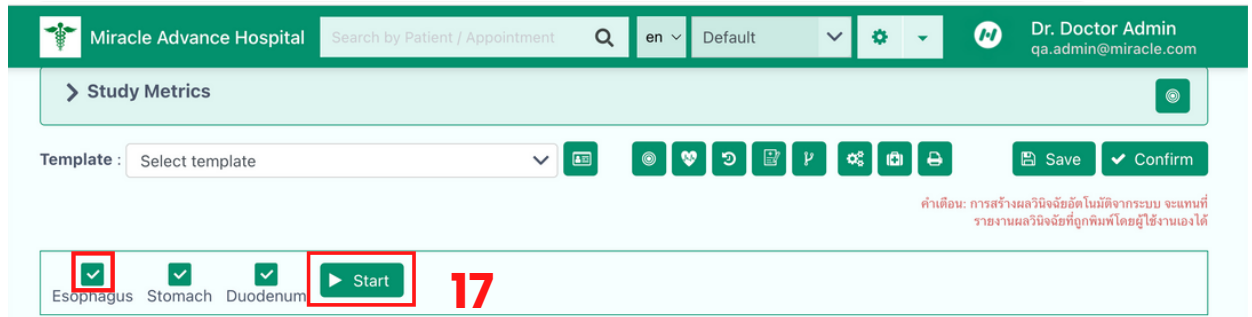

Miracle Advance Hospital Search by Patient / Appointment en Default

Dr. Doctor Admin qa.admin@miracle.com

> Study Metrics


Template : Select template

Save Confirm

Esophagus Stomach Duodenum Start

คำเตือน: การสร้างผลวินิจฉัยอัตโนมัติจากระบบ จะแทนที่ รายงานผลวินิจฉัยที่ถูกพิมพ์โดยผู้ใช้งานเองได้

5.4.1.2 After that, the user will find the Findings textbox of each organ where they have to input the finding details, intervention, impression and comment.

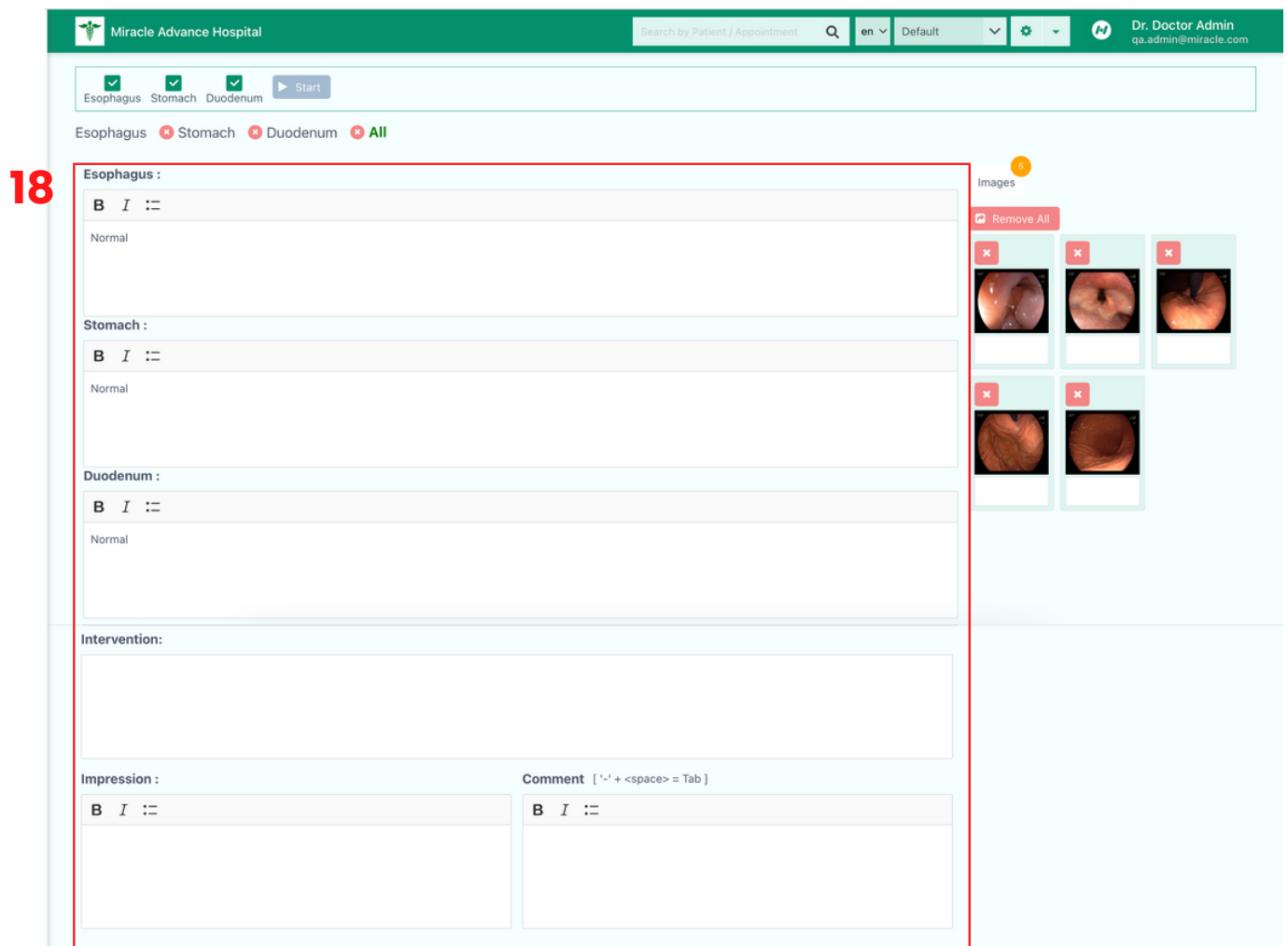

Miracle Advance Hospital Search by Patient / Appointment en Default

Dr. Doctor Admin qa.admin@miracle.com

Esophagus Stomach Duodenum Start

Esophagus Stomach Duodenum All

Esophagus :

B I ::

Normal

Stomach :

B I ::

Normal

Duodenum :

B I ::

Normal

Intervention:

Impression :

B I ::

Comment ['-' + <space> = Tab]

B I ::

Images

Remove All


18

14

5. REPORT

5.4.2 Simple reporting: Findings

5.4.2.1 the user will find the Findings textbox where they have to input the finding details, impression and comment.

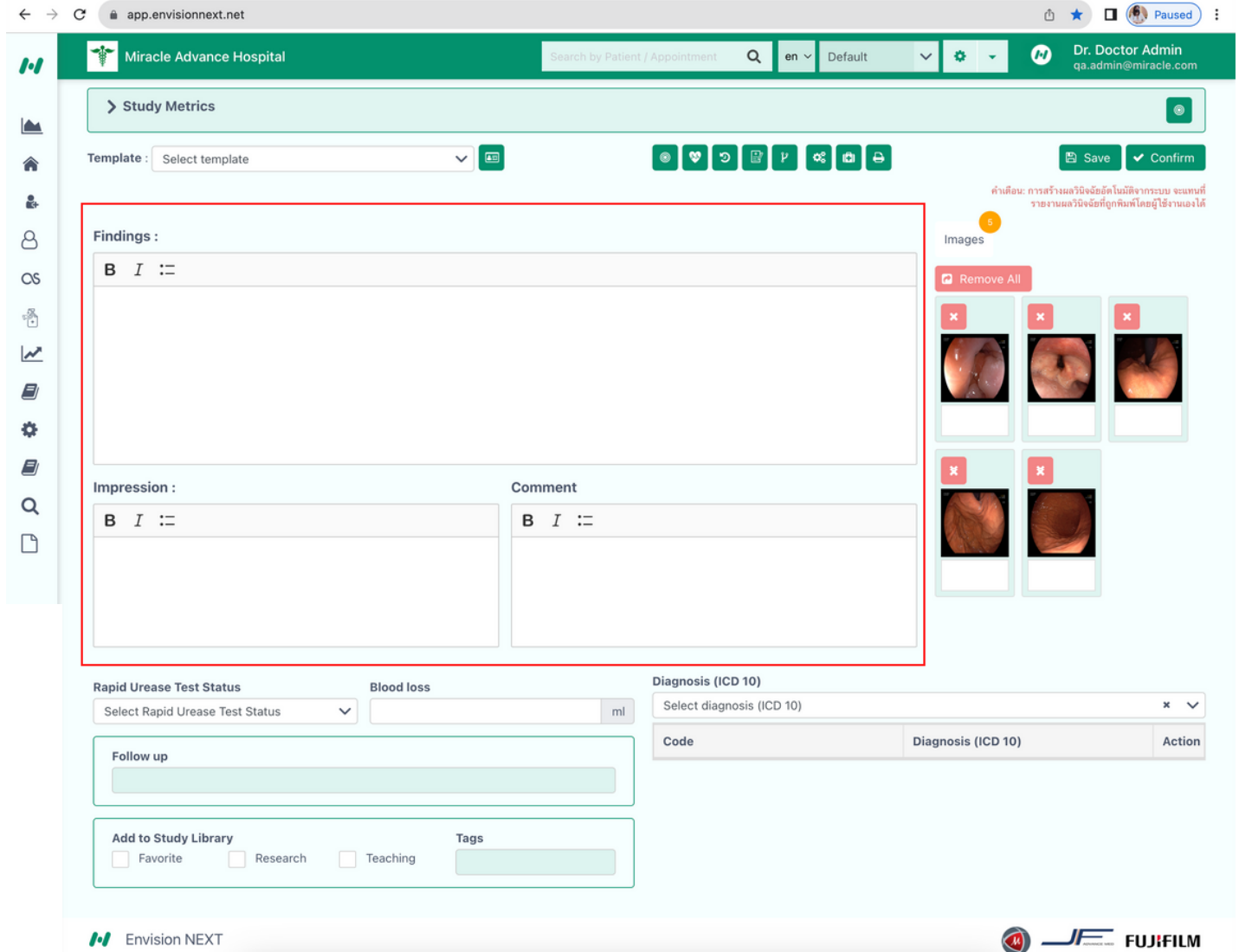

The screenshot shows the 'Miracle Advance Hospital' interface. The top navigation bar includes a search bar, language selector (en), and user profile (Dr. Doctor Admin). The main content area is titled 'Study Metrics' and contains a 'Template' dropdown, a 'Findings' section with a large text area, and 'Impression' and 'Comment' sections with smaller text areas. A red box highlights the 'Findings', 'Impression', and 'Comment' sections. To the right, there is an 'Images' section with a 'Remove All' button and five image thumbnails. Below the main sections, there are fields for 'Rapid Urease Test Status', 'Blood loss', and 'Diagnosis (ICD 10)'. At the bottom, there is a 'Follow up' section and an 'Add to Study Library' section with checkboxes for 'Favorite', 'Research', and 'Teaching', along with a 'Tags' input field.

Findings :

B I ::

Impression :

B I ::

Comment

B I ::

Rapid Urease Test Status

Select Rapid Urease Test Status

Blood loss

ml

Diagnosis (ICD 10)

Select diagnosis (ICD 10)

Code	Diagnosis (ICD 10)	Action

Follow up

Add to Study Library

☐ Favorite ☐ Research ☐ Teaching

Tags

5. REPORT

5.5 RAPID UREASE TEST/ BLOOD LOSS/ ICD.10/ FOLLOW UP PLAN/ STUDY LIBRARY

5.5.1 The user can also select the Rapid Urease Test Status, specify blood loss, select diagnosis(ICD.10)and select the follow up plan.

1 "Add to study library" option. There are 3 check boxes. they are favorite, research and teaching. the user can select the checkbox according to their importance.

5.5.2 The user will find different types of plans after clicking the comment plan. The user has to select the plan according to the patient situation and can add text.

5.6 Save, Preview, Print and Download the Report

5.6.1 Finally, The user has to click the **save** button for saving the report and click the **Confirm** button for completing the report.

5.6.2 To see the report, the user has to click the Preview Report button.

To preview the report.

5.5.3.3 You can also print and download the report by clicking the button as below:

Endoscopy Report

Gastroscopy

Procedure Date: 28-Nov-2023 , 19:30 to 19:45 (00:15)
 HN#: P2022032512290 Gender: Male
 DOB: 01 March, 2022 Age: 1 year-old
 Patient Name: Dok Phi

Accession# : 2023112815756 Room : Room 101
 Endoscopist : Miss. Dr.Thanyarat W. Scope :

