

Technical Manual

RadPortal

Developed By -

Miracle Advance Technology

Last Update Date -

12/22/2022

This Document Is for the RadPortal Application

Version – 1.1.1

Written By -

S. M. Al Shahaneaz Shaikat

Revision Sheet

Version	Date	Revision Description	Author Name
V 1.1.1	21/12/22	Draft status, Version Controlling, Study Status	S. M. Al Shahaneaz Shaikat
Rev. 2	DD/MM/YY	Short description of Major changes	Name
Rev. 3	DD/MM/YY	Short description of Major changes	Name
Rev. 4	DD/MM/YY	Short description of Major changes	Name

USER'S MANUAL OF RADPORTAL

Table of Contents

1. General Information about RadPortal.....	5
1.1 System Overview	5
1.2 Project References.....	6
1.3 Organization of the Manual	6
1.4 Acronyms and Abbreviations	6
1.5 Additional Help.....	7
2. System Summary	7
2.1 Minimum Hardware Requirements, Supported Browsers List, And Operating System List for RadPortal	7
2.1.1 Hardware Requirements Personal Computers and Laptops	7
2.1.2 Browsers.....	8
2.1.3 Operating System	8
2.2 User Roles and Permission Levels.....	8
2.2.1 Service Provider(Admin)	8
2.2.2 Customer	9
2.2.3 Hospitals(Healthcare Service Provider)	9
2.2.4 Radiologist(Subject Matter Experts)	9
2.2.5 Viewer.....	9
3. Getting Started.....	10
3.1 Logging In	11
3.2 System Menus.....	15
3.3 How to Change Password, Change Profile Picture, and Add Your Digital Signature.....	16
3.3.1 Changing Password	16
3.3.2 Changing Profile Photo	18
3.3.3 Adding Digital Signature	19
3.4 Software Version and Updates	20
3.5 Sign Out from The System.....	21
4. Exploring Worklist Menu	22
4.1 How to Search Studies.....	23
4.2 How to Sort Studies	24
4.3 How Study information is Organized In Study worklist.....	25

4.4 How to Review Images, Resend To PACS, Cancel, Edit, View Report and Attach Additional Documents to Studies	27
4.4.1 Review Images of a Specific Study.....	28
4.4.2 Resend to PACS (Picture archiving and communication system).....	30
4.4.3 Cancel (Remove) Studies	30
4.4.4 Reconcile Studies	31
4.4.5 How to Submit Study/Case	33
4.4.6 View Study Reports	34
4.4.7 Attach Additional Files To A Study.....	35
4.6 How to Submit Studies	37
4.7 How to Assign a Radiologist to a Case/Study.....	41
5. Registration Menu.....	43
6. User list.....	45
6.1 How to Add New Users	46
6.2 How to Approve and Update User Information	48
6.3 How to Organize Information of Users	49
7. Organization Setup	52
7.1 How to Add New Organization	52
7.2 How to Update Organization Information	55

1. General Information about RadPortal

Radportal is a Secured cloud-based radiology platform. The main problem it solves is that there are many hospitals in remote places worldwide with all the modern medical equipment for medical tests such as (MRI, CT scan, x-rays, nuclear medicine, ultrasound, and more). Still, hospitals may not have expert radiologists/doctors readily available to write the reports for these studies when patients come to these remote hospitals and complete their tests.

Sometimes the Patient has to physically go to hospitals/Other more developed cities where radiologists/Doctors are available, but this process is difficult and time-consuming.

Using the RadPortal System, the hospital can easily save patient test results. Then the hospitals can forward the test results to an expert doctor/radiologist who may not live nearby.

1.1 System Overview

RadPortal is a secured cloud-based radiology platform. It helps hospitals to send their Patient's radiology test images to doctors/radiologists worldwide and receive written reports from remote doctors/radiologists.

RadPortals Core functionalities are

- All hospitals/sites that have already configured modality devices can automatically send DICOM files to their routers, which are then forwarded to the RadPortal cloud database.
- RadPortal system reads the DICOM headers and collects information from the header, and saves it as an entity called "study".
- Hospitals also can upload single or multiple DICOM files manually.
- Hospitals can then access DICOMs from their RadPortal profiles and assign specific radiologists to write reports on that particular study or keep the study open for other radiologists.
- Suppose any open study is available on RadPortal. In that case, the independent radiologists can log in to their RadPortal profiles and independently select any open study. Then the individual Radiologist can write reports for that particular study.
- The Radiologist has a roaming profile, so they do not need specific devices to interact with the RadPortal system. The Radiologist can interact with the RadPortal system with laptops, desktops, or tabs and a reliable Internet connection.

- RadPortal also helps Government Institutes such as the National Health Security office to track the number of studies and reports the hospitals and the radiologists completed and disburse fees and payments accordingly.

1.2 Project References

References used in preparing this document are listed here. Project References may help you understand the complex topics and terms used in the RadPortal system.

<http://dicomlookup.com/modalities.asp>

1.3 Organization of the Manual

This document is patterned to familiarize you with the core functionalities of the RadPortal application. This document explains how the application is intended to use by properly explaining its functionalities and limitations.

You will also know how the RadPortal application saves information and how the information is shared across multiple user roles within the RadPortal system. We will discuss the most common elements so readers can later relate to the terms as they look deeper into the application's features and functionalities.

The functionalities will vary depending on user roles, but the role **Admin** will have access to all the functionalities of RadPortal. User roles are extensively explained in section [2.2, User Roles and Permission Levels](#).

1.4 Acronyms and Abbreviations

Common acronyms and abbreviations are introduced in this section so the reader can get comfortable with the acronyms and abbreviations used in this documentation and the RadPortal application.

DICOM Digital Imaging and Communications in Medicine.

PACS Picture archiving and communication system.

1.5 Additional Help

Common terms are introduced throughout the user manual to familiarize the reader with the common terms used in the application. Sections are indexed (numbered, bullet points, and roman numbers). Sections are also linked to common keywords so the reader can go directly to any relevant section to find more information. For easier navigation, keywords are marked in [blue](#).

Here are the styles used for lists

Bulleted lists:

- **Text 1**
- **Text 2**
- **Text 3**

Numerical numbered lists:

1. **Text 1**
2. **Text 2**
3. **Text 3**

Roman numbered lists:

- I. **Text 1**
- II. **Text 2**
- III. **Text 3**

Here are the styles for Notes and warnings.

| **Note:** Notes are marked in the color blue.

| **Warning:** Warning messages are marked in the color red.

2. System Summary

2.1 Minimum Hardware Requirements, Supported Browsers List, And Operating System List for RadPortal

2.1.1 Hardware Requirements Personal Computers and Laptops
Processor (CPU) with 2 gigahertz (GHz) frequency or above

- A minimum of 2 GB of RAM
- Monitor Resolution 1024 X 768 or higher
- A minimum of 20 GB of available space on the hard disk
- Internet Connection Broadband (high-speed) Internet connection with a speed of 4 Mbps or higher
- Keyboard and mouse or compatible pointing device

2.1.2 Browsers

- Chrome 36+
- Edge 20+
- Mozilla Firefox 31+
- Internet Explorer 11+ (Windows only)
- Safari 6+ (macOS only)

2.1.3 Operating System

- Windows 7, Windows 8, Windows 10, or Windows 11
- Mac OSX 10.8, 10.9, 10.10, or 10.11

2.2 User Roles and Permission Levels

Before we look deeper into the RadPortal system, you must understand the user roles properly. The user roles are listed below.

- I. Service Provider(Admin)
- II. Customer
- III. Hospitals(Healthcare Service Provider)
- IV. Radiologist(Subject Matter Experts)
- V. Viewer

2.2.1 Service Provider(Admin)

Miracle Advance Technology is the service provider, and Miracle Advance Technology plays the Administration role inside the RadPortal System. The main functionality of an administrator role is to

manage the system so that the system runs smoothly and without errors. For the RadPortal system, the Service Provider will also install, configure and manage DICOM routers for the healthcare service provider.

2.2.2 Customer

This role is very important in the RadPortal system. The only customer role can Configure Turnaround Time, Assign radiologists to a study and manage and verify Radiologists. They also have full control over the RadPortal system. The Customer is also in charge of disbursing Fees and payments to the Hospitals and Radiologists.

2.2.3 Hospitals(Healthcare Service Provider)

Hospitals are also called "**sites**" in the RadPortal system. Hospitals play an significant role in the RadPortal System. They have to perform medical tests such as (MRI, CT scan, x-rays, nuclear medicine, ultrasound, and more) and manually upload files and images from browsers or the modality devices can automatically send **DICOM** files to RadPortal servers to create **cases** which are also called **studies**. The hospital can access the files, images, and DICOMs uploaded to the RadPortal server.

2.2.4 Radiologist(Subject Matter Experts)

The Radiologist also plays a significant role in the RadPortal System. There are two types of radiologist roles in the RadPortal system.

- I. Normal Radiologists
- II. Supervisor Radiologist

The main difference between the Supervisor and a normal Radiologist is that a Supervisor Radiologist can view all cases or studies. A normal radiologist can only see the studies/cases they are assigned to or any studies that are kept open for radiologists. Only the verified normal Radiologist can write reports of submitted studies/cases.

2.2.5 Viewer


The Viewer role is a bit unique. The viewer does not interact with any aspects of the RadPortal system. The viewers only view how many cases/studies the hospital/sites completed and how many reports the radiologists submitted and finalized in the RadPortal system.

This role is given to people with authority to monitor and audit the RadPortal system. This role is also given to the government institutes that might want to monitor the whole RadPortal system.

3. Getting Started


Before using the RadPortal system, you have to log in to the RadPortal system.

To log in, you have to go to this URL – <https://app.radportal.net/#/login>. After the login page properly loads, you should see a page like this

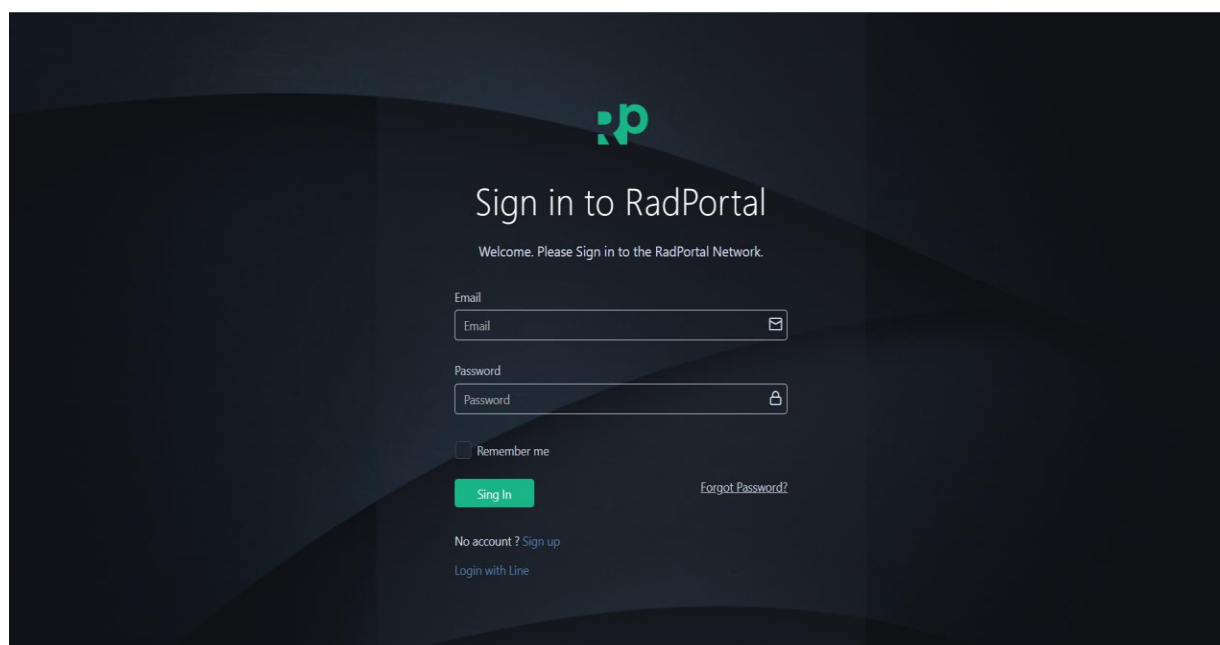


3.1 Logging In

This workflow diagram shows you the login process of the Service Provider (Admin)

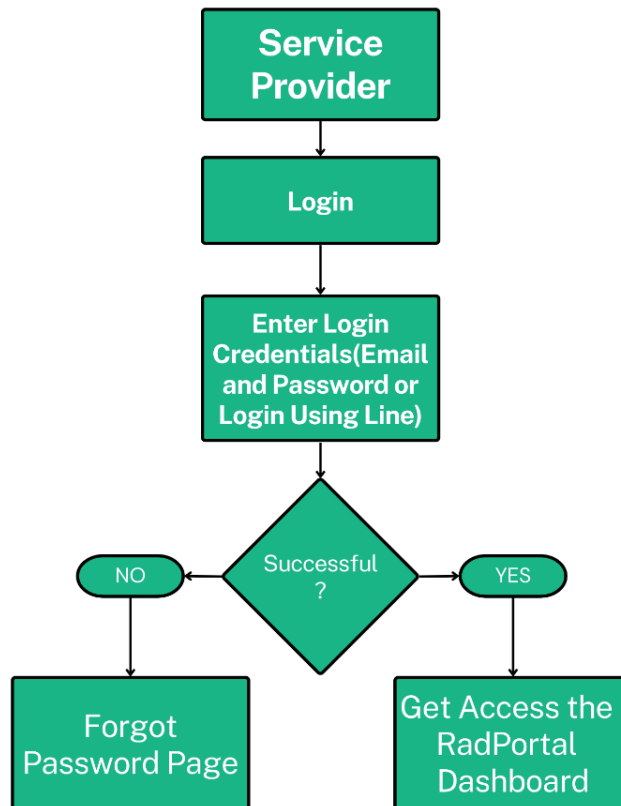

After you access the login page, you can use your **Email [1]** and **Password [2]** to log in to the RadPortal. You have to use the login credentials you used to create your account. Finally, click on the **Sign In [3]** button.

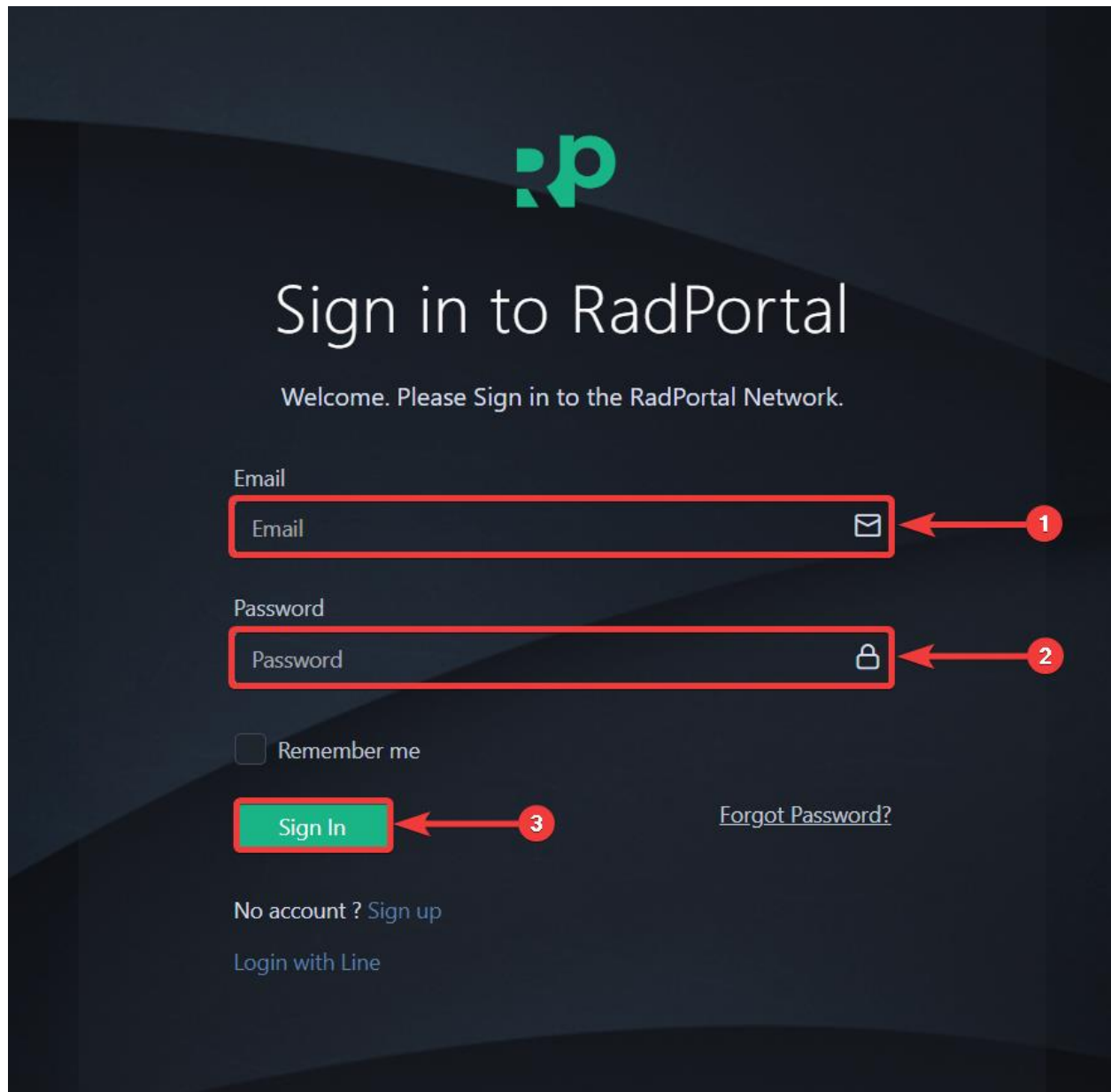

The image shows the RadPortal sign-in interface. At the top center is the RadPortal logo, a stylized 'rp' in green. Below the logo is the heading 'Sign in to RadPortal' in white, followed by the subtitle 'Welcome. Please Sign in to the RadPortal Network.' in a smaller white font. The main form area contains three input fields: 'Email' with a red box around it and a red arrow pointing to it labeled '1'; 'Password' with a red box around it and a red arrow pointing to it labeled '2'; and a 'Remember me' checkbox. Below the password field is a green 'Sign In' button with a red box around it and a red arrow pointing to it labeled '3'. To the right of the 'Sign In' button is a link that says 'Forgot Password?'. At the bottom left, there are two links: 'No account ? Sign up' and 'Login with Line'.

rp

Sign in to RadPortal

Welcome. Please Sign in to the RadPortal Network.

Email

Email

Password

Password

☐ Remember me

Sign In


[Forgot Password?](#)

No account ? [Sign up](#)

[Login with Line](#)

If you can not remember your Password, you can reset your Password.

To reset your Password, you have to click on the **Forgot Password?** Text.

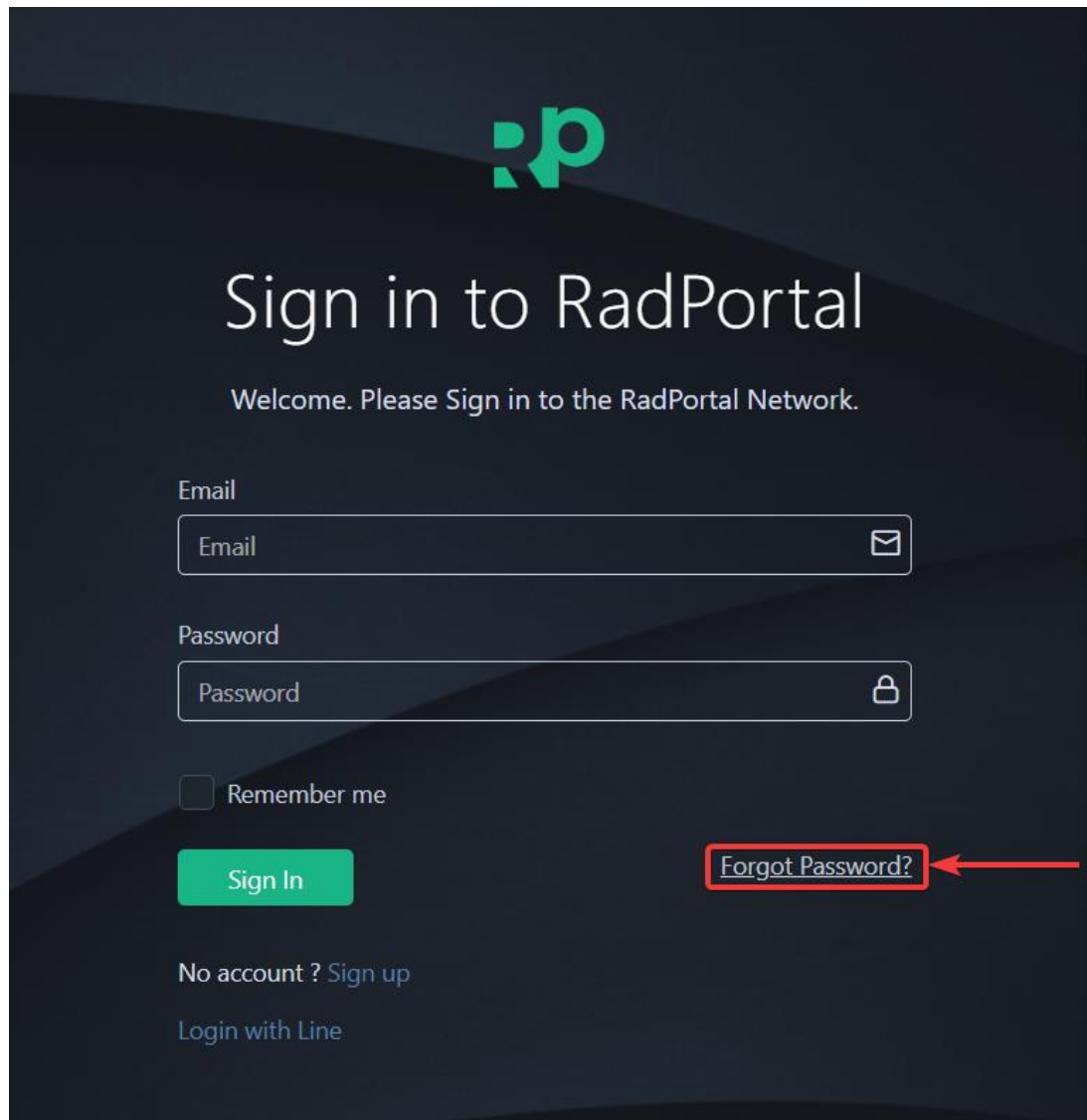

The image shows the RadPortal sign-in page. At the top center is the RadPortal logo, a stylized 'rp' in green. Below the logo is the heading 'Sign in to RadPortal' in large white text, followed by the subtitle 'Welcome. Please Sign in to the RadPortal Network.' in smaller white text. The form contains two input fields: 'Email' and 'Password', both with placeholder text and icons (an envelope for email and a lock for password). Below these fields is a checkbox labeled 'Remember me'. A green 'Sign In' button is positioned to the left of a red-bordered link that says 'Forgot Password?'. A red arrow points to this link. At the bottom left, there are links for 'No account ? Sign up' and 'Login with Line'.

rp

Sign in to RadPortal

Welcome. Please Sign in to the RadPortal Network.

Email

Email

Password

Password

☐ Remember me


Sign In

[Forgot Password?](#)

No account ? [Sign up](#)

[Login with Line](#)

After clicking, the system will take you to the forgot password page. The forgot password Page should look like this

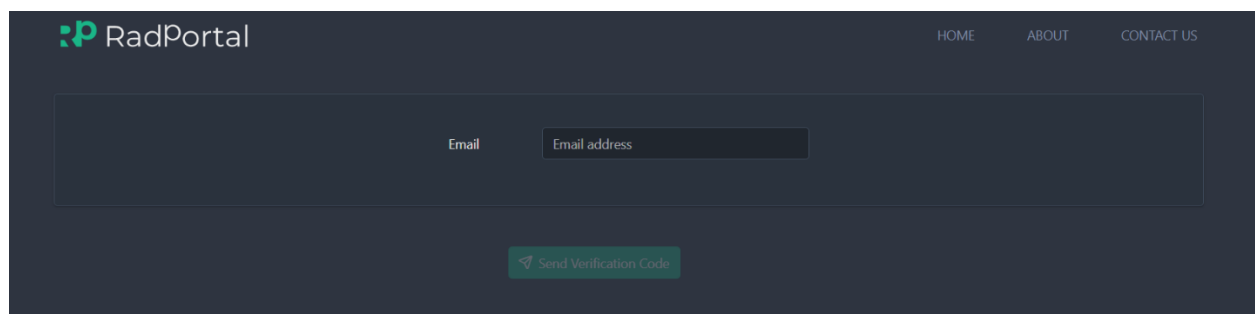

The image shows the RadPortal forgot password page. At the top left is the RadPortal logo. At the top right are links for 'HOME', 'ABOUT', and 'CONTACT US'. The main content area has a large input field with a placeholder 'Email' and a sub-label 'Email address'. Below this field is a green button with a checkmark icon and the text 'Send Verification Code'.

rp RadPortal


HOME ABOUT CONTACT US

Email

Email address


Send Verification Code

Now, insert your Email and click on **Send verification Code**.

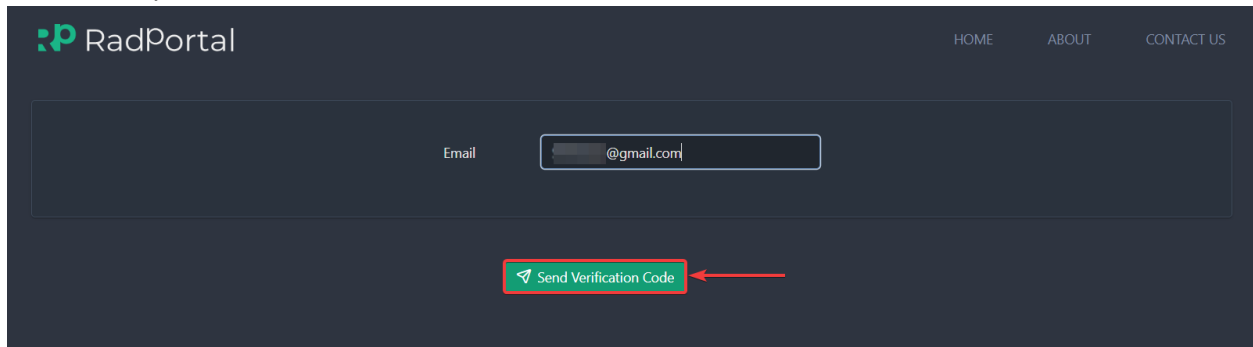

The image shows the RadPortal login interface. At the top left is the RadPortal logo. At the top right are links for HOME, ABOUT, and CONTACT US. In the center, there is an email input field with the placeholder text '@gmail.com'. Below the email field is a green button with a checkmark icon and the text 'Send Verification Code'. A red arrow points to this button.


After clicking, you should see a **Type code field [1]** popup under the email field where you have to insert the OTP (One Time Password). You should receive the OTP code in your mail. After inserting the OTP code, click on the **Confirm [2]** button.

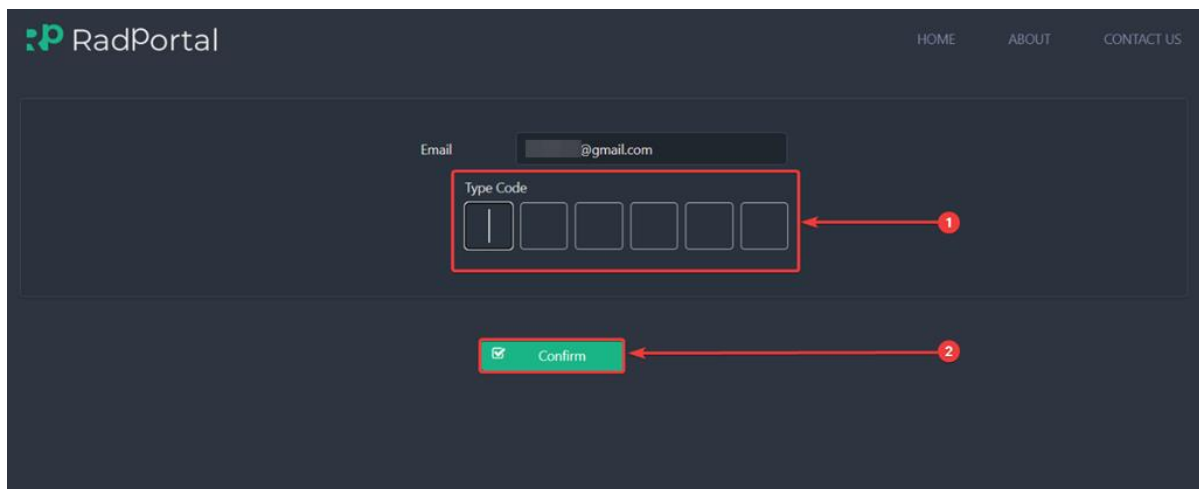


The image shows the RadPortal login interface after clicking the 'Send Verification Code' button. The email field remains at the top. Below it, a 'Type Code' field is highlighted with a red box and labeled with a red '1'. This field consists of six input boxes, with the first one containing a vertical line. Below the 'Type Code' field is a green button with a checkmark icon and the text 'Confirm', which is also highlighted with a red box and labeled with a red '2'. A red arrow points to this button.

After clicking, you will be presented with two new fields where you have to insert your new Password and click on the Conform button. Now you can log in to the RadPortal system with your new Password.

After successfully logging in, you should see the dashboard of the RadPortal system.

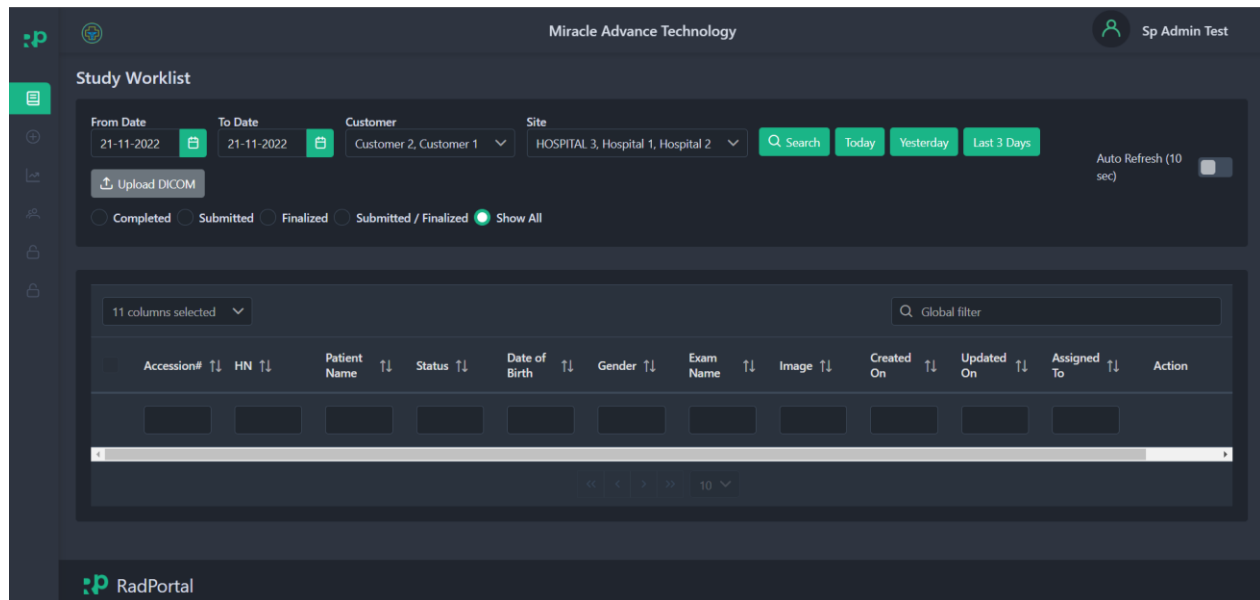


3.2 System Menus

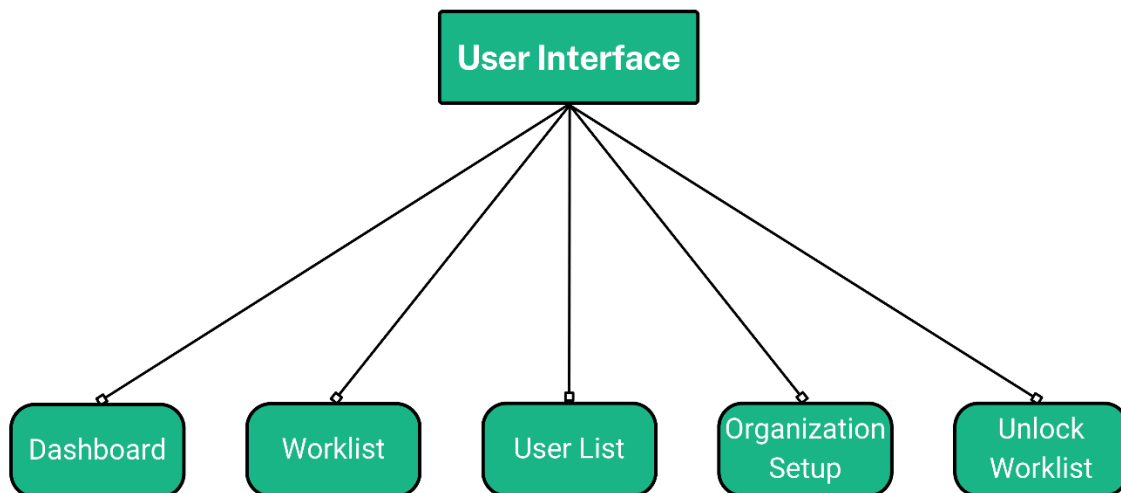

The RadPortal system has these menus currently available


- 1 **Dashboard**
- 2 **Worklist**
- 3 **User List**
- 4 **Organization Setup**
- 5 **Unlock Worklist**

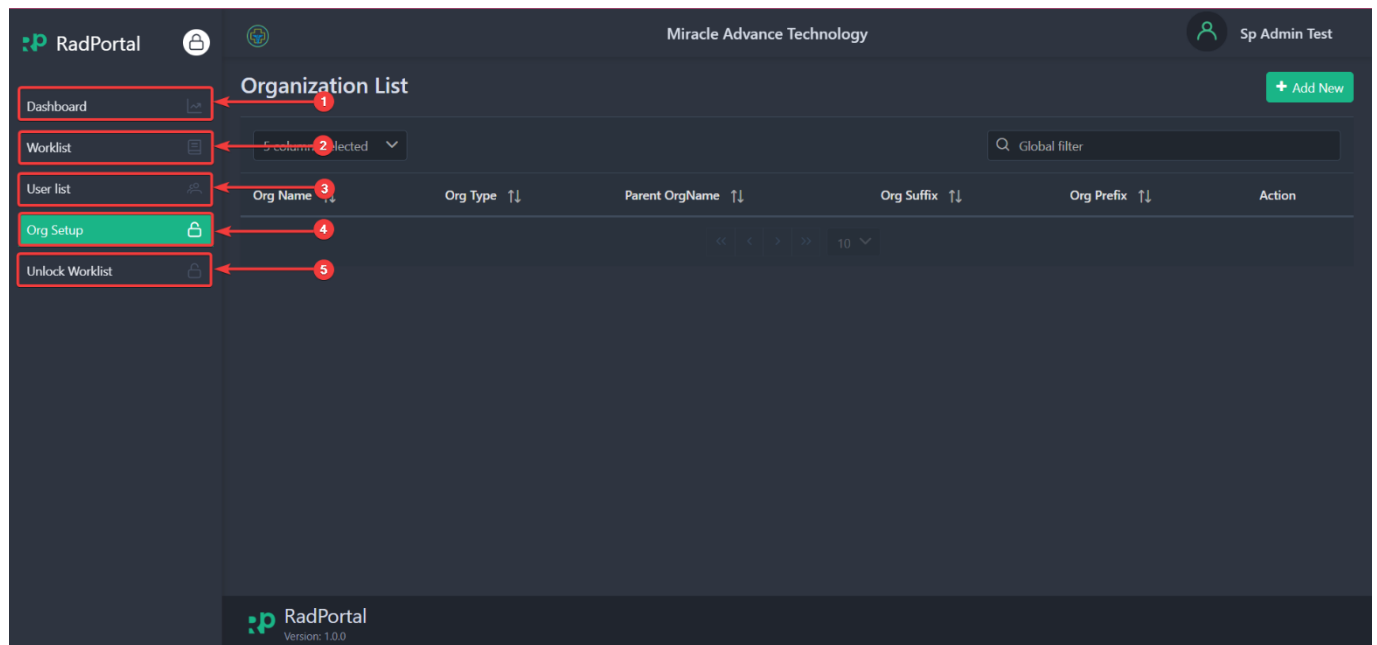


3.3 How to Change Password, Change Profile Picture, and Add Your Digital Signature

3.3.1 Changing Password

You can easily change your RadPortal login password from your profile. To access your profile, click on your **username [1]** in the top right corner and click on **Profile [2]**.

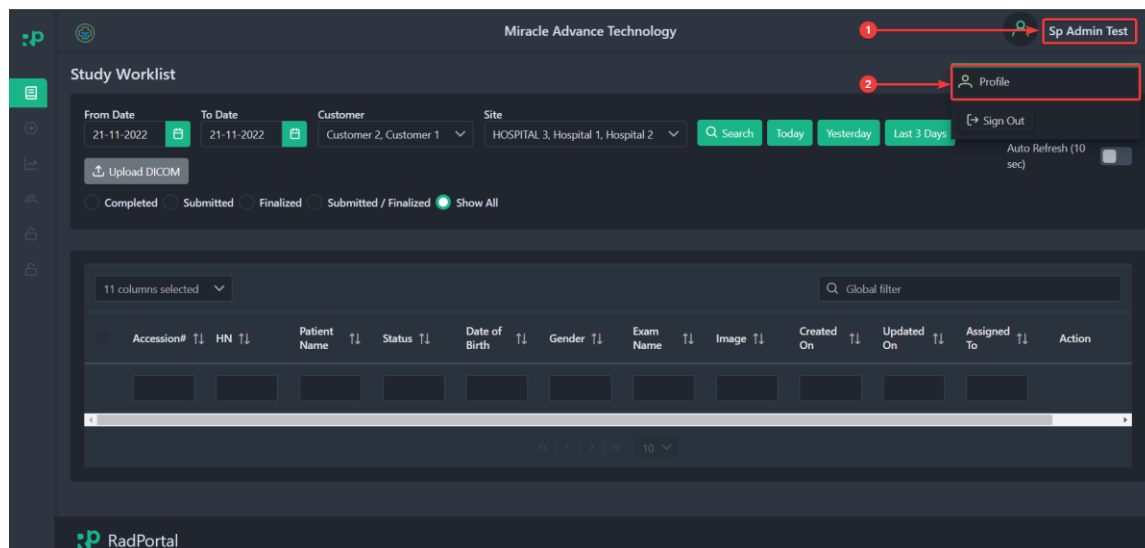


After clicking on The Profile, the system will redirect you to your profile page. On this page, you can change your Password. To change your user Password, click on the **Change User Password** button.

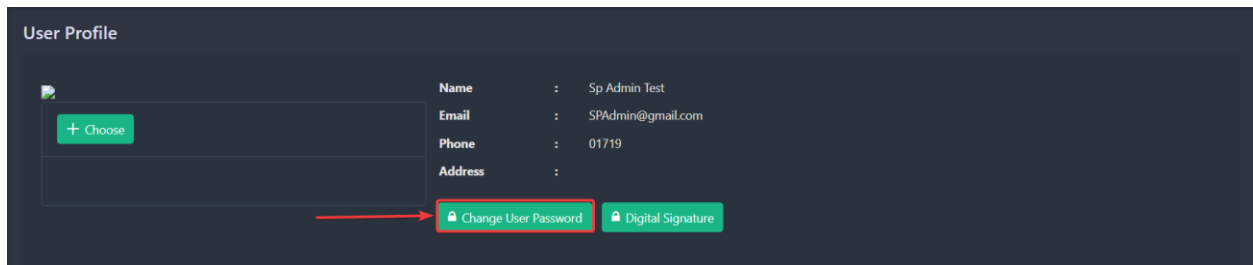

The screenshot shows the 'User Profile' page. On the left, there is a profile picture placeholder with a '+ Choose' button. To the right, user details are listed: Name (Sp Admin Test), Email (SPAdmin@gmail.com), Phone (01719), and Address. At the bottom, there are two buttons: 'Change User Password' (highlighted with a red box and an arrow) and 'Digital Signature'.


After clicking, you should see a new page where you have to confirm your **Old Password [1]**, enter your **New Password [2]**, and **confirm your new Password [3]**. Click on the **Update button [4]** to finalize your new password setup.

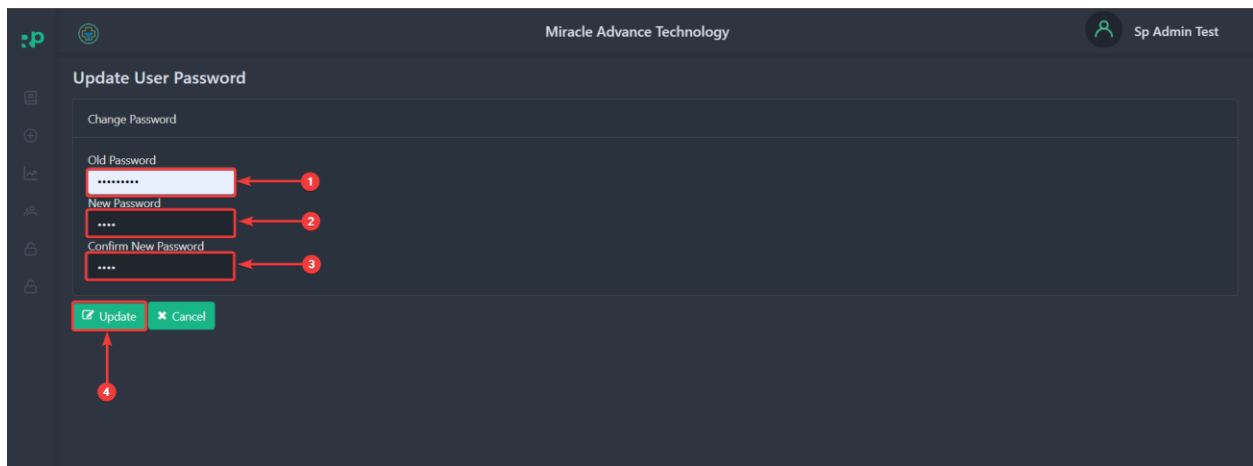


The screenshot shows the 'Update User Password' page. It has a header with the 'Miracle Advance Technology' logo and the user 'Sp Admin Test'. The main section is titled 'Change Password' and contains three input fields: 'Old Password' (marked with a red box and '1'), 'New Password' (marked with a red box and '2'), and 'Confirm New Password' (marked with a red box and '3'). At the bottom, there are two buttons: 'Update' (marked with a red box and '4') and 'Cancel'.

3.3.2 Changing Profile Photo


To change your profile photo, go back to your profile again and click on the Choose button.

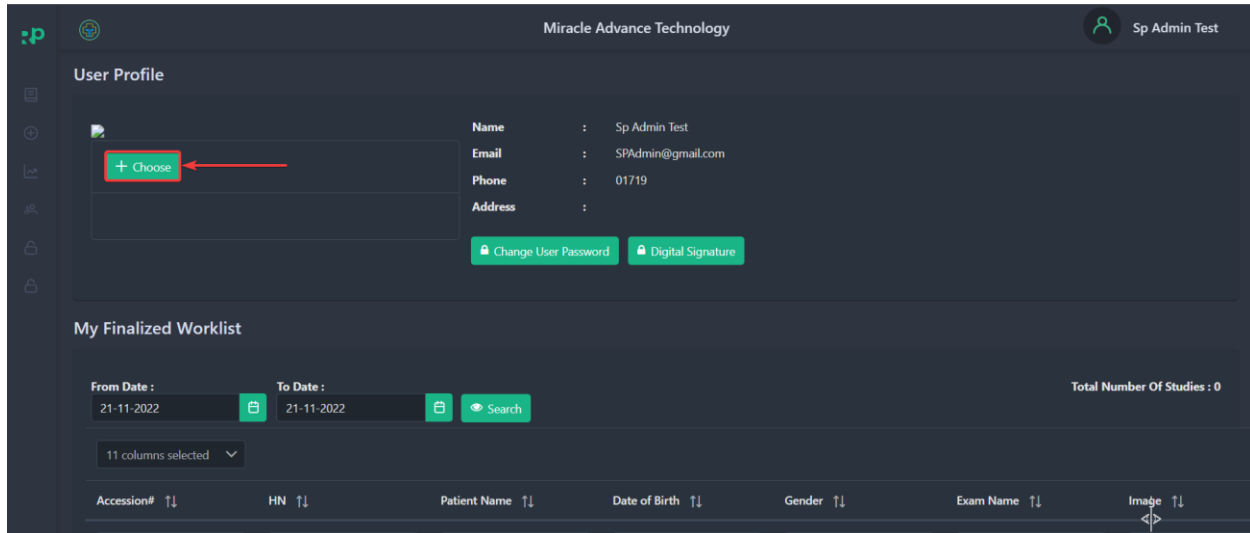

After clicking on the button, you can choose a photo for your profile. You can select your photos directly from your computer. After you upload your profile photo, you will see the preview in the section below.

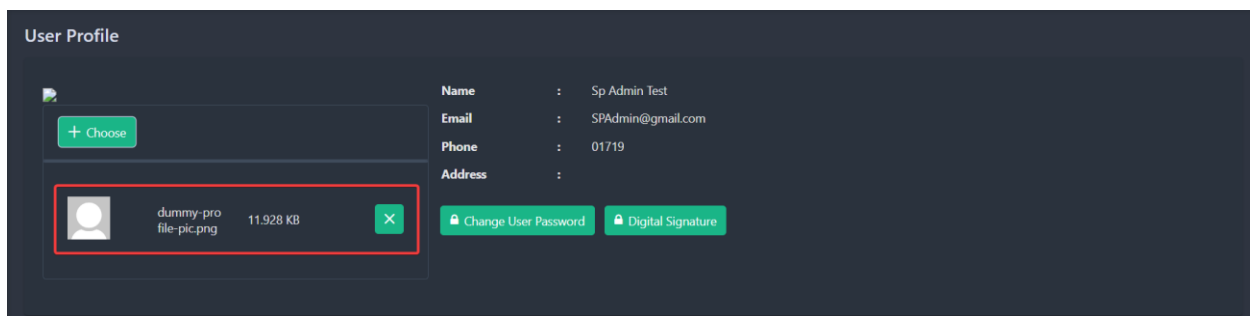


You can remove the photo by clicking on the green cross icon.

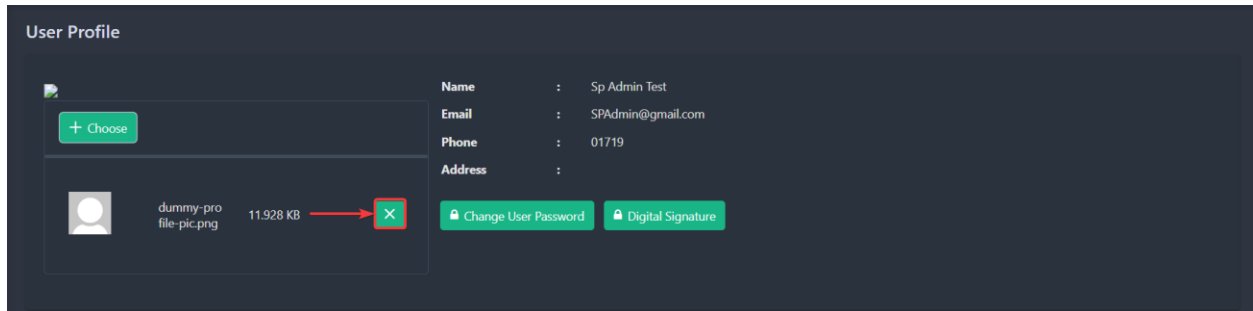


3.3.3 Adding Digital Signature


To add your Digital Signature, go back to your profile and click the **Digital Signature** button.

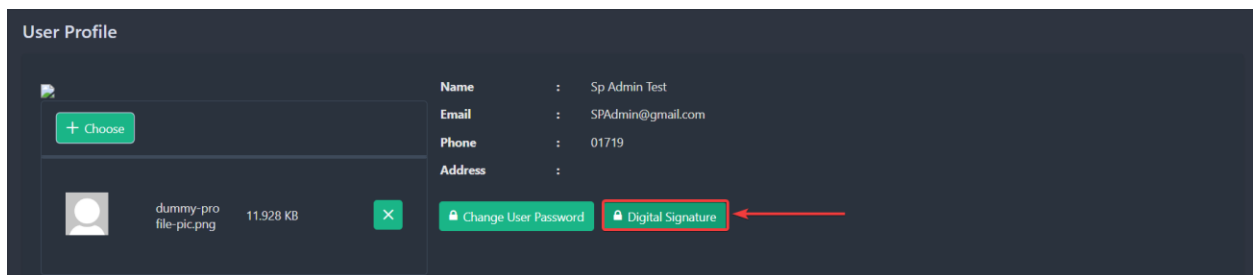


After clicking, you should see a new page where you will be able to upload your digital signature

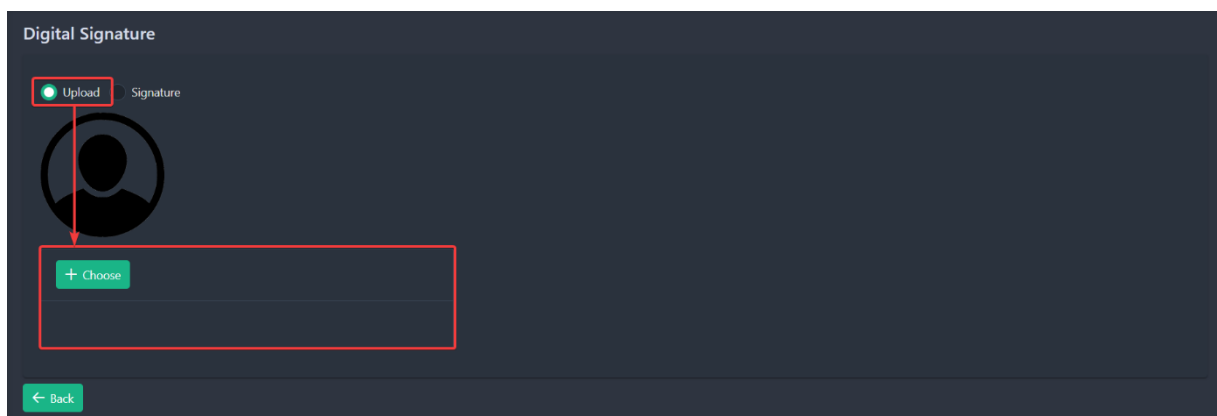

Or you can directly sign on the white blank section.

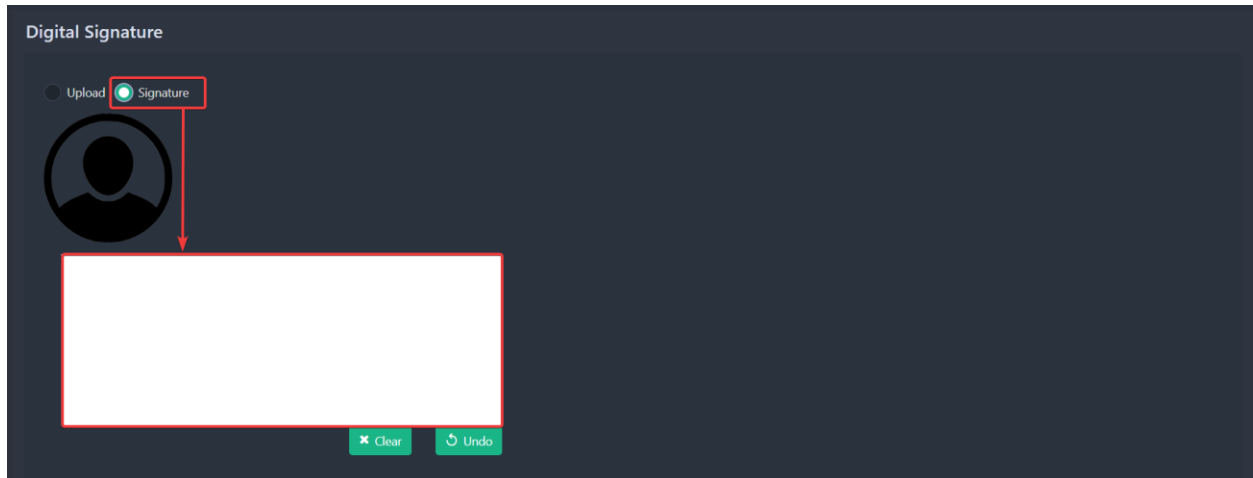


You can Clear your signature or undo a specific part of the signature by clicking on the **Clear** or **Undo** Button.

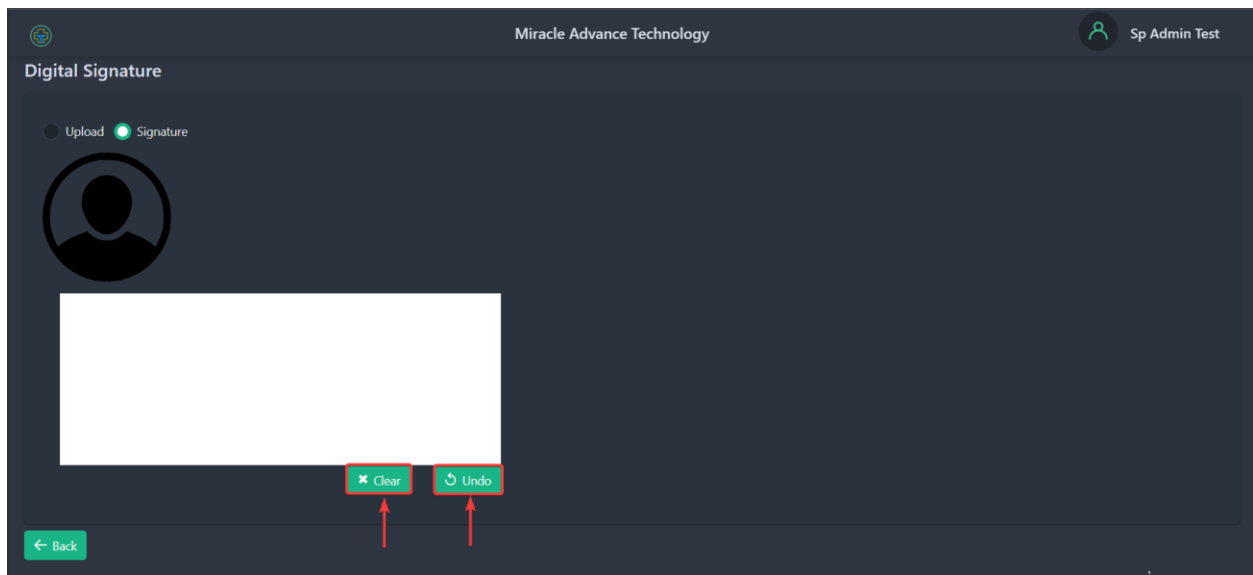


3.4 Software Version and Updates

You can check The RadPortal version in the footer of the RadPortal website. The latest Radportal version is 1.0.0

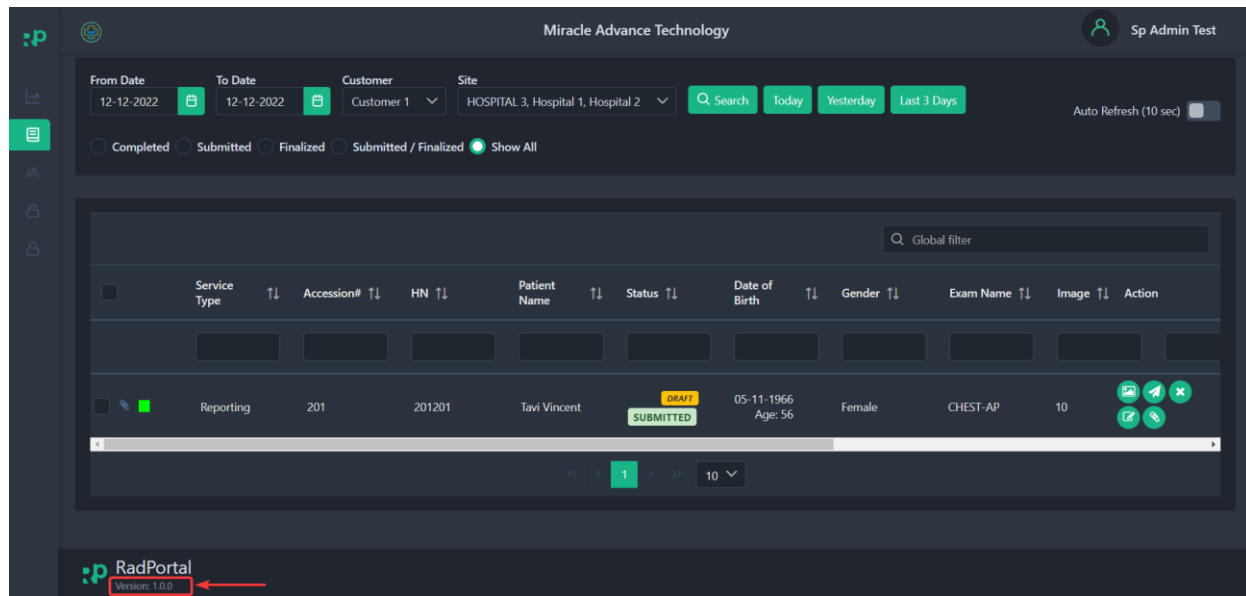


3.5 Sign Out from The System

To Sign Out from the RadPortal system, click on your **username** [1] in the top right corner and click on **Sign Out** [2].

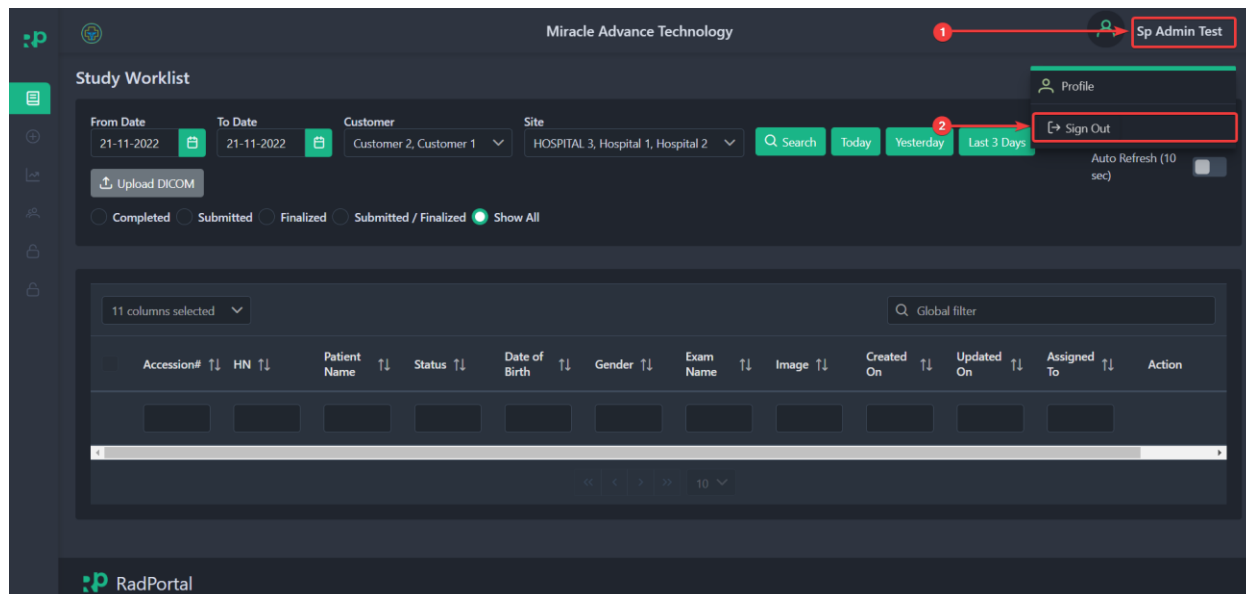

4. Exploring Worklist Menu

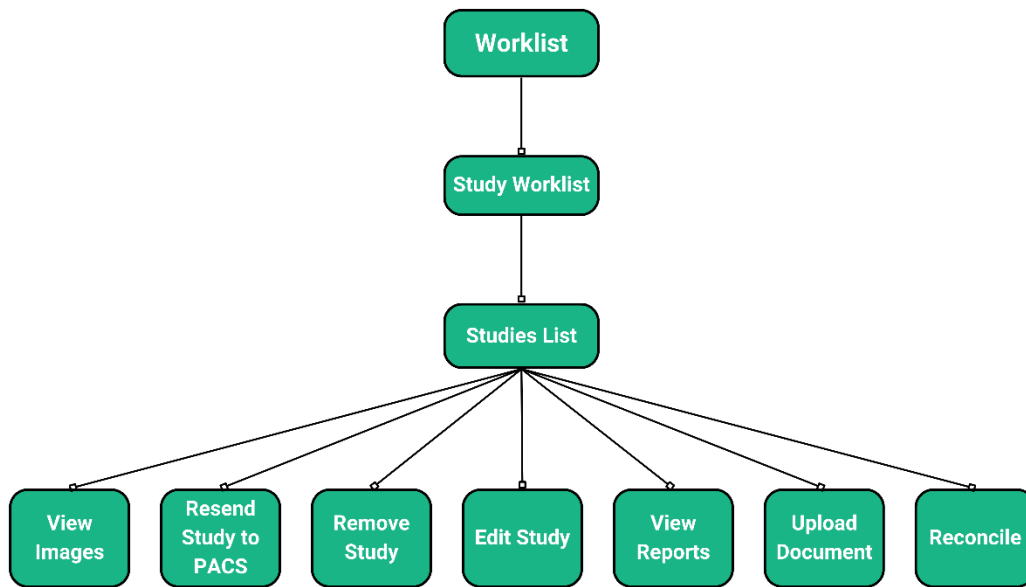

This menu in the RadPortal system allows you to keep track of all New, Submitted, and Finalized studies, or you can search for studies based on custom criteria. Using this Worklist Menu, the hospitals and the admin groups can Upload DICOM files manually. The worklist Menu page should look like the screenshot below.

Study Worksheet

From Date

22-11-2022

To Date

22-11-2022

Customer

Customer 2, Customer 1

Site

HOSPITAL 3, Hospital 1, Hospital 2

Search

Today

Yesterday

Auto Refresh (10 sec) ☐

Last 3 Days

Upload DICOM

☐ Completed

☐ Submitted

☐ Finalized

☐ Submitted / Finalized

☒ Show All

11 columns selected

Global filter

☐	Accession# ↑↓	HN ↑↓	Patient Name ↑↓	Status ↑↓	Date of Birth ↑↓	Gender ↑↓	Exam Name ↑↓	Image ↑↓	Created On ↑↓	Update On	Action
☐											

<

<<

>

>>

10

>

4.1 How to Search Studies

You can custom search based on criteria such as **From Date [1], To Date [2]**, or search based on the specific Customer by selecting a **customer [3]**, or you can also search by selecting a specific **Site or Hospital [4]**. Finally, click on the **Search [5]** button.

The screenshot shows the 'Study Worklist' search interface. At the top, there are five numbered red arrows pointing to specific search criteria: 1 points to the 'From Date' field (22-11-2022), 2 points to the 'To Date' field (22-11-2022), 3 points to the 'Customer' dropdown (Customer 2, Customer 1), 4 points to the 'Site' dropdown (HOSPITAL 3, Hospital 1, Hospital 2), and 5 points to the 'Search' button. Below these fields are buttons for 'Last 3 Days', 'Upload DICOM', and 'Today'/'Yesterday'. There is also an 'Auto Refresh (10 sec)' toggle. Below the search bar, there is a table with 11 columns: Accession#, HN, Patient Name, Status, Date of Birth, Gender, Exam Name, Image, Created On, Updat On, and Action. The table is currently empty, and the bottom shows pagination controls.


After clicking, you should see all studies/cases that match your custom criteria.

You can even see all studies/cases created **Today [1], Yesterday [2]**, or in the last 3 days [3].


This screenshot shows the 'Study Worklist' interface with a sidebar on the left containing links like 'Worklist', 'Registration', 'BI', 'User list', 'Org Setup', and 'Unlock Worklist'. The main search area has the same fields as the previous screenshot. In this view, three red arrows point to the 'Last 3 Days' button (labeled 3), the 'Today' button (labeled 1), and the 'Yesterday' button (labeled 2). The table below the search bar is empty, and the bottom shows pagination controls.

4.2 How to Sort Studies


You can easily sort study-list based on status, such as **Completed [1]**, **Submitted [2]**, **Finalized [3]**, or both **Submitted and Finalized [4]**. To sort studies, just click on the radio button you see before each status section.

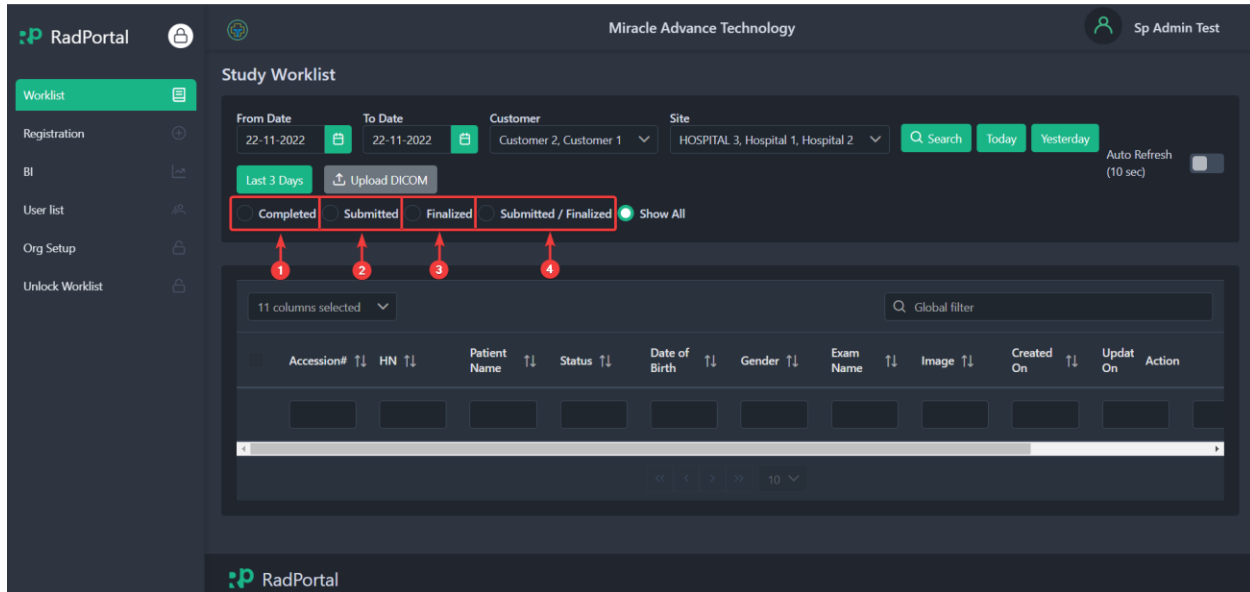

Or you can view all studies by clicking on the **Radio Button** you see before **Show All** text.

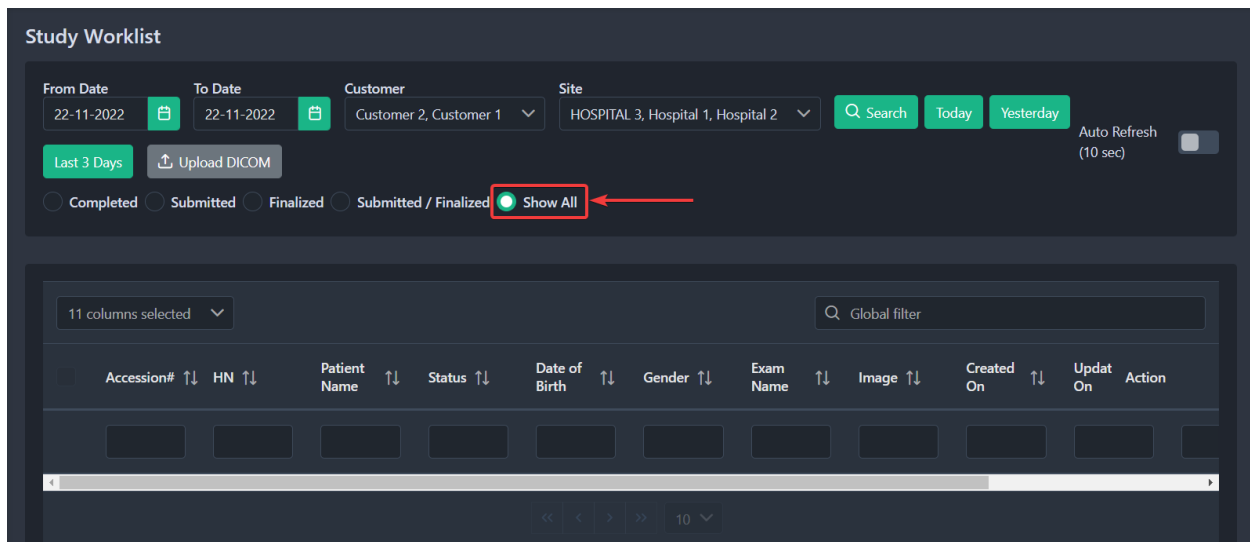

You can enable auto-refresh by clicking on the **auto-refresh slide button**. After enabling auto-refresh, the page will refresh continuously after every 10 seconds.

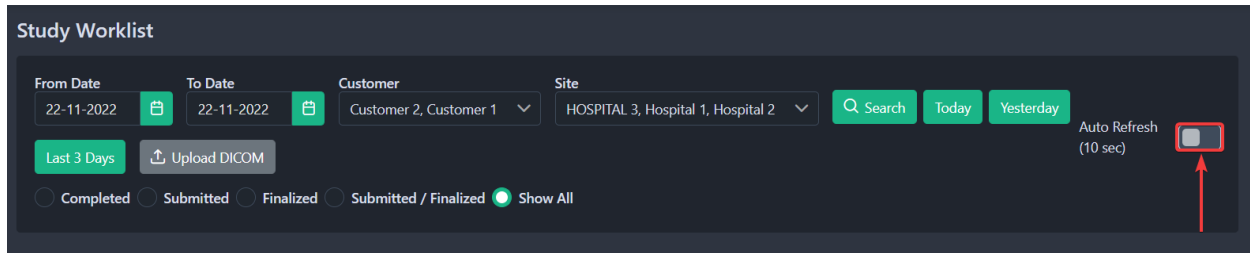


4.3 How Study information is Organized In Study worklist


Currently, there are 13 information types on the study worklist page, such as:

1. **Service Type:** It shows the service type of a specific study.
2. **Accession Number:** It is a unique identifier of a study for the RadPortal system.
3. **HN (Hospital Number):** is a unique identifier generated on the hospital end.
4. **Patient Name:** It shows the patient's Name.
5. **Status:** It shows the status of a particular case/study.
6. **Date of Birth:** Shows the patient's Date of birth.
7. **Gender:** Show the patient's gender.
8. **Exam Name:** It shows the exam name of a case/study
9. **Image:** It shows the total image count of a specific study.
10. **Created On:** It shows the creation date and time of a particular case/study.
11. **Updated On:** It shows the last update date of a particular case/study.
12. **Assign To:** It shows the Name of a radiologist assigned to a particular case/study.
13. **Modality:** It shows the test code (this code is collected from the DICOM Header) of a particular case/study.


To activate an information type, click on **Choose Columns [1]** and find the information type you want to activate. To activate, click on the **Checkbox [2]** before the information type.

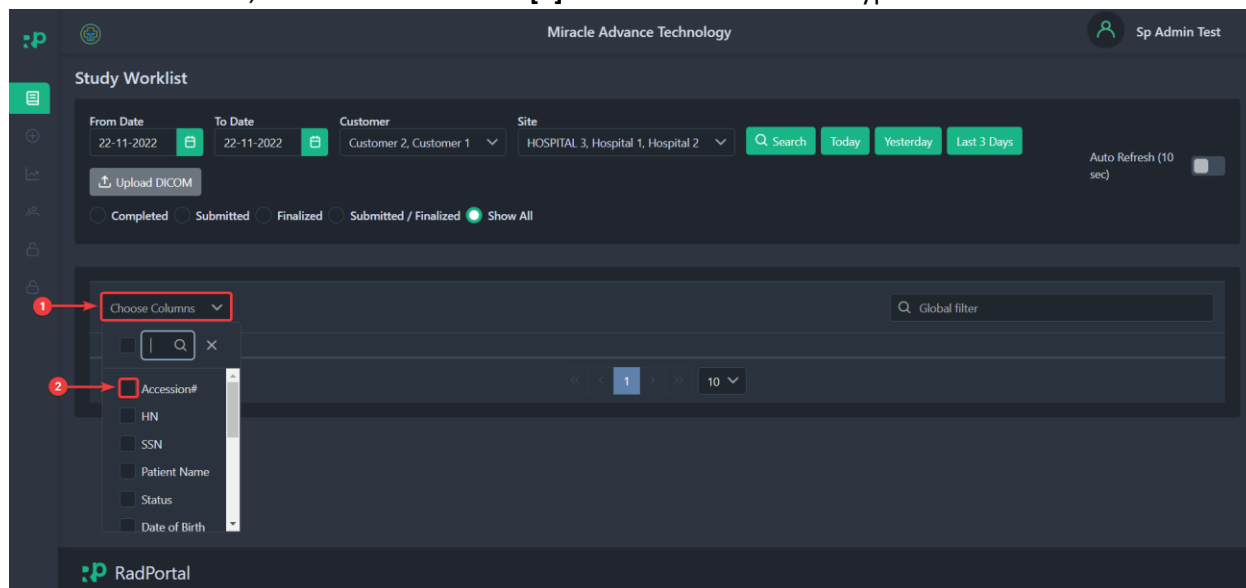


After selecting, you should see that the specific information type is now activated.

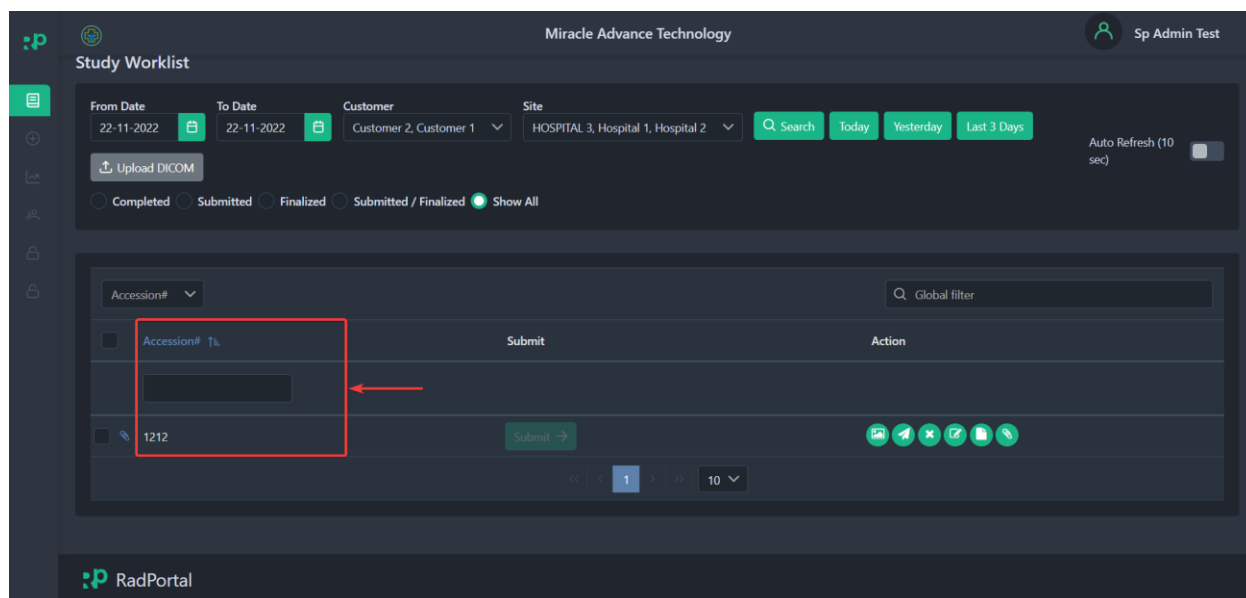

You can activate all information types by clicking on **Choose Columns [1]** and then clicking on **Checkbox [2]** at the top of the dropdown.

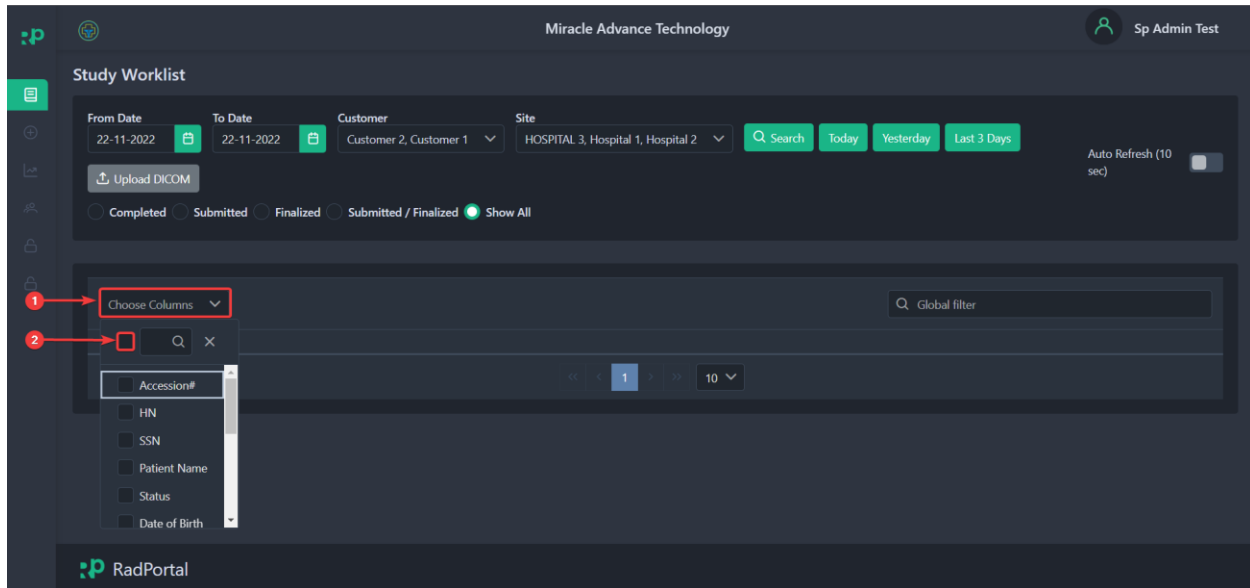

After clicking on the Checkbox, you should see that all information types are now activated.

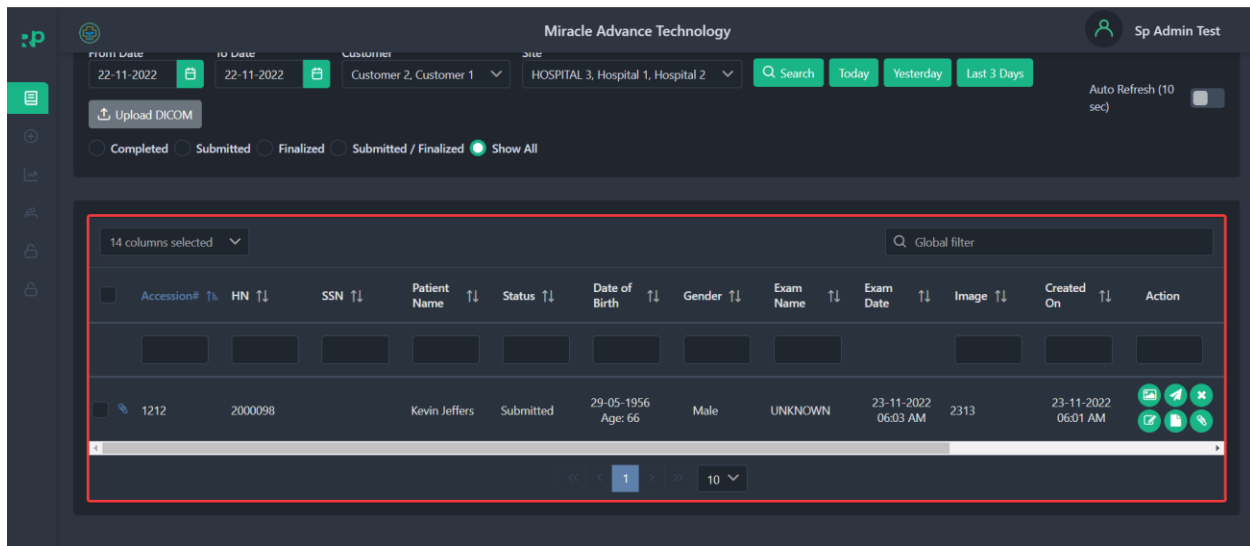


4.4 How to Review Images, Resend To PACS, Cancel, Edit, View Report and Attach Additional Documents to Studies

On the Study Worklist page, you can review Images of a Specific Study, resend them to PACS (Picture archiving and communication system), add, reconcile or delete studies, and attach necessary documents to existing studies. You can also assign radiologists to a specific study.

4.4.1 Review Images of a Specific Study

You can easily access all images of a specific study by clicking on the **Image Icon** in the screenshot below.

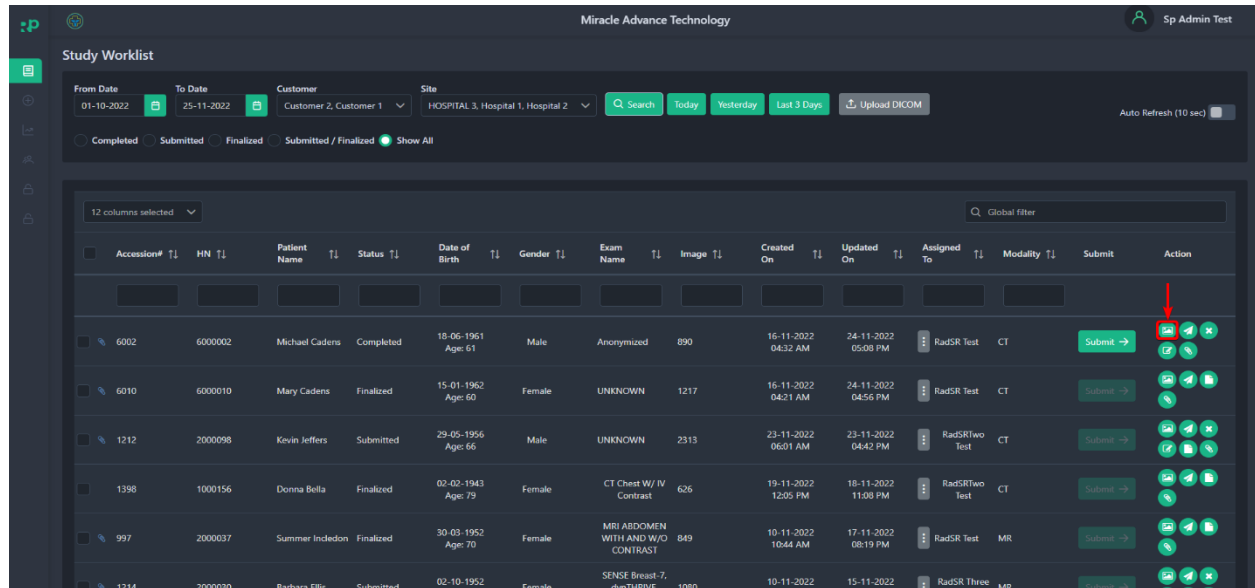


The screenshot shows the 'Study Worklist' page with a table of study records. The table has columns for Accession#, HN, Patient Name, Status, Date of Birth, Gender, Exam Name, Image, Created On, Updated On, Assigned To, Modality, Submit, and Action. A red arrow points to the 'Image Icon' in the 'Action' column for the first record.

Accession#	HN	Patient Name	Status	Date of Birth	Gender	Exam Name	Image	Created On	Updated On	Assigned To	Modality	Submit	Action
6002	6000002	Michael Cadens	Completed	18-06-1961 Age: 61	Male	Anonymized	890	16-11-2022 04:32 AM	24-11-2022 05:08 PM	RadSR Test	CT	Submit →	Image Icon
6010	6000010	Mary Cadens	Finalized	15-01-1962 Age: 60	Female	UNKNOWN	1217	16-11-2022 04:21 AM	24-11-2022 04:56 PM	RadSR Test	CT	Submit →	Image Icon
1212	2000098	Kevin Jeffers	Submitted	29-05-1956 Age: 66	Male	UNKNOWN	2313	23-11-2022 06:01 AM	23-11-2022 04:42 PM	RadSRtwo Test	CT	Submit →	Image Icon
1398	1000156	Donna Bella	Finalized	02-02-1943 Age: 79	Female	CT Chest W/ I/M Contrast	626	19-11-2022 12:05 PM	18-11-2022 11:08 PM	RadSRtwo Test	CT	Submit →	Image Icon
997	2000037	Summer Incledon	Finalized	30-03-1952 Age: 70	Female	MRI ABDOMEN WITH AND W/O CONTRAST	849	10-11-2022 10:44 AM	17-11-2022 08:19 PM	RadSR Test	MR	Submit →	Image Icon
1214	2000030	Barbara Ellis	Submitted	02-10-1952	Female	SENSE Breast-7, dynTHRIVE	1080	10-11-2022	15-11-2022	RadSR Three	MR	Submit →	Image Icon

After clicking, you should see a popup window where you must sign in with a specific **Username [1]** and **Password [2]**. This step ensures that only authorized personnel can view the medical test images.

Finally, click on the **Sign In[3]** button to sign in.

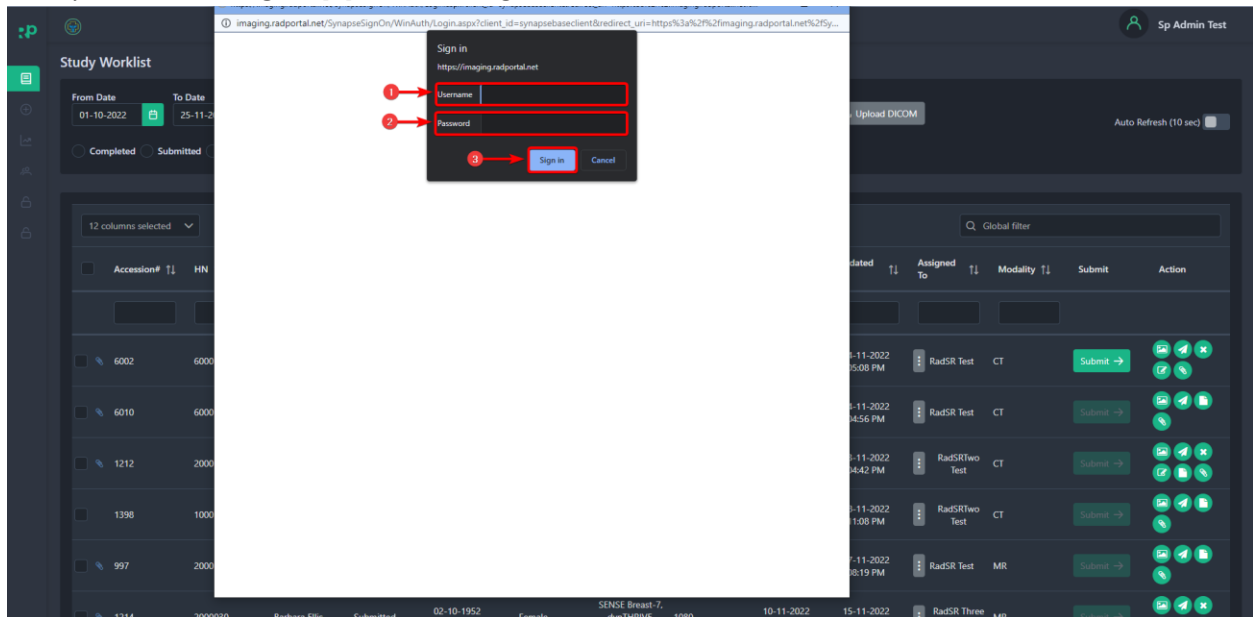

After successfully logging in to the system, you should be able to see all the test images of that specific Patient's medical test.



4.4.2 Resend to PACS (Picture archiving and communication system)

It's possible to resend studies from the local archive to the PACS if they were not correctly sent before.

To resend studies from the local archive to the destination PACS node, select the study (or multiple studies) and click on the **Resend to PACS** button.

The screenshot shows the 'Study Worklist' interface with a table of medical studies. The table has columns for Accession#, HN, Patient Name, Status, Date of Birth, Gender, Exam Name, Image, Created On, Updated On, Assigned To, Modality, Submit, and Action. A red arrow points to the 'Resend to PACS' button in the Action column for the first study (Accession# 6002).

Accession#	HN	Patient Name	Status	Date of Birth	Gender	Exam Name	Image	Created On	Updated On	Assigned To	Modality	Submit	Action
6002	6000002	Michael Cadens	Completed	18-06-1961 Age: 61	Male	Anonymized	890	16-11-2022 04:32 AM	24-11-2022 05:08 PM	RadSR Test	CT	Submit →	[Resend to PACS]
6010	6000010	Mary Cadens	Finalized	15-01-1962 Age: 60	Female	UNKNOWN	1217	16-11-2022 04:21 AM	24-11-2022 04:56 PM	RadSR Test	CT	Submit →	[Resend to PACS]
1212	2000098	Kevin Jeffers	Submitted	29-05-1956 Age: 66	Male	UNKNOWN	2313	23-11-2022 06:01 AM	23-11-2022 04:42 PM	RadSRtwo Test	CT	Submit →	[Resend to PACS]
1398	1000156	Donna Bella	Finalized	02-02-1943 Age: 79	Female	CT Chest W/ IV Contrast	626	19-11-2022 12:05 PM	18-11-2022 11:08 PM	RadSRtwo Test	CT	Submit →	[Resend to PACS]


4.4.3 Cancel (Remove) Studies

Studies can be easily deleted from the RadPortal System. To remove studies, click on the **Cancel button**.

The screenshot shows the 'Study Worklist' interface with a table of medical studies. The table has columns for Accession#, HN, SSN, Patient Name, Status, Date of Birth, Gender, Exam Name, Exam Date, Image, Created On, Updated On, Assigned To, and Action. A red arrow points to the 'Cancel' button in the Action column for the first study (Accession# 6002).

Accession#	HN	SSN	Patient Name	Status	Date of Birth	Gender	Exam Name	Exam Date	Image	Created On	Updated On	Assigned To	Action
6002	6000002		Michael Cadens	Completed	18-06-1961 Age: 61	Male	Anonymized	16-11-2022 04:34 AM	890	16-11-2022 04:32 AM	24-11-2022 05:08 PM	RadSR Test	[Cancel]
6010	6000010		Mary Cadens	Finalized	15-01-1962 Age: 60	Female	UNKNOWN	16-11-2022 04:23 AM	1217	16-11-2022 04:21 AM	24-11-2022 04:56 PM	RadSR Test	[Cancel]
1212	2000098		Kevin Jeffers	Submitted	29-05-1956 Age: 66	Male	UNKNOWN	23-11-2022 06:03 AM	2313	23-11-2022 06:01 AM	23-11-2022 04:42 PM	RadSRtwo Test	[Cancel]
1398	1000156		Donna Bella	Finalized	02-02-1943 Age: 79	Female	CT Chest W/ IV Contrast	19-11-2022 12:06 PM	626	19-11-2022 12:05 PM	18-11-2022 11:08 PM	RadSRtwo Test	[Cancel]
997	2000037		Summer Incledon	Finalized	30-03-1952 Age: 70	Female	MRI ABDOMEN WITH AND W/O CONTRAST	10-11-2022 10:45 AM	849	10-11-2022 10:44 AM	17-11-2022 08:19 PM	RadSR Test	[Cancel]
1214	2000030		Barbara Ellis	Submitted	02-10-1952 Age: 70	Female	SENSE Breast-7, dynTHRIVE, diffusion	10-11-2022 05:22 AM	1080	10-11-2022 05:21 AM	15-11-2022 05:03 PM	RadSR Three Test	[Cancel]

After clicking on the cancel button, you should see a popup window where you can specify the reason for deleting that particular study in the **Reason** [1] text field. To confirm the deletion of that study, click on the **Yes** [2] button, or if you do not want to delete the study, click on the **No** [3] button.

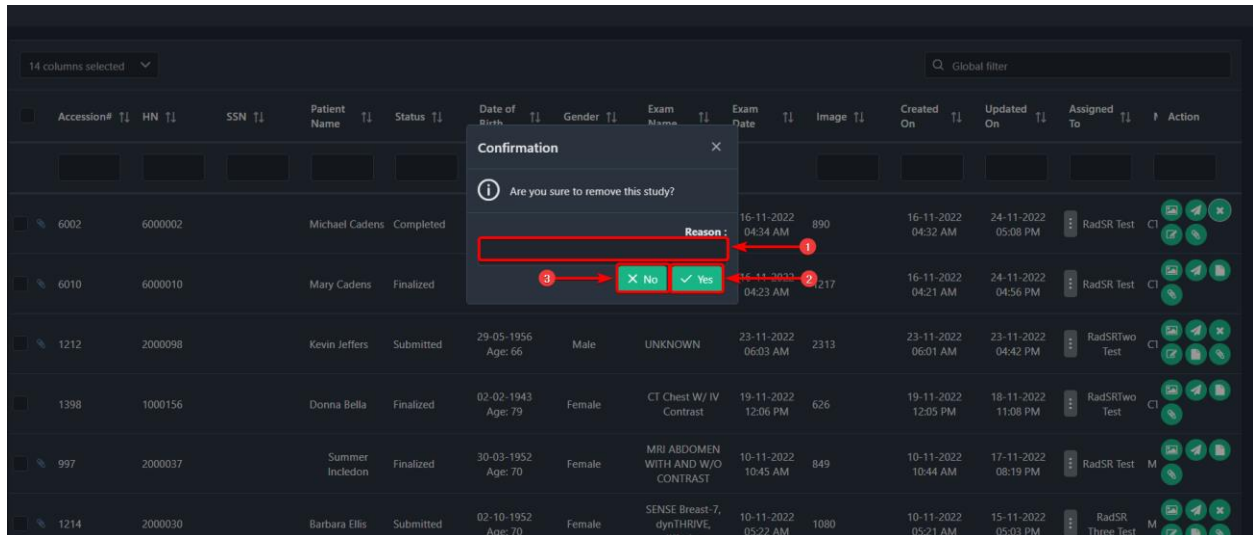

4.4.4 Reconcile Studies

You can easily reconcile studies in the RadPortal. By default, the RadPortal system collects information from the DICOM Files. So, you may need to reconcile/edit studies to attach additional information. You can edit/add information such as SSN, HN, First Name, Middle Name, Last Name, Gender, Age, DOB, Accession No, Study Date, and Study Govt Code.

To reconcile a particular study, click on the **Reconcile button**.

Miracle Advance Technology Sp Admin Test

T#	SSN	Patient Name	Status	Date of Birth	Gender	Exam Name	Exam Date	Image	Created On	Updated On	Assigned To	Modality	Submit	Action
I002		Michael Cadens	Completed	18-06-1961 Age: 61	Male	Anonymized	16-11-2022 04:34 AM	890	16-11-2022 04:32 AM	24-11-2022 05:08 PM	RadSR Test	CT	Submit →	
I010		Mary Cadens	Finalized	15-01-1962 Age: 60	Female	UNKNOWN	16-11-2022 04:23 AM	1217	16-11-2022 04:21 AM	24-11-2022 04:56 PM	RadSR Test	CT	Submit →	
I098		Kevin Jeffers	Submitted	29-05-1956 Age: 66	Male	UNKNOWN	23-11-2022 06:03 AM	2313	23-11-2022 06:01 AM	23-11-2022 04:42 PM	RadSRTwo Test	CT	Submit →	
I156		Donna Bella	Finalized	02-02-1943 Age: 79	Female	CT Chest W/ IV Contrast	19-11-2022 12:06 PM	626	19-11-2022 12:05 PM	18-11-2022 11:08 PM	RadSRTwo Test	CT	Submit →	
I037		Summer Incledon	Finalized	30-03-1952 Age: 70	Female	MRI ABDOMEN WITH AND W/O CONTRAST	10-11-2022 10:45 AM	849	10-11-2022 10:44 AM	17-11-2022 08:19 PM	RadSR Test	MR	Submit →	
I030		Barbara Ellis	Submitted	02-10-1952 Age: 70	Female	SENSE Breast-7, dynTHRIVE, diffusion	10-11-2022 05:22 AM	1080	10-11-2022 05:21 AM	15-11-2022 05:03 PM	RadSR Three Test	MR	Submit →	
I06		Carotids Vincent	Finalized	06-11-1966 Age: 56	Female	UNKNOWN	09-11-2022 09:26 AM	772	09-11-2022 09:25 AM	08-11-2022 03:24 PM	RadSR Test	CT	Submit →	
I01		Tavi Vincent	Finalized	05-11-1966 Age: 56	Female	UNKNOWN	08-11-2022 06:40 AM	1314	08-11-2022 06:39 AM	07-11-2022 04:58 PM	RadSRTwo Test	CT	Submit →	

After clicking, you should see a popup where you can edit or update patient and study information.

Miracle Advance Technology Sp Admin Test

T#	SSN	Patient Name	Status	Date of Birth	Gender	Exam Name	Exam Date	Image	Created On	Updated On	Assigned To	Modality	Submit	Action
I002		Michael Cadens	Completed	18-06-1961 Age: 61	Male	Anonymized	16-11-2022 04:34 AM	890	16-11-2022 04:32 AM	24-11-2022 05:08 PM	RadSR Test	CT	Submit →	
I010		Mary Cadens	Finalized	15-01-1962 Age: 60	Female	UNKNOWN	16-11-2022 04:23 AM	1217	16-11-2022 04:21 AM	24-11-2022 04:56 PM	RadSR Test	CT	Submit →	
I098		Kevin Jeffers	Submitted	29-05-1956 Age: 66	Male	UNKNOWN	23-11-2022 06:03 AM	2313	23-11-2022 06:01 AM	23-11-2022 04:42 PM	RadSRTwo Test	CT	Submit →	
I156		Donna Bella	Finalized	02-02-1943 Age: 79	Female	CT Chest W/ IV Contrast	19-11-2022 12:06 PM	626	19-11-2022 12:05 PM	18-11-2022 11:08 PM	RadSRTwo Test	CT	Submit →	
I037		Summer Incledon	Finalized	30-03-1952 Age: 70	Female	MRI ABDOMEN WITH AND W/O CONTRAST	10-11-2022 10:45 AM	849	10-11-2022 10:44 AM	17-11-2022 08:19 PM	RadSR Test	MR	Submit →	
I030		Barbara Ellis	Submitted	02-10-1952 Age: 70	Female	SENSE Breast-7, dynTHRIVE, diffusion	10-11-2022 05:22 AM	1080	10-11-2022 05:21 AM	15-11-2022 05:03 PM	RadSR Three Test	MR	Submit →	
I06		Carotids Vincent	Finalized	06-11-1966 Age: 56	Female	UNKNOWN	09-11-2022 09:26 AM	772	09-11-2022 09:25 AM	08-11-2022 03:24 PM	RadSR Test	CT	Submit →	
I01		Tavi Vincent	Finalized	05-11-1966 Age: 56	Female	UNKNOWN	08-11-2022 06:40 AM	1314	08-11-2022 06:39 AM	07-11-2022 04:58 PM	RadSRTwo Test	CT	Submit →	

Study Information

Patient Information

SSN: HN: First Name: Middle Name:

Last Name: Gender: DOB: Age:

Study Information

Accession No.: Study Date *: Study *: Govt Code *:

NOTE: When reconciling the study, you must fill in the Study Date, Study, and Govt Code

When you are done, click the **Save** button to save the updated information.

Study Information

Patient Information

SSN: 0011 HN: 6000002 First Name: Cadens Middle Name: Middle Name

Last Name: Michael Gender: Male DOB: 18-06-1961 Age: 61

Study Information

Accession No.: 6002 Study Date *: 16-11-2022 Study *: Chest Govt Code *: Chest Lung

Save Close

4.4.5 How to Submit Study/Case

You have to submit the Study/Case before you assign a radiologist. After you submit a study, it will automatically be sent to the PACS (Picture Archiving and Communication System). When a case has the **Completed** status, you should see the **Submit Button** Activated under the **Submit** field.

Now to submit the study, click on the **Submit** button.

Miracle Advance Technology

Sp Admin Test

13 columns selected


Global filter

	Patient Name	Status	Date of Birth	Gender	Exam Name	Image	Created On	Updated On	Assigned To	Modality	Submit	Action
124	Melissa Dayton	Finalized	13-12-1955 Age: 67	Female	MRI CARDIAC FUNCTION	12	26-11-2022 11:45 AM	30-11-2022 01:39 PM	RadSR Test	MR	Submit ->	
178	Photichanprasert Phumiphat	Completed	18-04-1989 Age: 33	Male	Chest	1	30-11-2022 08:35 AM	29-11-2022 06:31 PM	RadSR Test	DX	Submit ->	
156	Donna Bella	Submitted	02-02-1943 Age: 79	Female	BREAST*Routine	13	26-11-2022 11:45 AM	29-11-2022 05:56 PM	RadSR Test	MR	Submit ->	
16	Carotids Vincent	Finalized	06-11-1966 Age: 56	Female	UNKNOWN	10	26-11-2022 11:45 AM	28-11-2022 07:55 PM	RadSR Test	CT	Submit ->	
102	Michael Cadens	Completed	18-06-1961 Age: 61	Male	Anonymized	11	26-11-2022 11:45 AM	28-11-2022 07:55 PM	RadSR Test	CT	Submit ->	

NOTE: After the hospital submits a study, the hospital will not be able to cancel the study.

4.4.6 View Study Reports

When a radiologist submits a report for a specific case/study, you can view the report by clicking on the **View Report** button.

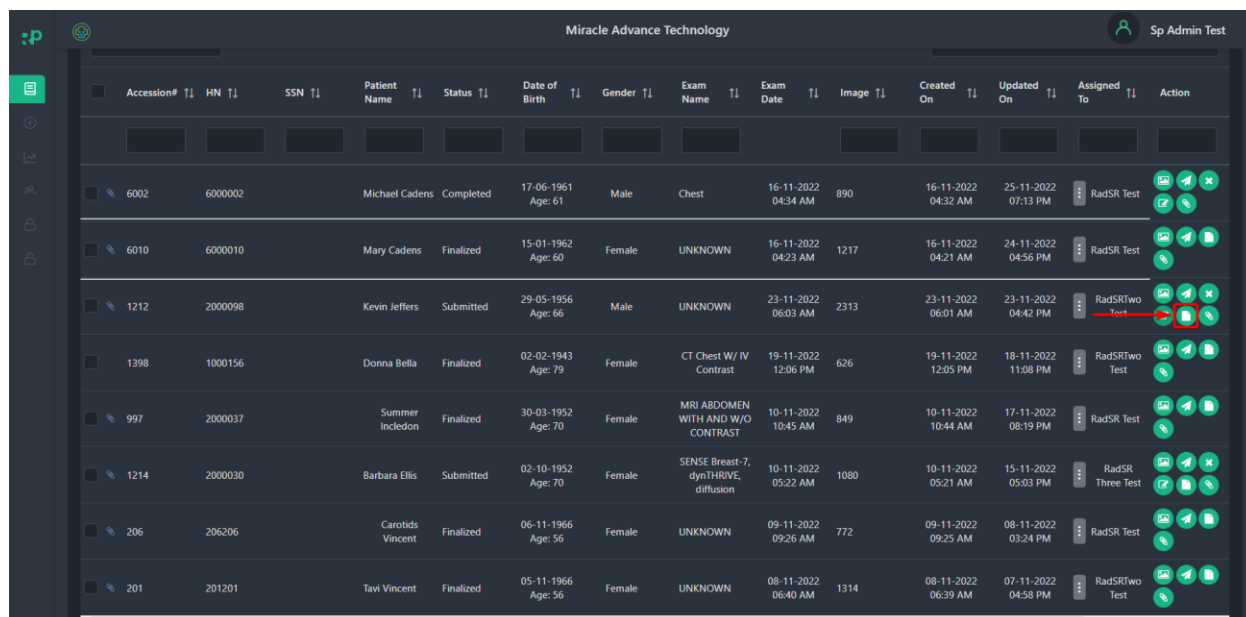

Accession#	HN	SSN	Patient Name	Status	Date of Birth	Gender	Exam Name	Exam Date	Image	Created On	Updated On	Assigned To	Action
6002	6000002		Michael Cadens	Completed	17-06-1961 Age: 61	Male	Chest	16-11-2022 04:34 AM	890	16-11-2022 04:32 AM	25-11-2022 07:13 PM	RadSR Test	
6010	6000010		Mary Cadens	Finalized	15-01-1962 Age: 60	Female	UNKNOWN	16-11-2022 04:23 AM	1217	16-11-2022 04:21 AM	24-11-2022 04:56 PM	RadSR Test	
1212	2000098		Kevin Jeffers	Submitted	29-05-1956 Age: 66	Male	UNKNOWN	23-11-2022 06:03 AM	2313	23-11-2022 06:01 AM	23-11-2022 04:42 PM	RadSRTwo Test	
1398	1000156		Donna Bella	Finalized	02-02-1943 Age: 79	Female	CT Chest W/ IV Contrast	19-11-2022 12:06 PM	626	19-11-2022 12:05 PM	18-11-2022 11:08 PM	RadSRTwo Test	
997	2000037		Summer Incedon	Finalized	30-03-1952 Age: 70	Female	MRI ABDOMEN WITH AND W/O CONTRAST	10-11-2022 10:45 AM	849	10-11-2022 10:44 AM	17-11-2022 08:19 PM	RadSR Test	
1214	2000030		Barbara Ellis	Submitted	02-10-1952 Age: 70	Female	SENSE Breast-7, dynTHRIVE, diffusion	10-11-2022 05:22 AM	1080	10-11-2022 05:21 AM	15-11-2022 05:03 PM	RadSR Three Test	
206	206206		Carotids Vincent	Finalized	06-11-1966 Age: 56	Female	UNKNOWN	09-11-2022 09:26 AM	772	09-11-2022 09:25 AM	08-11-2022 03:24 PM	RadSR Test	
201	201201		Tavi Vincent	Finalized	05-11-1966 Age: 56	Female	UNKNOWN	08-11-2022 06:40 AM	1314	08-11-2022 06:39 AM	07-11-2022 04:58 PM	RadSRTwo Test	

After clicking, you can see the Finalized report.

SENSE BREAST-7, DYNTHRIVE,

Site : Hospital 1
Code :
Assigned By : RadSR Three Test
Finalized By :
Accession# : 1214
Type : MR
Assigned On : 21/11/2022 11:47 AM
Finalized On :

Findings

test


NOTE: You will only see the **View Report** button once the Radiologist has finalized the report from their end.

4.4.7 Attach Additional Files To A Study


You can attach additional files to a study by clicking on the Document Upload button.

Miracle Advance Technology														Sp Admin Test	
Accession#	HN	SSN	Patient Name	Status	Date of Birth	Gender	Exam Name	Exam Date	Image	Created On	Updated On	Assigned To	Action		
6002	6000002		Michael Cadens	Completed	17-06-1961 Age: 61	Male	Chest	16-11-2022 04:34 AM	890	16-11-2022 04:32 AM	25-11-2022 07:13 PM	RadSR Test	  		
6010	6000010		Mary Cadens	Finalized	15-01-1962 Age: 60	Female	UNKNOWN	16-11-2022 04:23 AM	1217	16-11-2022 04:21 AM	24-11-2022 04:56 PM	RadSR Test	  		
1212	2000098		Kevin Jeffers	Submitted	29-05-1956 Age: 66	Male	UNKNOWN	23-11-2022 06:03 AM	2313	23-11-2022 06:01 AM	23-11-2022 04:42 PM	RadSR Two Test	  		
1398	1000156		Donna Bella	Finalized	02-02-1943 Age: 79	Female	CT Chest W/ IV Contrast	19-11-2022 12:06 PM	626	19-11-2022 12:05 PM	18-11-2022 11:08 PM	RadSR Two Test	  		
997	2000037		Summer Incledon	Finalized	30-03-1952 Age: 70	Female	MRI ABDOMEN WITH AND W/O CONTRAST	10-11-2022 10:45 AM	849	10-11-2022 10:44 AM	17-11-2022 08:19 PM	RadSR Test	  		
1214	2000030		Barbara Ellis	Submitted	02-10-1952 Age: 70	Female	SENSE Breast-7, dynTHRIVE, diffusion	10-11-2022 05:22 AM	1080	10-11-2022 05:21 AM	15-11-2022 05:03 PM	RadSR Three Test	  		
206	206206		Carotids Vincent	Finalized	06-11-1966 Age: 56	Female	UNKNOWN	09-11-2022 09:26 AM	772	09-11-2022 09:25 AM	08-11-2022 03:24 PM	RadSR Test	  		
201	201201		Tavi Vincent	Finalized	05-11-1966 Age: 56	Female	UNKNOWN	08-11-2022 06:40 AM	1314	08-11-2022 06:39 AM	07-11-2022 04:58 PM	RadSR Two Test	  		


After clicking, you should see a popup window where you can define document **Types [1]** such as **Report, Lab data, or Clinical diagnosis**. Select the document from the local storage by clicking the **Choose Files [2]** button. Add the Name of the document you just uploaded in the **Title Field [3]**, and finally, click on the **Upload button [4]**.

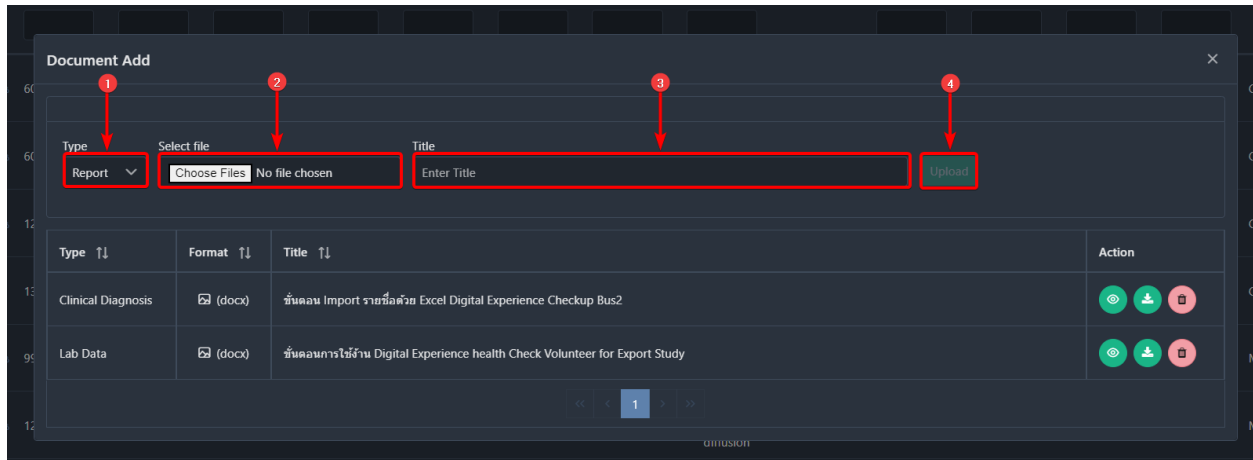

After completing the upload process, you should see the document in this list.

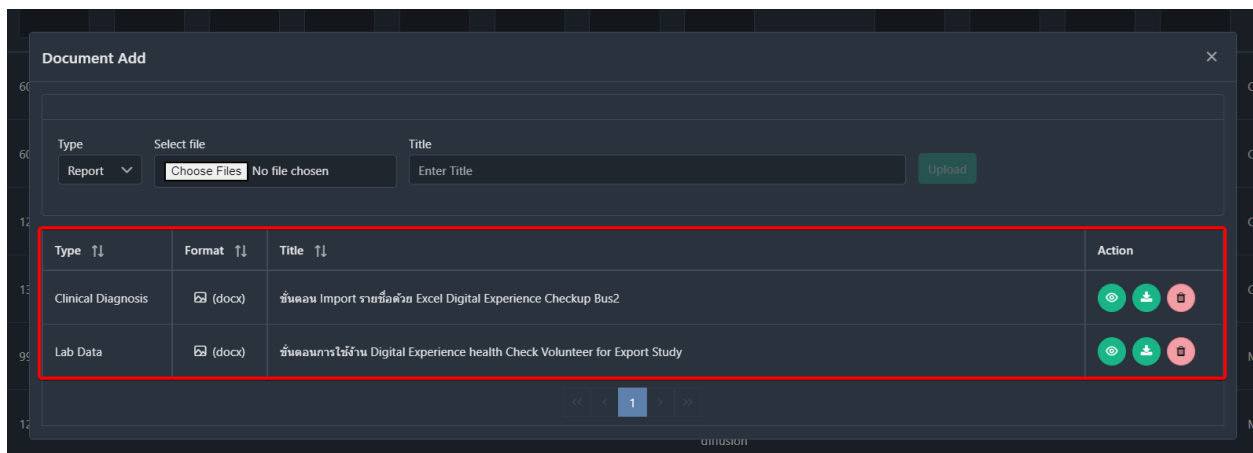


You can view, download or delete the document by clicking on the **View [1]**, **Download [2]**, and **Delete [3]** buttons.

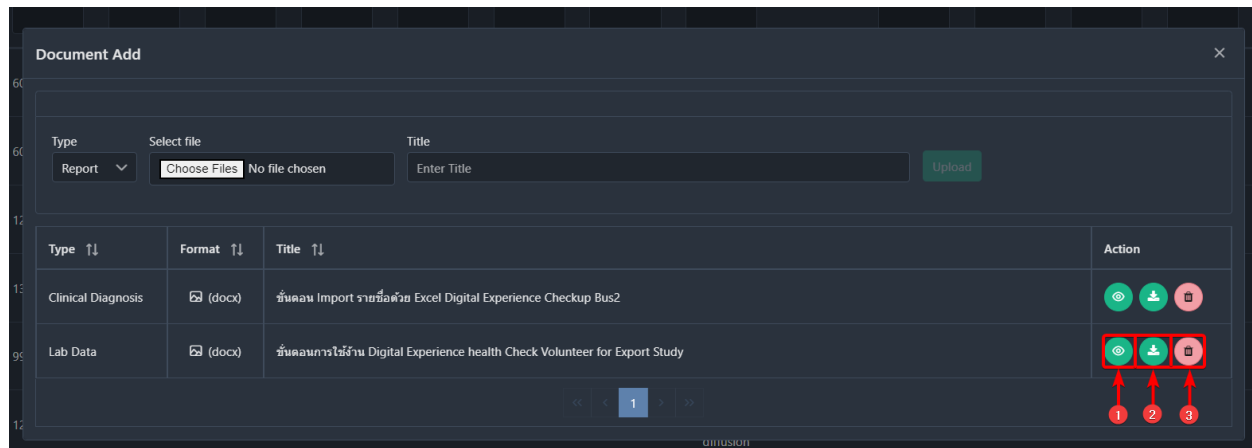

4.6 How to Submit Studies

When a Healthcare Service Provider uploads studies, you can see them in the worklist with the status of **Completed**.

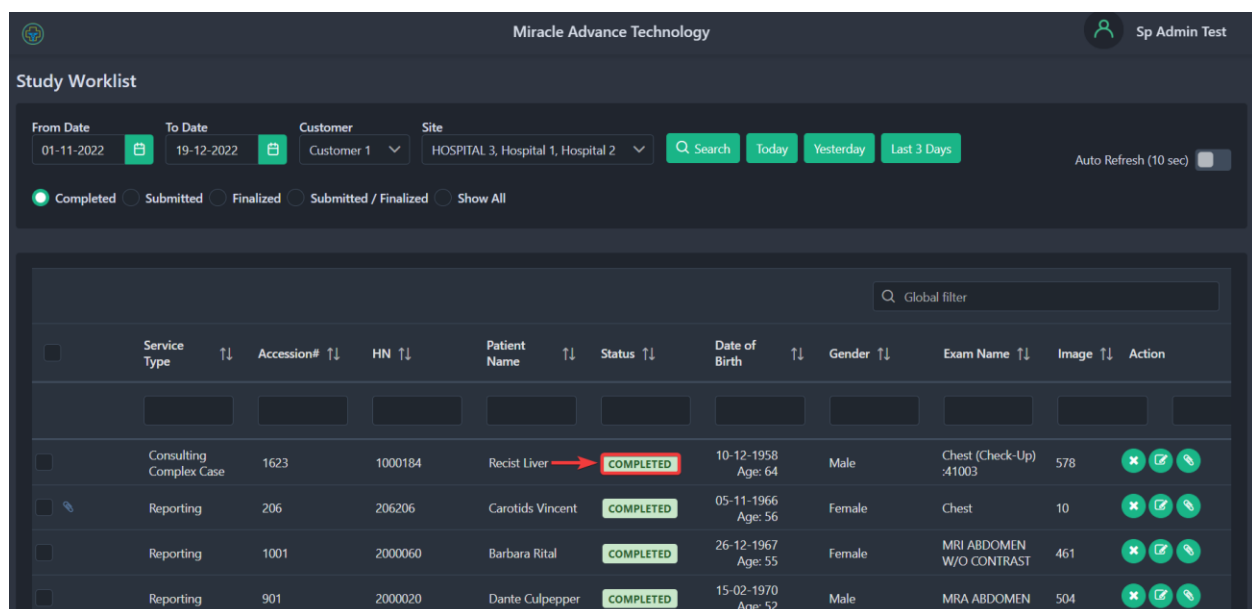

Now, you have to reconcile study data before you submit the study. To reconcile, click on the **Reconcile** button.

Miracle Advance Technology Sp Admin Test

Global filter

	Date of Birth	Gender	Exam Name	Image	Created On	Updated On	Assigned To	Modality	Submit	Action
PLETED	10-12-1958 Age: 64	Male	Chest (Check Up) :41003	578	06-12-2022 02:36 PM	14-12-2022 05:05 PM		CT	Submit →	
PLETED	05-11-1966 Age: 56	Female	Chest	10	09-12-2022 10:05 AM	14-12-2022 01:18 PM	RadSR Test	CT	Submit →	
PLETED	26-12-1967 Age: 55	Female	MRI ABDOMEN W/O CONTRAST	461	06-12-2022 02:40 PM	14-12-2022 01:18 PM	RadSR Test	MR	Submit →	
PLETED	15-02-1970 Age: 52	Male	MRA ABDOMEN	504	06-12-2022 02:43 PM	14-12-2022 01:18 PM	RadSRTwo Test	MR	Submit →	
PLETED	01-02-1943 Age: 79	Female	Cardiac^26_CoroCT/ (Adult)	13	09-12-2022 10:05 AM	14-12-2022 11:52 AM	RadSR Test	MR	Submit →	
PLETED	15-01-1962 Age: 60	Female	UNKNOWN	10	15-11-2022 02:21 PM			CT	Submit →	
PLETED	15-06-1961 Age: 61	Male	Chest	890	15-11-2022 02:32 PM			CT	Submit →	
PLETED	01-02-1943 Age: 79	Female	CT Chest W/ IV Contrast	596	06-12-2022 02:33 PM			CT	Submit →	

After clicking, you will see a popup window where you can update all patient information.

Study Information

Patient Information

Service Type* Consulting Complex Case Turnaround Time* GREEN (0 - 4) SSN* 221 HN* 1000184

First Name* Liver Middle Name Middle Name Last Name* Regist Gender* Male

DOB 10-12-1958 Age 64

Study Information

Accession No.* 1623 Study Date* 06-12-2022 Study* Select a Study Govt Code* Chest Lung

Save Close

To complete the reconciling process, you must assign **Study Code[1]** and **Govt Code[2]**.

Study Information [X]

Patient Information

Service Type*	Turnaround Time*	SSN*	HN*
Consulting Complex Case	GREEN (0 - 4)	221	1000184
First Name*	Middle Name	Last Name*	Gender*
Liver	Middle Name	Recist	Male
DOB	Age		
10-12-1958	64		

Study Information

Accession No.*	Study Date *	Study *	Govt Code *
1623	06-12-2022	Select a Study	Chest Lung

[Save] [Close]

After assigning study and govt codes, click on the **Save** button.

Study Information [X]

Patient Information

Service Type*	Turnaround Time*	SSN*	HN*
Consulting Complex Case	GREEN (0 - 4)	221	1000184
First Name*	Middle Name	Last Name*	Gender*
Liver	Middle Name	Recist	Male
DOB	Age		
10-12-1958	64		


Study Information

Accession No.*	Study Date *	Study *	Govt Code *
1623	06-12-2022	Chest	CHEST PA Upright

[Save] [Close]

After you complete the reconciliation process, you will be able to submit the studies. Only the Service Provider(admin) and the Customer are allowed to submit studies.


To submit the study, click on the **Submit** button.

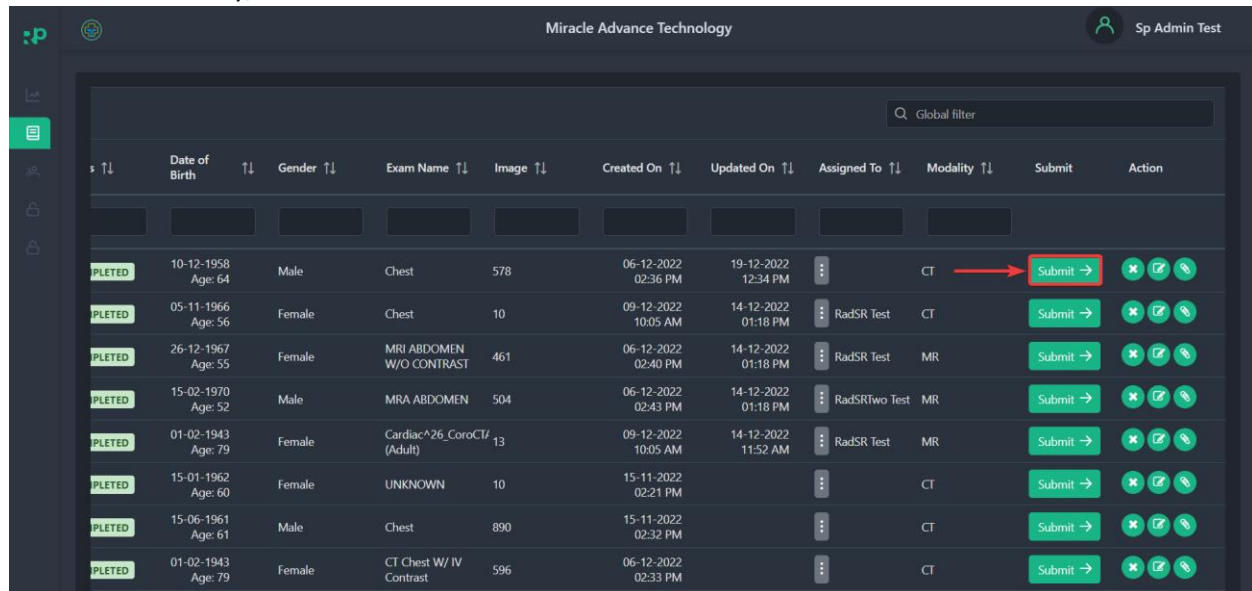

The screenshot shows a web application interface for 'Miracle Advance Technology'. At the top right, there is a user profile icon and the text 'Sp Admin Test'. Below the header is a table with columns: 's', 'Date of Birth', 'Gender', 'Exam Name', 'Image', 'Created On', 'Updated On', 'Assigned To', 'Modality', 'Submit', and 'Action'. The table contains several rows of study data. The first row is highlighted, and a red arrow points to the 'Submit' button in the 'Submit' column. The 'Action' column contains icons for edit, delete, and other actions.

s	Date of Birth	Gender	Exam Name	Image	Created On	Updated On	Assigned To	Modality	Submit	Action
PLETED	10-12-1958 Age: 64	Male	Chest	578	06-12-2022 02:36 PM	19-12-2022 12:34 PM		CT	Submit →	✕ ↻
PLETED	05-11-1966 Age: 56	Female	Chest	10	09-12-2022 10:05 AM	14-12-2022 01:18 PM	RadSR Test	CT	Submit →	✕ ↻
PLETED	26-12-1967 Age: 55	Female	MRI ABDOMEN W/O CONTRAST	461	06-12-2022 02:40 PM	14-12-2022 01:18 PM	RadSR Test	MR	Submit →	✕ ↻
PLETED	15-02-1970 Age: 52	Male	MRA ABDOMEN	504	06-12-2022 02:43 PM	14-12-2022 01:18 PM	RadSRTwo Test	MR	Submit →	✕ ↻
PLETED	01-02-1943 Age: 79	Female	Cardiac^26_CoroCTA (Adult)	13	09-12-2022 10:05 AM	14-12-2022 11:52 AM	RadSR Test	MR	Submit →	✕ ↻
PLETED	15-01-1962 Age: 60	Female	UNKNOWN	10	15-11-2022 02:21 PM			CT	Submit →	✕ ↻
PLETED	15-06-1961 Age: 61	Male	Chest	890	15-11-2022 02:32 PM			CT	Submit →	✕ ↻
PLETED	01-02-1943 Age: 79	Female	CT Chest W/ IV Contrast	596	06-12-2022 02:33 PM			CT	Submit →	✕ ↻

After clicking on Submit button, you should see the same popup window where now you should see the **Save & Submit** button. To submit, click on it.

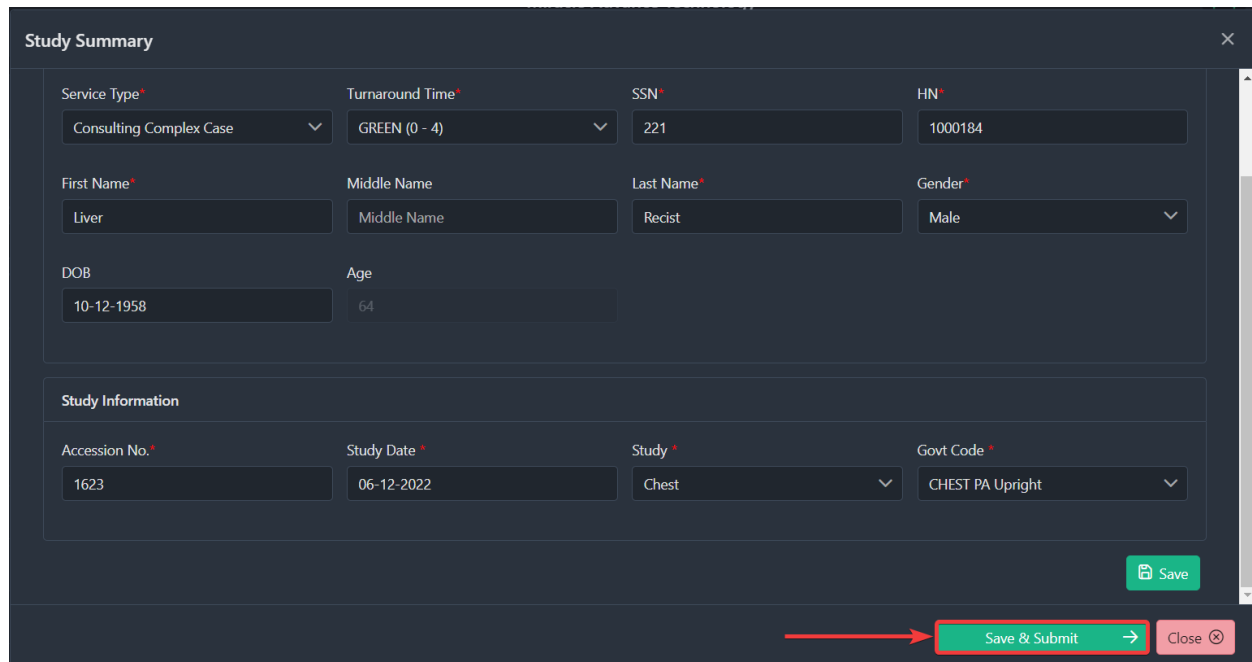

The screenshot shows a 'Study Summary' popup window. It contains several sections with form fields. The first section has 'Service Type' (Consulting Complex Case), 'Turnaround Time' (GREEN (0 - 4)), 'SSN' (221), and 'HN' (1000184). The second section has 'First Name' (Liver), 'Middle Name' (Middle Name), 'Last Name' (Recist), and 'Gender' (Male). The third section has 'DOB' (10-12-1958) and 'Age' (64). The fourth section, 'Study Information', has 'Accession No.' (1623), 'Study Date' (06-12-2022), 'Study' (Chest), and 'Govt Code' (CHEST PA Upright). At the bottom right, there is a 'Save' button and a 'Save & Submit' button, which is highlighted with a red arrow.

4.7 How to Assign a Radiologist to a Case/Study

Any Submitted studies can be assigned to a radiologist. You can easily assign available radiologists to a particular case/study. To assign a radiologist, you have to click on the **Kebab Menu Icon** under the **Assign To** section.

12 columns selected

Global filter

Accession#	HN	Patient Name	Status	Date of Birth	Gender	Exam Name	Image	Created On	Updated On	Assign To	Modality	Submit	Action
1212	2000098	Kevin Jeffers	Submitted	29-05-1956 Age: 66	Male	UNKNOWN	2313	23-11-2022 06:01 AM	23-11-2022 04:42 PM		CT	Submit →	
997	2000037	Summer Inclendon	Finalized	30-03-1952 Age: 70	Female	MRI ABDOMEN WITH AND W/O CONTRAST	849	10-11-2022 10:44 AM	17-11-2022 08:19 PM	RadSR Test	MR	Submit →	
1214	2000030	Barbara Ellis	Submitted	02-10-1952 Age: 70	Female	SENSE Breast-7, dynTHRIVE, diffusion	1080	10-11-2022 05:21 AM	15-11-2022 05:03 PM	RadSR Three Test	MR	Submit →	
206	206206	Carotids Vincent	Finalized	06-11-1966 Age: 56	Female	UNKNOWN	772	09-11-2022 09:25 AM	08-11-2022 03:24 PM	RadSR Test	CT	Submit →	
201	201201	Tavi Vincent	Finalized	05-11-1966 Age: 56	Female	UNKNOWN	1314	08-11-2022 06:39 AM	07-11-2022 04:58 PM	RadSRTwo Test	CT	Submit →	

1 10

After clicking, you should see a popup window.

Assignment

Study Metrics

Site : Hospital 1	Study : UNKNOWN	Accession# : 1212
Patient Name : Kevin Jeffers (2000098)	Age : Number of Image : 2313	Gender : Male
Assigned To :	Assigned By :	Assigned On :

Available Radiologists

4 columns selected

Global filter

Radiologist Name	Sub Specialty	Completed	Current Workload	Action
RadSR Test		0	3	
RadSRTwo Test		0	2	
RadSR Three Test		0	0	

1 10

There you will see the list of available radiologists.

The screenshot shows a dark-themed 'Assignment' window. At the top, 'Study Metrics' are displayed in a grid:

Site : Hospital 1	Study : UNKNOWN	Accession# : 1212
Patient Name : Kevin Jeffers (2000098)	Age : Number of Image : 2313	Gender : Male
Assigned To :	Assigned By :	Assigned On :

Below this is the 'Available Radiologists' section, which is highlighted with a red border. It contains a table with the following data:


Radiologist Name ↑↓	Sub Specialty ↑↓	Completed ↑↓	Current Workload ↑↓	Action
RadSR Test		0	3	
RadSRTwo Test		0	2	
RadSR Three Test		0	0	

At the bottom of the table is a pagination bar showing '1' of 10 items.

To assign a radiologist, click on the **Plus Icon**.

This screenshot is identical to the previous one, but with a red arrow pointing to the 'Plus icon' in the 'Action' column for the 'RadSR Test' row, indicating the click action to assign a radiologist.

After clicking, you should see another popup window asking for your confirmation. Click on **Yes** to confirm or **no** to return to the Assignment page.

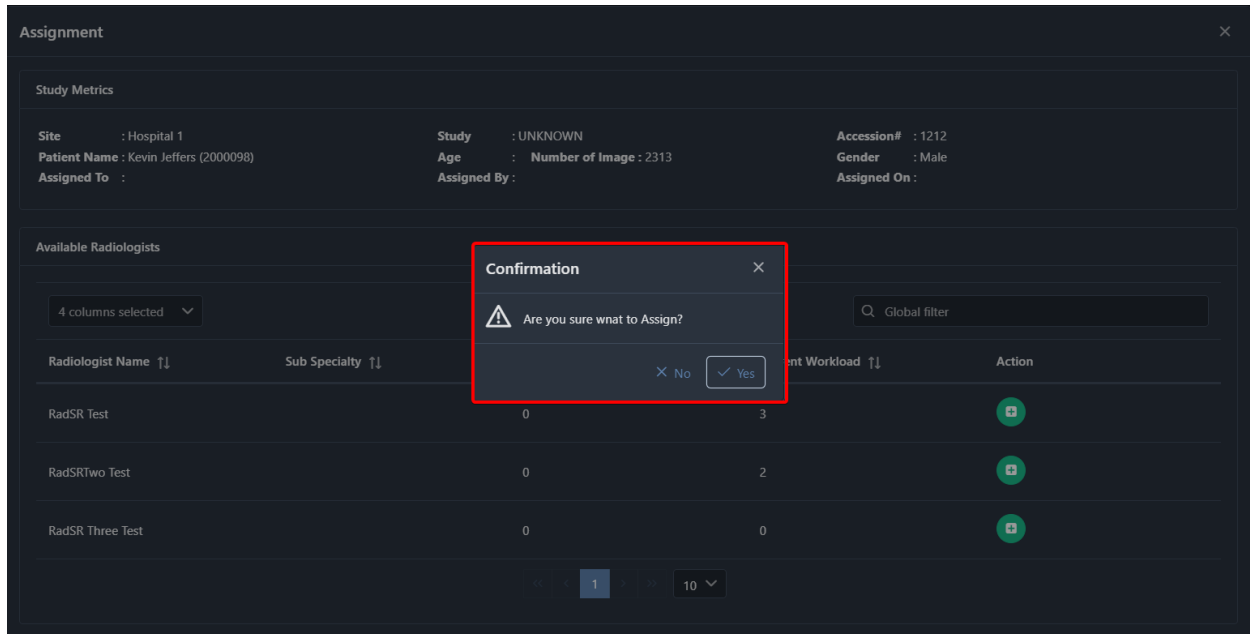

5. Registration Menu

This menu in the RadPortal system allows you to manually register new studies with patient information. There are 2 sections you have to fill in before you try to register a new patient.

You have to fill in

- Patient Information
- Study Information

RadPortal Miracle Advance Technology Sp Admin Test

Study Registration

Patient Information

First Name: Test Last Name: Patient

DOB: 01-01-2020 Gender: Male

Age: 23 HN: HN ☒ Auto Generate

SSN: SSN ☒ Auto Generate

Study Information

Accession No: Accession# ☒ Auto Generate Selected Group: Select Hospital

Queue#: 1 Select Study:

Patient Information:

To register a new Case/Study, you have to submit information such as the Patient's **First Name [1]**, **last Name [2]**, **DOB (Date of Birth) [3]**, **Gender [4]**, **Age [5]**, **HN (Hospital Number) [6]**, and **SSN (Social Security Number) [7]**.

Study Registration

Patient Information

1 → First Name: Test 2 → Last Name: Patient

3 → DOB: 01-01-2020 4 → Gender: Male

5 → Age: 23 6 → HN: HN ☒ Auto Generate

7 → SSN: SSN ☒ Auto Generate

Study Information

In the next Field, you have to fill in study information such as Accession No. [1], Select Groups [2], Queue Number [3], and Select Studies [4].

The screenshot shows a 'Study Information' form. It contains the following fields and controls:

- Accession No:** A text input field with a placeholder 'Accession#'. A red circle with the number '1' and an arrow points to this field.
- Auto Generate:** A green checkmark icon followed by the text 'Auto Generate'.
- Selected Group:** A dropdown menu with 'Select Hospital' as the selected option. A red circle with the number '2' and an arrow points to this dropdown.
- Queue#:** A text input field with the value '1'. A red circle with the number '3' and an arrow points to this field.
- Select Study:** A dropdown menu with a red circle with the number '4' and an arrow pointing to it.

Finally, click on the **Save [1]** button to save or click on the **Save & Preview [2]** to see a preview after saving.


The screenshot shows a full form titled 'Miracle Advance Technology' with a user 'Sp Admin Test'. It is divided into two main sections:

- Patient Information:** Contains fields for First Name (Test), Last Name (Patient), DOB (01-01-2020), Gender (Male), Age (23), HN (HN), and SSN* (SSN). There are 'Auto Generate' checkboxes for HN and SSN.
- Study Information:** Contains fields for Accession No (Accession#), Selected Group (Select Hospital), Queue# (1), and Select Study*. There is an 'Auto Generate' checkbox for Accession No.

At the bottom right, there are two buttons: 'Save' and 'Save & Preview'. Red circles with numbers '1' and '2' and arrows point to these buttons respectively.


6. User list

This menu in the RadPortal system allows you to review all users in the RadPortal system. You also can Add, Modify or Remove users. The user list menu page should look like the screenshot below.

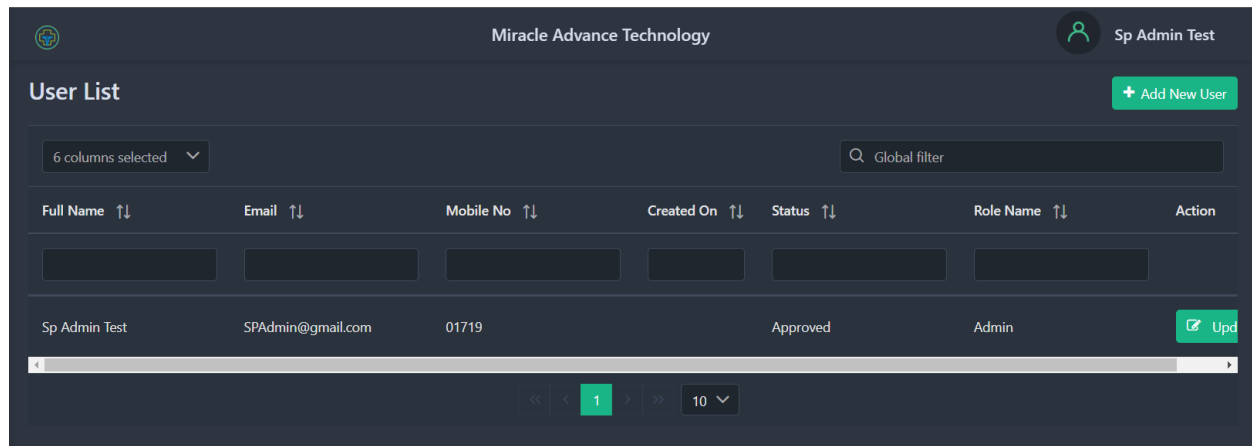

6.1 How to Add New Users

You can easily add users to the RadPortal system. To add a new user, click the **Add New User** button on the top right side of the page.

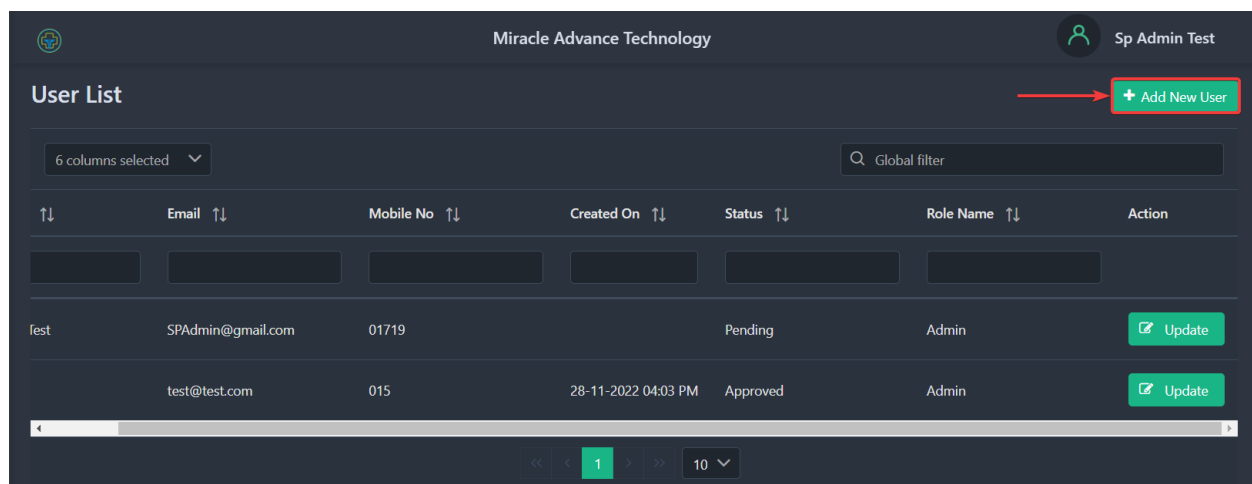

After clicking, you should see a new popup page where you can add information about the new user. You can add information such as First Name [1], Last Name [2], Phone No [3], Email [4], Password [5], Confirm Password [6], and User Role [7].

User Information Entry

1 → First Name (First Name)

2 → Last Name (Last Name)

3 → Phone No (Phone No)

4 → Email (Email)

5 → Password (Password)

6 → Confirm Password (Confirm Password)

7 → User Role (Select Role)

Buttons: Create, Create & Approve, Cancel

In the user role drop-down, you can choose user roles such as **Admin[I]**, **Radiologist [II]**, **Clinician [III]**, and **Help Desk [IV]**.

RadPortal | Miracle Advance Technology | Sp Admin Test

User List | 6 columns selected | Global filter

+ Add New User

User Information Entry

First Name: First Name | Last Name: Last Name

Phone No: Phone No | Email: Email

Password: Password | Confirm Password: Confirm Password

User Role: Select Role

Dropdown Options:

- I → Admin
- II → Radiologist
- III → Clinician
- IV → Help Desk

Buttons: Create, Create & Approve, Cancel

Actions: Update, Update

After you fill in all the information, you can click on the **Create** button to create the new user, and you can click on **Save & Approve** to create and approve the new user and click on **Cancel** to go back to the User list page.

The screenshot shows the 'User List' page in RadPortal. A modal titled 'User Information Entry' is open, allowing the creation of a new user. The fields are as follows:

Field	Value
First Name	Test
Last Name	Admin
Phone No	01
Email	test@test1.com
Password	*****
Confirm Password	*****
User Role	Admin

At the bottom of the modal, there are three buttons: 'Create' (labeled 1), 'Create & Approve' (labeled 2), and 'Cancel' (labeled 3).

If you have clicked on the create button, you have to approve the new user before they try to log in to the RadPortal system. The new user will not be able to log in if their account is not approved.

6.2 How to Approve and Update User Information

You can approve the users you just created. You can also edit or update your user's information by clicking on the **Update** button.

Full Name	Email	Mobile No	Created On	Status	Role Name	Action
Test Admin	test@test1.com	01	28-11-2022 04:36 PM	Pending		Update
Sp Admin Test	SPAdmin@gmail.com	01719		Pending	Admin	Update
Test Role	test@test.com	015	28-11-2022 04:03 PM	Approved	Admin	Update

After clicking on the Update, you should see a new popup window to edit user information.

Update User

First Name: Test

Last Name: Admin

Phone No: 01

Email: test@test1.com

New password:
visibility_off

Status: Deny

User Role: Admin

Update

To approve a user, you have to click on the Status **Dropdown [1]** and select **Approve [2]** and click on **Update [3]**.

Update User

First Name: Test

Last Name: Admin

Phone No: 01

Email: test@test1.com

New password:
visibility_off

User Role: Admin

Status: Deny

Approve

Deny

Delete

Update


6.3 How to Organize Information of Users

Currently, you can view 8 types of information that you can show in the study worklist page, such as:


1. **Organization ID:** Show's the Organization unique Identity number.
2. **Organization Type:** Show the organization type.
3. **Full Name:** Shows the user's full Name.
4. **Email:** Shows the users login email address

5. **Mobile No:** shows the user's mobile number
6. **Created On:** Show the Date the user was created.
7. **Status:** It shows the status of a particular user role.
8. **Role Name:** Shows the user's role name.


To activate an information type, click on **Choose Columns [1]** and find the information type you want to activate. Click on the **Checkbox [2]** before the information type to activate.

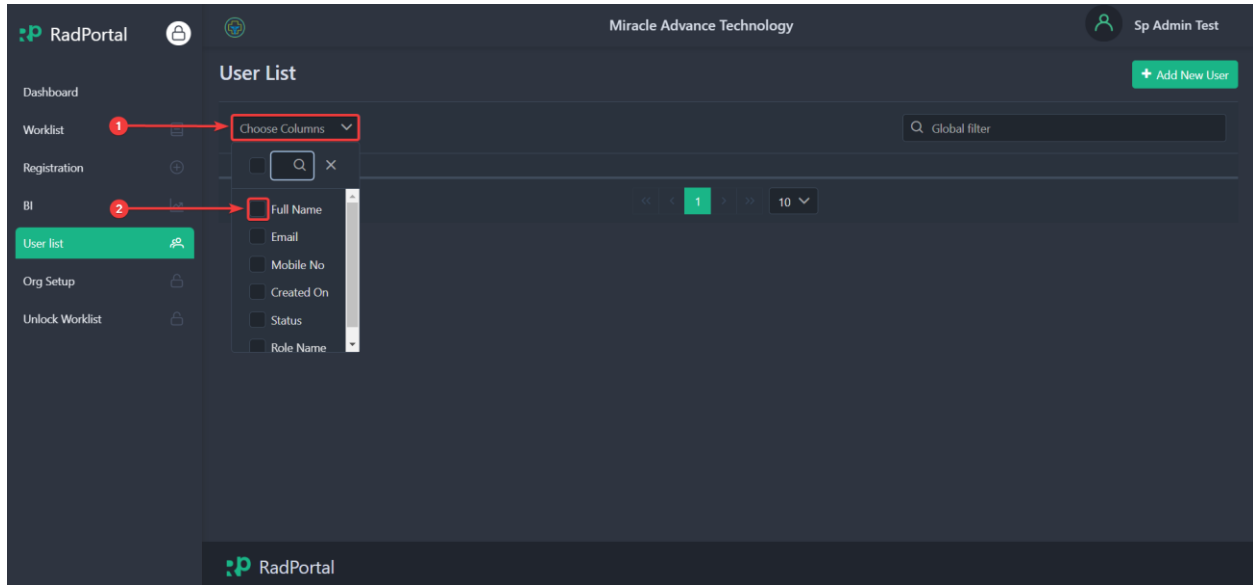


After selecting, you should see that the specific information type is activated.

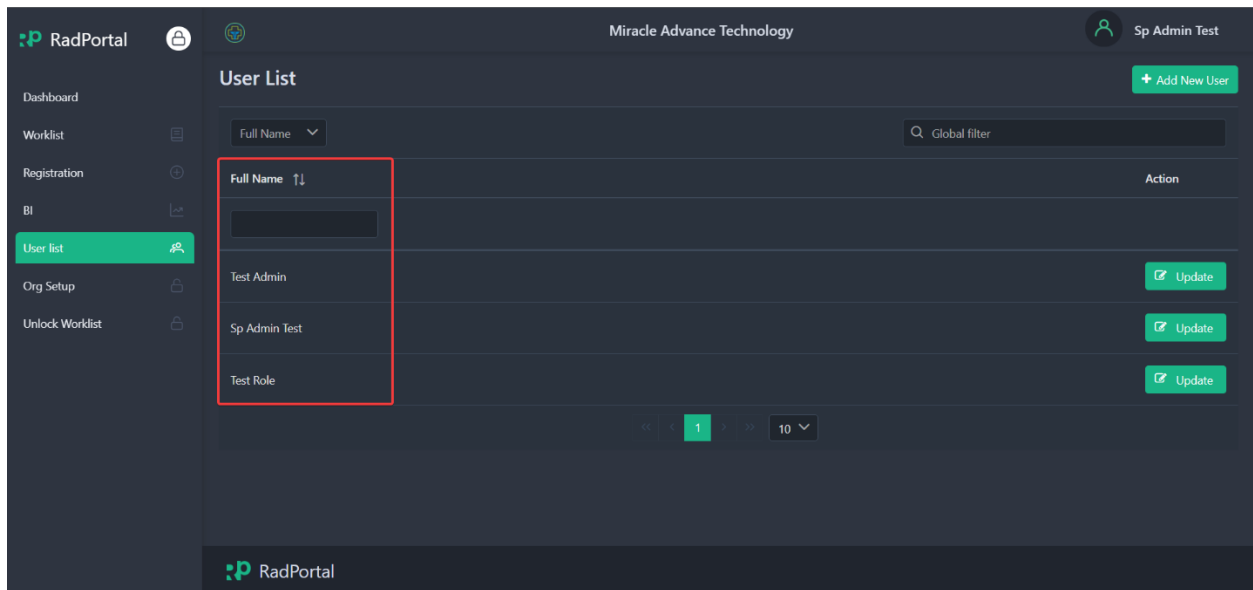

You can activate all information types by clicking on **Choose Columns [1]** and then clicking on **Checkbox [2]** at the top of the dropdown.

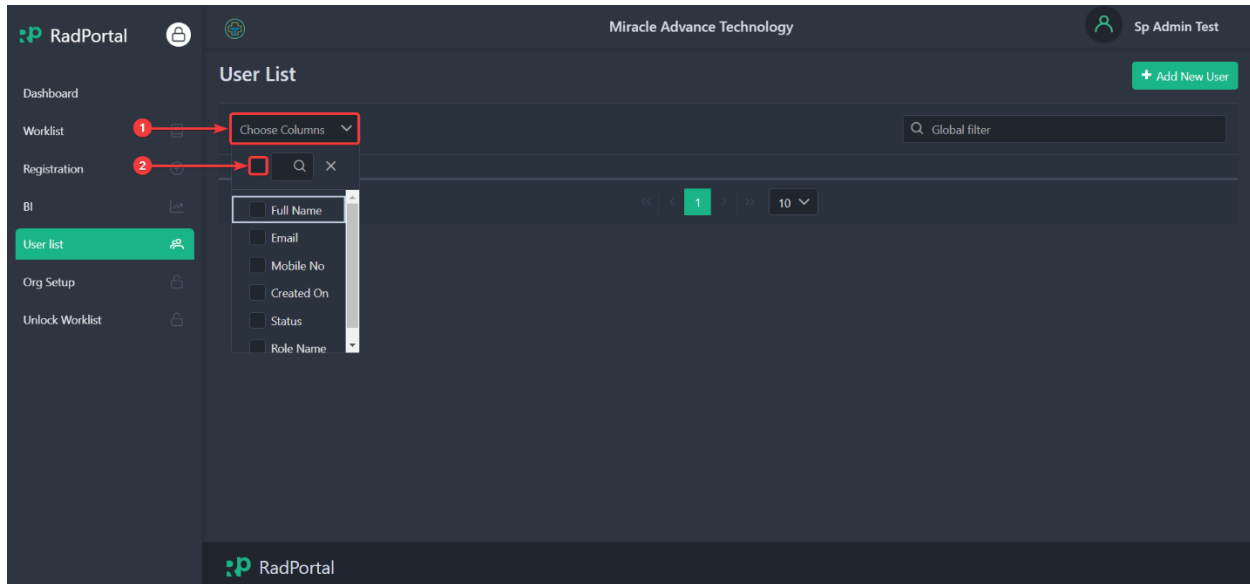


After clicking on the Checkbox, you should see that all information types are activated.

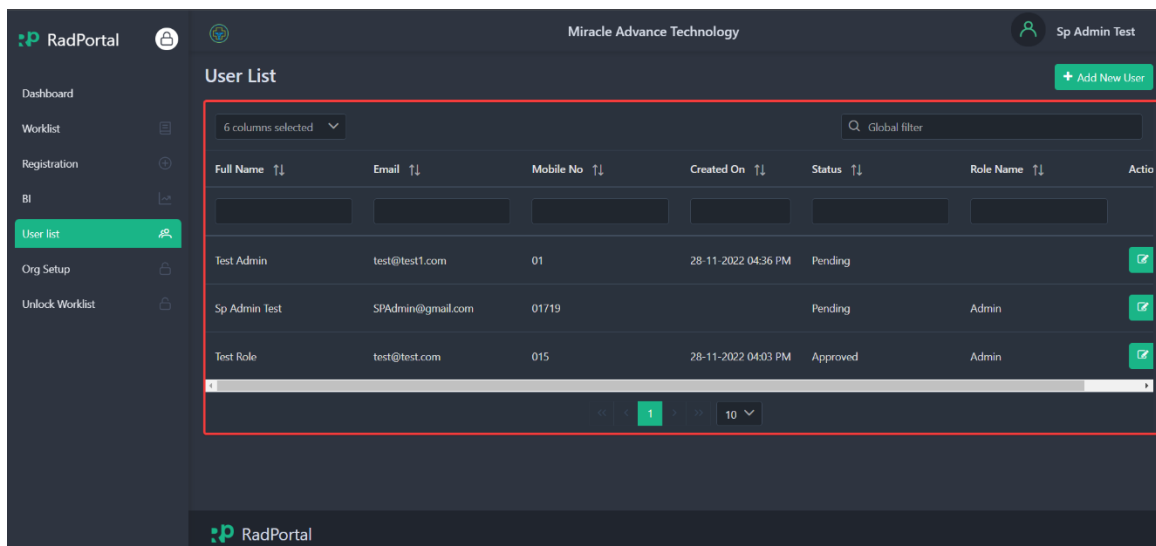

7. Organization Setup

This menu in the RadPortal system allows you to review all organizations configured in the RadPortal system. You also can Add, Modify or Remove organizations. The Organization Setup page should look like the screenshot below.

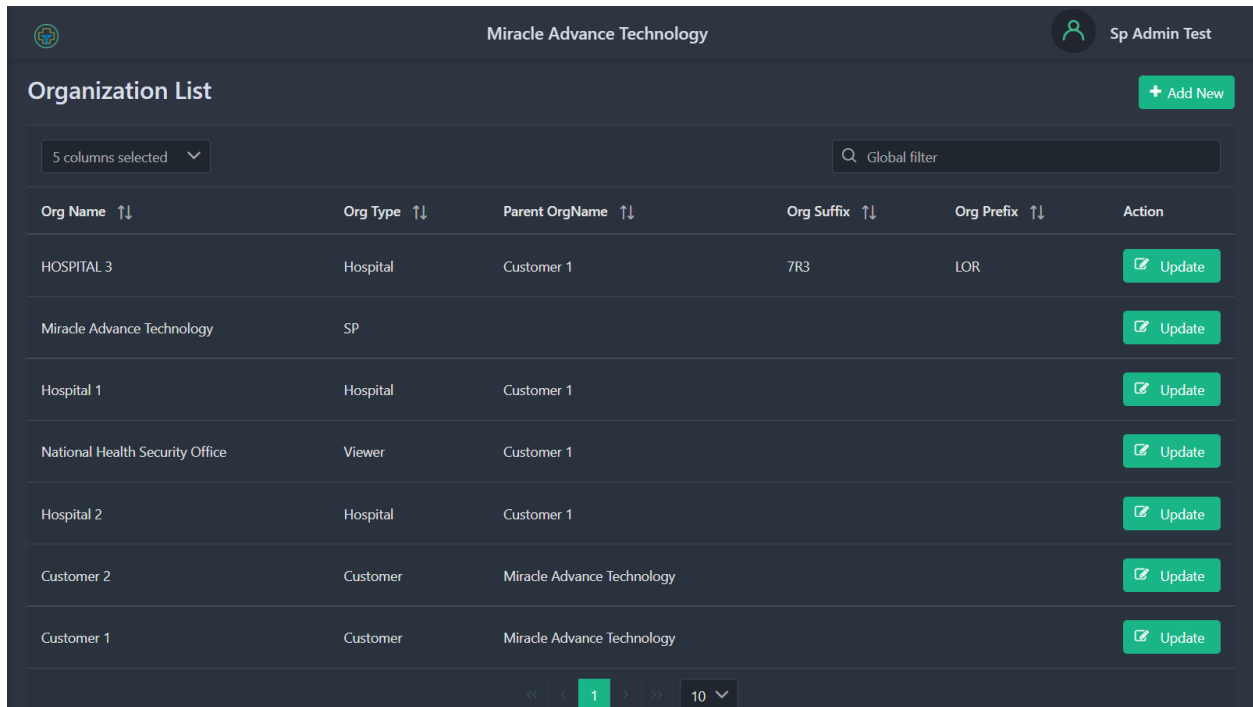

Org Name ↑↓	Org Type ↑↓	Parent OrgName ↑↓	Org Suffix ↑↓	Org Prefix ↑↓	Action
HOSPITAL 3	Hospital	Customer 1	7R3	LOR	Update
Miracle Advance Technology	SP				Update
Hospital 1	Hospital	Customer 1			Update
National Health Security Office	Viewer	Customer 1			Update
Hospital 2	Hospital	Customer 1			Update
Customer 2	Customer	Miracle Advance Technology			Update
Customer 1	Customer	Miracle Advance Technology			Update

7.1 How to Add New Organization

You can easily add Organization to the RadPortal system. To add a new organization, click the **Add New** button on the top right side of the page.

Organization List + Add New

5 columns selected Global filter

Org Name ↑↓	Org Type ↑↓	Parent OrgName ↑↓	Org Suffix ↑↓	Org Prefix ↑↓	Action
HOSPITAL 3	Hospital	Customer 1	7R3	LOR	Update
Miracle Advance Technology	SP				Update
Hospital 1	Hospital	Customer 1			Update
National Health Security Office	Viewer	Customer 1			Update
Hospital 2	Hospital	Customer 1			Update
Customer 2	Customer	Miracle Advance Technology			Update
Customer 1	Customer	Miracle Advance Technology			Update

<<
<
1
>
>>
10

To add a new organization, you have to fill in 2 sections:

- General Information
- Organization Address

Add Organization Information

General Information

Org Type: Parent Org:

Org Name: Web Url:

Tel. No.: Email:

Line ID: Org Prefix:

Org Suffix:

Organization Address

Country: Province:

District: Sub District:

Zip Code:

General Information:


To register an organization, you have to submit information such as the Organization Type [1], Parent Organization [2], Organization Name [3], Website URL [4], Telephone Number [5], Organization Email [6], Line ID [7], Organization Prefix [8], and Organization Suffix [9].

In the Organization Type role drop-down, you can choose organization types such as **Customer [I]**, **Hospital [II]**, or **Viewer [III]**.

Organization Address

In the following Field, you have to fill in the Organization Address, such as Country Name [1], Province [2], District [3], Sub District [4], and Zip Code [5].

Finally, click on the **Save [1]** button to save or click on **Cancel [2]** to go back to the Organization management page.

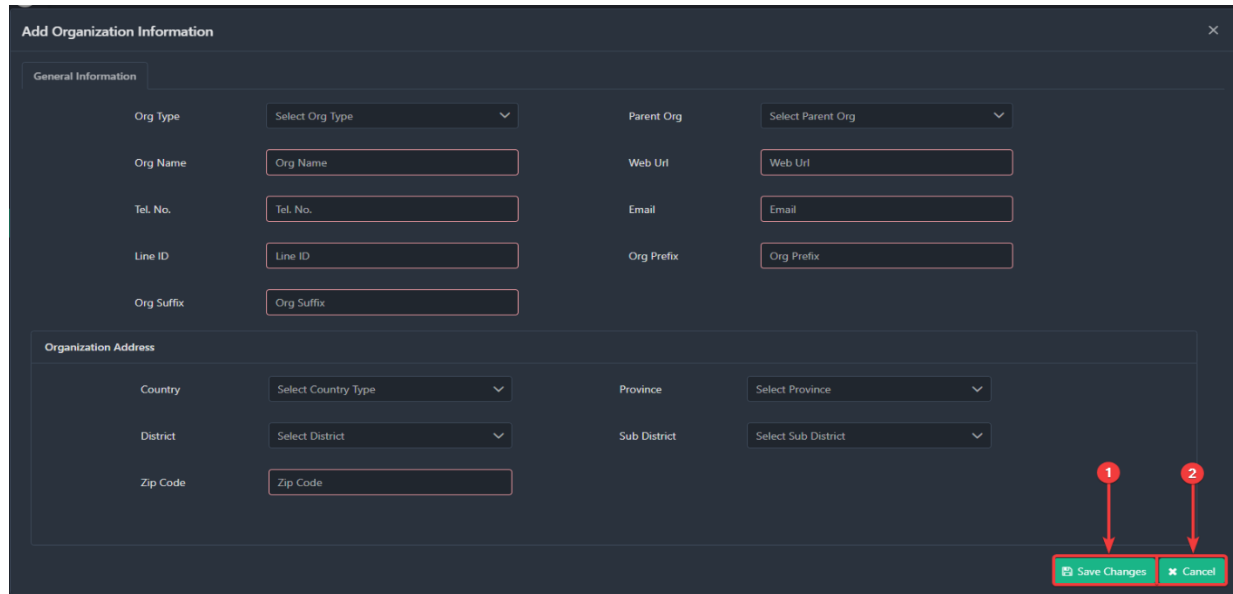

Add Organization Information

General Information

Org Type: Parent Org:

Org Name: Web Url:

Tel. No.: Email:

Line ID: Org Prefix:

Org Suffix:

Organization Address

Country: Province:


District: Sub District:

Zip Code:

1 **2**

7.2 How to Update Organization Information

You can approve the users you just created. You can also edit or update your user's information by clicking on the **Update** button.

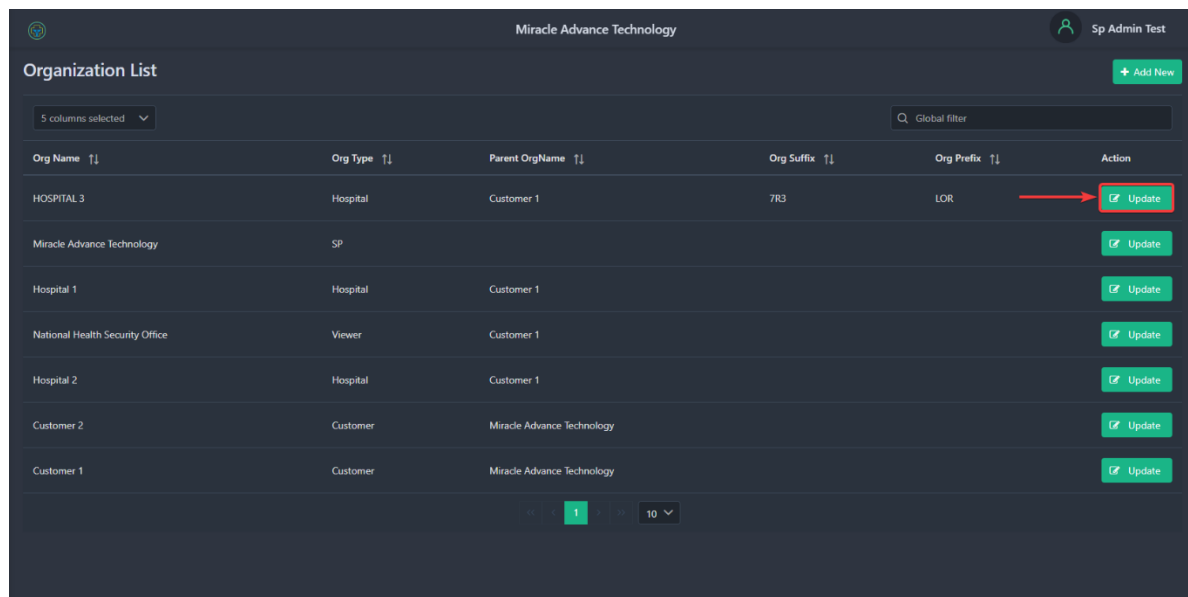

Miracle Advance Technology

Sp Admin Test

Organization List

5 columns selected

Org Name	Org Type	Parent OrgName	Org Suffix	Org Prefix	Action
HOSPITAL 3	Hospital	Customer 1	7R3	LOR	Update
Miracle Advance Technology	SP				Update
Hospital 1	Hospital	Customer 1			Update
National Health Security Office	Viewer	Customer 1			Update
Hospital 2	Hospital	Customer 1			Update
Customer 2	Customer	Miracle Advance Technology			Update
Customer 1	Customer	Miracle Advance Technology			Update

1 10

After clicking on the Update, you should see a new popup window where you can update organization information.

Miracle Advance Technology

Update Organization Information

General Information

Org Type	Customer	Parent Org	Miracle Advance Technology
Org Name	Customer 1	Web Url	Web Url
Tel. No.	Tel. No.	Email	Email
Line ID	Line ID	Org Prefix	Org Prefix
Org Suffix	Org Suffix		

Organization Address

Country	Select Country Type	Province	Select Province
District	Select District	Sub District	Select Sub District
Zip Code	Zip Code		

Save Changes Cancel

RadPortal