



Envision NEXT

USER MANUAL

Contents



3

Signing Up and Signing In 3

 5

Navigation 5

 13

Patient Registration & Appointment 13

Patient Dashboard: 16

Generating ID Card for Patient: 16

Creating an Appointment for the Patient: 17

 21

Whiteboard 21

Whiteboard Features: 22

Search Features: 22

Column Features: 24

Reservation Worklist: 26

Preparation Worklist: 28

Vital Signs 35

Imaging: 36

 49

Nursing 49

 62

Device Management 62

 78

Disinfection Worklist 78

 78

BI 78

 78

Administration 78

Document Version No- 1.1.2


Last update Datetime- 14/03/2022

Last update by-Wasiur Rahman

Applicable Software no-0.8.0


Signing Up and Signing In

Before you can login to Envision Next, you will need to gather your username and password from an administrator. An **administrator will have to register** you from within the system, and provide you with your proper credentials. The **Sign In** page looks like this:-

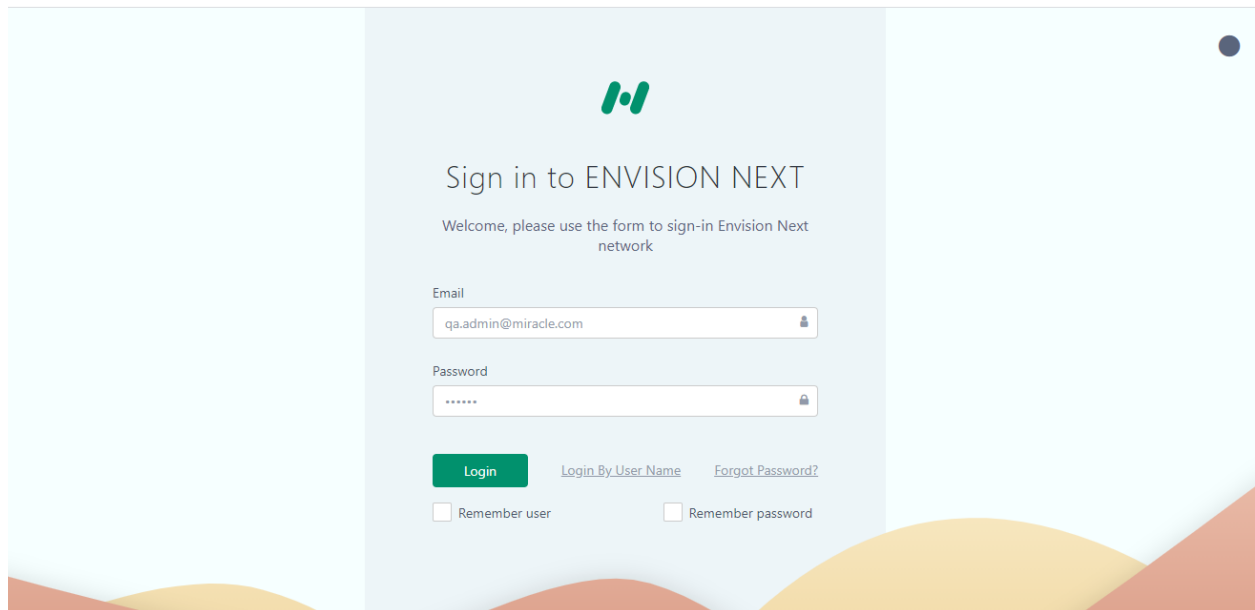

You will have a number of options in the sign in page that you can use to customize or facilitate your next logins.


Sign in Page Theme: You can for example, **click on the dark grey dot on the top right** corner of the login page to change the theme of the login portal to a dark one.

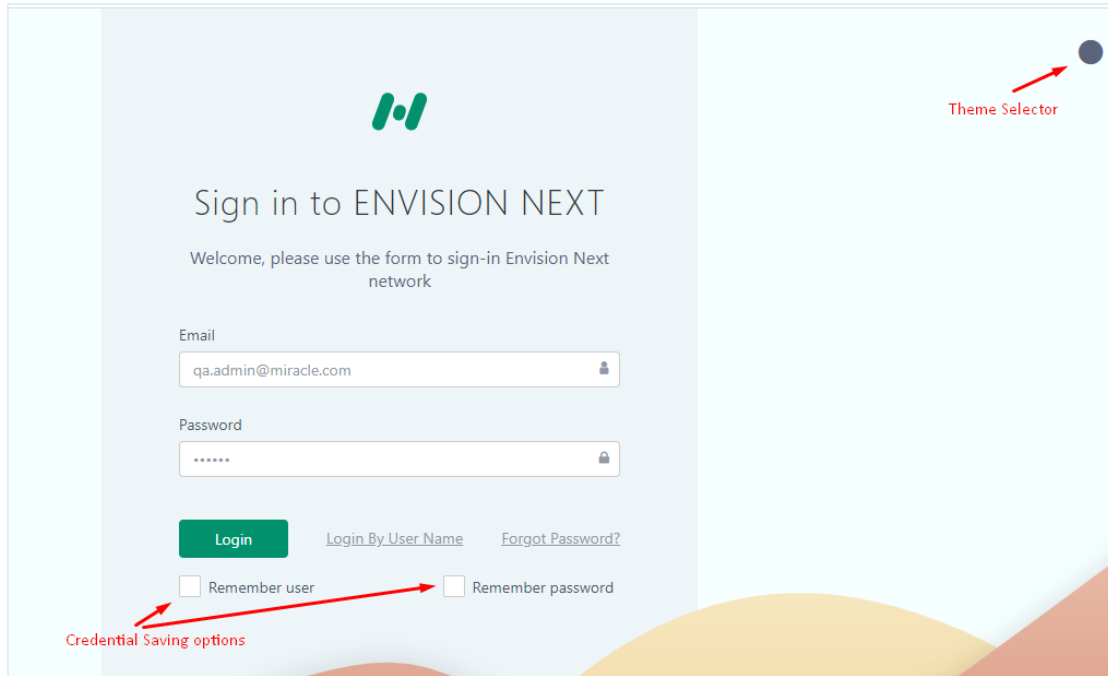


You can also:

- Login with your **username** instead of your email
- **Choose which information (user/pass)** you want the application to remember


In case you're curious, here's what the dark login theme looks like!

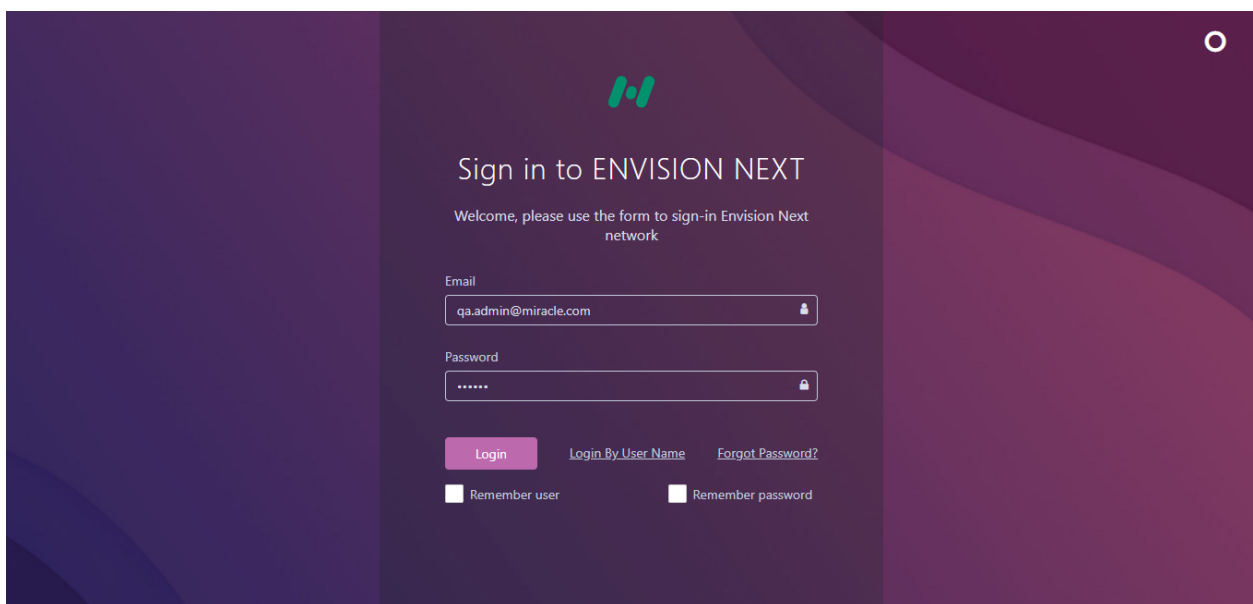

After clicking the Login button, the user will get a **Confirmation** message for replacing another session.


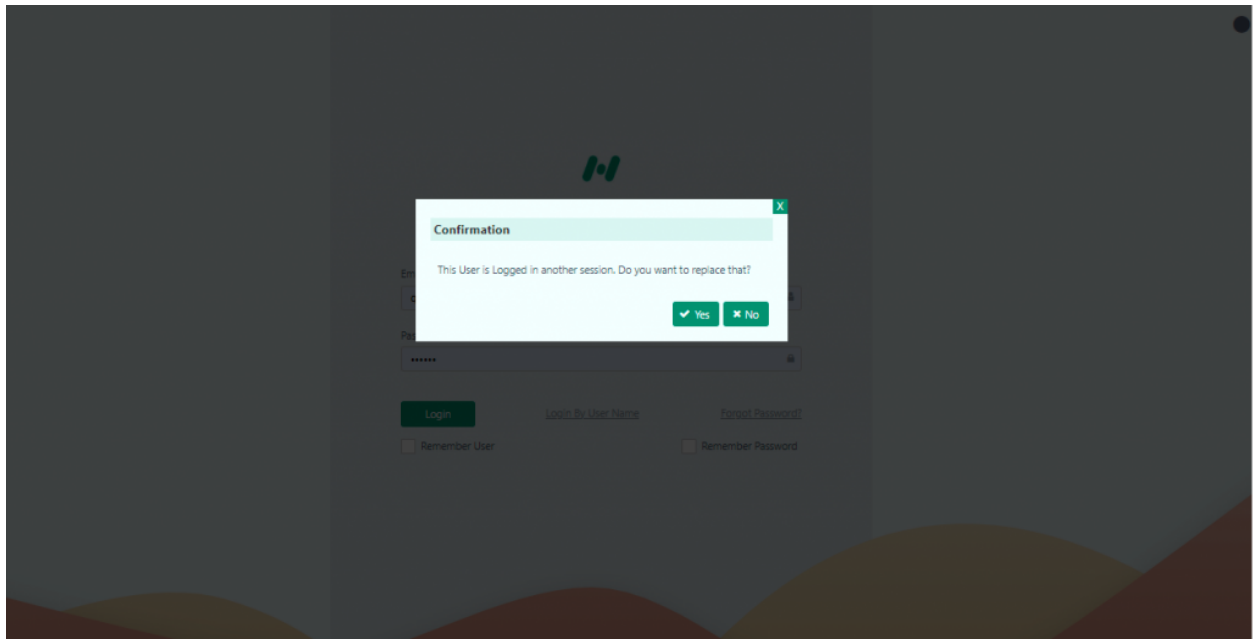


After clicking the 'Yes' button the user will successfully login to the system.

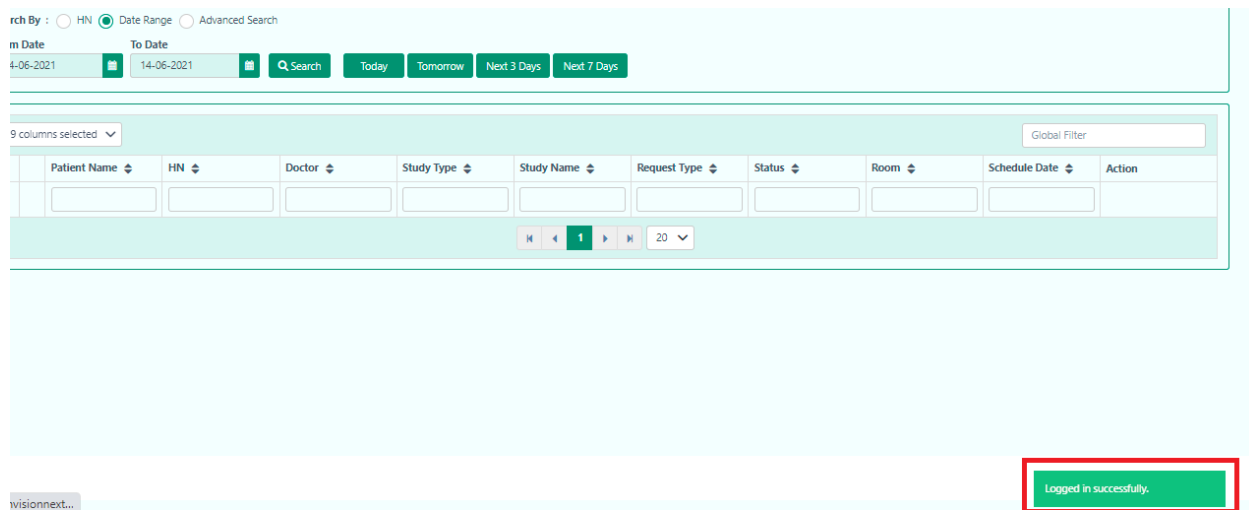

Navigation



Here's what the application will look like once you're logged in. The page you see may be different depending on what you want to see as your homepage.

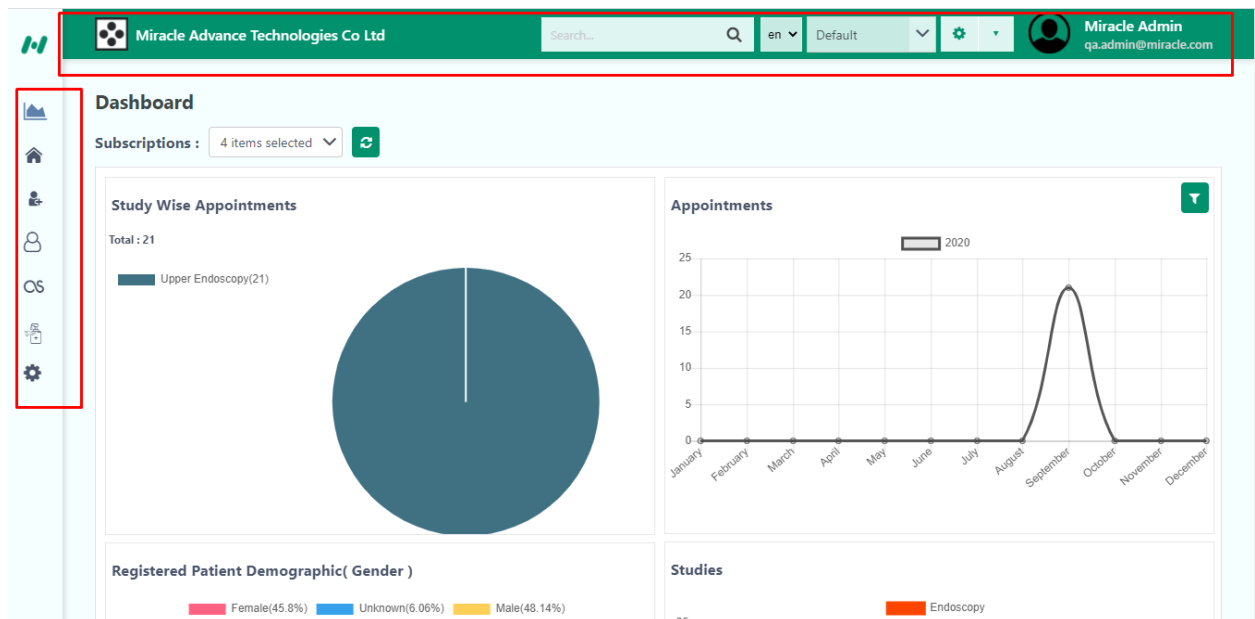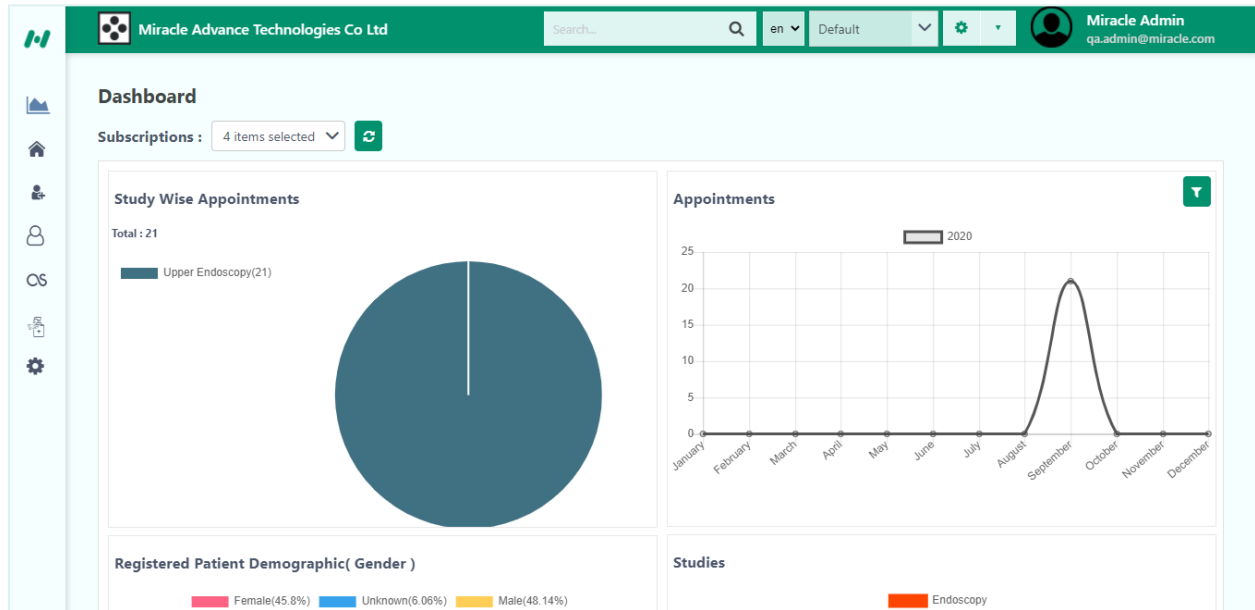


The left-hand **Main Menu** contains all the sections that you will need to use throughout the application.

The **Top Bar/Menu** contains the following features:


- **Search bar:** Here you can **search anything regarding** any patients or appointments, or both at the same time (see the 2nd screenshot below) just by selecting or deselecting the option.

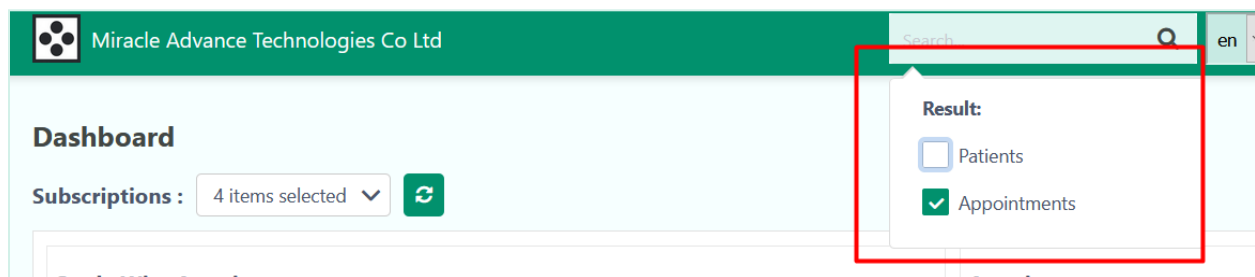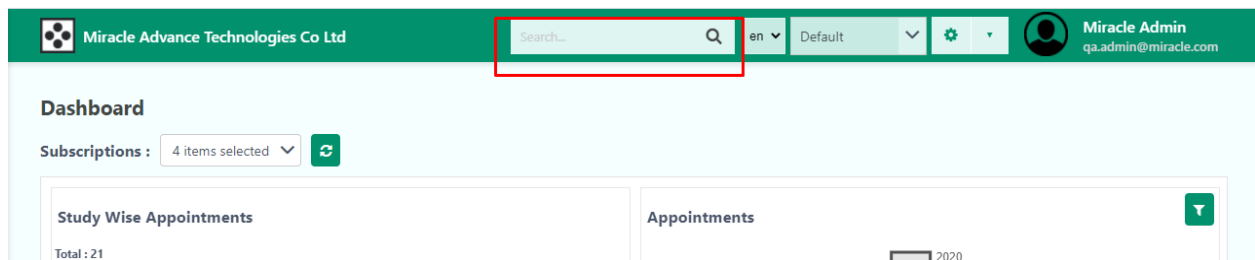

- **Language Selector:** You can select what language you want to use.

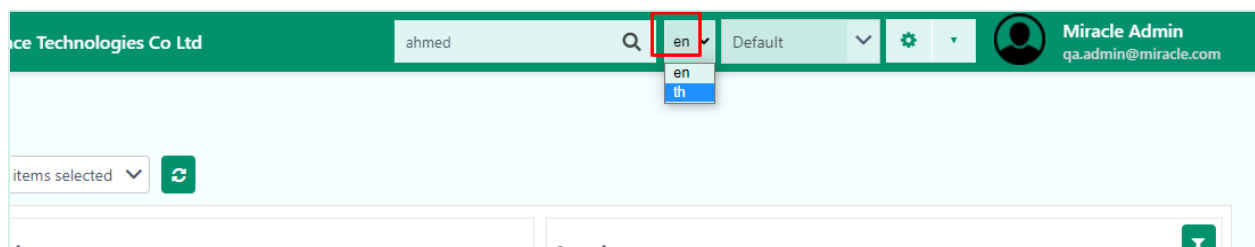


- **Theme Selector:** You can choose what type of background you can use , for instance, ‘ Default’, ‘ Blue’ or ‘Dark’.

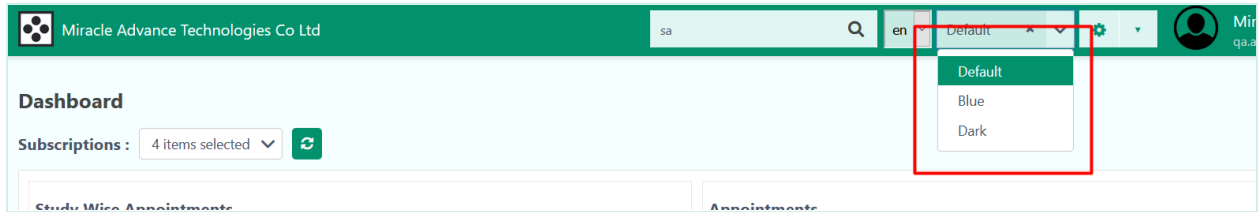

→ Settings

- **Set As Home:** By pressing this, you can set any page as the home page when you log in, so for example if you want to set “Whiteboard” as the home page, you can just go to WhiteBoard and set it as home screen.

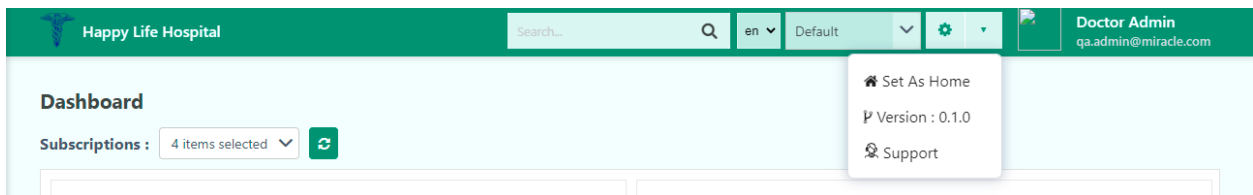


- **Version Number:** This shows what version of the software you are using, what date it was published on and any new features /updates that version brings to the software.
- **Support:** This helps to get any kinds of complaints from the users.

 A screenshot of the 'Support Ticket' form. The form has a light blue header with the title 'Support Ticket'. Below the header, there are three input fields: 'Subject*' (with a placeholder 'Subject'), 'Priority*' (with a placeholder 'Select priority' and a dropdown arrow), and 'Support Type*' (with a placeholder 'Select support type' and a dropdown arrow). Below these fields is a 'Notes*' section with a rich text editor. The editor has a toolbar with various icons for text formatting (bold, italic, underline, link, unlink, list, indent, outdent, undo, redo) and a 'Save' button. At the bottom of the form, there is a 'Choose' button and a 'Save' button. A small text at the bottom left states: 'Available file types are: .jpeg, .jpg, .png, .gif, .mp4, .mpg, .flv, .avi, .mov, .pdf, .doc, .docx, .txt, .xls, .xlsx, .csv, .zip, .rar'.


The user when they find any difficulties using this application, he/she simply uses this support option. The support option provides a form. It contains-

- **Subject:** The user writes the issue in this field when they are really facing the problems.

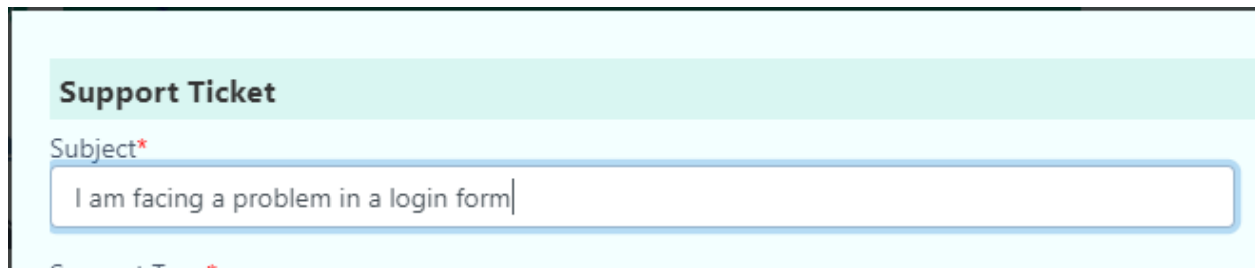

The screenshot shows a form titled "Support Ticket" with a light blue header. Below the header, the label "Subject*" is followed by a text input field containing the text "I am facing a problem in a login form".

- **Priority:** The user can select the priority when they are facing the problems with the help of High, Medium and Normal. According to the priority the developer will try to solve the problems.

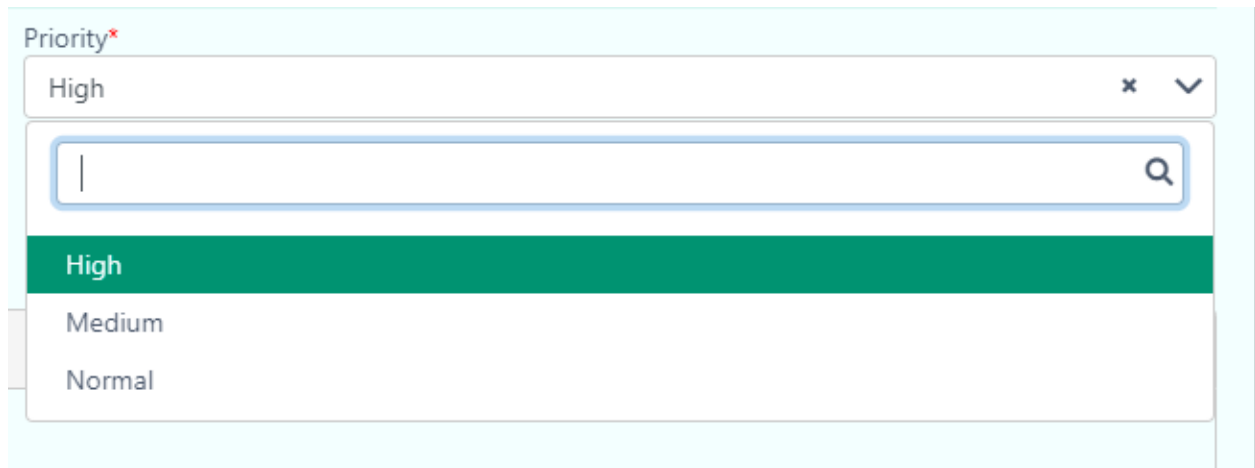

The screenshot shows a dropdown menu for "Priority*". The current selection is "High". Below the dropdown, a search bar is visible. The dropdown list is open, showing three options: "High" (highlighted in green), "Medium", and "Normal".

- **Select support type:** The user can select Error, Bug, Query according to the issues while using this application.

Support Type*

Select support type

|

Error

Bug

Query

- **Notes:** The user can explain details about the issues while using the application through this section. They can include files in any format. Example- .jpeg, .jpg, .png, .gif, .mp4, .mpg, .flv, .avi, .mov, .pdf, .doc and so on.

Notes*

Normal Sans Serif B I U A [icon] [icon] [icon] [icon] [icon] [icon] [icon] [icon] [icon] [icon]

Available file types are: .jpeg, .jpg, .png, .gif, .mp4, .mpg, .flv, .avi, .mov, .pdf, .doc, .docx, .txt, .xls, .xlsx, .csv, .zip, .rar

+ Choose

Save

When all the procedure is over the user has to press the 'Save' button in the below right corner of the Support page.

Support Ticket

Subject* Priority*

Support Type*

Notes*

Normal Sans Serif B I U A [color picker] [background color picker] [list icon] [link icon] [unlink icon] [code icon]

Available file types are: .jpeg, .jpg, .png, .gif, .mp4, .mpg, .flv, .avi, .mov, .pdf, .doc, .docx, .txt, .xls, .xlsx, .csv, .zip, .rar

→ Profile Options:


Search... en Default [gear icon] [user icon]

Miracle Admin
qa.admin@miracle.com

-
-
-

Appointments


- **Profile:** Here you can view your Profile and change your photo, edit user, provide signature and your Password.

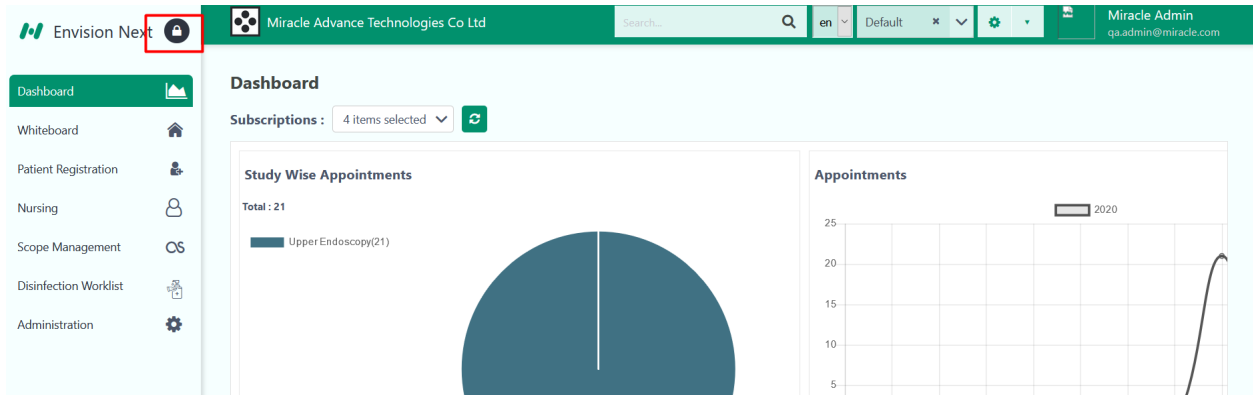

You can lock the main menu, so that it always stays there even if you move the mouse away from it, by **selecting the lock button** on the top right side of the main menu.

Patient Registration & Appointment

From the **Main Menu** on the left, you can start a **Patient**: -

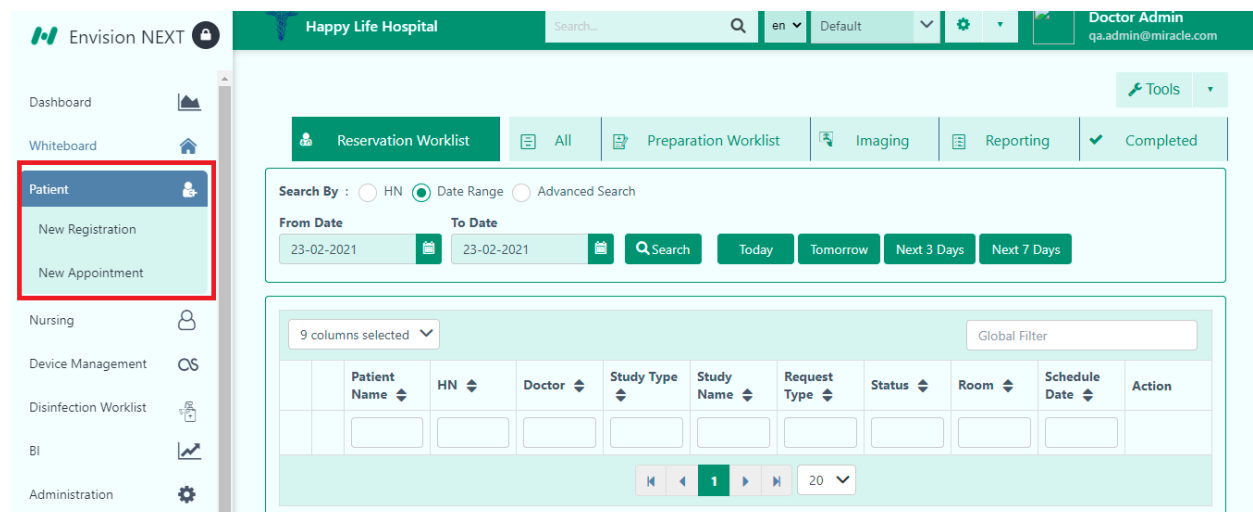

This takes you to sub menus- New Registration and New Appointment.

The **New Registration** form has 3 sections that you can fill up -

- Patient Information
- Additional
- Emergency contact
- Patient Allergy

Note: You will need to fill up the (minimum) mandatory fields in the form before you can successfully complete a patient registration. Mandatory fields to be filled are marked with red asterisks (*)

The screenshot shows a 'Patient Information' form with the following fields and options:

- HN:** Input field with an 'Auto Generate' checkbox.
- Title:** Dropdown menu with 'Select Title'.
- First Name*:** Input field (mandatory).
- DOB *:** Input field with a calendar icon (mandatory).
- Phone:** Input field.
- Gender:** Dropdown menu with 'Select Gender'.
- Country:** Dropdown menu with 'Thailand' selected.
- Address:** Text area.
- Known Disease:** Dropdown menu with 'Choose'.
- Special Instruction:** Dropdown menu with 'Choose'.
- Is VIP:** Checkbox.
- Auto Generate:** Checked checkbox.
- Profile Picture:** Section with a person icon.
- SSN:** Input field.
- Last Name*:** Input field (mandatory).
- Age:** Input fields for 'Years', 'Months', and 'Days'.
- Email:** Input field.
- Blood Group:** Dropdown menu with 'Select'.
- Province:** Dropdown menu with 'กรุงเทพมหานคร' selected.
- Zip Code:** Input field.
- Known Disorder:** Dropdown menu with 'Choose'.
- Patient Type:** Dropdown menu with 'Select'.

Mandatory fields are highlighted in red*

Note that you can add multiple options for fields like 'Known Disease' and 'Known Disorder'. You can search for the diseases at your convenience.

Miracle Advance Technologies Co Ltd

Search...

en

Default

Miracle Admin
qa.admin@miracle

Country: Thailand

Province: กรุงเทพมหานคร

Address:

Zip Code:

Known Disease: Lung disease, Bronchial asthma

Known Disorder: Vision impairment flag, Hearing impairment flag, Language impairment flag

Special Instruction:

Emergency Contact:

Title: Select Title

Relation:

Patient Type: Select

Search diseases here

Liver dysfunction

Renal dysfunction

Bronchial asthma

Lung disease

Heart disease

Vision impairment flag, Hearing impairment flag, Language impairment flag.

Note: Is VIP option protects the user data. Only the admin can read the data.

Phone *

Email:

Gender: Select Gender

Blood Group: Select

Country: Thailand

Province: กรุงเทพมหานคร

Address:

Zip Code:

Known Disease: Bronchial asthma, Lung disease

Known Disorder: Choose

Special Instruction: Choose

Patient Type: Select

Is VIP: ☒

Bronchial asthma, Lung disease.

Emergency Contact

Title

Mr. x v

First Name*

Edmund

Phone *

7867786

Address

Relation

Frenemy

Last Name*

Tully

Email

Patient Allergy

Allergy Name

Food Allergy x v

Comments

Nut allergy

+ ADD

Allergy Name	Comments	
Dust Allergy	Allergic to deserts	

Back

Save

Mandatory fields are highlighted in red*

Note: When adding allergies, please note that you will have to click on the ‘ADD’ button to be able to list the allergies successfully.

Once you’ve filled all the sections, you can click on **Save**. You will now be taken to the **Patient Dashboard**: -

Patient Dashboard

EMR

Appointment

Tools

Back

HN

:

P2020091611956

Name

:

Mr. Kencho Chopel

DOB

:

02-03-1966

Phone

:

985454665

Gender

:

Male

Address

:

Registration Date

:

16-09-2020

SSN

:

57867TN-5657

Age

:

54 Years 7 Months 1 Days

Email

:

kencho@chopel.com

Blood Group

:

O-

Zip Code

:

Edit

History

Pathology

Invoice

Appointment

Allergy


History

Patient Dashboard:

The Patient Dashboard is where you can access information about the patient, edit the information and view useful details such as Invoices and History.

Generating ID Card for Patient:

From Patient Dashboard, go to Tools on the top right corner:

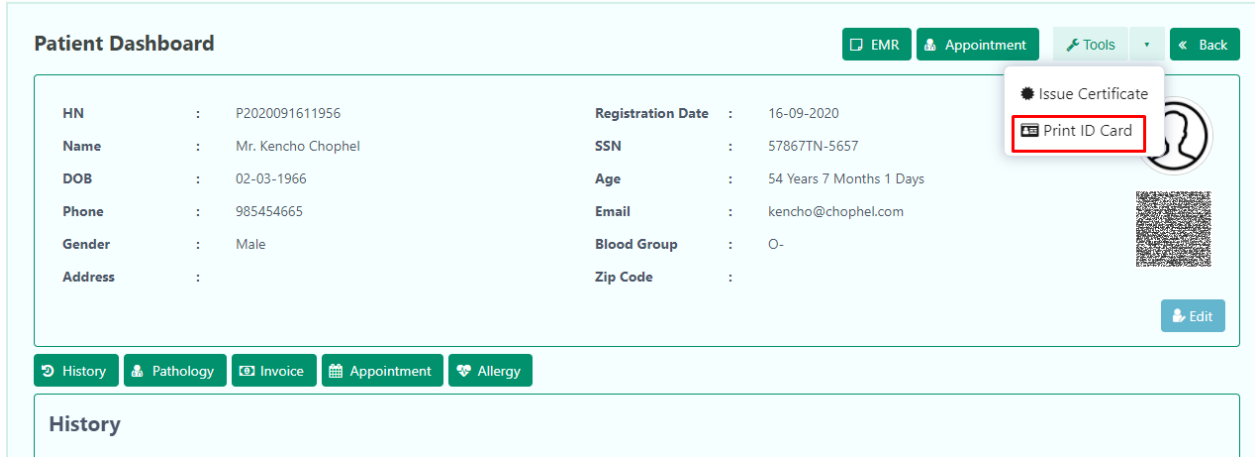


The screenshot shows the 'Patient Dashboard' interface. At the top right, there are buttons for 'EMR', 'Appointment', 'Tools', and 'Back'. The 'Tools' button is highlighted with a red box, and a dropdown menu is open showing 'Issue Certificate' and 'Print ID Card'. The 'Print ID Card' option is also highlighted with a red box. Below the buttons, there is a patient information card for 'Mr. Kencho Chopel' with details like HN, Name, DOB, Phone, Gender, Address, Registration Date, SSN, Age, Email, Blood Group, and Zip Code. A QR code and an 'Edit' button are also visible. At the bottom, there are buttons for 'History', 'Pathology', 'Invoice', 'Appointment', and 'Allergy'.

HN	:	P2020091611956	Registration Date	:	16-09-2020
Name	:	Mr. Kencho Chopel	SSN	:	57867TN-5657
DOB	:	02-03-1966	Age	:	54 Years 7 Months 1 Days
Phone	:	985454665	Email	:	kencho@chopel.com
Gender	:	Male	Blood Group	:	O-
Address	:		Zip Code	:	

An ID card print is generated!

Creating an Appointment for the Patient:

Right before **Tools** is the **Appointment** button. It takes you to the Patient Appointment.

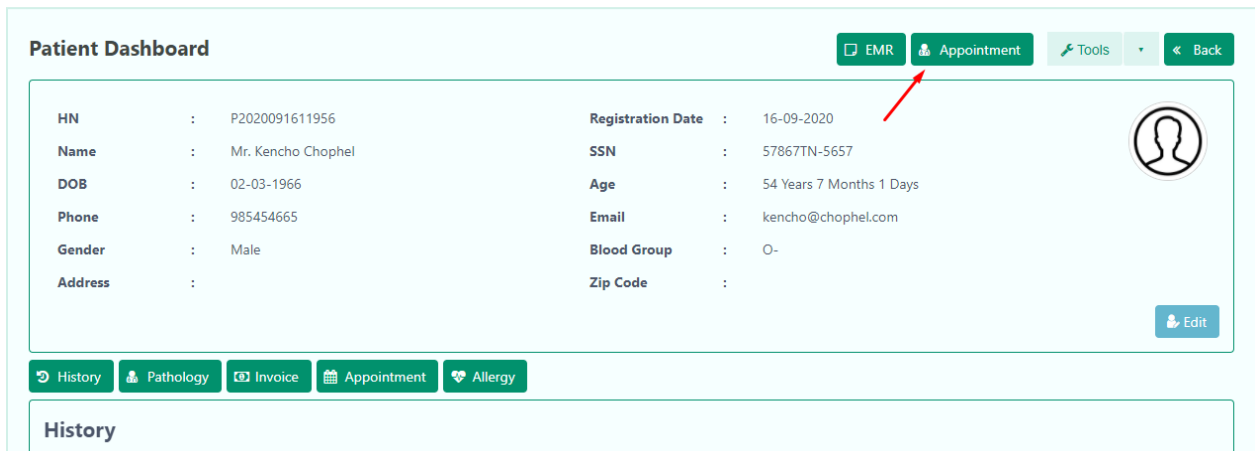

The screenshot shows the 'Patient Dashboard' interface. At the top right, there are buttons for 'EMR', 'Appointment', 'Tools', and 'Back'. A red arrow points to the 'Appointment' button. The patient information card for 'Mr. Kencho Chopel' is visible below the buttons. At the bottom, there are buttons for 'History', 'Pathology', 'Invoice', 'Appointment', and 'Allergy'.

HN	:	P2020091611956	Registration Date	:	16-09-2020
Name	:	Mr. Kencho Chopel	SSN	:	57867TN-5657
DOB	:	02-03-1966	Age	:	54 Years 7 Months 1 Days
Phone	:	985454665	Email	:	kencho@chopel.com
Gender	:	Male	Blood Group	:	O-
Address	:		Zip Code	:	

Patient Appointment ← Back

Clinic Type: Elective × ▼
 Doctor: Select doctor ▼
 Room: Choose ▼
 Study Type: Select exam type ▼
 Study: Select procedure Type ▼

Study Dependency
Study Dependency
No Dependency
Arrived

June 15, 2021 day Week Month

	Room 101	Room 102	Room 105	Room 106	Room 107	Test Room	Test Room 2
09:00							
10:00							
11:00							
12:00							
13:00							
14:00							
15:00							
16:00	16:00 - 17:00 EGD/ERCP			16:00 - 17:00 - EGD			
17:00	17:00 - 18:00 EGD/Colonoscopy						
18:00			18:00 - 19:00 - EGD				
19:00							

The Appointment calendar is where you can select a slot [example: **red area**] (by clicking on it) and set an appointment for the patient from the Patient Dashboard that you just came from. You can select an appointment Room Wise.

67ba8ed6a1ef.ngrok.io/#/dashboard/patient-appointment-by-room/20849 Doctor Admin

Patient Appointment ← Back

Clinic Type: Elective × ▼
 Doctor: Select doctor ▼
 Room: Room 101 ▼
 Study Type: Select exam type ▼
 Study: Select procedure Type ▼

Study Dependency
Study Dependency
No Dependency
Arrived

June 15, 2021 day Week Month

Room 101

	Room 101	Room 102	Room 105	Room 106	Room 107
09:00					
10:00					
11:00					
12:00					
13:00					
14:00					
15:00					

Clicking a slot takes you to the Appointment card where you can input details-

Patient Appointment

Close Save

Patient Information

Name : Mr. Kencho Chophel

Phone : 985454665

Address :

Doctor Information

Name :

Phone :

Department :

> Study

> Additional

Close Save

Note- When you click **Mark As Arrived** it will take you to the **Preparation Worklist**. Otherwise, the patient will go to the **Reservation Worklist**.

Fill in the rest of the information in the two tabs.

Please note that you will not be able to Save the information until you have filled out all the required fields.

The fields are highlighted in the screenshots below-

▼ Study

Date* 15-06-2021

Duration (in Minutes)* 30

Doctor* Select doctor

Study Type* Bronchoscopy

Clinic Type* Elective

Chief Complaint Choose

Start Time* 22:00

End Time* 22:30

Room* Room 107

Study* Select procedure

Request Type Routine

Diagnosis

Mandatory fields are highlighted in red*

Note- You can also search patients in the search bar.

▼ Additional

Department

Internal Medicine x ▼

Floor

Report

Please select report type ▼

Consent

Please Select ▼

Nursing Class

Please select nursing class ▼

Indication

Choose ▼

B I U

Patient Status

Please select patient status ▼

Phone

Finding

Choose ▼

Walking

Please select walking ▼


Purpose

Choose ▼



Comment

B I U


Once you have saved the information, you will be taken back to the calendar to show the appointment that was being made.

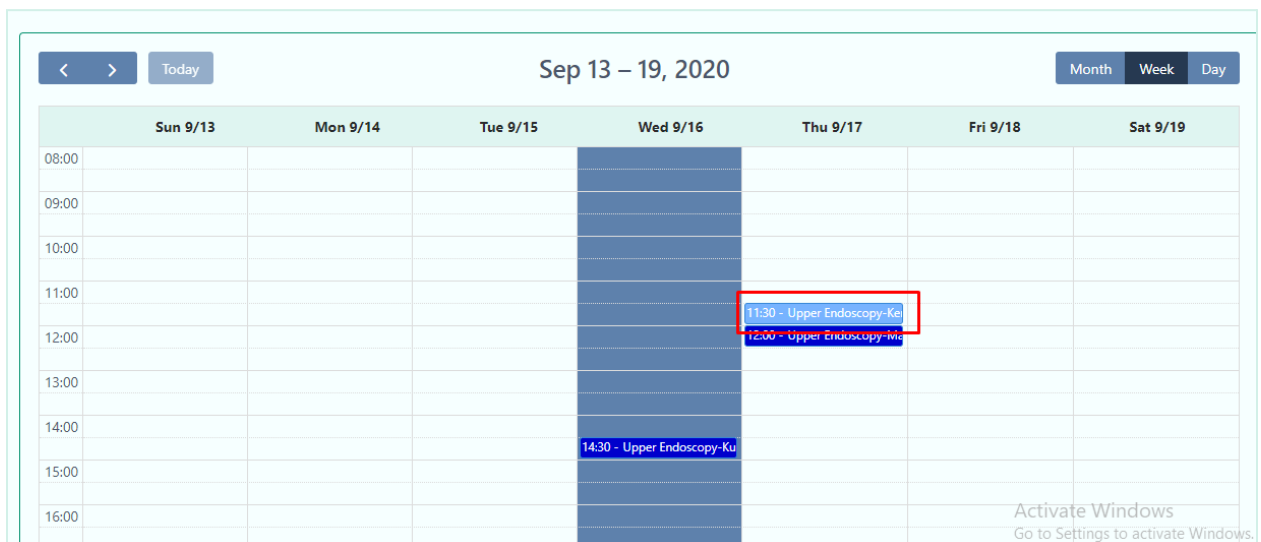

You can extend the slot duration by simply dragging the bottom of the coloured box that reserves the slot. You will consecutively be asked for a reason for the time change-

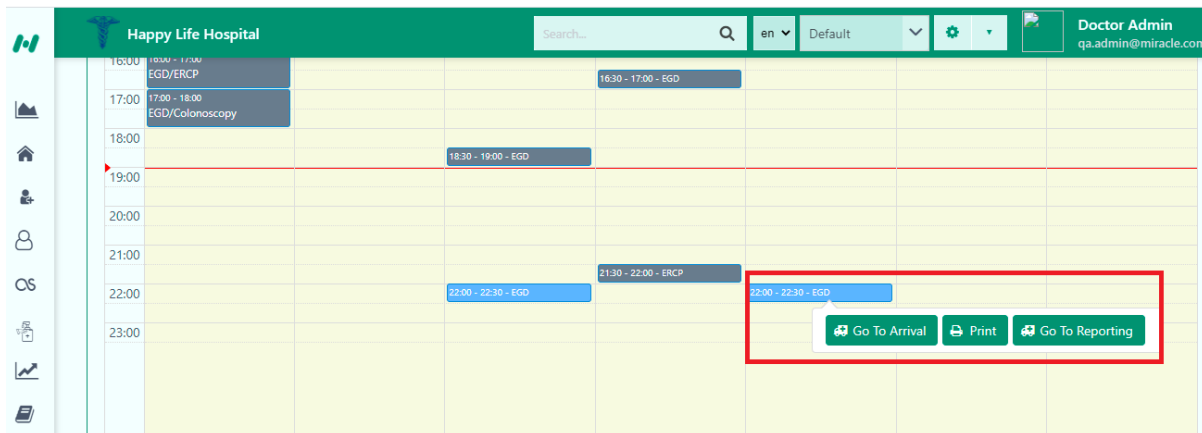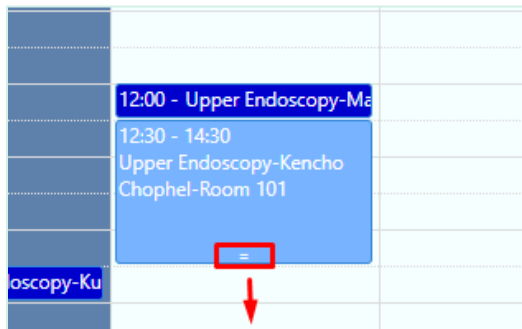


You can also register a patient and make an appointment from the **Tools** section in the **Whiteboard**.

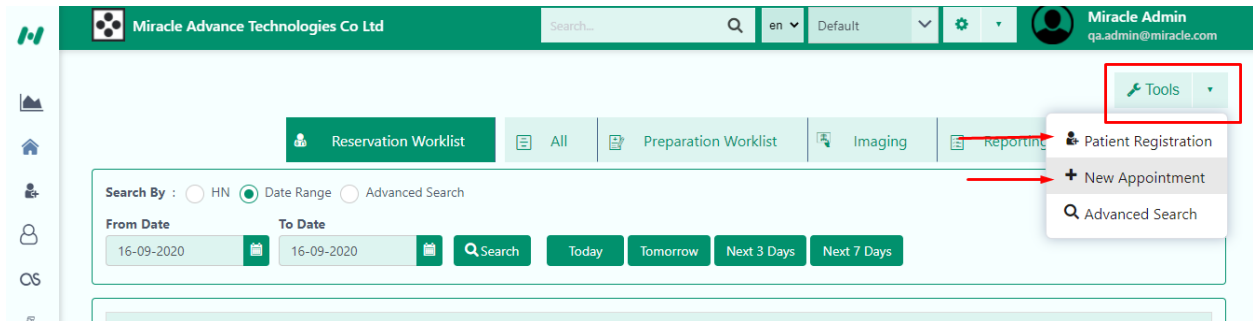


Whiteboard

Inside the Whiteboard, you can **view statuses** of a patient, what **stage** they are in, and in general, **the whole process/cycle** that a patient will go through from the moment s/he registers until the completion of the billing.


Whiteboard is organized with tabs that will carry a patient record from one tab to the next, until the client has completed the payments.

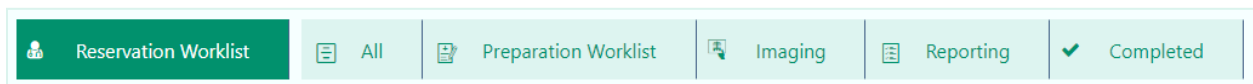

Once you have made an appointment (either from Tools or from Patient Dashboard), it will show in the **Reservation Worklist** under the respective date inside the Whiteboard.

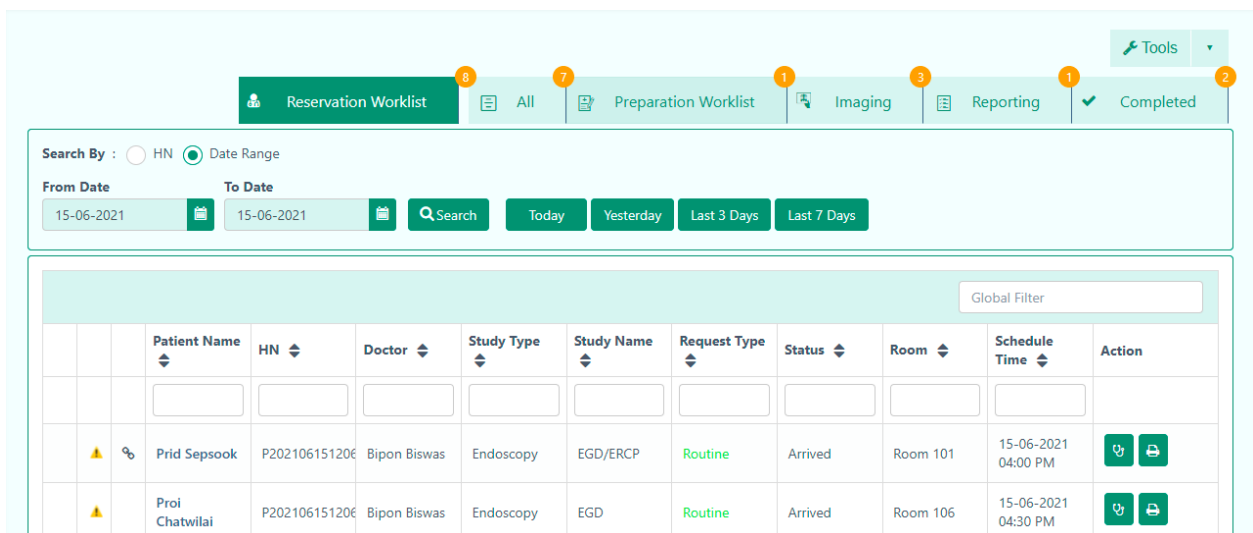

Tip: You can directly go to **Whiteboard** from the **Main Menu**, and even set it as your default home page.

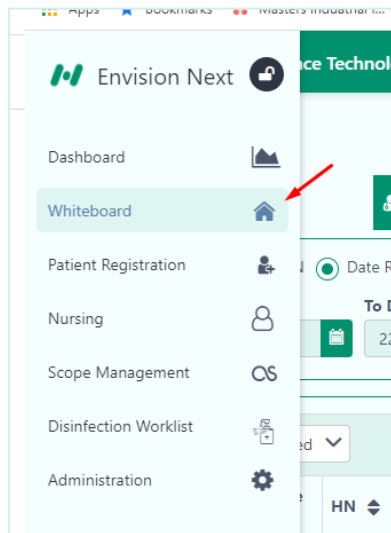

Whiteboard Features:

Search Features:

The application offers 3 types of search, namely **HN**, **Data Range** and **Advanced Search**. **Date Range** is selected as the default search type.

Date Range: The *From* and *To* dates are selected to be the present day by default. It means you will only see reservations made for the day in the **Reservation Worklist**, until you change the date range.

Search By : ☐ HN ☒ Date Range ☐ Advanced Search

From Date: 09-09-2020 To Date: 17-09-2020 Q Search Today Tomorrow Next 3 Days Next 7 Days

< September 2020 >

Su	Mo	Tu	We	Th	Fr	Sa
30	31	1	2	3	4	5
6	7	8	9	10	11	12
13	14	15	16	17	18	19
20	21	22	23	24	25	26
27	28	29	30	1	2	3

Physician	Study Type	Request Type	Procedure Type	Status	Room
1956 Miracle Admin	Endoscopy	Stat	Upper Endoscopy	Scheduled	Room

1 20

Click on **Search** to show the results from the date range.

To reduce your work, we have also set buttons that can directly show results for the past one, two and next three and next seven days.

Search By : ☐ HN ☒ Date Range ☐ Advanced Search

From Date: 17-09-2020 To Date: 17-09-2020 Q Search Today Tomorrow Next 3 Days Next 7 Days

9 columns selected

Patient Name	HN	Physician	Study Type	Request Type	Procedure Type	Status	Room
Kencho Chophel	P2020091611956	Miracle Admin	Endoscopy	Stat	Upper Endoscopy	Scheduled	Room 101

HN (Hospital Number):

Let's you search according to a specific hospital number-

Search By : ☒ HN ☐ Date Range ☐ Advanced Search

HN

Advanced Search:

In any case, the results may require you to make your filters as specific as possible, requiring many fields. Advanced Search gives you that option, and includes both HN and date range searches also.

Search By : ☐ HN ☐ Date Range ☒ Advanced Search ☒ Include Date Range in Search

From Date: 17-09-2020 To Date: 17-09-2020 HN: First Name: Last Name: Phone:

Requesting Dept: Study Type: Study: Arrival Status: Room: Clinical Diagnosis:

Clinical Finding: Study Purpose: Study Region: Chief Complaint: Request Type: Nursing Class:

Ambulatory Class:

Column Features:


- Selecting fields to show: select the specific fields you want/don't want to see

Search By : ☐ HN ☒ Date Range ☐ Advanced Search


From Date: 17-09-2020 To Date: 17-09-2020

9 columns selected


Patient Name	HN	Physician	Study Type	Request Type	Procedure Type
Kencho Chopel	P2020091611956	Miracle Admin	Endoscopy	Stat	Upper Endoscopy

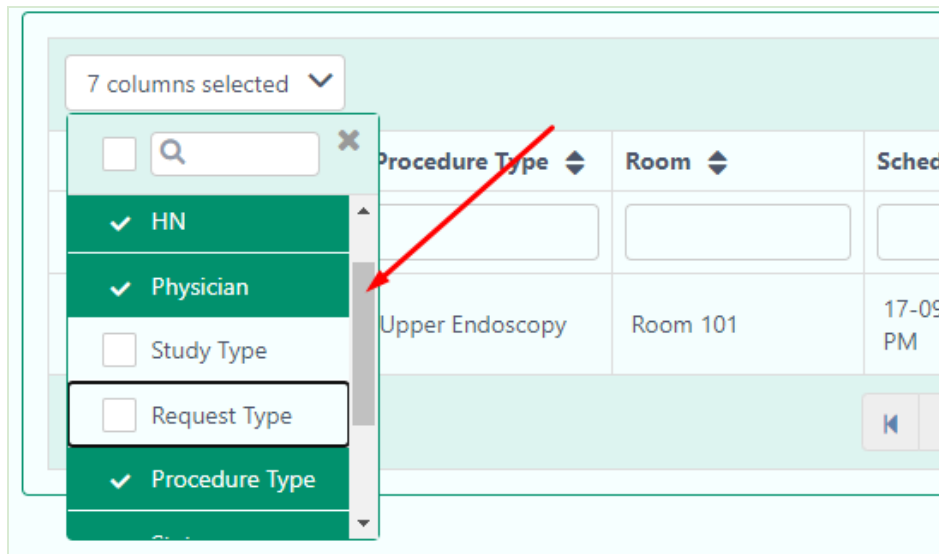

b. Type Search in every column: Use this feature to look up a specific value in a column

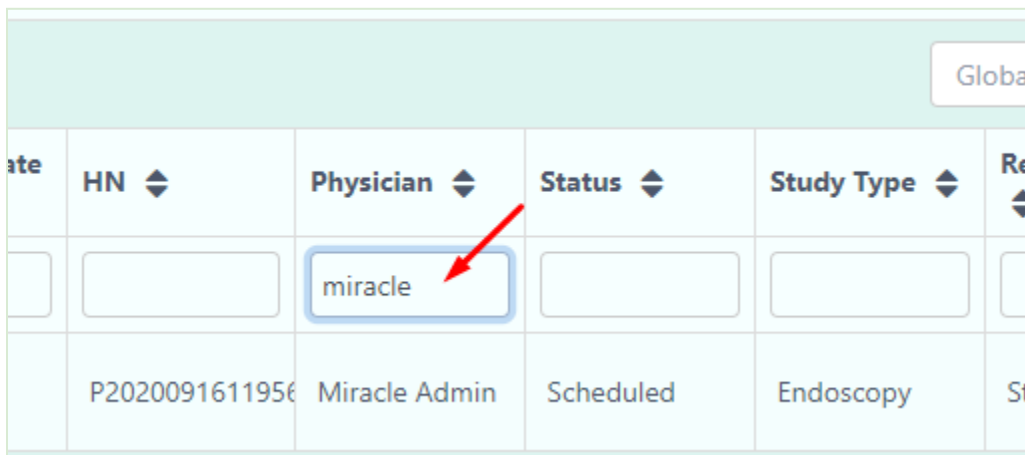

c. **Global Filter:** keywords you type are run in a search that covers data from all columns.

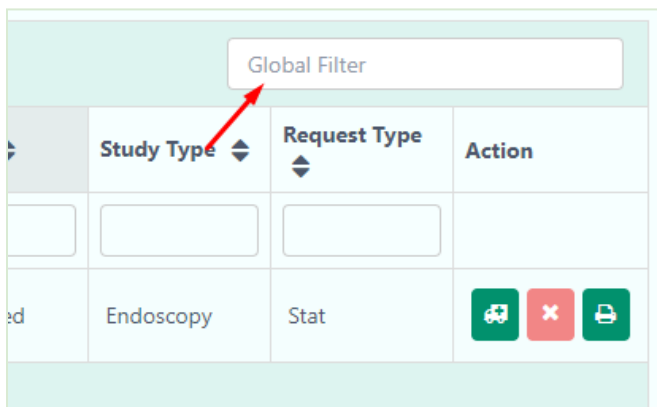

Reservation Worklist:

This worklist shows all the appointments for the day. From this worklist, ***you will have to mark arrival for the patient*** before he/she can be sent for preparation. This is done from the Action column:

	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Schedule Date	Action
	*****	P202102231202	Doctor Admin	Endoscopy	Colonoscopy	Routine	Scheduled	Room 102	24-02-2021 08:30 AM	   

Note: You can also arrive a patient from the Appointment Calendar

This takes the patient information to the **Patient Arrival** page-

Patient Arrival

Reservation Worklist

All

Preparation Worklist

Imaging

Reporting

Completed

Update

Arrived

Back

Patient Information

Name : Mr. Kencho Chopel

Phone : 985454665

Address :



Doctor Information

Name : Miracle Admin

Phone : +9999999999

Department : Physical Therapy Clinic



> Study

> Additional

Update

Arrived

Study and Advanced sections are made “collapsible”-

Study

Date*17-09-2020

Start Time*18:30

Duration (In Minutes)*30

End Time*20:00

Doctor*Miracle Admin

Room*Room 101

Study Type*Endoscopy

Study*Select procedure

Service Type	Study Name	Start Time	End Time	Action
Endoscopy	Upper Endoscopy	12:30 PM	02:00 PM	

Additional

Department
Internal Medicine
x
v

Patient Status
Out Patient
x
v

Floor

Phone

Request Type
Stat
x
v

Report

Content

Region
Choose
v

Finding
Choose
v

Intraoral bleeding, Abdominal rigidity.

Diagnosis
Choose
v

Gastric ulcer, Mitral insufficiency.

Consent
yes

Walking
Transportation with a stretcher
x
v

Chief Complain
Choose
v

Heartburn.

Nursing Class
Requires constant nursing
x
v

Purpose
Choose
v

Suspected tuberculosis.

Indication
B I U

Comment
B I U

Activate Windows
Go to Settings to activate Windows.

You can change information from the Patient Arrival page, and **Update** the information before you mark the patient “arrived”.

You will be asked to select the Clinic Type after you click on **Arrive**:

Confirm Clinic Type

Do You Want To Change Clinic Type?

Select Clinic Type :

Regular
v

OT

Emergency

Special

Regular

No

Yes

The patient information is now moved to the next stage: **Preparation Worklist**

Preparation Worklist:

The **Preparation** is the next step after a patient has arrived. It shows up in the **Whiteboard** like this-

Tools

Reservation Worklist All Preparation Worklist Imaging Reporting Completed

Search By : ☐ HN ☒ Date Range

From Date 23-02-2021 To Date 23-02-2021 Search Today Yesterday Last 3 Days Last 7 Days

10 columns selected Global Filter

Accession No	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Date	Action
20210223138	*****	P2021022312	Doctor Admin	Endoscopy	Colonoscopy	Routine	Arrived	Room 102	23-02-2021 02:21 PM	     

You can start preparing (and perform a number of other actions) the patient from this section.

Accession No	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Date	Action
20210223138	*****	P2021022312	Doctor Admin	Endoscopy	Colonoscopy	Routine	Arrived	Room 102	23-02-2021 02:21 PM	     

The button indicates **Resource Allocation**. After clicking the button, the user will get the **Endoscopy Preparation** page. Here they will find-

- **Team settings**
- **Study**
- **Scope List**
- **Endoscopy Machine**
- **Patient Information**

In the team settings, the user can set their desired team from the dropdown. They will get confirmation notification when adding the team.

HN : 123456
Name : นาย สมศรี สาระชาติ
DOB : September 5, 2016
Gender : Male
Study Name : EGD
Accession No : 2021103014381
Doctor Name : Dr. Doctor Admin



Primary Endoscopist

Dr. Doctor Admin

Recorder

Mr. Bipon Biswas

Team

Primary Endoscopist

Interventional Endoscopist

Supervising Endoscopist

Senior Fellow

Sedating Physician

Anesthesiologist

Nurse

Assistant Nurse

Study

EGD

Scope List

Scope List ☐ Include all Scopes

5 columns selected					
Global Filter					
Barcode	SerialNo	Scope Name	Study Name	Status	Action
20					

Endoscopy Machine

Endoscopy Machine List ☐ Include all Endoscopy Machine

Barcode, Machine Name		
Global Filter		
Barcode	Machine Name	Action
20		

In the **Examination Team**, the user can add team members to unlock some features of the Action and also as a helping hand to operate the patient efficiently.

Global Filter													
		Patient Name ▾	HN ▾	Doctor ▾	Study Type ▾	Study Name ▾	Request Type ▾	Status ▾	Room ▾	Examination Team ▾	Arrival Time ▾	Waiting Time	Action
		Prottoy Roy	P202106171	Panida Piyachaturav	Endoscopy	EGD / Colonoscopy	Routine	Arrived	Room 107	Panida Piyachaturav	17-06-2021 10:49 AM	Less than an hour	   
⏪ ⏩ 1 ⏪ ⏩ 20													

The team member who has no access, they cannot operate to-

- **Transfer to Imaging**
- **Vital Sign**
- **Patient History**

When the team member has the access, the user will find the interface like-

Global Filter													
		Patient Name ▾	HN ▾	Doctor ▾	Study Type ▾	Study Name ▾	Request Type ▾	Status ▾	Room ▾	Examination Team ▾	Arrival Time ▾	Waiting Time	Action
		Prottoy Roy	P202106171	Panida Piyachaturav	Endoscopy	EGD / Colonoscopy	Routine	Arrived	Room 107	Doctor Admin, Panida Piyachaturav	17-06-2021 10:49 AM	1 hour ago	      
⏪ ⏩ 1 ⏪ ⏩ 20													

the user can include all Endoscopy Machine, Study, and Scope List.

Study

EGD

Scope List

Scope List

☐ Include all Scopes

5 columns selected

Global Filter

Barcode	SerialNo	Scope Name	Study Name	Status	Action

1

20

To add **Vital Signs**, scroll, simply select the Vital sign icon from the **Action** column-

Global Filter

Room	Arrival Date	Action
Room 103	22-09-2020 10:07 PM	

This opens a panel from the left, where you can put in the information, **Save**, or **Close** the panel. This panel can also be helpful in providing you the patient's basic details (name, DoB, Gender, HN) at a glance.

Vital Sign

Patient Name : Mrs Carolina Statham

HN : P2020092211958

Gender : Female

DOB : March 5, 1985



Save

Height (cm)

Weight (kg)

BMI

BP (mmHg)

Systolic / Diastolic

SPO₂

Heart Rate (BPM)

Temperature (F)

Temperature (C)

Respiration Rate (RPM)

Pulse rate (BPM)

Capillary Refill Time (Sec)

BCS

Pain Score

Hydration Status

We can now move on to **Preparation**.

The first option is ***Prepare PreScope***, that preps the patient for the endoscopy-

10 columns selected												Global Filter	
		Accession No	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Date	Action	
		20210223138	*****	P2021022312	Doctor Admin	Endoscopy	Colonoscopy	Routine	Arrived	Room 102	23-02-2021 02:21 PM	     	

This takes you to the **Prepare Prescope** section:

(not full section shown here)

Pre-Preparation

Clear All
 Save
 Save & Transfer

Patient Information

Name : Carolina Statham
Phone : 9635474857
Address :

Doctor Information

Name : Miracle Admin
Phone : +9999999999
Department : Physical Therapy Clinic

Patient Information

Procedure/Study

PEG

Status

OPD

Diagnosis

Mitral insufficiency, Gastric ulcer

Indication

+

Purpose

Suspected tuberculosis

Finding

Intraoral bleeding, Abdominal rigidity

This section contains a number of forms that are not collapsible. If you keep scrolling, you will find all the options that you can select/fill in. The sections include:

- Patient Information

Patient Information

Clear

Procedure/Study

PEG

Status

OPD

Diagnosis

Mitral insufficiency, Gastric ulcer

Indication

+

Purpose

Suspected tuberculosis

Finding

Intraoral bleeding, Abdominal rigidity

Nursing Class

Requires occasional nursing

Region

+

Chief Complain

Dysphagia

Clicking on the plus sign lets you enter text information:

Select Indication

Indication

B I U

Fever, Nausea

Cancel

Save

- Check List

Check List				Clear	
Confirm Patient Name	☒ Yes	Confirm Procedure	☒ No	Consent Form Sign	☒ Yes
NOP	☒ No	Antiplatelet/ Anticoagulant	☒ Yes	NO PPI	☒ Yes
Companion	☒ No	Remove False Teeth/Acc	☒ No		
Attachment	☐ OPD Card ☐ Other ☐ Film				
Premed	☒ Domicum ☐ Pethidine				

The default choice for all items in the list is set to 'No'. Clicking on it will set it to 'Yes'

● Patient History

Patient History				Clear	
Drug Allergy	☒ No	Egg Soybean Allergy	☒ No	Previous Surgery	☒ No
Pregnancy	☒ No	Pacemaker/ICD	☒ No	Special Equipment	☒ Yes
Prophylactic Antibiotic	☒ Yes	Difficult Airway/ Aspiration	☒ No	Anxiety	☒ No
Smoking	☒ Yes				
Alcohol drinking	☐ Regular ☐ Social ☐ No				
Comorbidities					
None	☒ No	DM	☒ Yes	HT	☒ No
Heart Disease	☒ Yes	Cirrhosis	☒ No	CRF	☒ Yes
Bleed Disorder	☒ No	Other	☒ No		
Current Medication					
None	☒ Yes	Aspirin	☒ No	NSAIDs	☒ No
Plavix	☒ No	Warfin	☒ No	Icon	☒ No

(not full section shown here)

This is the same as **Check List**, the default choice for all items in the list is set to 'No'. Clicking on it will set it to 'Yes'.

● Patient History

Patient History				Clear			
Drug Allergy	☒ No			Egg Soybean Allergy	☒ No	Previous Surgery	☒ No
Pregnancy	☒ No			Pacemaker/ICD	☒ No	Special Equipment	☒ No
Prophylactic Antibiotic	☒ No			Difficult Airway/ Aspiration	☒ No	Anxiety	☒ No
Smoking	☒ No						
Alcohol drinking	Regular	Social	No				
Comorbidities							
None	☒ No			DM	☒ No	HT	☒ No
Heart Disease	☒ No			Cirrhosis	☒ No	CRF	☒ No
Bleed Disorder	☒ No			Other	☒ No		
Current Medication							
None	☒ No			Aspirin	☒ No	NSAIDs	☒ No
Plavix	☒ No			Warfin	☒ No	Icon	☒ No
NOACs	☒ No			Other	☒ No		
LOC							
Alert	☒ No			Confuse	☒ No	Drowsy	☒ No
Stuporous	☒ No			Semicoma	☒ No	Coma	☒ No

● Breathing

Breathing			Clear
Spontaneous	☒ Yes	Room Air	☒ Yes
O2 Mask C Bag	☒ No	Artificial	☒ No
ET Tube	☒ No	Ambu Bag	☒ Yes
		O2 Cannula	☒ No
		T Tube	☒ Yes
		Ventilator	☒ No

● Blood Component

Blood Component				Clear	
Red Cell(Unit)	5	☒ Yes	Plasma(Unit)	12	☒ Yes
CRYO(Unit)		☒ No	Platelet(Unit)		☒ No

This is also the same as **CheckList**, the default choice for all items in the list is set to **‘No’**. Clicking on it will set it to **‘Yes’**. When you click **‘Yes’**, the user can give input. Otherwise, the user can’t give input if it is **‘No’**.

● Complete Blood Count

Complete Blood Count				Clear
Hct	7	Plt	7	
PT	4	INR	4	
PTT	6.8	CBC Date	08-09-2020	

● Blood Chemistry

Blood Chemistry				Clear
Bun	3	Cr	76.4	
Na	5.8	K	5	
Cl	10.7	CO2	5.4	
Chem Date	02-09-2020			

● Liver Function Test

Liver Function Test				Clear
TB	67.5	DB	23.9	
SGOT	0	SGPT	457	
ALP	12.9	ALB	98.56	
LFT Date	16-09-2020	Other Labs		
Labs Date	11-09-2020			

● Check List

Check List				Clear	
Confirm Name Procedure	☒ Yes	Confirm Correct HN	☒ Yes	Patient Signed Consent	☒ Yes
Antibiotic before procedure	☒ No	No Antibiotic before procedure	☒ Yes	Documents Check	☒ Yes
OPD Card	☒ Yes	Other Check	☒ No	Fim	☒ No
Critical Event	☒ No				

● Nurse Info

Nurse Info	Alert Message
Nurse Sign Out: Please Select Nurse Nurse User	Set Alert Message: No

The nurse can set an Alert Message for the doctor if needed.

Alert Message

Set Alert Message: ✓ Yes

The user has to enable **Yes** to write the alert message for the patient.

The last option is Chief Complaint and Diagnosis, that prepares the patient for the endoscopy-

Chief Complaint and Diagnosis

Chief Complaint

Heartburn, Constipation

Heartburn, Constipation.

Diagnosis

- Anastomotic stricture
- Angioectasia
- Atrophic gastritis
- Barretts esophagus

Qualifier

- ☐ Sampling of
- ☐ Screening of
- ☒ Suspected
- ☐ Treatment of

+ Add

Diagnosis	Qualifier	Action
Atrophic gastritis		☐ ☐

✕ Close
Save

After selecting the Diagnosis and Qualifier, the user will get a confirmation message by clicking on the **Save** button.

		20210201138	Mahfuz05Pet Mahfuz05Pet	P20200505-11485	Doctor Admin	Endoscopy	Upper Endoscopy	Routine		Room 102	01-02-2021 11:02 AM	  
		20210201138	MAX MAX	P20200427-10723	Doctor User	Endoscopy	Upper Endoscopy	Routine		Room 104	01-02-2021 04:17 PM	  
		20210201138	MAX MAX	P20200427-10723	Doctor Admin	Endoscopy	ERCP	Routine		Room 105	01-02-2021 04:26 PM	  
		20210202138	*****	P2021012212	Bipon Biswas	Endoscopy	Upper Endoscopy	Routine		Room 101	02-02-2021 01:46 PM	  
		20210202138	MAX MAX	P20200427-10723	Doctor Admin	Endoscopy	Upper Endoscopy	Routine		Room 102	02-02-2021 02:36 PM	  
		20210202138	MAX MAX	P20200427-10723	Doctor Admin	Endoscopy	Colonoscopy	Routine		Room 102	02-02-2021 02:36 PM	  
		20210202138	MAX MAX	P20200427-10723	Doctor Admin	Endoscopy	ERCP	Routine		Room 102	02-02-2021 02:36 PM	  

Chief Complaint and Diagnosis
Updated Successfully

Vital Signs

You can add a patient's Vital signs from the **Preparation Worklist** or from the **Pre-preparation** section. To add from the latter, click on the **Vital Signs** icon at the end of the form:

This opens the Vital Sign form as a panel from the left-

Patient Name : Mrs Carolina Statham

HN : P2020092211958

Gender : Female

DOB : March 5, 1985



Update

Height (cm)

153

Weight (kg)

56

BMI

23.92

BP (mmHg)

Systolic / Diastolic

SPO₂

Heart Rate (BPM)

Temperature (F)

189

Temperature (C)

87

Respiration Rate (RPM)

Pulse rate (BPM)

Capillary Refill Time (Sec)

BCS

Pain Score

4

Hydration Status

67

Vital signs saved successfully.

You can now save the details and transfer the patient for **Imaging**, using the **Save & Transfer** button. If you wish to just save the information and not **Transfer** yet, you can select the **Save** button, and **Transfer** later.

Nurse Info

Nurse Sign Out:

Nurse User x v

Clear All

Save

Save & Transfer

Note that the *Vital Signs button* is now red, indicating that it has information in it. Once you transfer the patient, the patient information moves to the **Imaging** tab.

Imaging:

From the Imaging section, you can-

Tools

Back

Reservation Worklist

All 171

Preparation Worklist 4

Imaging 40

Reporting 50

Completed 77

Search By :

☐ HN
☒ Date Range
☐ Only My Cases

From Date

To Date

01-09-2021

16-02-2022

Search

Today

Yesterday

Last 3 Days

Last 7 Days

☐ Show Canceled

12 columns selected v

Global Filter

		Patient Image	Patient Name	HN	Doctor	Study Type	Study Name	Accession	Request Type	Status	Room	Examination Team	Arrival Time	Waiting Time	Action
			นาย สมศรี สดอาด	123456	Bipon Biswas	Endoscop	Colonoscopy / EGD	20210927 / 20210927	Routine	Imaging In Progress P(1/2)	Room 105	Bipon Biswas, David Attenborc	27-09-2021 09:24 AM	168 day(s) ago	

- **Start imaging**
- **Patient PreScope**
- **Patient History**
- **Chief Complaint and Diagnosis**
- **Resource Allocation**
- **Add Consumable Item**

Patient History: As a doctor, you might want to take a look at the patient's **History**. Click on the icon that takes you to the patient's past records:




Updated:23/2/2021


Start imaging: When the patient is ready, the user can start imaging by clicking the start imaging button.

Patient Image	Patient Name	HN	Doctor	Study Type	Study Name	Accession	Request Type	Status	Room	Examination Team	Arrival Time	Waiting Time	Action
	นางสาว อรุณพร อรุณพร	123456	Bipon Biswas	Endoscopy	Colonoscopy / EGD	20210927 / 20210927	Routine	Imaging in Progress P(1/2)	Room 105	Bipon Biswas, David Attenbore	27-09-2021 09:24 AM	168 day(s) ago	



Note: If you do not select the scope in the previous step , it will continuously give notification for selecting the scope.

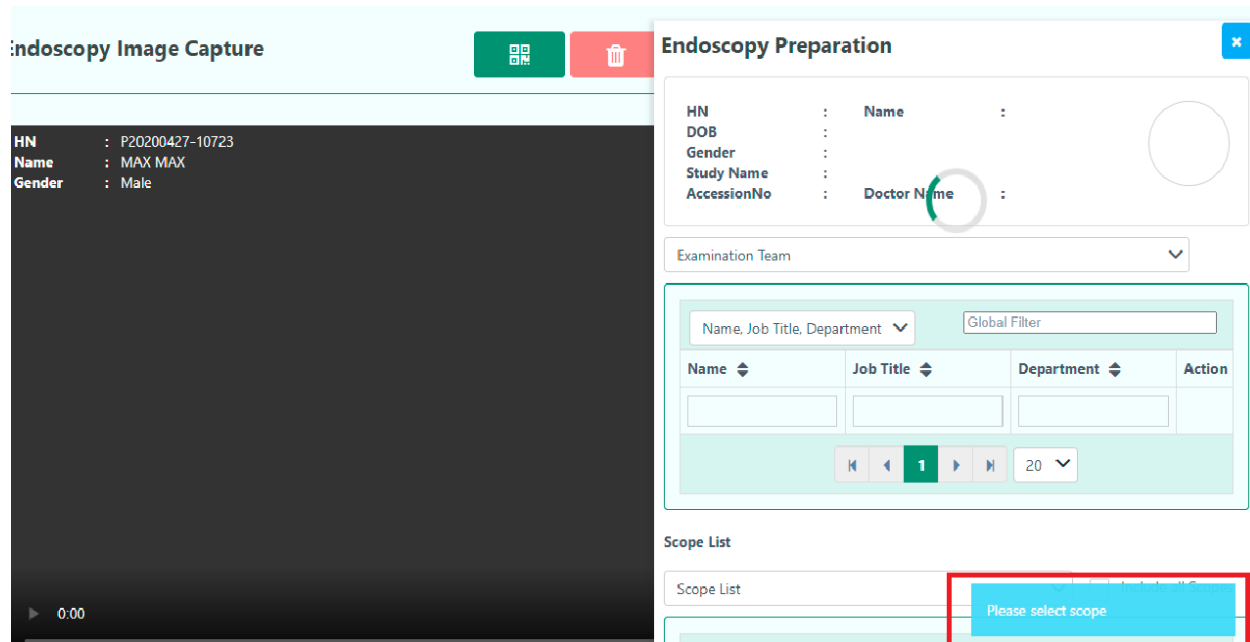


The user has to select the Resource Allocation for selecting the scope.

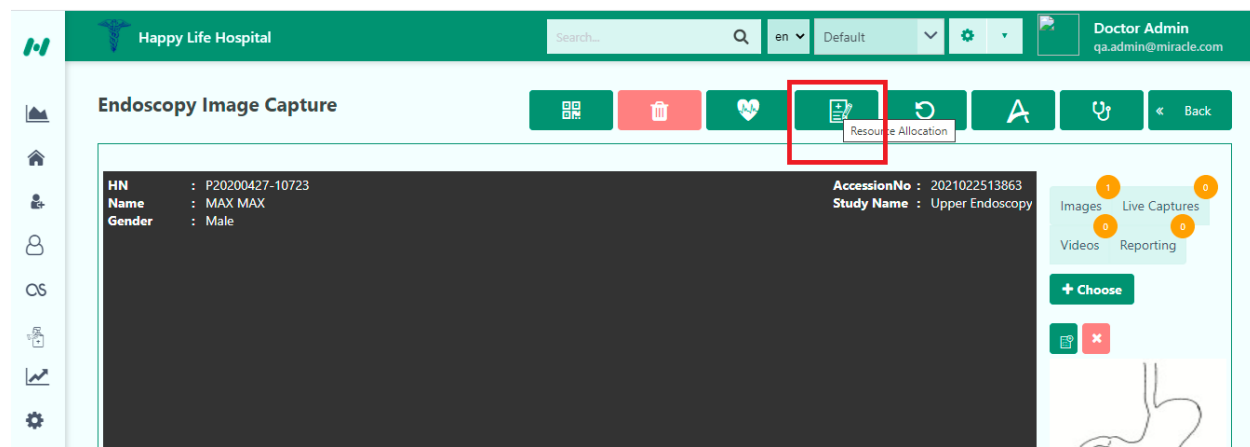


After selecting the **Resource Allocation**, the user has to select **Include all Scopes**. Then, select Scopes from the dropdown.

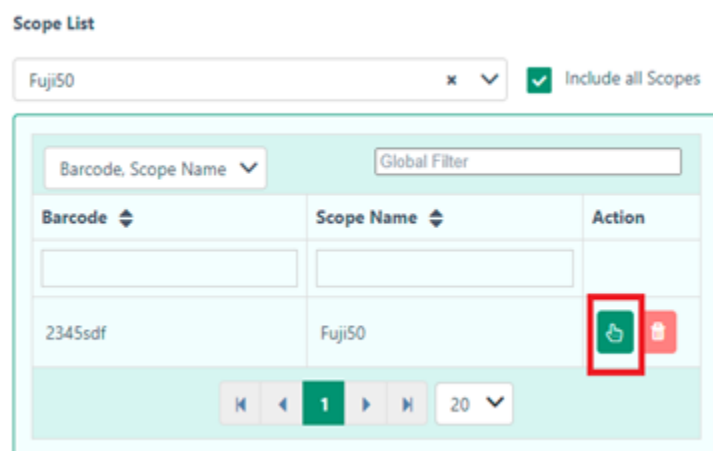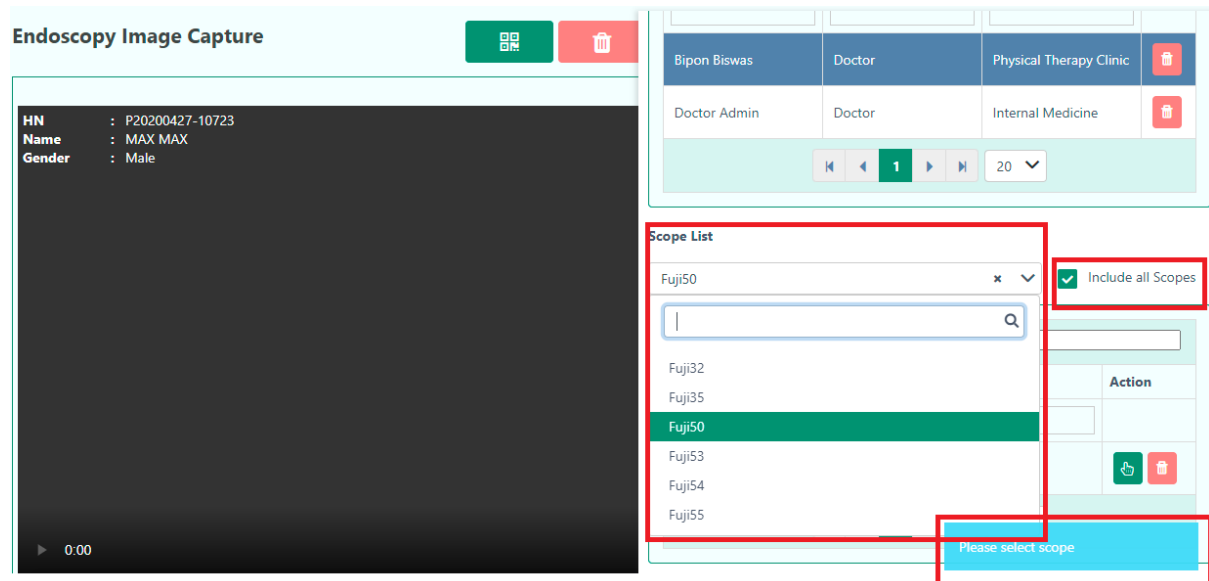



After pressing the select button, the user then selects the **Scope in** button to start the procedure.

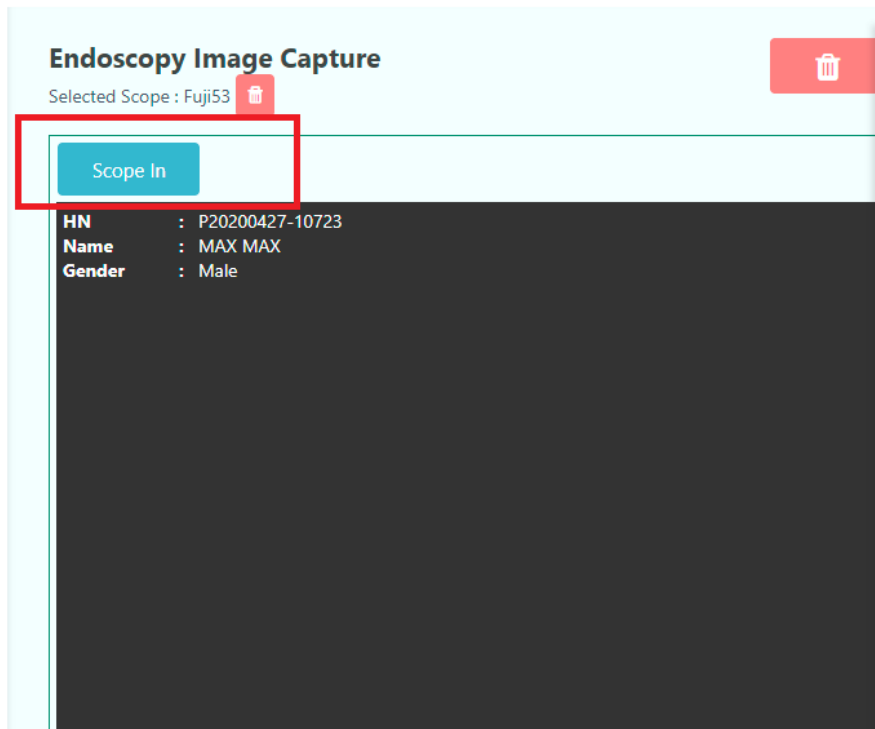


Note: When the user presses the **Scope in** button, the Action list in the right corner becomes disable. The user should disinfect the scope later to use it again.

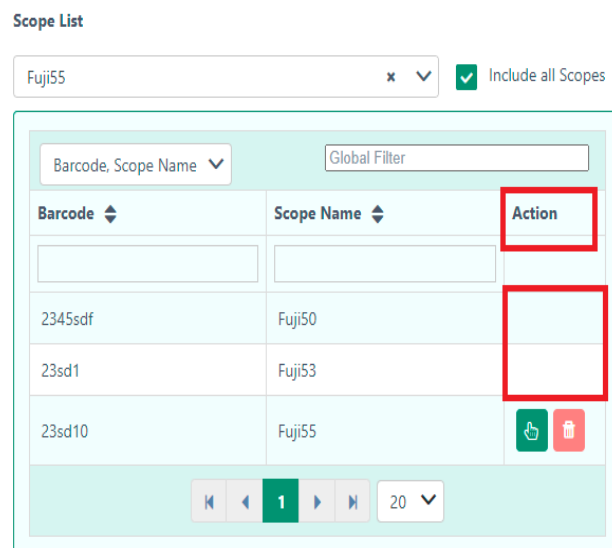


The user then starts the imaging process by selecting the video source.

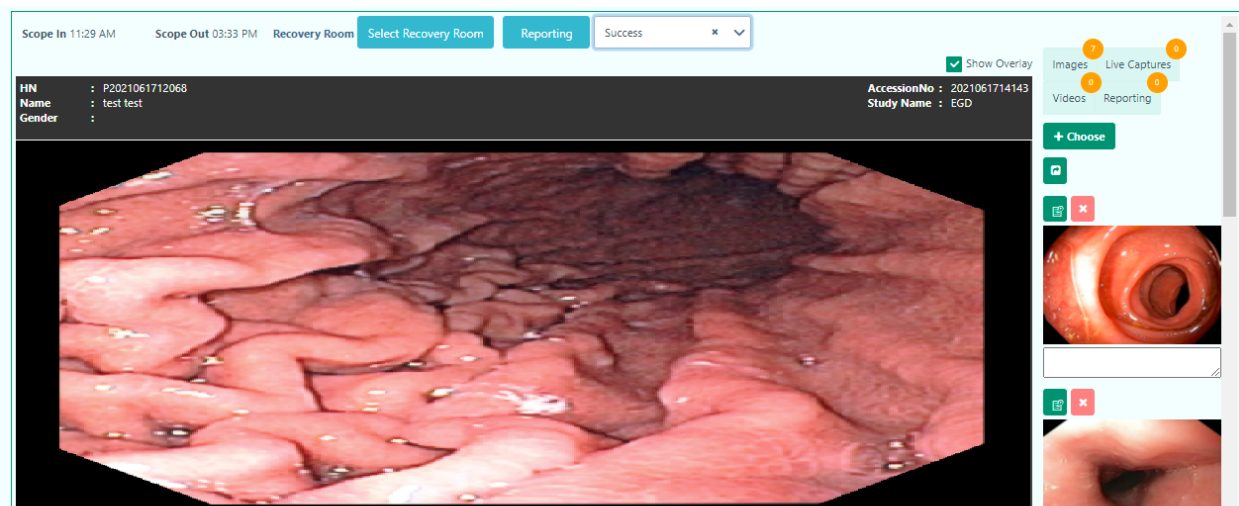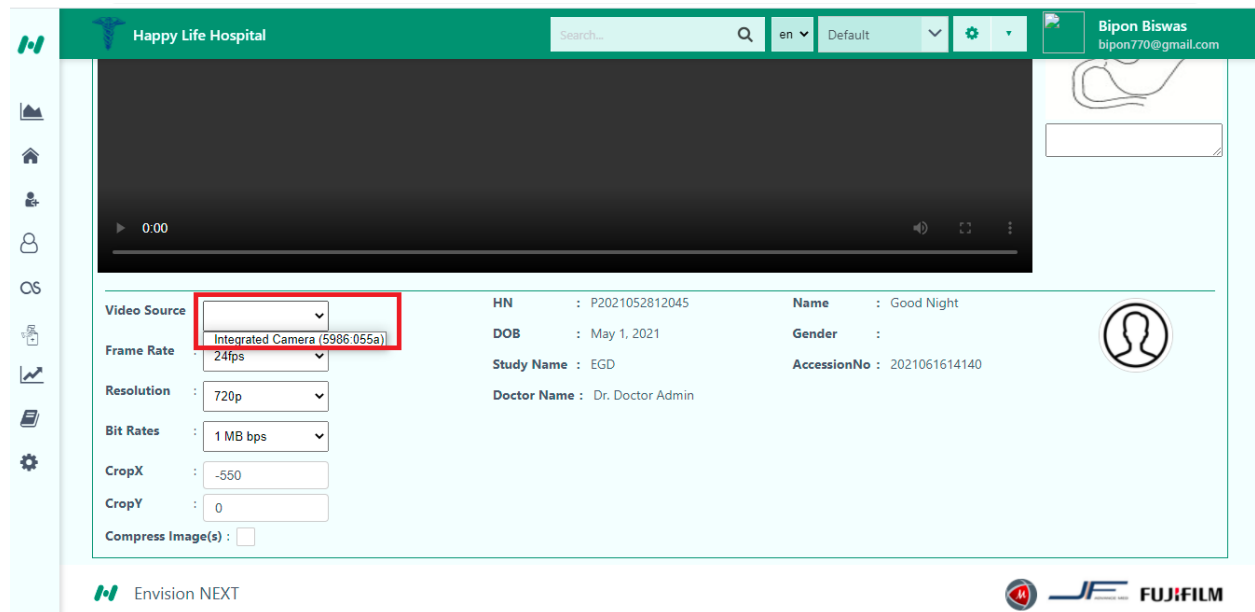

The user can notice the number of images and Live captures from the above.

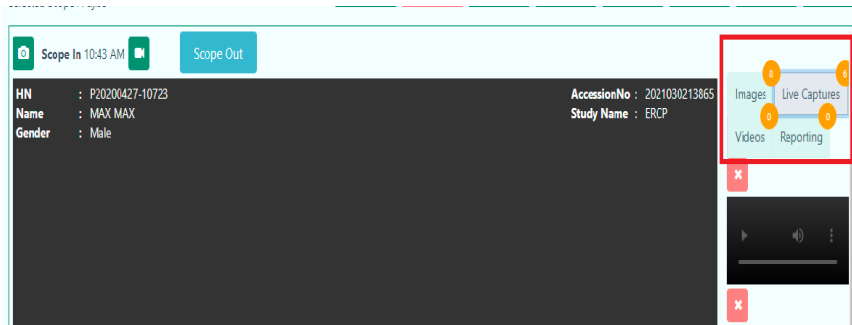



The user can see the video where the image is captured by simply clicking this button.

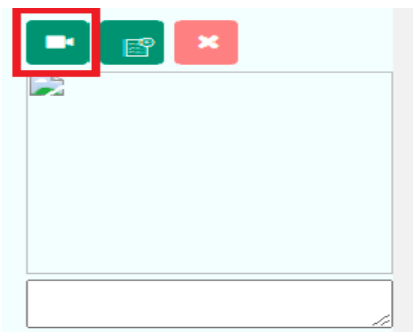


Note- Keyboard Shortcut-

- Alt+ c= Take Picture.
- Alt+ r= Start Recording.
- Alt+ s= Save.
- - = Capture Image

Moreover, the user can crop video by giving start and end value.

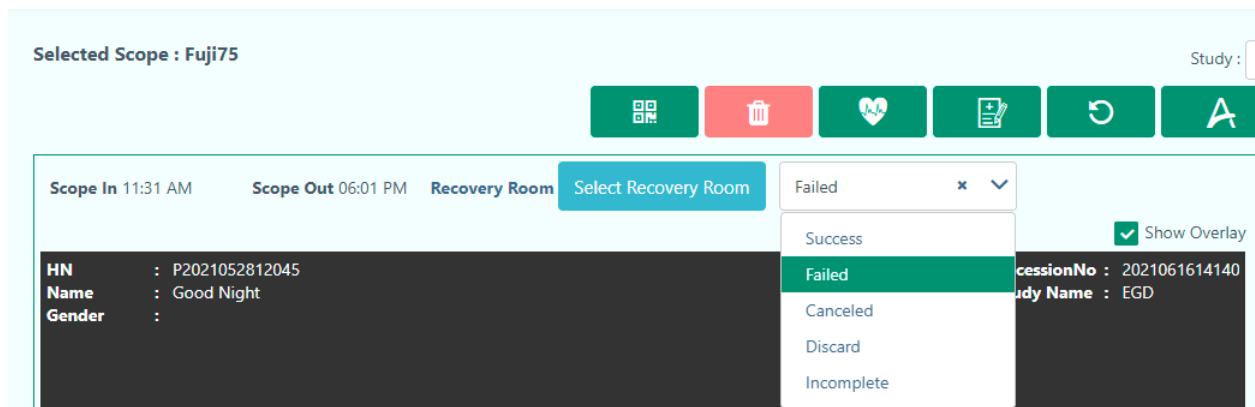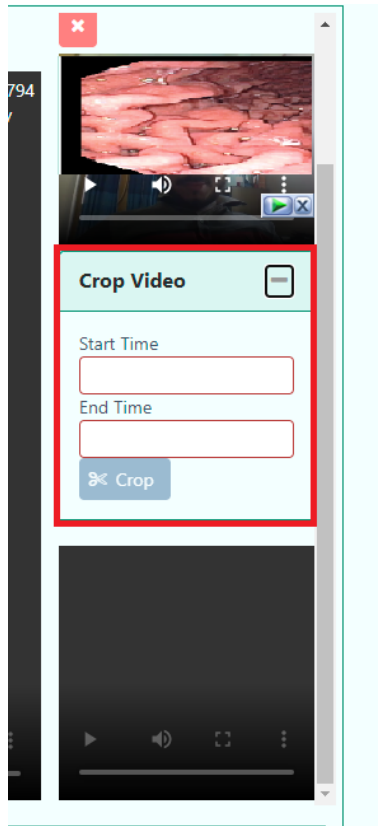


The user has to select procedure status to enable the Reporting button. The button will enable when the status is **Success**.

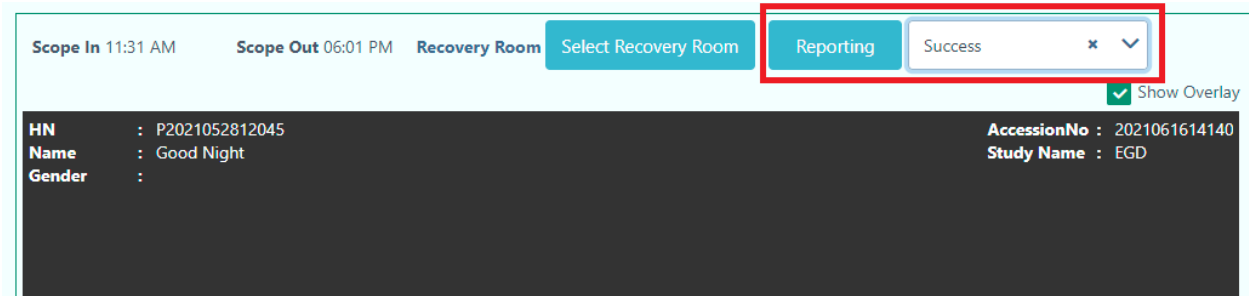

There are five procedure states.

- Success
- Failed
- Canceled
- Discard
- Incomplete

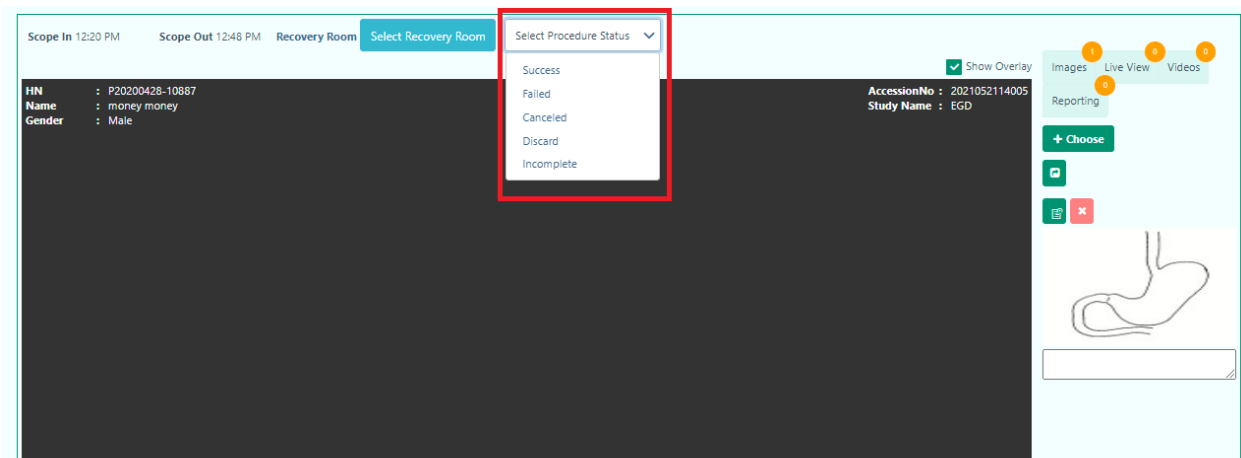

The procedure depends on the study. When you finish the study successfully, you can select Success. Otherwise, you can select according to your result.

The user can also control the overlay of the screen by simply clicking the “Show Overlay” option.

☒ Show Overlay

HN	: P2021052812045	AccessionNo	: 2021061614140
Name	: Good Night	Study Name	: EGD
Gender	:		

☐ Show Overlay

In the imaging section, the user can also check the **Vital Sign** of the patient.

Endoscopy Image Capture

Scope In 12:32 PM

Scope Out

Choose

HN	: P20200427-10723	AccessionNo	: 2021022513863
Name	: MAX MAX	Study Name	: Upper Endoscopy
Gender	: Male		


Patient Name : Mrs Carolina Statham HN : P2020092211958 Gender : Female DOB : March 5, 1985		
Update		

Height (cm)	153	Weight (kg)	56
BMI	23.92	BP (mmHg)	Systolic / Diastolic
SPO ₂		Heart Rate (BPM)	
Temperature (F)	189	Temperature (C)	87
Respiration Rate (RPM)		Pulse rate (BPM)	
Capillary Refill Time (Sec)		BCS	
Pain Score	4	Hydration Status	67

Vital signs saved successfully.

The user can also check Patient History and Chief Complaint and Diagnosis accordingly.


The user can add all Images to **Reporting** and also they can select one by one to **Reporting**.

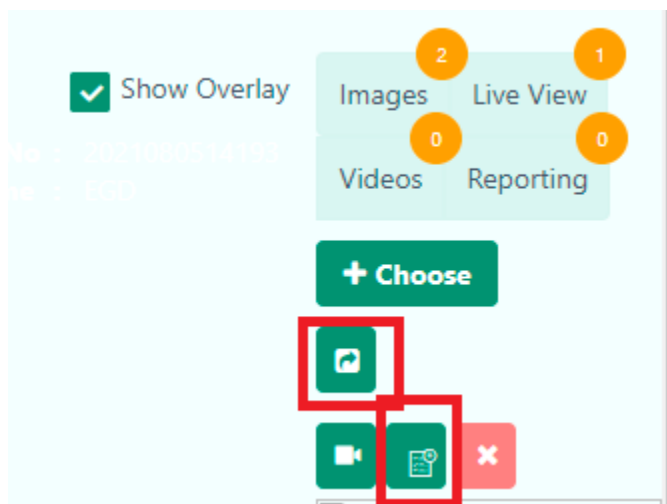

Updated:25/2/2021

Reporting:

The user can go to Reporting from Imaging by selecting the **Procedure status** to **Success**.

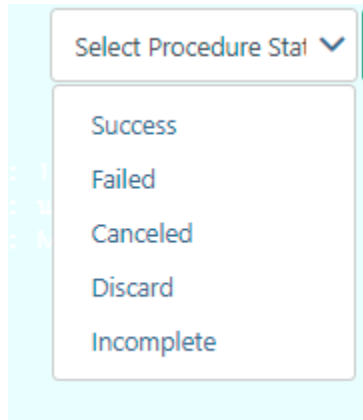


A dropdown menu titled "Select Procedure Status" with a downward arrow. The menu is open, showing five options: "Success", "Failed", "Canceled", "Discard", and "Incomplete".



A button labeled "Success" with a close icon (X) and a dropdown arrow, next to a green button labeled "Reporting".

Also the user can easily generate reports for the patient from this section.

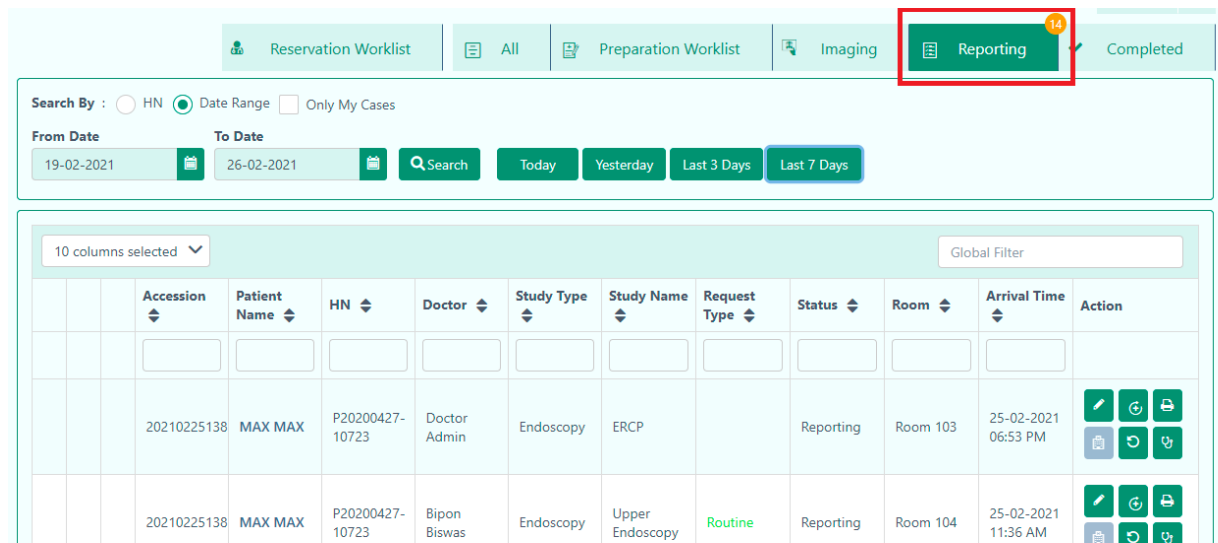

The screenshot shows the "Reporting" section of the application. The top navigation bar includes "Reservation Worklist", "All", "Preparation Worklist", "Imaging", and "Reporting" (highlighted with a red box and a notification badge). Below the navigation bar, there is a search section with "Search By" options (HN, Date Range, Only My Cases) and date range filters (From Date: 19-02-2021, To Date: 26-02-2021). A "Search" button and time filters (Today, Yesterday, Last 3 Days, Last 7 Days) are also present. The main area displays a table with 10 columns selected. The table has headers: Accession, Patient Name, HN, Doctor, Study Type, Study Name, Request Type, Status, Room, Arrival Time, and Action. Two rows of data are visible, both with "Reporting" status.


Accession	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Time	Action
20210225138	MAX MAX	P20200427-10723	Doctor Admin	Endoscopy	ERCP		Reporting	Room 103	25-02-2021 06:53 PM	[Edit] [Refresh] [Print] [Share] [Download]
20210225138	MAX MAX	P20200427-10723	Bipon Biswas	Endoscopy	Upper Endoscopy	Routine	Reporting	Room 104	25-02-2021 11:36 AM	[Edit] [Refresh] [Print] [Share] [Download]

From the **Reporting** the user can-

- Write Reporting.

- Go to PostScope.
- Discharge.
- Patient History.
- Chief Complaint and Diagnosis.
- Report.
- Consumable Item.

Write Reporting: When the imaging part is done, then the user can start this operation by clicking this button.

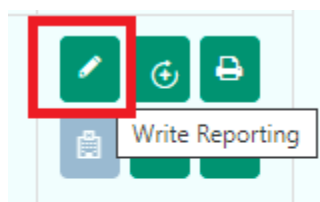


In this section, the user will find two sections.

- The Information of the Patient
- Captured Files.


In the first dropdown, the user will find the patient information.


Finally, the user can add the images to the reporting.

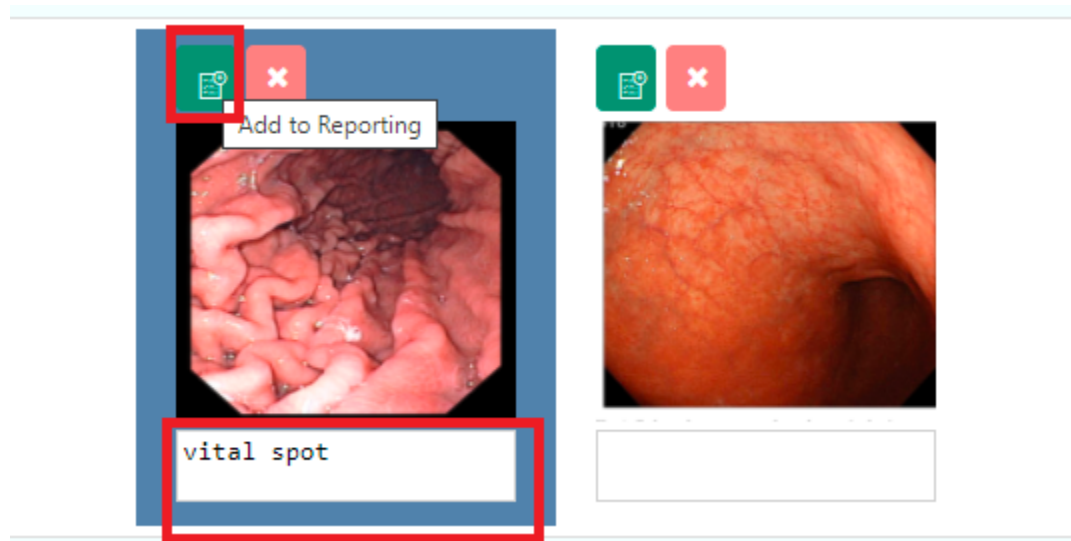

The user can easily delete the picture with the above cross sign button.

Note: If the user deletes the picture mistakenly, they can simply restore the photos.

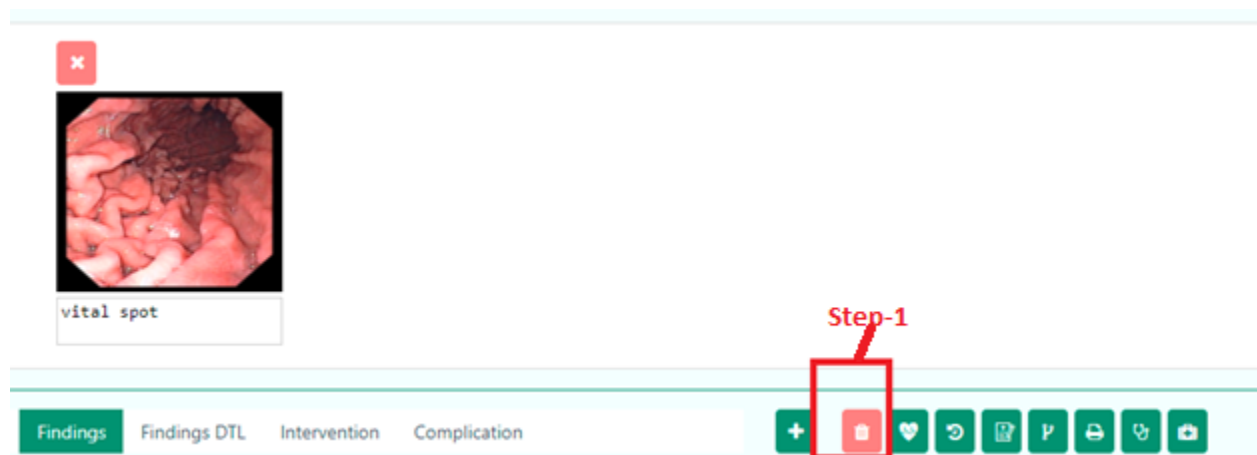

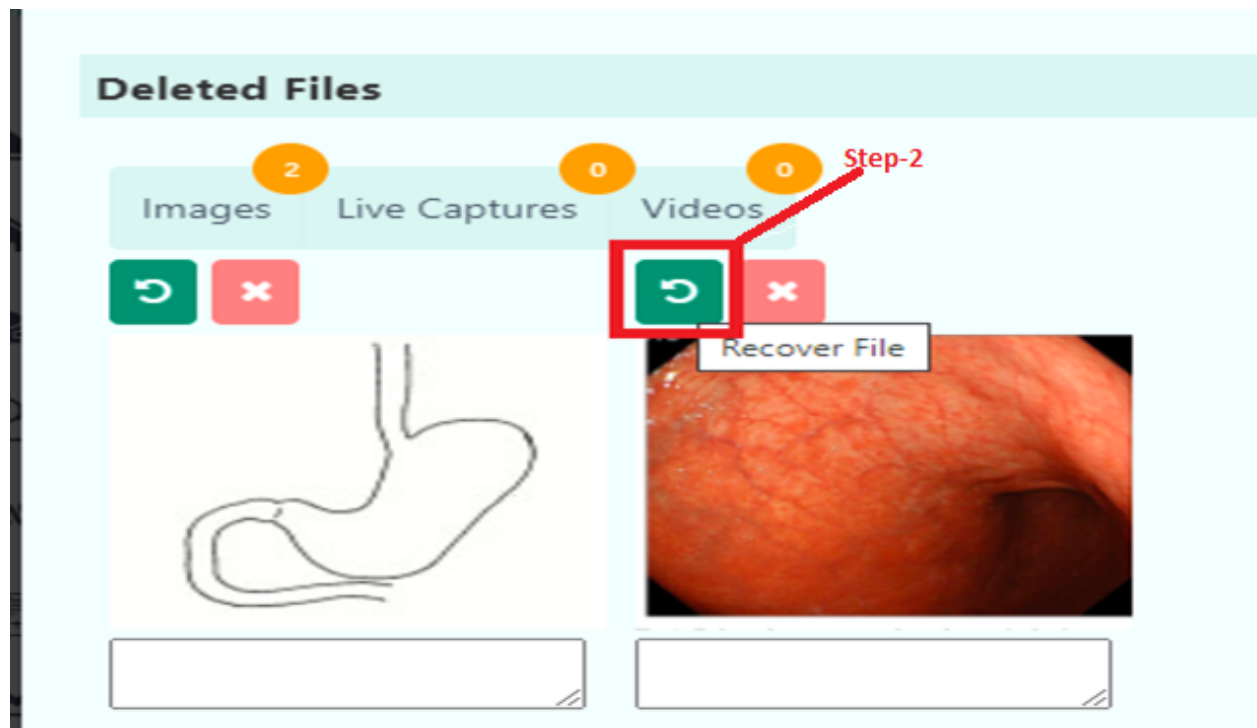

Study Metrics

The screenshot shows the 'Study Metrics' interface. At the top, there are tabs for 'Indication/Diagnosis', 'Sedation', 'Bowel Preparation', 'Timing', and 'Team'. The 'Timing' tab is selected. Below the tabs, there is a section titled 'Chief Complaint and Diagnosis'. This section contains three columns: 'Chief Complaint', 'Diagnosis', and 'Qualifier'. Each column has a search bar and a list of items. The 'Chief Complaint' column lists 'asd', 'Cancer screening', 'Constipation', and 'Defecation disorder'. The 'Diagnosis' column lists 'Bile duct stone', 'Bile duct stricture', 'Caroli disease', and 'Caustic ingestion'. The 'Qualifier' column lists 'Established', 'Follow-up of', 'Sampling of', and 'Screening of'. There are 'Add' and 'Save' buttons at the bottom right.

Chief Complaint	Diagnosis	Qualifier
☐ asd	☐ Bile duct stone	☐ Established
☐ Cancer screening	☐ Bile duct stricture	☐ Follow-up of
☐ Constipation	☐ Caroli disease	☐ Sampling of
☐ Defecation disorder	☐ Caustic ingestion	☐ Screening of

The user will find five sections in **Study Metrics**.

- **Indication/Diagnosis**
- **Sedation**

- **Bowel Preparation**
- **Timing**
- **Team**

Indication/Diagnosis-

Indication/Diagnosis
Sedation
Bowel Preparation
Timing
Team

Chief Complaint and Diagnosis

Chief Complaint	Diagnosis	Qualifier
☐ ☐ asd ☐ Cancer screening ☐ Constipation ☐ Defecation disorder	☐ Bile duct stone Bile duct stricture Caroli disease Caustic ingestion	☐ ☐ Established ☐ Follow-up of ☐ Sampling of ☐ Screening of

+ Add

Add Sedation

Indication/Diagnosis
Sedation
Bowel Preparation
Timing
Team






Add Sedation

Local xylocaine throat spray ☒ Yes ☐ No

Sedation ☒ Yes ☐ No

Save

Bowel Preparation

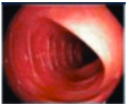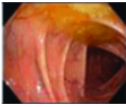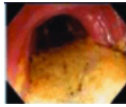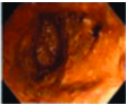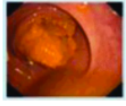Bowel Preparation					
BBPS	3	2	1	0	
3=Excellent 2=Good 1=Poor 0=Inadequate					
					
LC	☐	☐	☐	☐	☐
TC	☐	☐	☐	☐	☐
RC	☐	☐	☐	☐	☐
BBPS = 0					
Agent:	☐ Niftec ☐ PicoPrep ☐ PEG ☐ Swift ☐ Unison Enema ☐ Cleansing Enema ☐ Bisacodyl Suppository				
Bowel Preparation:					
					Reset Save

Timing

The user can view the start and end procedural time of the scope. It also provides the duration of the procedural time.

🔍 Indication/Diagnosis 💊 Sedation 🚽 Bowel Preparation ⌚ Timing 👤 Team				
Colonoscopy				
Study Time				
Start Procedural Time	End Procedural Time	Duration (hh:mm:ss)		
15:49	16:04	00:15:00		
Scope List v		☐ Include all Scopes		
Save Changes				

The user can also select scope from Scope list. Moreover, they can select **Include all Scopes** checkbox where they will find all study scopes.

Study Time

Start Procedural Time

15:58

End Procedural Time

16:13

Duration (hh:mm:ss)

00:15:00

Scope List

▼

☐ Include all Scopes

Scope Time

Scope Name / Serial No	Cecal Intubation Time	Withdrawal Time
Fuji702 / 702	0	0

Save Changes

Note: The user can also edit the start and end procedural time.

Indication/Diagnosis

Sedation

Bowel Preparation

Timing

Team

Colonoscopy

Study Time

Start Procedural Time

16:04

End Procedural Time

16:19

Duration (hh:mm:ss)

00:15:00

16 : 04

▼

☐ Include all Scopes

Save Changes

Team

The user can add more teams from this section. They can select the team member from each of the dropdowns.

The user can choose the findings of the patient by

- **Lesion**
- **Findings**
- **Findings DTL**
- **Intervention**
- **Complication**

The user will find detailed options in every finding. It will help to make a detailed report for the patient.

The user can add a multi-region and repeat the findings of the patient. After selecting the **Lesion**, the user can add multiple findings by clicking this  button. It will help the user to add more findings of the patient if needed.

+ Lesion Findings Findings DTL Intervention Complication

Stomach

1) Atrophy 2) Bezoar 3) Bezoar

Template : Select Template

Finding 1	Number	Size	Region 1	Region 2
- Angioectasia - Atrophy - Bezoar - Cautic injury - Dieulafoy's lesion - Dilation			☐ ☐ Anastomosis site ☐ Anterior wall of body ☒ Antrum ☐ Body ☐ Cardia	

Findings :
Stomach:
 1) Atrophy at Antrum stomach.
 2) Bezoar
 3) One Bezoar size 1.5 cm

Images

The user can add this report as a template by simply clicking  this button. After that, the user will find a pop up page where they have to input template name, access type, and description.

X

Add Template

Study Name

EGD

Template Name

p1

Access Type

Private

▼

Description

multiple findings-1

×

Close

Save

Finally, the user have to click the **Save** button.

Template :

Select Template

▼

Region 2

Q

Q

Custom Tem

×

EE

×

Tem-110101

×

p1

×

Lesion Findings Findings DTL Intervention Complication

Ileocecal valve Colon

Mass

Study: Colonoscopy Template: Select Template

Finding 1	Number	Size	Region 1	Region 2
Mass	1	Custom	☐ Appendiceal orifice	
Melanosis	2	0.5 cm	☐ Cecum	
Metal clip	3	1.5 cm	☒ Ascending colon	
Nodular mucosa	>3	2 cm	☐ Transverse colon	
Nodule		2.5 cm	☐ Descending colon	
Parasite		3 cm	☐ Sigmoid colon	
Perforation		3.5 cm	☐ Rectosigmoid junction	
Polyp		4 cm	☐ Rectum	
Previous surgery		4.5 cm	☐ Diffused	

Findings :

[Ileocecal valve]:

1) More Than Three pedunculated polyp(s) size 1.5 cm at Diffused Colon, Bleeding -Yes.

Intervention :

ESD/Result: En-bloc resection .

EMR/Technique: Banding . Result: complete piece-meal resection . Defect closure: Yes .

Biopsy forceps .

Type: Biopsy .

[Colon]:

1) More Than Three Mass size 1.5 cm at Ascending colon .

Complications: Yes:

Complication Type: Bleeding, Location: Sphincterotomy Site.

Intervention: Yes:

Type: Injection, Result: Success.

Type: Hemostasis, Result: Success.

Moreover, the user can save the template of the reporting for similar types of patients.

Add Template

Study Name: Colonoscopy

Template Name: Temp1

Access Type: Private

Description:

Close Save

The user will get a confirmation message after clicking the **Save** button.

The user will find this template in the Template dropdown.

Lesion Findings Findings DTL Intervention Complication

Terminal ileum Ileocecal valve Colon Anus

Mass

Study: Colonoscopy

Template: Select Template

Finding 1	Number	Size	Region 1
Mass	1	Custom	☐ Appendiceal orifice
Melanosis	2	0.5 cm	☐ Cecum
Metal clip	3	1.5 cm	☒ Ascending colon
Nodular mucosa	>3	2 cm	☐ Transverse colon
Nodule		2.5 cm	☐ Descending colon
Parasite		3 cm	☐ Sigmoid colon
Perforation		3.5 cm	☐ Rectosigmoid junction
Polyp		4 cm	☐ Rectum
Previous surgery		4.5 cm	☐ Diffused

Comment

Findings :

ileum]: Normal

ervention :
ampling(random) at; Ileocecal valve: 2, Transverse colon:

valve]: Normal

ervention :
ampling(random) at; Terminal Ileum: 2, Ascending colon:

3.

[Anus]: Normal

Intervention :
Sampling(random) at; Terminal Ileum: 3, Ileocecal valve: 4

Complications: Yes;
Complication Type: Bleeding; Location: Sphincterotomy Site;
Intervention: Yes;
Type: Injection; Result: Success;
Type: Hemosty: Result: Success;

Moreover, the user can select Sedation and Bowel Preparation for Reporting. The user can choose Bowel Preparation only in Colonoscopy.

Detected

Lesion Findings Findings DTL Intervention Complication

Ileocecal valve

Study: Colonoscopy

Findings :

[Ileocecal valve]:

1) More Than Three pedunculated polyp(s) size 1.5 cm at Diffused Colon. Bleeding -Yes.

Intervention :
EMR; Technique: Banding , Result: complete piece-meal resection. Defect closure: Yes .
Biopsy forceps .
Type: Biopsy .
ESD; Result: En-bloc resection .

Complications: Yes;
Complication Type: Bleeding; Location: Sphincterotomy Site;
Intervention: Yes;
Type: Injection; Result: Success;
Type: Hemosty: Result: Success;

Images

Detected

Add Sedation

Local xylocaine throat spray ☒ Yes


Sedation ☒ Yes

Type of Sedation Conscious Sedation


Sedative Agent Pethidine Amount 2 mg


+ Add

Close Update

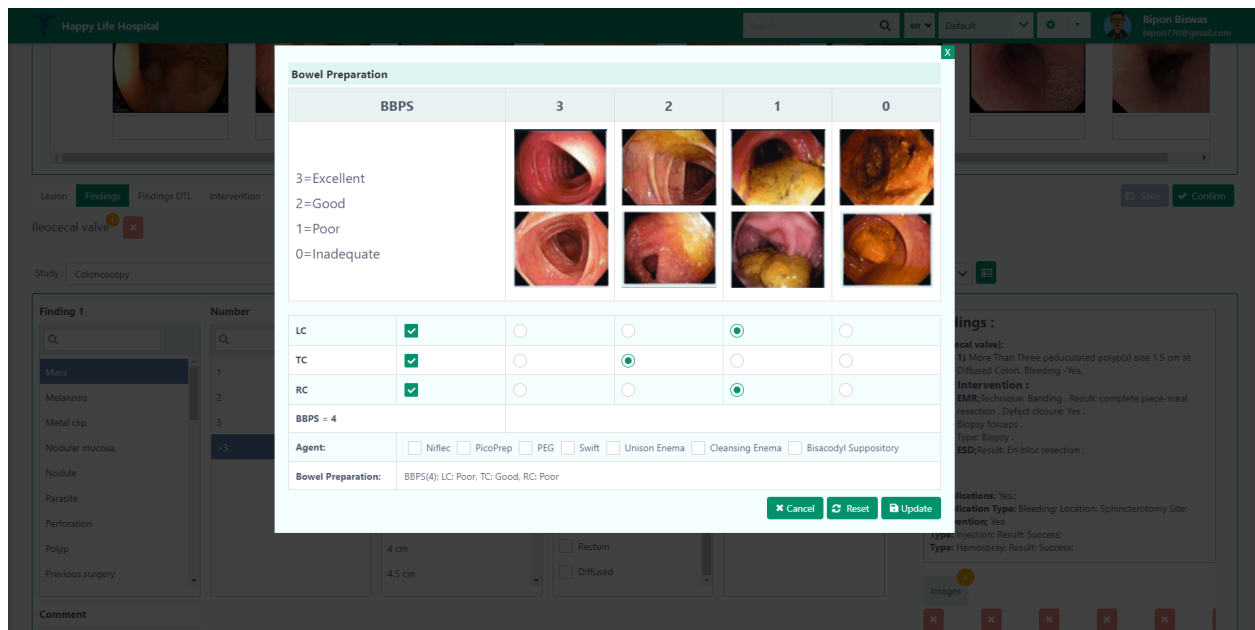



Finally, the user has to select the **Confirm** button, after completing all the procedure.

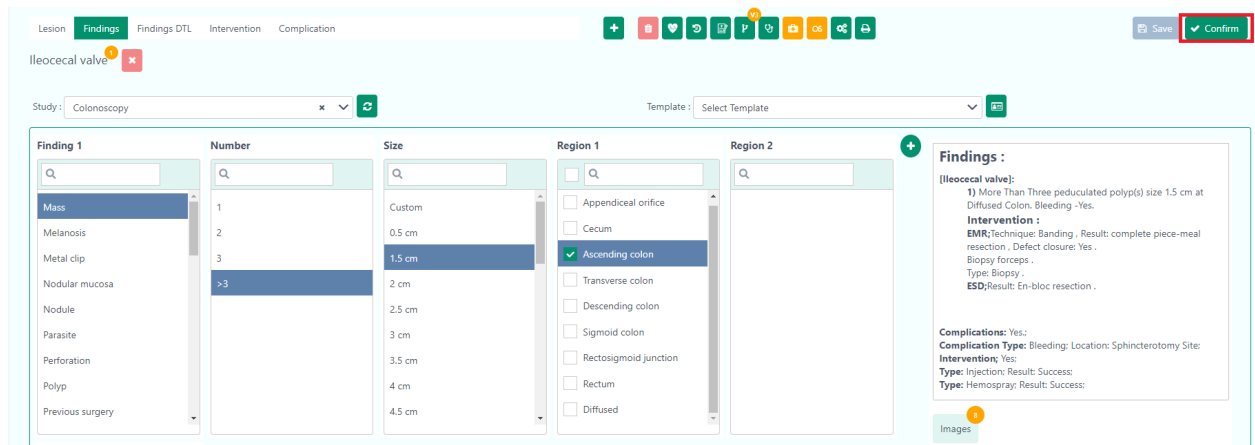

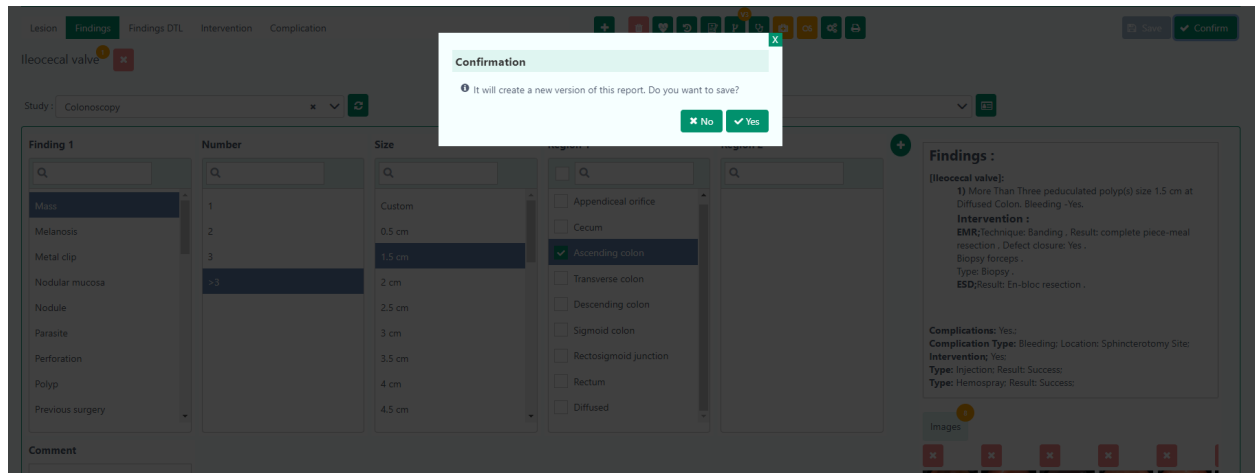

To see the report, the user has to click the Preview Report button.

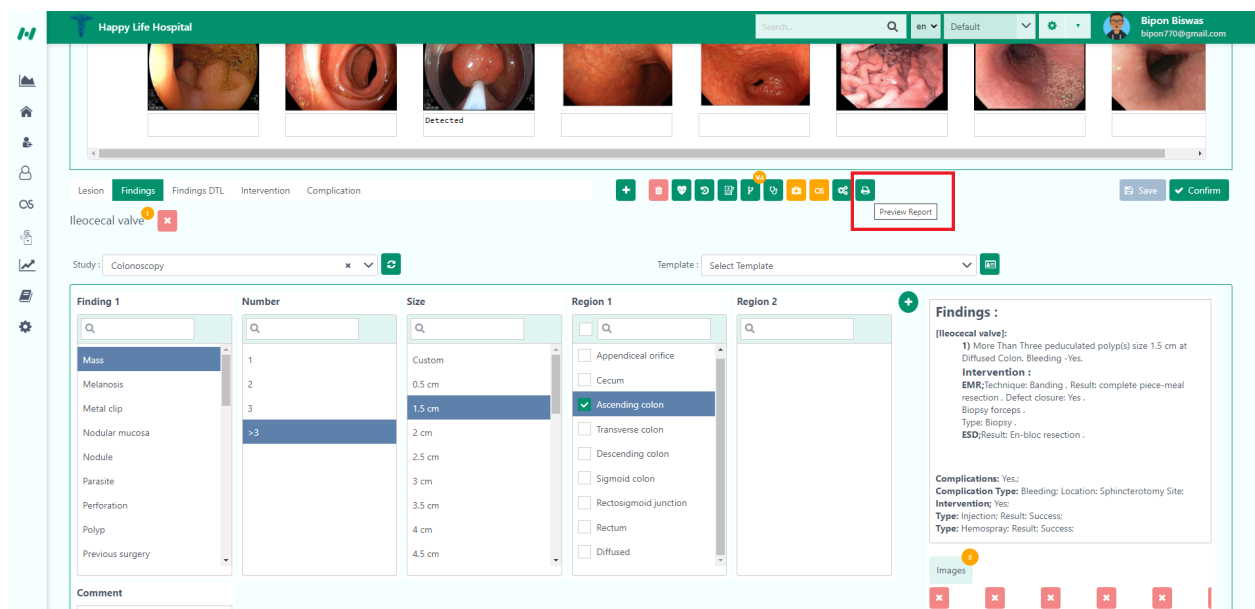

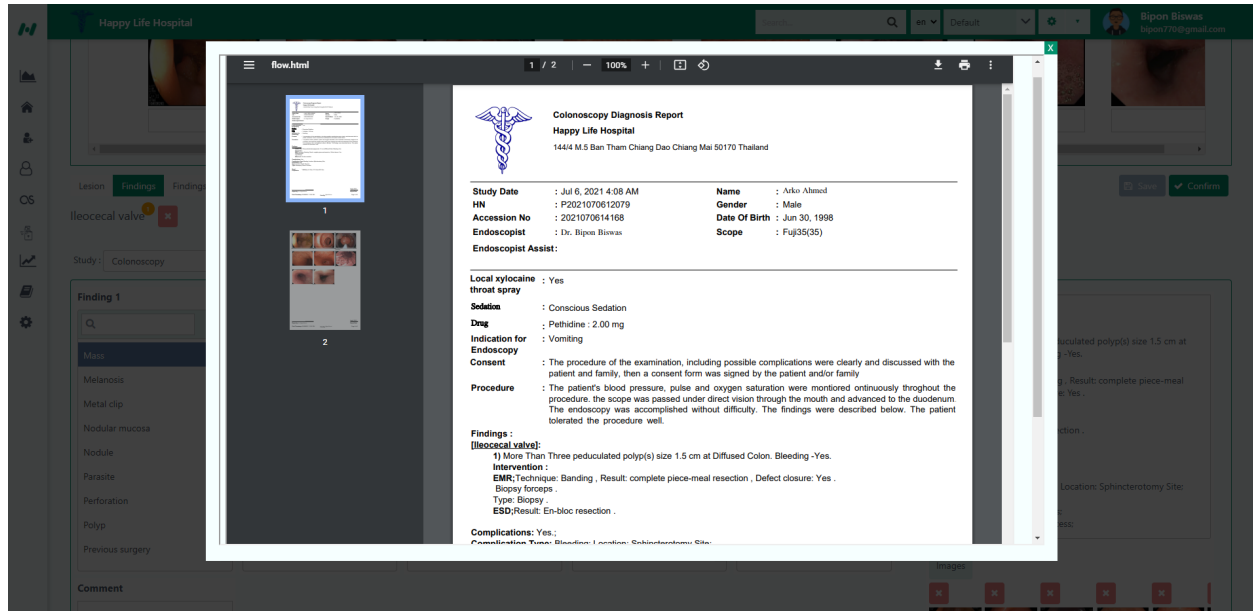

Updated: 3/2/2021

Completed:

Nursing

Inside the Nursing, you can **view statuses** of a patient, what **stage** they are in, and in general, **the whole process/cycle** that a patient will go through from the moment s/he discharged until the completion of the nurses.

Nursing is organized with two tabs that will carry a patient record form one tab to the next, until the client has completely discharged.

Pre Scope Worklist

Pre-Scope Worklist

Post-Scope Worklist

Search By :

HN

Date Range

Advanced Search

From Date

To Date

02-03-2021

02-03-2021

Search

Today

Tomorrow

Next 3 Days

Next 7 Days

10 columns selected

Global Filter

		Accession No	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Date	Action
		20210302138	MAX MAX	P20200427-10723	Bipon Biswas	Endoscopy	ERCP	Stat	Imaging	Room 104	02-03-2021 10:28 AM	

1

20

Updated-2/3/2021

Pre-Scope Worklist:

When the user clicks the **Nursing** button, this page is shown as a default page. When the patient registered, the user can see the data of that patient in this field.

Pre Scope Worklist

Pre-Scope Worklist

28

Post-Scope Worklist

Search By : ☐ HN ☒ Date Range ☐ Advanced Search

From Date

To Date

01-02-2021

08-03-2021

Search

Today

Tomorrow

Next 3 Days

Next 7 Days

10 columns selected

Global Filter

	Accession No	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Date	Action
	20210201138	Mahfuz05Pet Mahfuz05Pet	P20200505-11485	Doctor Admin	Endoscopy	Upper Endoscopy	Routine		Room 102	01-02-2021 11:02 AM	
	20210201138	MAX MAX	P20200427-10723	Doctor User	Endoscopy	Upper Endoscopy	Routine		Room 104	01-02-2021 04:17 PM	

From the **Pre-Scope** the user can-

- Prepare PreScope.
- Transfer to Image
- Chief Complaint and Diagnosis.

	Accession No	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Date	Action
	20210201138	Mahfuz05Pet Mahfuz05Pet	P20200505-11485	Doctor Admin	Endoscopy	Upper Endoscopy	Routine		Room 102	01-02-2021 11:02 AM	

Prepare PreScope

The first option is **Prepare PreScope**, that preps the patient for the endoscopy-

10 columns selected

Global Filter

	Accession No	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Date	Action
	20210201138	Mahfuz05Pet Mahfuz05Pet	P20200505-11485	Doctor Admin	Endoscopy	Upper Endoscopy	Routine		Room 102	01-02-2021 11:02 AM	
	20210201138	MAX MAX	P20200427-10723	Doctor User	Endoscopy	Upper Endoscopy	Routine		Room 104	01-02-2021 04:17 PM	

This takes you to the **Pre-Preparation** section:

Pre-Preparation



Clear All

Save

Save & Transfer

Patient Information

Name : Carolina Statham

Phone : 9635474857

Address :



Doctor Information

Name : Miracle Admin

Phone : +9999999999

Department : Physical Therapy Clinic



Patient Information

Clear

Procedure/Study	PEG	Status	OPD
Diagnosis	Mitral insufficiency, Gastric ulcer	Indication	+
Purpose	Suspected tuberculosis	Finding	Intraoral bleeding, Abdominal rigidity

This section contains a number of forms that are not collapsible. If you keep scrolling, you will find all the options that you can select/fill in. The sections include:

- Patient Information

Patient Information

Clear

Procedure/Study	PEG	Status	OPD
Diagnosis	Mitral insufficiency, Gastric ulcer	Indication	+
Purpose	Suspected tuberculosis	Finding	Intraoral bleeding, Abdominal rigidity
Nursing Class	Requires occasional nursing	Region	+
Chief Complain	Dysphagia		

Clicking on the plus sign lets you enter text information:

Select Indication

X

Indication

B **I** **U**

Fever, Nausea

Cancel

Save

- Check List

Check List				Clear	
Confirm Patient Name	☒ Yes	Confirm Procedure	☒ No	Consent Form Sign	☒ Yes
NOP	☒ No	Antiplatelet/Anticoagulant	☒ Yes	NO PPI	☒ Yes
Companion	☒ No	Remove False Teeth/Acc	☒ No		
Attachment	☐ OPD Card ☐ Other ☐ Film				
Premed	☒ Domicum ☐ Pethidine				

The default choice for all items in the list is set to 'No'. Clicking on it will set it to 'Yes'.

- Patient History

Patient History				Clear	
Drug Allergy	☒ No	Egg Soybean Allergy	☒ No	Previous Surgery	☒ Yes
Pregnancy	☒ No	Pacemaker/ICD	☒ No	Special Equipment	☒ No
Prophylactic Antibiotic	☒ Yes	Difficult Airway/Aspiration	☒ Yes	Anxiety	☒ No
Smoking	☒ No				
Alcohol drinking	☐ Regular ☐ Social ☐ No				
Comorbidities					
None	☒ No	DM	☒ Yes	HT	☒ No
Heart Disease	☒ No	Cirrhosis	☒ No	CRF	☒ No
Bleed Disorder	☒ No	Other	☒ No		
Current Medication					
None	☒ No	Aspirin	☒ No	NSAIDs	☒ No
Plavix	☒ No	Warfin	☒ Yes	Icon	☒ No
NOACs	☒ No	Other	☒ No		
LOC					
Alert	☒ Yes	Confuse	☒ No	Drowsy	☒ No
Stuporous	☒ Yes	Semicoma	☒ No	Coma	☒ No

This is the same as **CheckList**, the default choice for all items in the list is set to 'No'. Clicking on it will set it to 'Yes'.

- Breathing

Breathing						Clear
Spontaneous	☒ Yes	Room Air	☒ Yes	O2 Cannula	☒ No	
O2 Mask C Bag	☒ No	Artificial	☒ No	T Tube	☒ Yes	
ET Tube	☒ No	Ambu Bag	☒ Yes	Ventilator	☒ No	

- Blood Component

Blood Component						Clear
Red Cell(Unit)	5	☒ Yes	Plasma(Unit)	12	☒ Yes	
CRYO(Unit)		☒ No			☒ No	

This is also the same as **Check List**, the default choice for all items in the list is set to '**No**'. Clicking on it will set it to '**Yes**'. When you click '**Yes**', the user can give input. Otherwise, the user can't give input if it is '**No**'.

- Complete Blood Count

Complete Blood Count						Clear
Hct	7	Plt	7			
PT	4	INR	4			
PTT	6.8	CBC Date	08-09-2020			

- Blood Chemistry

Blood Chemistry						Clear
Bun	3	Cr	76.4			
Na	5.8	K	5			
Cl	10.7	CO2	5.4			
Chem Date	02-09-2020					

- Liver Function Test

Liver Function Test

Clear

TB	67.5	DB	23.9
SGOT	0	SGPT	457
ALP	12.9	ALB	98.56
LFT Date	16-09-2020	Other Labs	
Labs Date	11-09-2020		

- Check List

Check List

Clear

Confirm Name Procedure	☒ Yes	Confirm Correct HN	☒ Yes	Patient Signed Consent	☒ Yes
Antibiotic before procedure	☒ No	No Antibiotic before procedure	☒ Yes	Documents Check	☒ Yes
OPD Card	☒ Yes	Other Check	☒ No	Film	☒ No
Critical Event	☒ No				

- Nurse Info

Nurse Info

Alert Message

Nurse Sign Out:

Please Select Nurse

Nurse User

Set Alert Message:

☒ No

The nurse can Set Alert Message for the doctor if needed.

Alert Message


Set Alert Message:

☒ Yes

The user has to enable **Yes** to write the alert message for the patient.

When the nurse provides the Alert message, the user can see it in the left corner.

10 columns selected												Global Filter	
	Accession No	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Date	Action		
	20210201138	Mahfuz05Pet Mahfuz05Pet	P20200505-11485	Doctor Admin	Endoscopy	Upper Endoscopy	Routine		Room 102	01-02-2021 11:02 AM			

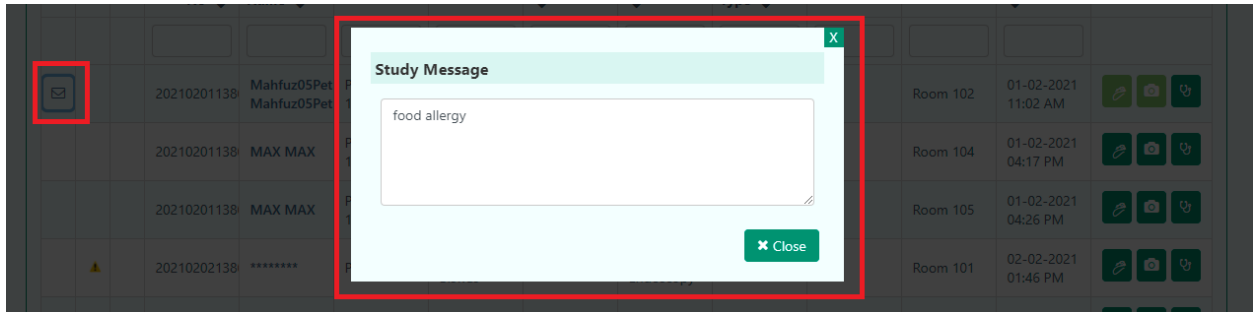

Transfer to Image

The second option is ***Transfer to Image***, that takes the patient to **Imaging-**


10 columns selected

Global Filter

		Accession No	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Date	Action
		20210201138	Mahfuz05Pet Mahfuz05Pet	P20200505-11485	Doctor Admin	Endoscopy	Upper Endoscopy	Routine		Room 102	01-02-2021 11:02 AM	
		20210201138	MAX MAX	P20200427-10723	Doctor User	Endoscopy	Upper Endoscopy	Routine		Room 104	01-02-2021 04:17 PM	

Transfer to Image

Note: When the user presses the ‘Yes’ button, this will take to **Imaging** section. Above you will find the details of **Imaging**.

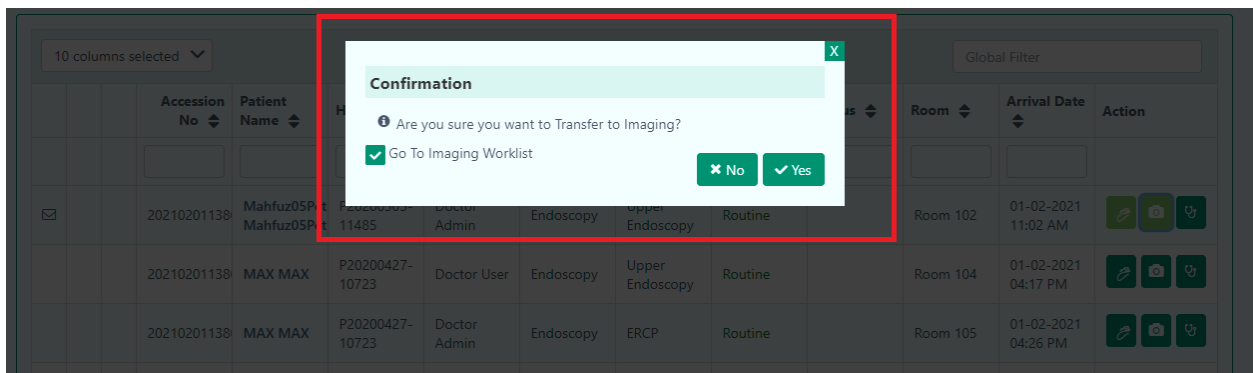

When the Imaging part is completed, the user can see that patient details in the **Post-Scope Worklist**.

Updated: 8/3/2021

Post-Scope Worklist

After completing the **Imaging** of the patient, the user can start doing **Post-Scope Worklist**.

In Post-Scope Worklist, the user can

- Go to Post Scope
- Discharge

Go to Post Scope

By clicking the button, the user will find the **Post Scope** page.

10 columns selected												Global Filter
		Accession	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Time	Action
		20210305138	Patient 101	P2021020812	Sureeporn Jangsirikul	Endoscopy	Colonoscopy	Routine	Reporting	Room 103	05-03-2021 09:07 AM	
		20210305138	Patient 102	P2021020812	Sureeporn Jangsirikul	Endoscopy	Upper Endoscopy	Routine	Reporting	Room 102	05-03-2021 08:48 AM	

Post ScopeBack

Patient : Patient 101**Age** : 39yrs

HN : P2021020812014**SSN** :

Allergy : No

Study Information

Code : 21172

Name : Colonoscopy

Datetime : 05-03-2021

Doctor Information

Name : Sureeporn Jangsirikul

Phone :

Department :

Time in Recovery Room

Recorded Nurse
Dr. Doctor Admin

Recorded Date
09/03/2021

Recorded Time
11:49

Then, the user can check the patient information in the **Vital Sign** option where they can **Record** basic information of the patient for giving discharge or further procedure.

Vital Sign

Time(HH:mm | 24 hour format) Free Text Record

Global Filter

Action	Time	Free Text	Respiration	Consciousness	SBP	O2Sat

When the patient enters to Post Scope, the timing is also captured in the **Vital Sign**.

Patient Name : Patient 101
HN : P2021020812014
Gender :
DOB : February 9, 1962

Save

Respiration	86	Consciousness	no
SBP	32	O2Sat	
Height (cm)	168	Weight (kg)	86
BMI	30.47	BP (mmHg)	Systolic / Diastolic
SpO ₂	96	Heart Rate (BPM)	86
Temperature (F)	98	Temperature (C)	
Respiration Rate (RPM)	90	Pulse rate (BPM)	
Capillary Refill Time (Sec)	5	BCS	
Pain Score	11	Hydration Status	

Close Save

After clicking the **Save** button the user will find the information in the **Vital Sign**.

Note: The user can update the patient record or delete the existing one from the **Vital Sign**.

Vital Sign

Time(HH:mm | 24 hour format) Free Text Record

Global Filter

Action	Time	Free Text	Respiration	Consciousness	SBP	O2Sat	Height
✓ ✕	11:49		86	no	32		168

1 20

After successfully completing the **Vital Sign**, the user will find another section Total score for discharge. The patient will be discharge when the score is above or equal 9.

The point will count by checking the below requirements.

- Motor
- Respiration
- Consciousness
- SBP
- O2Sat

Motor

The screenshot shows the 'Motor' assessment section. At the top, it says 'Total Score ≥ 9 for Discharge'. Below this, there's a 'Time' field with '12:40:07' and a 'Point' field with '2'. The 'Motor' dropdown menu is open, showing options: 'Able to move arms and legs', 'Able to move arms and legs but quite weak', and 'Unable to move'. The 'Respiration' dropdown is set to 'Select Respiration', 'Consciousness' to 'Select Conscious', 'SBP' to 'Select SBP', and 'O2Sat' to 'Select O2Sat'. There is an '+ ADD' button. At the bottom right, there are buttons for 'Back', 'Save', and 'Save & Discharge'. The Fujifilm logo is in the bottom right corner.

When the user selects the options according to the condition of the patient, the point will be calculated automatically.

Respiration

The screenshot shows the 'Respiration' assessment section. At the top, it says 'Total Score ≥ 9 for Discharge'. Below this, there's a 'Time' field with '12:45:13' and a 'Point' field with '3'. The 'Motor' dropdown is set to 'Able to move arms and legs'. The 'Respiration' dropdown menu is open, showing options: 'Able to breath', 'Able to breath and cough', 'Able to breath', and 'Need device to maintain airway'. The 'Envision NEXT' logo is in the bottom left corner.

Consciousness

The screenshot shows the 'Consciousness' assessment section. At the top, it says 'Total Score ≥ 9 for Discharge'. Below this, there's a 'Time' field with '12:47:10' and a 'Point' field with '5'. The 'Motor' dropdown is set to 'Able to move arms and legs', and the 'Respiration' dropdown is set to 'Able to breath'. The 'Consciousness' dropdown menu is open, showing options: 'Fully Conscious', 'Need to be waken to be Conscious', and 'Not Conscious'. The 'Envision NEXT' logo is in the bottom left corner.

SBP

Total Score ≥ 9 for Discharge

Time: 12:49:36

Point: 7

Motor: Able to move arms and legs

Respiration: Able to breath

Consciousness: Fully Conscious

SBP: Between ± 20

Envision NEXT

O2Sat

Total Score ≥ 9 for Discharge

Time: 12:51:18

Point: 9

Motor: Able to move arms and legs

Respiration: Able to breath

Consciousness: Fully Conscious

SBP: Between ± 20

O2Sat: $\geq 92\%$ Room air

Envision NEXT

After successfully selecting all the sections, the user has to press **Add** for further procedure.

Total Score ≥ 9 for Discharge

Time: 12:55:34

Point: 0

Motor: Select Motor

Respiration: Select Respiration

Consciousness: Select Conscious

SBP: Select SBP

O2Sat: Select O2Sat

+ ADD

Post Scope Time	Point	Motor	Respiration	Consciousness	SBP	O2Sat	
12:55	9	Able to move arms and legs	Able to breath	Fully Conscious	Between ± 20	$\geq 92\%$ Room air	

Back Save Save & Discharge

If the user wants to observe the patient more, they can simply press **Save** button. Otherwise, the user can press the **Save & Discharge** button.

Note: When the discharge point is below 9, the user can give conditional discharge.

Confirmation

Total point is ≤ 9 . Are you sure you want to discharge the patient?

No Yes

Total Score ≥ 9 for Discharge

Time: 14:31:43

Point: 0

Motor: Select Motor

Respiration: Select Respiration

Consciousness: Select Conscious

SBP: Select SBP

O2Sat: Select O2Sat

+ ADD

Post Scope Time	Point	Motor	Respiration	Consciousness	SBP	O2Sat	
14:25	7	Able to move arms and legs but quite weak	Able to breath and cough	Not Conscious	Between ± 20	$\geq 92\%$ Room air	

Back Save Save & Discharge

The user has to press ‘Yes’ to open the conditional discharge form where the user provides the instructions for discharging the patient.

Discharge

Transfer To* Home Ward Critical Care Unit

Ambulatory Status* Ambulate WheelChair Stretcher Other

Note

Have to take full bed rest for 7days, Antibiotic course will continue till this days.

Cancel Discharge

After pressing the discharge button, the user will get confirmation box to proceed.

Confirmation

Are you sure to discharge the patient?

No Yes

Post Scope Worklist Pre-Scope Worklist Post-Scope Worklist

Search By : ☐ HN ☒ Date Range ☒ All

From Date: 01-02-2021 To Date: 09-03-2021 Search Today Yesterday Last 3 Days Last 7 Days

10 columns selected Global Filter

	Accession	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Time	Action
	2021030513882	Patient 101	P2021020812014	Sureeporn Jangsirikul	Endoscopy	Colonoscopy	Routine	Completed	Room 103	05-03-2021 09:07 AM	
	2021030513881	Patient 102	P2021020812015	Sureeporn Jangsirikul	Endoscopy	Upper Endoscopy	Routine	Reporting	Room 102	05-03-2021 08:48 AM	
	2021030413880	Patient 101	P2021020812014	Santi Kulpacharapong	Endoscopy	ERCP	Routine	Reporting	Room 101	04-03-2021 09:27 PM	
	2021030413878	MAX MAX	P20200427-10723	Bipon Biswas	Endoscopy	Sigmoidoscopy		Completed	Room 104	04-03-2021 09:16 PM	
	2021030413877	MAX MAX	P20200427-10723	Bipon Biswas	Endoscopy	Sigmoidoscopy	Routine	Reporting	Room 103	04-03-2021 07:47 PM	
	2021030413875	Dipu Siam	P2021020412011	Doctor Admin	Endoscopy	Colonoscopy	Routine	Reporting	Room 113	04-03-2021 07:47 PM	

Discharged Successfully

Finally, the user will get successful notification after discharging the patient. The button of discharge become disable when the patient successfully discharged by the user.

10 columns selected												Global Filter	
			Accession	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Time	Action
			20210305138	Patient 101	P2021020812	Sureeporn Jangsirikul	Endoscopy	Colonoscopy	Routine	Completed	Room 103	05-03-2021 09:07 AM	 

Updated:9/3/2021

Device Management

This section will help the user to maintain the device. The Device Management contains-

- Scope Management
- Endoscopy Management
- Washing Machine Management

Envision NEXT

WRITE YOUR OWN

Patient

Nursing

Device Management

Scope Management

Endoscopy Management

Washing Machine Management

Happy Life Hospital

Search...

en ▼

Default ▼

⚙️ ▼

👤

Doctor Admin

qa.admin@miracle.com

Accession	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Status	Room	Arrival Time	Action
202103051	Patient 101	P20210206	Sureeporn Jangsirikul	Endoscopy	Colonosco	Routine	Completed	Room 103	05-03-2021 09:07 AM	
202103051	Patient 102	P20210206	Sureeporn Jangsirikul	Endoscopy	Upper Endoscopy	Routine	Reporting	Room 102	05-03-2021 08:48 AM	

Scope Management

In the **Scope Management**, the user will find the Scope lists.

Scope Lists

+ Add Scope

< Back

Scope Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Fuji 001	12345	Coloplast		10-09-2020 12:00 AM	04-03-2021 09:20 PM	Dirty		
Fuji55	55	Fujifilm		04-02-2021 12:00 AM	04-03-2021 06:18 PM	Ready		
Fuji35	35	Coloplast		13-01-2021 12:00 AM	25-02-2021 02:52 PM	Ready		
Fuji32	32	Olympus		10-01-2021 12:00 AM	25-02-2021 11:51 AM	In Use		

In the Scope List, the user can know the barcode, details of the scope and can active or inactive the scope from this list.

Scope Lists

+ Add Scope

< Back

Scope Barcode

fsfsfsfsf

Scope Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Fuji 001	12345	Coloplast		10-09-2020 12:00 AM	04-03-2021 09:20 PM	Dirty		
Fuji55	55	Fujifilm		04-02-2021 12:00 AM	04-03-2021 06:18 PM	Ready		

Fuji55	55	Fujifilm		04-02-2021 12:00 AM	04-03-2021 06:18 PM	Ready	Details
--------	----	----------	--	---------------------	---------------------	-------	--------------------

The user can easily see the scope details where the user can **Edit** Scope Dashboard.

Scope Dashboard

< Back

Scope Name : Fuji55

Purchase Date : 05-01-2021

First Use Date : 04-02-2021

Manufacturer : Fujifilm

Total Used : 2 times

Serial No : 55

Status : Ready

Last Use Date : 04-03-2021

Notes :

Total Disinfected : 2 times

Edit

Scope Update

Scope Picture

Close

Save

Scope Update

Scope Name*

Fuji55

Serial No*

55

Manufacturer

Fujifilm

✕

▼

Purchase Date

05-01-2021

📅

First Use

04-02-2021

📅

Notes

Notes

Barcode

23sd10

23sd10

Print

Close

Save

The user also finds **Study Information** and can **Add** new Study Information in Scope Dashboard.

Study Information

+ ADD

Study Code	Study Name	Action
Upper Endoscopy 2	Upper Endoscopy	
COLONOSCOPY	Colonoscopy	
ERCP 2	ERCP	

+ ADD

Study Registration

Study Name

Upper Endoscopy

▼

Study Code	Study Name	Action
Upper Endoscopy 2	Upper Endoscopy	

Close

Save

In the **History** area, the user can find Usage, Disinfection and Maintenance.

History

Usage Disinfection Maintenance

Usage History

Global Filter

Patient Name	Accession No	Study	Device Status	Doctor Name	Room	Start Time	End Time	Duration
Dipu Siam	2021030413875	Colonoscopy	Used	Doctor Admin	Room 113	04-03-2021 06:18 PM	04-03-2021 06:18 PM	00:00:08
Sud Ratana	2021021013838	Upper Endoscopy	Used	Mostafa Chowdhury	Room 106	10-02-2021 03:59 PM	10-02-2021 04:00 PM	00:00:01

1 20

The user will find the used scope in this **Usage** field and disinfection scope information in **Disinfection** part.

History

Usage Disinfection Maintenance

Disinfection History

Global Filter

Machine	Study	Device Status	Started By	Ended By	Start Time	End Time	Duration	Comments
Washing 102	Colonoscopy	Disinfected	Doctor Admin	Doctor Admin	04-03-2021 09:19 PM	04-03-2021 09:19 PM	00:00:03	


The user can simply **Active** or **Inactive** the Scope by clicking this button. While clicking the button, the user will get confirmation messages.

Scope Lists

+ Add Scope < Back

Global Filter

Scope Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Fuji 001	12345	Coloplast	Barcode	10-09-2020 12:00 AM	04-03-2021 09:20 PM	Inactive		Active

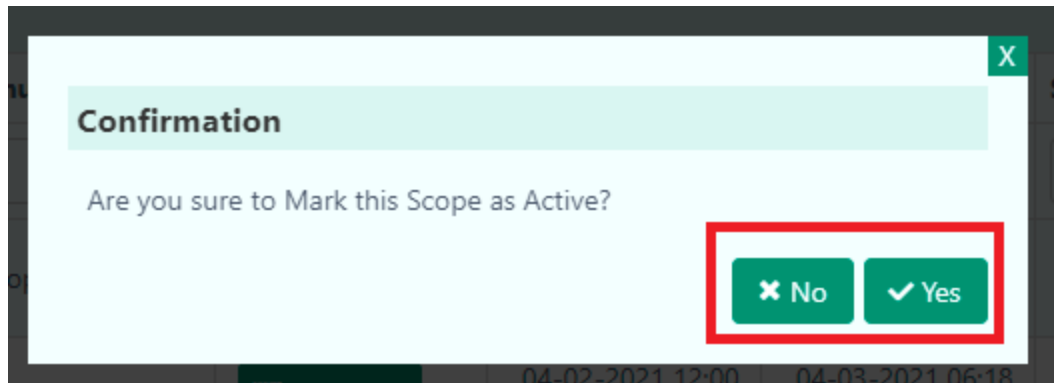

Scope Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Fuji 001	12345	Coloplast		10-09-2020 12:00 AM	04-03-2021 09:20 PM	Dirty		
Fuji55	55	Fujifilm		04-02-2021 12:00 AM	04-03-2021 06:18 PM	Ready		
Fuji35	35	Coloplast		13-01-2021 12:00 AM	25-02-2021 02:52 PM	Ready		
Fuji32	32	Olympus		10-01-2021 12:00 AM	25-02-2021 11:51 AM	In Use		
Fuji50	50	Fujifilm		01-02-2021 12:00 AM	24-02-2021 12:15 PM	Ready		
Fuji30	30	Boston Scientific		04-01-2021 12:00 AM	23-02-2021 01:08 PM	Ready		

Scope is activate successfully.

Note: When the user active the scope, the scope is ready for use in Imaging. However, when the user inactive the scope, they can't use this machine.

Note: If the user wants to active or inactive the machine, they can simply click the (Tick) button. When the button is blur, it will indicate as Inactivate. On the other hand, when the button is green, it indicates active mode.

Fuji 10001	1212121213312	R. Bard		09-03-2021 12:00 AM	10-03-2021 05:22 PM	Inactive		
Fuji50	50	Fujifilm		01-02-2021 12:00 AM	09-03-2021 03:52 PM	Ready		

The user can add new scope through Scope lists. After clicking the **Add Scope**, it will take them to the **Scope Registration** page.

Scope Lists

+ Add Scope

< Back

Global Filter								
Scope Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Fuji 001	12345	Coloplast	 Barcode	10-09-2020 12:00 AM	04-03-2021 09:20 PM	Dirty		 
Fuji55	55	Fujifilm	 Barcode	04-02-2021 12:00 AM	04-03-2021 06:18 PM	Ready		 

Scope Registration

Scope Picture



Close Save

Scope Information

Scope Name*

Scope Name

Manufacturer

Select manufacturer

First Use

Barcode

Barcode

Serial No*

Serial No

Purchase Date

Notes

Notes

Print

Study Information

Study Name

Choose

Study Type

Study Code

Study Name

Action

Close Save

Mandatory fields are highlighted in red*

The user has to press **‘Save’** after completing the registration. Finally, the user will get a confirmation message.

Scope Information

Scope Name*

Fuji3102021

Manufacturer

Fujifilm

First Use

Barcode

fuji3102021

Serial No*

3102021

Purchase Date

10-03-2021

Notes

New Equipment

Print

Q

Upper Endoscopy

Colonoscopy

ERCP

Sigmoidoscopy

Study Information

Study Name

Study Type

Endoscopy

Study Code

Upper Endoscopy 2

Study Name

Upper Endoscopy

Action

Close

Save

Happy Life Hospital

Search...

en

Default

Doctor Admin

qa.admin@miracle.com

Scope Lists

+ Add Scope

< Back

Scope Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Fuji 001	12345	Coloplast	Barcode	10-09-2020 12:00 AM	04-03-2021 09:20 PM	Dirty		
Fuji55	55	Fujifilm	Barcode	04-02-2021 12:00 AM	04-03-2021 06:18 PM	Ready		
Fuji35	35	Coloplast		13-01-2021 12:00 AM	25-02-2021 02:52 PM	Ready		
Fuji32	32	Olympus		10-01-2021 12:00 AM	25-02-2021 11:51 AM	In Use		
Fuji50	50	Fujifilm	Barcode	01-02-2021 12:00 AM	24-02-2021 12:15 PM	Ready		
Fuji30	30	Boston Scientific		04-01-2021 12:00 AM	23-02-2021 01:08 PM	Ready		
Fuji36	36	R, Bard		02-01-2021 12:00 AM	19-02-2021 06:52 PM	Assigned		
Fuji53	53	Fujifilm	Barcode	04-02-2021 12:00 AM	18-02-2021 07:08 PM	Assigned		
Fuji54	54	Fujifilm	Barcode	04-02-2021 12:00 AM	10-02-2021 05:32 PM	Ready		
Fuji 002	stsfst	Johnson & Johnson	Barcode	28-08-2020 12:00 AM	10-02-2021 03:19 PM	In Use		
dgdsdgdsd	5252525	Johnson & Johnson	Barcode	02-10-2020 12:00 AM	09-02-2021 10:48 AM	inactive		

Scope information is created successfully.

The user can use this new scope in the **Imaging** section.

Endoscopy Image Capture

Selected Scope : Fuj32

Scope In

HN : P20200427-10723

Name : MAX MAX

Gender : Male

LABORATORY WEB

Name, Job Title, Department

Name Job Title Department Action

Scope List

Fuj3102021

Barcode, Scope Name

Barcode Scope Name Action

Fuj32

fuj3102021 Fuj3102021

Endoscopy Machine

Endoscopy Machine List

Barcode, Machine Name

Endoscopy Management

In the **Endoscopy Management**, the user will find the Endoscopy Machine Lists.

Endoscopy Machine Lists								
Global Filter								
Machine Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Endo Machine 1002	12121323312		Barcode	07-09-2020 05:54 PM	07-09-2020 05:55 PM	Inactive		
Test Machine	0989090898					Inactive		
Test Machine 2	010101010101	Johnson & Johnson	Barcode	10-09-2020 12:00 AM		Ready	My sample notes	

In the Endoscopy Management Lists, the user can see the Barcode, details of the Endoscopy and active or inactive the machine for use.

Endoscopy Machine Lists								
+ Add Endoscopy Machine ◀ Back								
Global Filter								
Machine Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Endo Machine 1002	12121323312		 Barcode	07-09-2020 05:54 PM	07-09-2020 05:55 PM	Ready		 

Updated-10/3/2021



Note: When the user activates the Endoscopy machine, the machine is ready for use in Imaging. However, when the user inactive the Endoscopy machine, they can't use this machine.

Note: If the user wants to activate or inactive the machine, they can simply click the (Tick) button. When the button is blurred, it will indicate as Inactive. On the other hand, when the button is green, it indicates active mode.

Endo Machine 1002	12121323312		 Barcode	07-09-2020 05:54 PM	07-09-2020 05:55 PM	Ready		 
Test Machine	0989090898					Inactive		 

The user can easily see the Endoscopy machine details where the user can **Edit** Endoscopy Machine Dashboard.

Endoscopy Machine Dashboard

[< Back](#)

Machine Name : Endo Machine 1002
Purchase Date :
First Use Date : 07-09-2020
Manufacturer :
Total Used : 3 times

Serial No : 12121323312
Status : Ready
Last Use Date : 07-09-2020
Notes :

[Edit](#)

The user can also see the **History** of usages and maintenance record.

History

[Usage](#)[Maintenance](#)

Usages History

Scope Name	Study Name	Used By	Start Time	End Time	Use Duration	Comments

The user can add a new Endoscopy machine through Endoscopy Machine lists. After clicking the **Add Endoscopy Machine**, it will take them to the **Endoscopy Machine Registration** page.

Endoscopy Machine Lists

[+ Add Endoscopy Machine](#)[< Back](#)

Machine Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Endo Machine 1002	12121323312		Barcode	07-09-2020 05:54 PM	07-09-2020 05:55 PM	Ready		

Endoscopy Machine Registration

Machine Picture

Close

Save

Endoscopy Machine Information

Endoscopy Machine Name*

Johnson2

Serial No*

1132021

Manufacturer

Johnson & Johnson

Purchase Date

11-03-2021

First Use

Notes

Notes

Barcode

john11320201

john11320201

Print

Mandatory fields are highlighted in red*

The user has to press ‘Save’ after completing the registration. Finally, the user will get a confirmation message.

Endoscopy Machine Lists


+ Add Endoscopy Machine

< Back

Machine Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Endo Machine 1002	12121323312			07-09-2020 05:54 PM	07-09-2020 05:55 PM	Ready		
Test Machine	0989090898					Inactive		
Test Machine 2	010101010101	Johnson & Johnson		10-09-2020 12:00 AM		Ready	My sample notes	
Test endoscopy	11111			17-09-2020 12:00 AM		Ready	Testing notes	
Test Endoscopy - 3	5353535	Olympus		10-09-2020 12:00 AM		Ready	My notes is for	
Test Endoscopy - 4	242424242	Coloplast		10-09-2020 12:00 AM		Ready		

Endoscopy machine information is created successfully.

The user can use this new Endoscopy Machine in the **Imaging** section.

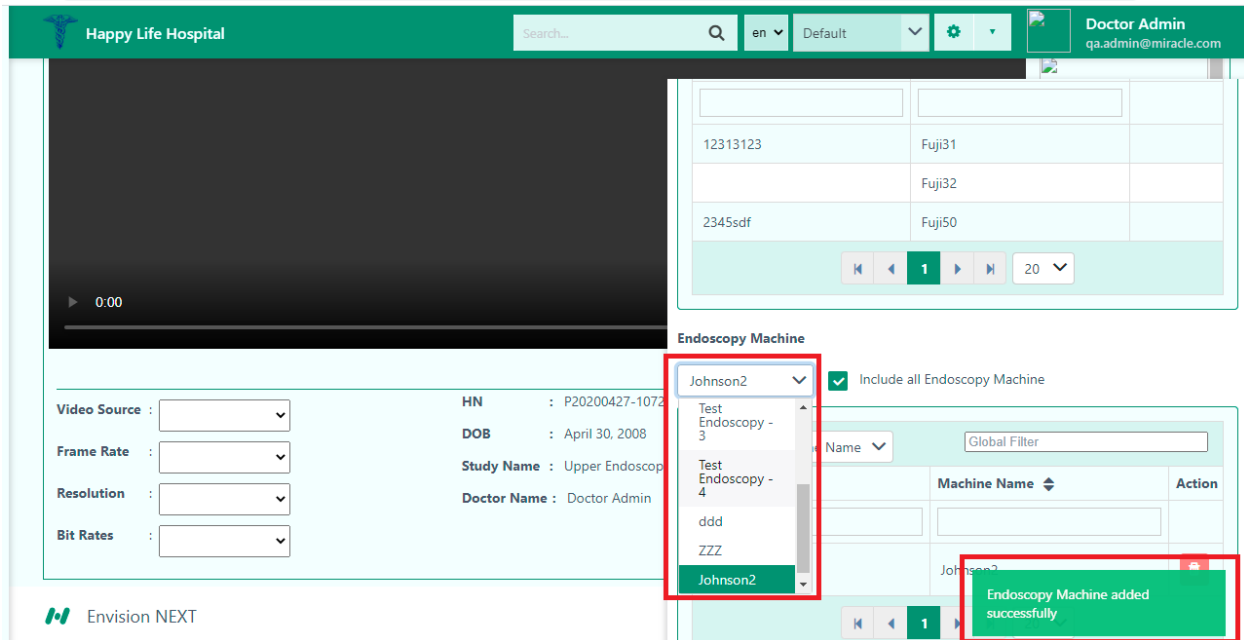

Washing Machine Lists

In the **Washing Machine Lists**, the user will find the lists of washing machines.

Washing Machine Lists								
+ Add Washing Machine Back								
Machine Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Washing 101	123456789	Johnson & Johnson	Barcode	01-09-2020 12:00 AM		Assigned	This is sample notes	
Washing 102	000121	Boston Scientific	Barcode	02-09-2020 12:00 AM		Ready	My sample notes	

In the Washing Machine Lists, the user can see the Barcode, details of the Washing Machine and active or inactive the machine for use.

Machine Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Washing 101	123456789	Johnson & Johnson	Barcode	01-09-2020 12:00 AM		Assigned	This is sample notes	



Note: When the user activates the Washing machine, the machine is ready for use in Imaging. However, when the user inactive the Endoscopy machine, they can't use this machine.

Note: If the user wants to activate or inactive the machine, they can simply click the (Tick) button. When the button is blurred, it will indicate as Inactive. On the other hand, when the button is green, it indicates active mode.

Washing 101	123456789	Johnson & Johnson		01-09-2020 12:00 AM	Inactive	This is sample notes	
Washing 102	000121	Boston Scientific		02-09-2020 12:00 AM	Ready	My sample notes	

The user can easily see the Washing machine details where the user can **Edit** Washing Machine Dashboard.

Washing Machine Dashboard

Machine Name : Washing 101

Purchase Date : 01-08-2020

First Use Date : 01-09-2020

Manufacturer : Johnson & Johnson

Total Used : 0 time

Serial No : 123456789

Status : Assigned

Last Use Date :

Notes : This is sample notes

Edit

The user can also see the **History** of usages and maintenance records.

History

Usage
Maintenance

Usages History

Global Filter

Scope Name	Started By	Ended By	Start Time	End Time	Use Duration	Comments
Fuji32	Doctor Admin		09-03-2021 03:52 PM		0	
Fuji31	Dilruba Tariq	Dilruba Tariq	09-03-2021 03:43 PM	09-03-2021 03:45 PM	00:02:04	
Fuji31	Doctor Admin	Doctor Admin	04-03-2021 09:21 PM	04-03-2021 09:21 PM	00:00:05	

The user can add a new Endoscopy machine through Endoscopy Machine lists. After clicking the **Add Washing Machine**, it will take them to the **Washing Machine Registration** page.

Washing Machine Lists

+ Add Washing Machine
Back

Global Filter

Machine Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Washing 101	123456789	Johnson & Johnson	Barcode	01-09-2020 12:00 AM		Inactive	This is sample notes	

Washing Machine Registration

Machine Picture

Close
Save

Washing Machine Information

Washing Machine Name*
Washing Wolf

Serial No*
3112020

Manufacturer
Johnson & Johnson

Purchase Date

First Use

Notes
Notes

Barcode
John3112021

John3112021
Print

Close
Save

Mandatory fields are highlighted in red*

The user has to press ‘Save’ after completing the registration. Finally, the user will get a confirmation message.

Happy Life Hospital

Search... en Default

Doctor Admin qa.admin@miracle.com

Washing Machine Lists

+ Add Washing Machine < Back

Global Filter

Machine Name	Serial No	Manufacturer	Barcode	First Use	Last Use	Status	Notes	Action
Washing 101	123456789	Johnson & Johnson	Barcode	01-09-2020 12:00 AM		Inactive	This is sample notes	
Washing 102	000121	Boston Scientific	Barcode	02-09-2020 12:00 AM		Ready	My sample notes	
Machine 1	10001	Johnson & Johnson	Barcode	14-09-2020 12:00 AM		Ready	My Notes111	
Test Machine 2	14141414141	Boston Scientific	Barcode	10-09-2020 12:00 AM		Ready	My testing is going on.	
ABCD Machine	0799797987	Johnson & Johnson	Barcode	10-09-2020 12:00 AM		Ready	Testing is done.	
gulgulgul	44353	Johnson & Johnson	Barcode	02-10-2020 12:00 AM		Ready		

Wash machine information is created successfully.

The user can use this new Washing Machine in **Washing and Disinfection** section.

Washing and Disinfection

Department : Internal Medicine

Maker Name : Olympus

Washing

Select Operator

Brushing : No

Cap Check : No

Leakage Check : Normal Abnormal Not Checked

Disinfection

Washing Wolf

Washing 102

Machine 1

Test Machine 2

ABCD Machine

gulgulgul

Washing Wolf

Close Save

After clicking the **Save** button, the user will get a confirmation message.

Happy Life Hospital Search... en Default ⚙️ ⌵ Doctor Admin qa.admin@miracle.com									
Fuji 001	12345	09-03-2021 03:43 PM	09-03-2021 03:45 PM	00:01:04	Coloplast		Disinfection Completed		
Fuji31	31	04-03-2021 09:21 PM	04-03-2021 09:21 PM	00:00:05	Olympus		Disinfection Completed		
Fuji33	33	04-03-2021 09:20 PM	04-03-2021 09:20 PM	00:00:03	Boston Scientific		Disinfection Completed		
Fuji34	34	04-03-2021 09:20 PM	04-03-2021 09:20 PM	00:00:03	Fujifilm		Disinfection Completed		
Fuji 004	123456789	04-03-2021 09:20 PM	04-03-2021 09:20 PM	00:00:03	Medtronic		Disinfection Completed		
Fuji 001	12345	04-03-2021 09:20 PM	04-03-2021 09:20 PM	00:00:03	Coloplast		Disinfection Completed		
Fuji32	32	04-03-2021 09:19 PM	04-03-2021 09:19 PM	00:00:03	Olympus		Disinfection Completed		
Fuji35	35	04-03-2021 09:19 PM	04-03-2021 09:19 PM	00:00:04	Coloplast		Disinfection Completed		
Fuji55	55	04-03-2021 09:19 PM	04-03-2021 09:19 PM	00:00:03	Fujifilm		Disinfection Completed		
Fuji 002	sfsfsf				Johnson & Johnson		Disinfection Saved successfully.		

Updated-3/12/2021

Disinfection Worklist

In the Disinfection Worklist, the user can disinfect scope when it is used for the patient in the **Imaging**. The user can't use the scope when it is already used. They have to confirm the disinfection to use the scope again.

Happy Life Hospital

en

Default

⚙️

⌵

Doctor Admin
qa.admin@miracle.com

0:00

Video Source :

HN : P20200427-1072

Frame Rate :

DOB : April 30, 2008

Scope List

✕

⌵

☒ Include all Scopes

Barcode, Scope Name

Global Filter

Barcode	Scope Name	Action
23sd12	Fuji54	
21213131313144121	Fuji 10001	

⏮

⏪

1

⏩



⏭

20

⌵

Select Scope

When the user selects the scope for the **Imaging**, the scope will be ready to use. After clicking the **Scope In** button, the user can't use this scope for another patient until the scope is disinfected in the **Disinfection Worklist**.

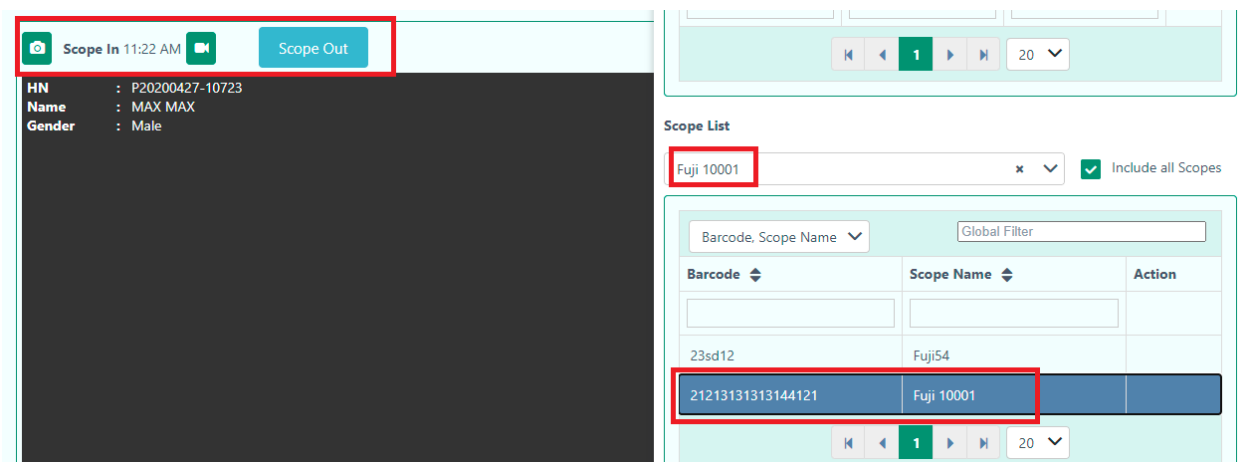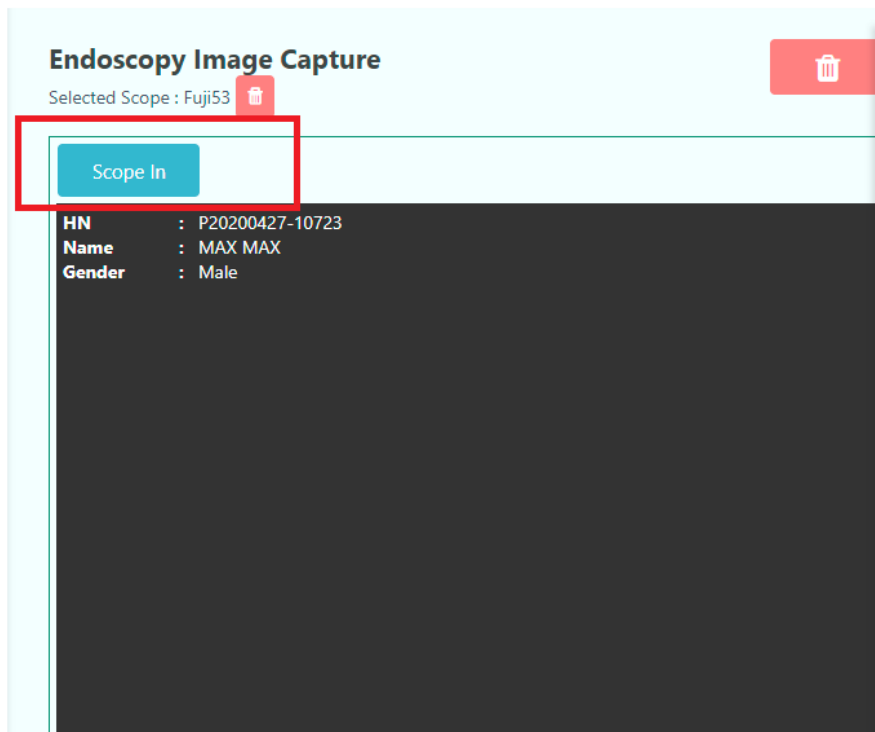

Note: When the user presses the **Scope in** button, the Action list in the right corner becomes disable. The user should disinfect the scope later to use it again.

Scope List

Fuji55 x ☒ Include all Scopes

Barcode, Scope Name Global Filter

Barcode	Scope Name	Action
2345sdf	Fuji50	
23sd1	Fuji53	
23sd10	Fuji55	👍 🗑️

The user will find this scope in the **Disinfection Worklist**. There are two parts in this section.

- **Washing**
- **Disinfection**

Washing Worklist

Back

From Date
To Date
Barcode
Serial Number

☐ Show All Dirty

9 columns selected

Scope Name	Serial No	Barcode	Start Time	End Time	Duration	Maker Name	Notes	Status	Action
EG-580UT	4U402K326	Barcode 112445568	07-03-2022 06:54 PM	07-03-2022 06:54 PM	00:00:00	Fujifilm		Washing Completed	👁
EG-580UT	4U402K326	Barcode 112445568	07-03-2022 06:54 PM	07-03-2022 06:54 PM	00:00:00	Fujifilm		Washing Completed	👁

Note: The user can search any user scope with **Barcode** and **Serial Number**. The user will see the dirty scope by searching with **Date**. It is a by default searching option.

Washing Worklist

« Back

Washing

Disinfection

From Date

To Date

Barcode

Serial Number

Q Search

Show All Dirty

9 columns selected

Global Filter

Scope Name	Serial No	Barcode	Start Time	End Time	Duration	Maker Name	Notes	Status	Action
EG-580UT	4U402K326	Barcode 112445568	07-03-2022 06:54 PM	07-03-2022 06:54 PM	00:00:00	Fujifilm		Washing Completed	
EG-580UT	4U402K326	Barcode 112445568	07-03-2022 06:54 PM	07-03-2022 06:54 PM	00:00:00	Fujifilm		Washing Completed	

When the user searches the dirty scope with **Barcode** or **Serial Number**, the searching option with Date will disappear.

Washing and Disinfection

« Back

Barcode

Serial Number

Q Search

8 columns selected

Global Filter

Scope Name	Serial No	Start Time	End Time	Duration	Maker Name	Notes	Status	Action
Fuji 10001	1212121213312				R. Bard			

1

20

The user presses the **Scope Out** button in the **Imaging** section, the Status will show **Dirty** in the **Working Worklist**.

From Date

To Date

Barcode

Serial Number

Search

Show All Dirty

9 columns selected

Global Filter

Scope Name	Serial No	Barcode	Start Time	End Time	Duration	Maker Name	Notes	Status	Action
scp	1616							Dirty	
sc	141414							Dirty	
zz	1414							Dirty	

Note: In the above right corner, the user will find the **Show All Dirty** checkbox where they can see only dirty scopes in the lists.

Envision NEXT

Happy Life Hospital

Search by Patient / Appointment

en

Default

Bipon Biswas

bipon770@gmail.com

Washing

Disinfection

From Date

To Date

Barcode

Serial Number

Search

Show All Dirty

9 columns selected

Global Filter

Scope Name	Serial No	Barcode	Start Time	End Time	Duration	Maker Name	Notes	Status	Action
scp	1616							Dirty	
sc	141414							Dirty	
zz	1414							Dirty	

Washing

In the **Washing** part, the user will find both completed and dirty scope. But if the user check the **Show All Dirty** then they will only get dirty status.

Envision NEXT | Happy Life Hospital | Search by Patient / Appointment | en | Default | Bipon Biswas | bipon770@gmail.com

Washing | Disinfection

From Date: 01-02-2022 | To Date: 14-03-2022 | Barcode: | Serial Number: | Search | **Show All Dirty**

9 columns selected | Global Filter

Scope Name	Serial No	Barcode	Start Time	End Time	Duration	Maker Name	Notes	Status	Action
scp	1616							Dirty	
sc	141414							Dirty	
zz	1414							Dirty	

By clicking the view washing  button, the user will find the details of the Washing.

Washing

Operator

Operator Name : Bipon Biswas

Gender :

Department :

Device

Device Name : EG-580UT

Serial Number : 4U402K326

Maker Name : Fujifilm

Washing

Bipon Biswas

Leakage Check * : Normal Abnormal **Not Checked**

Close

The user will find the dirty scope and can wash by clicking this  button.

In the **Washing** part, the user has to select the operator from the dropdown for washing.

Washing

Operator

Operator Name : Bipon Biswas

Gender :

Department :

Device

Device Name : scp

Serial Number : 1616

Maker Name :

Washing

Select Operator

นางสกลกัน ขันดีน (Nurse)

นางยุวดี พลอุทัย (Nurse)

Bipon Biswas (Administration)

Miss Kasama Pengsomboon (Nurse)

Miss Jariya Nurat (Nurse)

Mr. Kanokpon Thongpanchang (Nurse)

Close Save

After selecting the operator and leakage check status, the user has to click the **Save** button. The user will get a confirmation message after that.

Confirmation

☒ Do you sure you want to finish the Washing?

No Yes

Note: Leakage check is mendatory *.

Disinfection

In the disinfection part, the user must select the washing machine.

Disinfection

Washing 102

Washing 102

Machine 1

Test Machine 2

ABCD Machine

gdgdgdg

▶ Start

■ Finish

✕ Close

💾 Save

The user can check some fields and provide notes if necessary, in the **Disinfection** section. Then, the user can press the **Start** button to begin the disinfection procedure of the scope.

Disinfection

Washing 102

Alcohol Used : ✓ Yes

Concentration Check : Normal Abnormal Not Checked

Washing Status : Normal Abnormal Not Checked

Comment : Disinfection Fini Automatic Termi Not Checked

Notes : Try to disinfect for 50 second extra

▶ Start

■ Finish

✕ Close

💾 Save

Disinfection

Washing 102

Alcohol Used

:

✓

Yes

Concentration Check

:

Normal

Abnormal

Not Checked

Washing Status

:

Normal

Abnormal

Not Checked

Comment

:

Disinfection Finis

Automatic Termi

Not Checked

Notes

:

Try to disinfect for 50 second extra

▶ Start

00:00:07

■

Finish

It also shows the time scale when the user presses the **Start** button. When the user ensures the disinfection of the scope, they have to press **Finish** to stop the procedure. They will get a confirmation message for finishing the procedure. After clicking the ‘**Yes**’ button it will give a successful message to the user. This disinfected scope can be ready to use for another patient.

Confirmation

☒ Are you sure you want to finish the disinfection?

Comment : Disinfection Finis Automatic Termi Not Checked

Notes : Try to disinfect for 50 second extra

00:23:34

Washing and Disinfection

Department : Maker Name : Fujifilm

Washing

Doctor Admin (Doctor)

Brushing :

Cap Check :

Leakage Check :

Disinfection

Washing 102

Alcohol Used :

Concentration Check :

Washing Status :

Comment : Disinfection Finis Automatic Termina Not Checked

Notes : Try to disinfect for 50 second extra

00:25:21

Disinfection finished successfully.

The user will get detailed information of this disinfected scope in **Washing and Disinfection** lists.

Washing and Disinfection

Back

From Date

To Date

Barcode

Serial Number

Search

01-02-2021

15-03-2021

Search

8 columns selected

Global Filter

Scope Name	Serial No	Start Time	End Time	Duration	Maker Name	Notes	Status	Action
Fuji54	54	15-03-2021 03:56 PM	15-03-2021 04:21 PM	00:25:22	Fujifilm	Try to disinfect for 50 second extra	Disinfection Completed	View

The user can know the start and end time of disinfection procedure, duration of the disinfection, status and notes of the scope.

The user can also view the washing and disinfection page again by simply click the right corner button.

8 columns selected

Global Filter

Scope Name	Serial No	Start Time	End Time	Duration	Maker Name	Notes	Status	Action
Fuji54	54	15-03-2021 03:56 PM	15-03-2021 04:21 PM	00:25:22	Fujifilm	Try to disinfect for 50 second extra	Disinfection Completed	View

Updated-15/3/2021



BI

In this Business Intelligence section, the user will get raw data and statistical data for better visualization.

Happy Life Hospital

Search...

en

Default

Settings

Doctor Admin

Dashboard
Whiteboard
Patient
Nursing
Device Management
Disinfection Worklist
BI
Administration

Business Intelligence

From Date

To Date

16-03-2021

23-03-2021

Search

Include Chart

Study Report

Connect

Open

Save

Export

Format

Options

Fields

Fullscreen

	1	2	3	4	5	6	
1	STUDYNAME	Total Count of StudyName					
2	Colonoscopy	2					
3	ERCP	1					
4	Sigmoidoscopy	1					
5	Upper Endoscopy	2					
6	Grand Total	6					
7							
8							
9							
10							
11							
12							
13							

In the first part, the user can search data, **From Date -To Date**. They can also select whether they want to see **Study Report** or **Scope Report** from the dropdown.

Business Intelligence ← Back

From Date

16-03-2021

To Date

23-03-2021

Search

☐ Include Chart

Study Report

Study Report

Scope Report

Connect

Open


Save

Export

Study Report

Scope Report

If the user wants to visualize the data, they have to click the **Include Chart**.

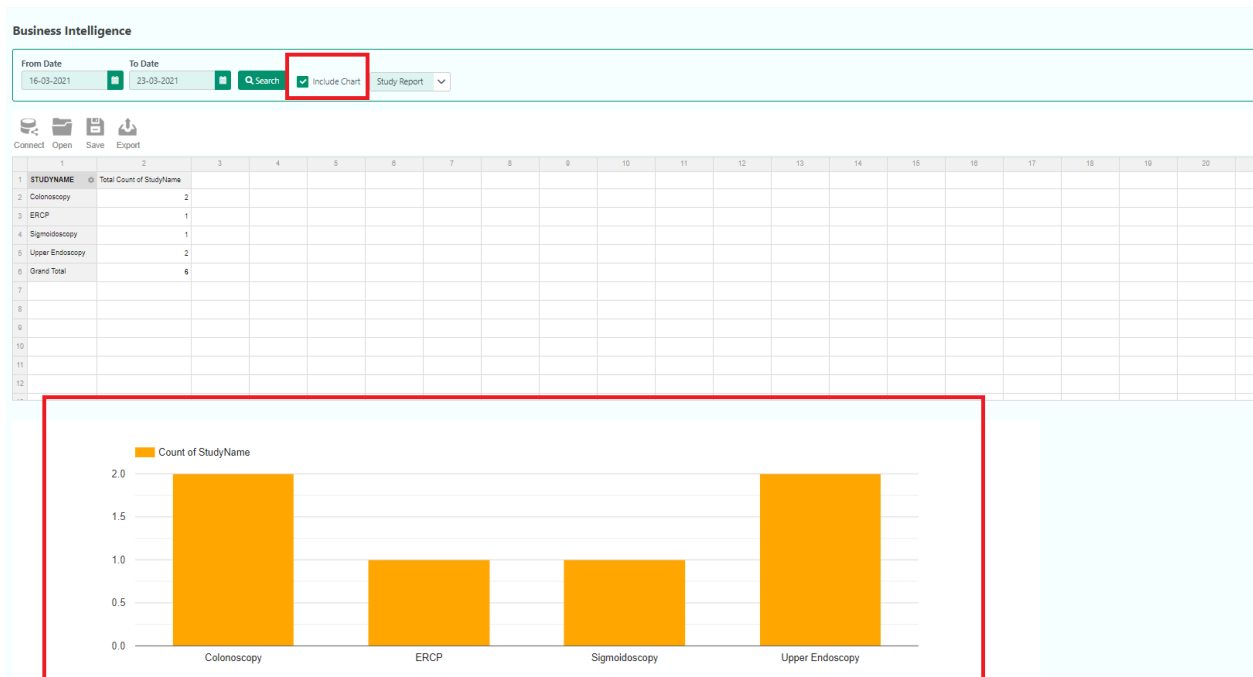


In the left corner users will find some functions like- **Connect, Open, Save, and Export**.


Connect

In this option the user can connect any data from-


- To Local CSV
- To Local JSON
- To remote CSV
- To remote JSON

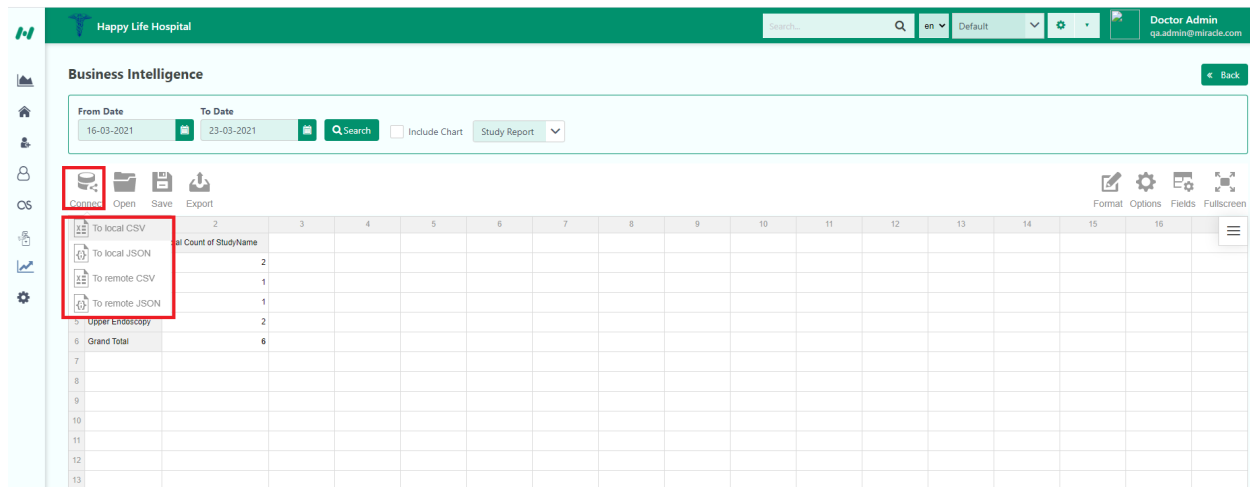

They can also open Local Report or Remote Report from **Open** option.

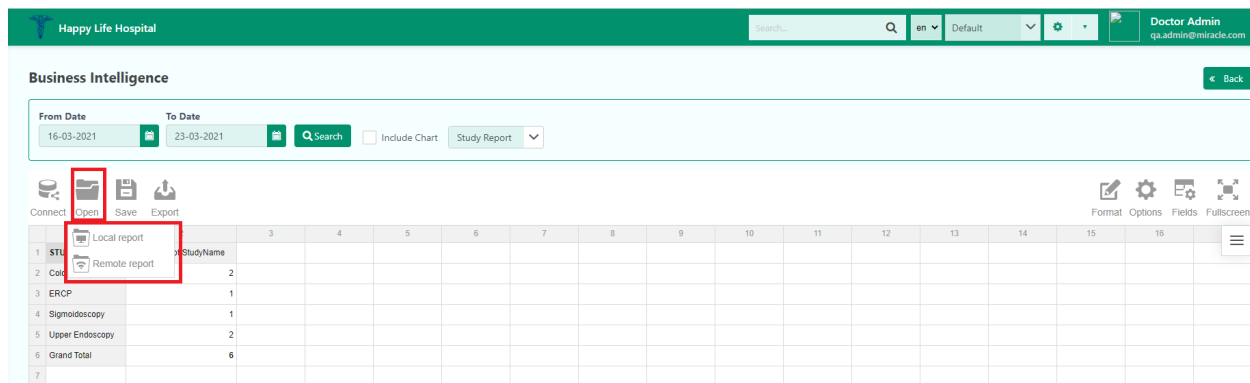


The user can also Print, export to HTML, Excel and PDF format by clicking in the **Export** option.

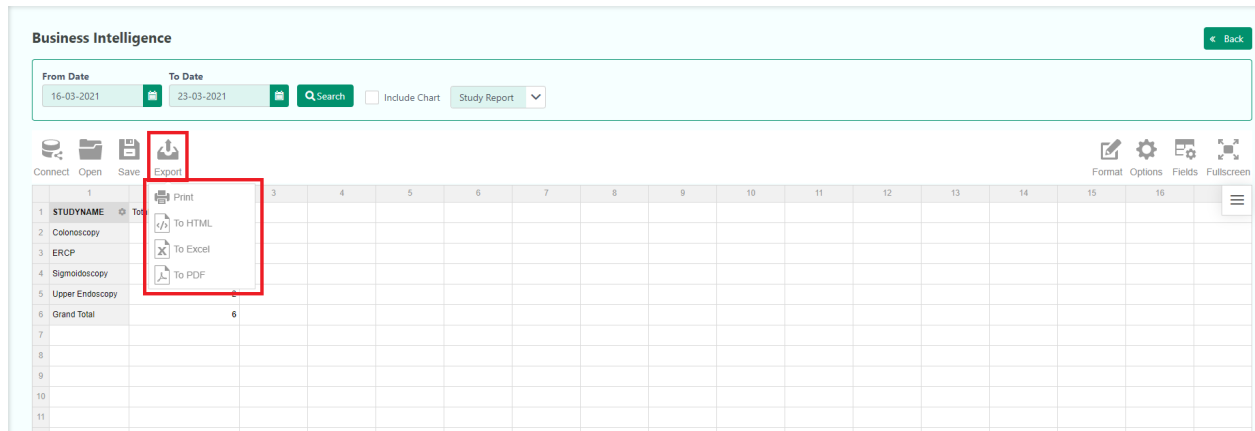

In the right corner, the user will find **Format**, **Options**, **Fields**, and **Full Screen** functions.

The user can easily format the data by clicking in the **Format**. The user will find **Format cells** and **Conditional formatting** inside the **Format** dropdown.

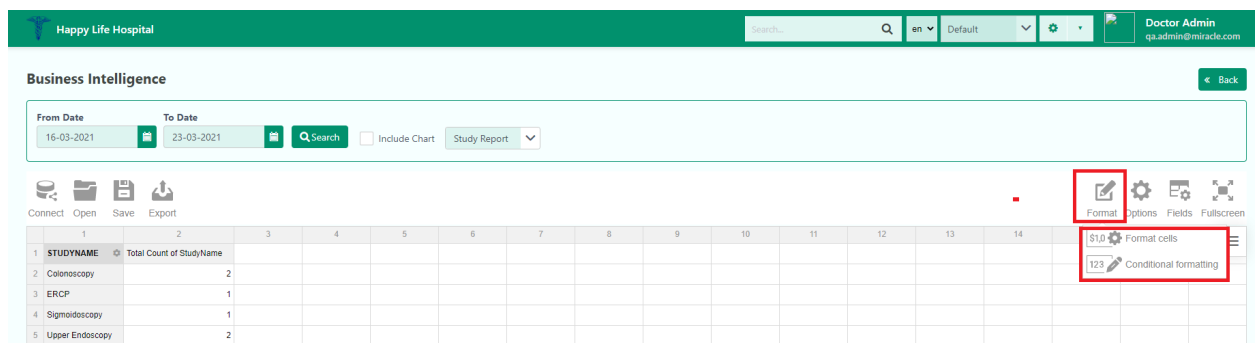


In the **Format cells**, the user will find some basic format changes options.

Format cells **APPLY** **CANCEL**

CHOOSE VALUE	Choose value	▼
Text align	right	▼
Thousand separator	(Space)	▼
Decimal separator	.	▼
Decimal places	None	▼
Currency symbol		
Currency align	left	▼
Null value		
Format as percent	false	▼

In the **Conditional formatting**, the user will first find a pop-up window where they have to choose '+' button for further proceeding.

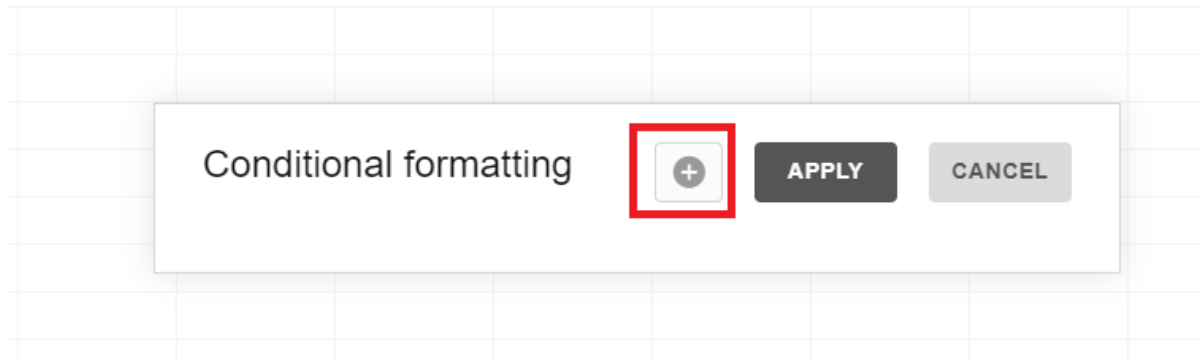

There they can choose the format either less than, less than or equal and so on. Moreover, they can also choose colour to those format changes.

Business Intelligence

From Date: 16-03-2021 To Date: 23-03-2021 Search Include Chart Study Report

Connect Open Save Export

STUDYNAME	Total Count of StudyName
Colonoscopy	2
ERCP	1
Sigmoidoscopy	1
Upper Endoscopy	2
Grand Total	6

Conditional formatting dialog box showing options for Value (All values), Format (Arial), and a list of conditions: Less than, Less than or equal to, Greater than, Greater than or equal to, Equal to, Not equal to, Between, Empty. The 'Less than' option is highlighted.

1	2	3	4	5	6	7	8	9
STUDYNAME	DOCTOR							
ROOM	Doctor Admin	Laura Doc	Mostafa Chowdhury	Santi Kulpacharapong	Sureeporn Jangsirikul	karim das	Total Count of StudyName	
> Colonoscopy	2,0000	1,0000	0,0000	1,0000	2,0000	0,0000	6,0000	
> EGD	10,0000	2,0000	0,0000	0,0000	0,0000	1,0000	13,0000	
> ERCP	7,0000	0,0000	0,0000	1,0000	3,0000	1,0000	12,0000	
> Sigmoidoscopy	2,0000	0,0000	1,0000	0,0000	0,0000	0,0000	3,0000	
Grand Total	21,0000	3,0000	1,0000	2,0000	5,0000	2,0000	34,0000	

In the Option, the user can choose Layout options of the data.

Layout options **APPLY** **CANCEL**

GRAND TOTALS

- ☐ Do not show grand totals
- ☒ **Show grand totals**
- ☐ Show for rows only
- ☐ Show for columns only

SUBTOTALS

- ☐ Do not show subtotals
- ☒ **Show subtotals**
- ☐ Show subtotal rows only
- ☐ Show subtotal columns only

LAYOUT

- ☒ **Compact form**
- ☐ Classic form
- ☐ Flat form

The user can drag and drop fields to arrange the data inside the **Fields** option.

Fields

Drag and drop fields to arrange

All Fields
Collapse All

OrderDate

☒ OrderDate.Day

☒ OrderDate.Month

☒ OrderDate.Year

☒ Doctor

☒ Room

☒ StudyName

☒ StudyType

Add calculated value

APPLY

CANCEL

Report Filters

OrderDate.Month

OrderDate.Day

Drop field here

Columns

Room

Σ Values

Drop field here

Rows

Doctor

OrderDate.Year

StudyType

Values

Count of Doctor
Σ ∨

Count of Room
Σ ∨

Count of StudyName
Σ ∨

Note: The user can **Add calculated value** from the **Fields**.

Calculated Value

APPLY

CANCEL

Drag values to formula

Value name

OrderDate.Month (Count)

OrderDate.Year (Count)

Room (Count)

StudyName (Count)

StudyType (Count)

☐ Calculate individual values

+	-	x	÷	^	=	<	>	≤
≥	==	!=	OR	AND	IF	ABS	MIN	MAX

Drop values and edit formula here

Finally, they can maximize the screen with the **Fullscreen** option.

Updated-23/3/2021

Administration

This section will help the user to set up the whole settings. The user can customize-

© 2021 Miracle Advance Technologies Co. Ltd

- **Video Settings**
- **Update Organization**
- **Room Management**
- **User Management**
- **Study Library**
- **Capture Settings**
- **Study Management**
- **Support Ticket Lists**
- **Locked Reports**
- **Examination Team**
- **PACS Worklist**
- **Consumable Items**
- **Patient List**
- **Dynamic Dashboard**
- **Add Post**
- **Calendar Settings**
- **Audit Trails**

Envision NEXT Happy Life Hospital

Search... en Default

Bipon Biswas bipon370@gmail.com

Video Settings + Add Device Back

8 columns selected Global Filter

Room Name	Device Name	Device Code	Resolution	Bit Rates	Frame Rate	Maker Name	Purchase Date	Action
Room 106	Device 001	121212123	720p	1 MB bps	24fps	Coleplast	02-09-2020	
Room 102	Integrated		1080p	100 MB bps	30fps		26-08-2021	
Room 101	Blue Scope		1080p	1 MB bps	30fps	Johnson & Johnson	09-10-2020	

Envision NEXT FUJIFILM

Video Settings

In this section, the user can Add Device for capturing video.

Video Settings + Add Device Back

8 columns selected Global Filter

Room Name	Device Name	Device Code	Resolution	Bit Rates	Frame Rate	Maker Name	Purchase Date	Action

By clicking in the **Add Device**, the user will get a pop-up form. They have to fill up the form for adding a new device for the Video.

The screenshot shows a web form titled "Add Device" with a sub-header "Device Information". The form contains two columns of input fields. The left column includes "Device Name*" (text), "Room*" (dropdown), "Frame Rate" (dropdown), "Maker" (dropdown), and "Price" (text). The right column includes "Device Code" (text), "Resolution" (dropdown), "Bit Rate" (dropdown), and "Purchase Date" (text with a calendar icon). The asterisks on "Device Name" and "Room" indicate they are mandatory fields. At the bottom right, there are "Close" and "Save" buttons.

Mandatory fields are highlighted in red*

After filling the form the user has to save the device by clicking the **Save** button.

This screenshot shows the same "Add Device" form, but now it is filled with data. The "Device Name" field contains "Fuji123456", "Room" is set to "Room 103", "Frame Rate" is "15fps", "Maker" is "Fujifilm", and "Price" is "200000". On the right, "Device Code" is "fuji123456", "Resolution" is "1080p", "Bit Rate" is "1 GB tps", and "Purchase Date" is "29-03-2021". The "Save" button at the bottom right is highlighted with a red rectangle, indicating the next step in the process.

The user will get successful notification after clicking the **Save** button.

The screenshot displays the 'Video Settings' page in the Envision NEXT application. The page header includes the Envision NEXT logo, 'Happy Life Hospital', a search bar, and user information for 'Doctor Admin'. The left sidebar lists various modules like Dashboard, Whiteboard, Patient, Nursing, etc. The main content area features a table with 8 columns selected. The table data is as follows:

Room Name	Device Name	Device Code	Resolution	Bit Rates	Frame Rate	Maker Name	Purchase Date	Action
Room 103	Fuji123456	Fuji123456	1080p	1 GB bps	15fps	Fujifilm	29-03-2021	[Add] [Delete]

A red box highlights the row for 'Room 103' and its corresponding action buttons. Another red box highlights a green notification message at the bottom right: 'Video setting saved successfully.'

Note: The user can delete or edit the video settings with the help of the above red marked button.

Update Organization

In this section, the user can update organizational information.

The screenshot shows the 'Update Organization' form. The form is titled 'Update Organization' and contains a section for 'Organization Information'. The form is divided into two columns of input fields:

Field	Value
Name	Happy Life Hospital
Address	144/4 M.5 Ban Tham Chiang Dao Chiang Mai 50170 Thailand
Email	admin@miracle.com
Client Logo	[Choose]
Alias	Happy
Contact	+8888888888
Website	https://happyliife.med

At the bottom right of the form, there are two buttons: 'Back' and 'Save'.

The user can also change some organizational settings by clicking the right corner button “Org Settings”.

After clicking the button, the user will get the Update **Organization** and **Custom Form** section.

Update Organization

Workflow Settings	DICOM Settings
Disinfection Type: Simple	Send to PACS?: ☐
Washing Start/Stop: No	PACS Send Status: Imaging
Report Layout: Layout One	PACS IP: 101 Ping
Auto Arrived: ☒	PACS Port: 104
Skip Pre Scope: ☒	PACS AE Title: ROUTER_AE Echo
Skip Capture Time: ☒	Router IP: 192.168.7.103
Skip Post Scope: ☒	Router Port: 110
Study Completed On: Imaging	Router AE Title: ROUTER_AE
Enable Quick Registration: ☒	

Back Save

Inside the update organization, there are two parts.

- Workflow settings.
- DICOM settings.

In the Workflow settings, the user can configure Disinfection type, Washing Start/Stop, Report Layout, and Study Completed on.

Moreover, the user can also configure DICOM settings.

The user has to click the **Save** button after configuring the Update Organization.

Workflow Settings	DICOM Settings
Disinfection Type: Simple	Send to PACS?: ☒
Washing Start/Stop: No	PACS Send Status: Imaging
Report Layout: Layout Three RH 👁	PACS IP: 101
Auto Arrived: ☒	PACS Port: 104
Skip Pre Scope: ☒	PACS AE Title: ROUTER_AE
Skip Capture Time: ☐	Router IP: 192.168.7.103
Skip Post Scope: ☒	Router Port: 110
Required Sedation: ☒	Router AE Title: ROUTER_AE
Make Sound On Image Capture: ☒	
Prevent Double Click In Imaging: ☒	
Make Sound On Start Video Recording: ☒	
Make Sound On Stop Video Recording: ☒	
Required Report Approval: ☒	
Skip Team Check: ☒	
Auto Select Intervention: ☐	
Skip Default Image: ☒	
Study Completed On: Reporting	
Enable Quick Registration: ☒	
Patient Name Format: (Title) (First Name) (Last Name)	
Display Patient Name(Lang): ☒ English ☐ Thai	
Doctor Name Format: (Title) (First Name) (Last Name)	
Display Doctor Name(Lang): ☒ English ☐ Thai	
Tab Refresh Times(sec): 2	
Assigned Scope Status: Ready	

Back Save

The user will get a notification after successfully updating the organization.

Room Management

To add a new room or know about the details and make the room active/inactive, the user can manage from this section. The user can see lists of the room from this part.

Room Lists								
+ Add Room ← Back								
Global Filter								
Room Name	Room Type	Building	Floor	Ward	Last Use	Total Study	Status	Action
Room 101		Test Building	ชั้น 2	Non-infectious cat ward	05-03-2021 01:14 PM	38	In Use	✓
Room 102		Test Building	ชั้น 1	Blood parasite dog ward	09-03-2021 04:41 PM	12	Available	✓
Room 103		Test Building	ชั้น 2	Infectious dog ward	29-03-2021 03:23 PM	10	In Use	✓
Room 104		Test Building	ชั้น 2	Non-infectious cat ward	04-03-2021 10:16 PM	7	Available	✓
Room 105		Test Building	ชั้น 3	Infectious cat ward	11-03-2021 11:06 AM	4	In Use	✓
Room 106		Test Building	ชั้น 3	CCU	11-03-2021 11:30 AM	7	In Use	✓

The user can add a new room by clicking in the **Add Room**. It will take them to **Room Registration**.

Room Registration

Room Picture

Close

Save

Room Information

Room Name*

Room 343

Room Type*

RECOVERY

Building

Test Building

Floor

ชั้น 3

Ward*

Select ward

Close

Save

Mandatory fields are highlighted in red*

User Management

The user can edit, add a new user, and map the role of the users in this section.

User Information

Upload Users

Add User

Role Mapping

Back

8 columns selected

Global Filter

First Name	Last Name	User Name	Email	Time Zone	JobType	Created On	Department Name	Action
somsak								
Somsak	Phophong	Somsak	somsak@miracle.co	(UTC +07:00:00)Asia/Ban	Doctor	April 5, 2022	Surgery department	

1

20

Add User

User Demographic

Title

Select Title

First Name*

Last Name*

Gender

Select Gender

Date of Birth

Phone

Line ID

User Information

Email*

Skip Verification

Password*

Confirm Password*

Role*

Choose

Job Type*

Select Job Type

Time Zone*

Select Time zone

Department Name*

Select Department

Organization

Happy Life Hospital

Back

Save

Envision NEXT

FUJIFILM

Mandatory fields are highlighted in red*

In the **Role Mapping**, the user can add roles and map roles for other users. They can also edit and assign the role of the user with the help of **mapping**.

Role Mapping			
			+ Add Role ← Back
Global Filter			
Role Code	Role Name	Description	Action
ADM001	Administrator	admin roles@@@	 
DOC001	Doctor	Doctor	 
END001	Endoscopist	Endoscopist	 
CNT001	Counter	Counter	 

Menu Assign to Role				
Role Assign to Menu	Submenu User Name	Router Link	Description	Action
Administration				
SUBM006	Administration	administration/user-list	Administration	✓
SUB008	Video Settings	administration/video-setting	Video Settings	✓
SUB009	Update Organization	update-org	Update Organization	✓
SUB0010	Room Management	administration/room-listings	Room Management	✓
SUB0011	User Management	administration/user-list	User Management	✓
SUB0011	Study Library	/dashboard/administration/study-library	Study Library	✓
SUB0013	Capture Settings	/dashboard/capture-setting	Capture Settings	✓
SUB0015	Study Management	/dashboard/administration/study-listings	Study Management	✓
SUB0016	Support Ticket Lists	/dashboard/support-ticket-listings	Support Ticket Lists	✓
SUB0015	Locked Reports	/dashboard/administration/report-locked	Locked Reports	✓

In the user information there are four **Action** buttons.

- **Edit Team Role**
- **Edit**
- **Reset Password**
- **Delete**

Edit Team Role

The user can update team role by clicking  this button.

X

Update Team Role

User Name : Somsak

Team Role :

☐

☐

Primary Endoscopist

☒

Interventional Endoscopist

☐

Supervising Endoscopist

☒

Senior Fellow

☒

Sedating Physician

☐

Anesthesiologist

☐

Nurse

☐

Assistant Nurse

☐

Recorder

✕ Cancel

✓ Update

After selecting the role, the user has to click the Update button.

User Information

Upload Users Add User Role Mapping Back

8 columns selected Global Filter

First Name	Last Name	User Name	Email	Time Zone	JobType	Created On	Department Name	Action
somsak								
Somsak	Phophong	Somsak	somsak@miracle.co	(UTC +07:00:00)Asia/Ban	Doctor	April 5, 2022	Surgery department	

Envision NEXT

Team Role(s) updated successfully.

Finally, they will get a confirmation message after clicking the update button.

Study Library

When the doctor saves the report for study, those data are stored in the study library.

Study Library

All Favorite Research Teaching

From Date To Date Tags

01-04-2021 01-04-2021 Search Add

Own Study Shared Study

13 columns selected Global Filter

Study Date	Accession No	Patient Name	Age	Gender	Study Code	Study Name	Doctor Name	Favorite	Research	Teaching	Tags	Shared	Action
25-03-2021	202103251392	jone doe	1 Year	Male	20171	Upper Endoscopy	Doctor Admin						

In this section, the user can

- View Report
- View Schedule
- View Preparation Data
- View Recovery Data
- PACS(PICTURE ARCHIVING & COMMUNICATION SYSTEM)
- Edit
- Share

- Resource Allocation
- Vital Sign
- Study Files

Capture Settings

The user can set the capture settings in the **Capture Settings Entry**. The user has to select Study, Time, Capture Type. When the user clicks **None**, it won't take any captured file. The **Pre** capture image from the starting portion of the video. It will capture the first 2 seconds of the video image. The **Post** captures an image from the ending portion of the video. It will capture the first 2 seconds of the video image. Both capture images from both portions of the video. It will capture the first 6 seconds of the video image.

The user can also update or reset the data of the capture settings.

Capture Settings

Capture Settings Entry

Study*

Select procedure

Capture Type

None Pre Post Both

Time*

2

Allowed Video

No

Update

Reset

Study	CaptureType	Time (in Sec)	Allowed Video	Action
ERCP	Pre	2	No	
PEG	Both	2	Yes	
Upper Endoscopy	Pre	2	Yes	
Colonoscopy	Both	4	Yes	

Study Management

The user will get study lists from the Study Management. They can see the **details** of the Study management and can **Active/Inactive** the study management from the **Study Lists**.

Study Lists								+ Add Study	← Back
Global Filter									
Study Code	Study UID	Study Name	Procedure Type	Completed	Scheduled	Imaging Done	Reporting Done	Action	
2512	ERCP	ERCP	Endoscopy	1	0	1	2		
2513	PEG	PEG	Endoscopy	0	0	0	0		
2514	Upper Endoscopy	Upper Endoscopy	Endoscopy	23	0	23	24		

To add a new study, the user has to click the **Add Study** button. This action will take to **Study Registration** pop-up form.

Study Registration

Study Information

Study UID*

Study UID

Study Name*

Study Name

Procedure Type*

Select Procedure Type

Close

Save

Mandatory fields are highlighted in red*

Support Ticket Lists

When the user faces problems, they can take help from the Support.

Happy Life Hospital

Search...

en

Default

Doctor Admin
qa.admin@miracle.com

Support Ticket Lists

Set As Home

Version : 0.1.0

Support

← Back

Ticket Id	Subject	Support Type	Error Code	Priority	Page Name	Ticket Status	Create Date	Action
-----------	---------	--------------	------------	----------	-----------	---------------	-------------	--------

They will find those support at the **Support Ticket Lists**.

Support Ticket Lists

Back

Global Filter								
Ticket Id	Subject	Support Type	Error Code	Priority	Page Name	Ticket Status	Create Date	Action
	1008 Reporting Intervention no tcomplete	Error		High	/dashboard/whiteboard	New	05-03-2021 10:31 AM	
	1007 ss	Bug		High	/dashboard/administrat library	Closed	23-12-2020 03:50 PM	
	1006 Support with ZIP file	Error		Normal	/dashboard/support-ticket-listings	New	21-12-2020 07:43 PM	
	1005 Support with RAR File	Query		Medium	/dashboard/support-ticket-listings	New	21-12-2020 07:42 PM	

The user can see the details of the tickets. They can also **Active/Inactive** the tickets from the **Support ticket Lists**.

Support Ticket

Notes

Write the details

Attachment

Emulation Next

AppWhisperer (Helpful)

AppWhisperer

+

Not secure | app.emulationnext.net/vs/dashboard/whisperboard/2

History | Life Preserver

Support Ticket

Subject*

Reporting Intervention not complete

Priority*

High

Support Type*

Error

Error Code

Error Code

Normal

Source Serif

B

I

U

A

Link

Table

Code

Quote

Fullscreen

More

Enter text domain

Available file types are: .png, .jpg, .gif, .mp4, .mpeg, .flv, .avi, .mkv, .pdf, .doc, .docx, .txt, .xls, .xlsx, .zip, .rar

Choose

Locked Reports

The user can not access other user's reports. However, the admin can unlock the report by clicking on the **Unlock** button. It will help the other user to access that report.

Search By : From Date
To Date
01-04-2021
01-04-2021
Search

Global Filter

	Accession No	Patient Name	HN	Doctor	Study Type	Study Name	Request Type	Locked Status	Locked On	Locked By	Requested By	Action
▲	2021020713833	Ayan Lee	P2021020712013	Santi Kulpacharapong	Endoscopy	Upper Endoscopy	Routine	Locked	10-02-2021 09:58 AM	Dr. Doctor Admin		Unlock
▲	2021021913849	MAX MAX	P20200427-10723	Blipon Biswas	Endoscopy	Colonoscopy		Locked	23-02-2021 08:43 AM	Dr. Doctor Admin		Unlock
▲	2021022513860	MAX MAX	P20200427-10723	Blipon Biswas	Endoscopy	Upper Endoscopy	Routine	Locked	02-03-2021 12:54 PM	Dr. Doctor Admin		Unlock
▲	2021022513866	Fahim7 Hasan	P2021022312021	Doctor Admin	Endoscopy	Sigmoidoscopy	Routine	Locked	04-03-2021 06:52 AM	Dr. Doctor Admin		Unlock

Examination Team

Here, the user can add an examination team. The user can add a team, edit, inactive, and delete an examination team from this option.

Examination Team

Add Team
Back

Team Name, Team Members, Status

Global Filter

Team Name	Team Members	Status	Action
Group 1	Dr. Nayem Siddique, Miss Kasama Pengromboon	Active	Add Edit Inactive Delete
Group 2	Dr. Nayem Siddique, Miss Kasama Pengromboon, David Attenborough	Active	Add Edit Inactive Delete

1
20

The user can Add a team by clicking on the **Add Team** button.

Add Team

Team Name*

Team Name

Add Member*

Select team member

mahfuz shazol (Doctor)
Mr. Md. Akter Hossain (Administration)
Fujfilm Salesrep (Doctor)
Dr. Ostafa Howidhu (Doctor)
Mr. Chevaranon Pitakjaruporn (Doctor)
Sitti Borisuit (Administration)
Miss. Warunee Kithibongkarn (Doctor)

Close
Save

Mandatory fields are highlighted in red*

PACS Worklist

The user can configure PACS(PICTURE ARCHIVING & COMMUNICATION SYSTEM) and see the details of the patient regarding PACS settings.

PACS Configuration

PACS IP: 101
PACS Port: 104
PACS AE Title: ROUTER_AE

Search By : From Date: 18-09-2021 To Date: 19-09-2021

Table Columns: Patient Name, Study Date, Study Name, Accession Number, Number Of Image, PACS Send Status, Action

Footer: Envision NEXT, FUJIFILM

Consumable Item

The user can add consumable item from this settings. They can also add new consumable item by clicking **Add Item**.

Consumable Item

Buttons: + Add Item, < Back

Table Columns: Item Uid, Item Name, Item Type, Brand, Vendor, Unit, Price, Action

Item Uid	Item Name	Item Type	Brand	Vendor	Unit	Price	Action
CMBL202107091002	Savlon	CONSUMABLE	Boston Scientific	Thai Medical Development Ltd. Part.	ml	220.5555	
CMBL202107091001	Disinfectant	CONSUMABLE	Boston Scientific	Medtop Co. Ltd.	Ps	70	

Add Consumable



Consumable Item

Item Code*: Item Code
Item Type*: Select Type
Brand: Select Brand
Item Unit: Select Unit
Auto Generate: ☐
Item Name*: Item Name
Vendor: Select Vendor
Item Price: Item Price

Buttons: Close, Save

Mandatory fields are highlighted in red*

The user can also add consumable item from Imaging, and Reporting.

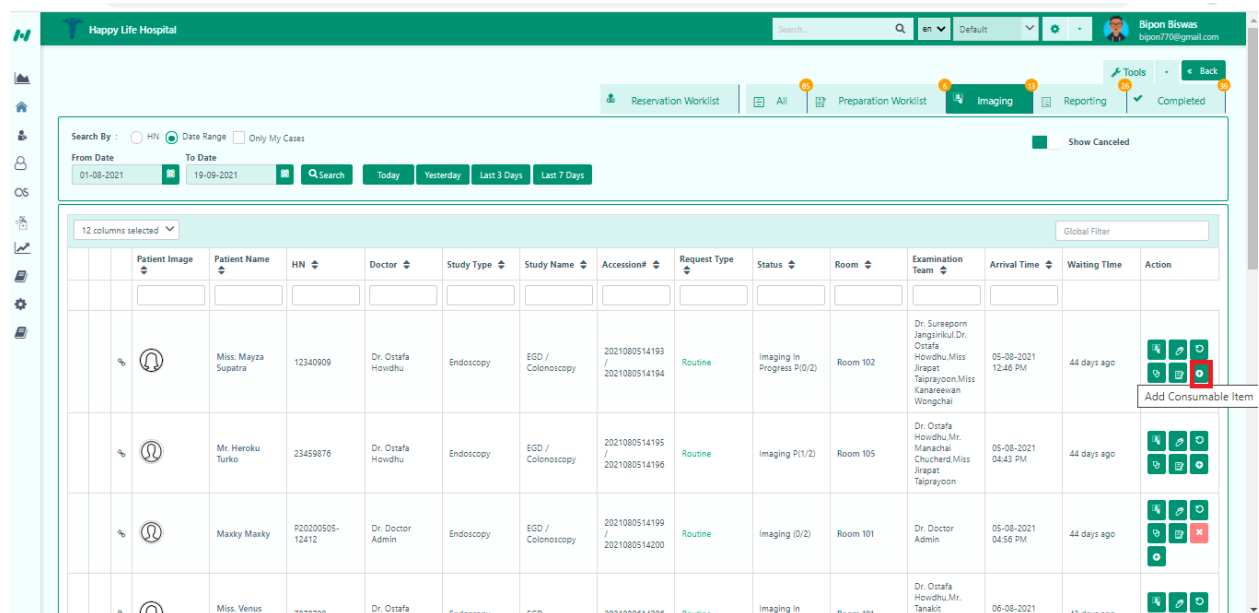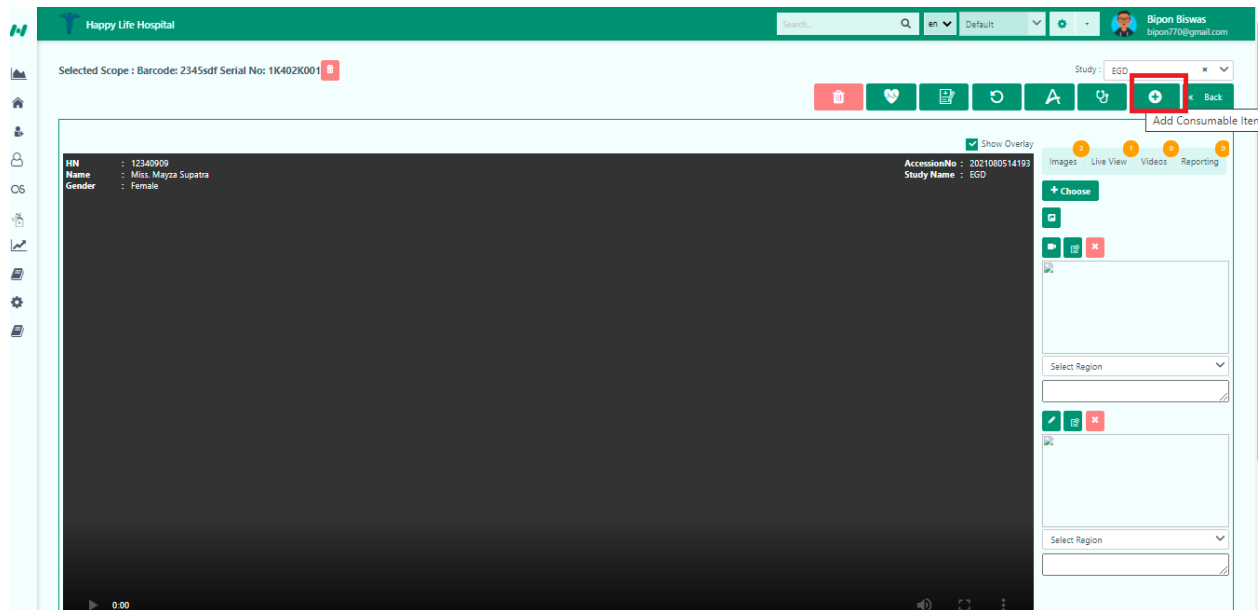

Patient List

The user can check the patient lists from here. They can Active/Inactive, Delete and Add to blacklist from this patient list.

Patient List Show Deleted Patient Back

7 columns selected Global Filter

Patient Image	Hn	Patient Name	Gender	Blood Group	Email	Status	Action
	P2021091712111	นายแพทย์ Jone Doe				Active	
	P2021091512110	Mahfuz Chan				Active	
	P2021091512109	Mr. Nathapong Sirisomboon	Male			Active	
	09864328	นายแพทย์ Sataporn Boongard				Active	
	96782819	นายแพทย์ Sataporn Boongard				Active	
	098737593	นายแพทย์ Sataporn Boongard	Male			Active	
	P2021090812108	Mr. Harry Williams				Active	
	P2021090312107	Mr. Sanan Supasawat	Male			Active	
	P2021090212106	Miss. Lalia Kundu				Active	
	123456	นายแพทย์ Sataporn Boongard	Male			Active	

Moreover, they can also see the deleted patient by simply clicking in the **Show Deleted Patient**.

Patient List Show Deleted Patient Back

7 columns selected Global Filter

Patient Image	Hn	Patient Name	Gender	Blood Group	Email	Status	Action
	P2021091712111	นายแพทย์ Jone Doe				Active	
	P2021091512110	Mahfuz Chan				Active	

Dynamic Dashboard

The user can see the daily update from this dashboard. They can know about the Recovery Room, Arrival Worklist, Imaging Worklist.

Add Post

Here, the user can post content regarding health.

Add Post

Thumbnail Image : + Choose


Title :

Subtitle :

Content :

Post


The user can add a thumbnail image, title of the content, subtitle and main content.

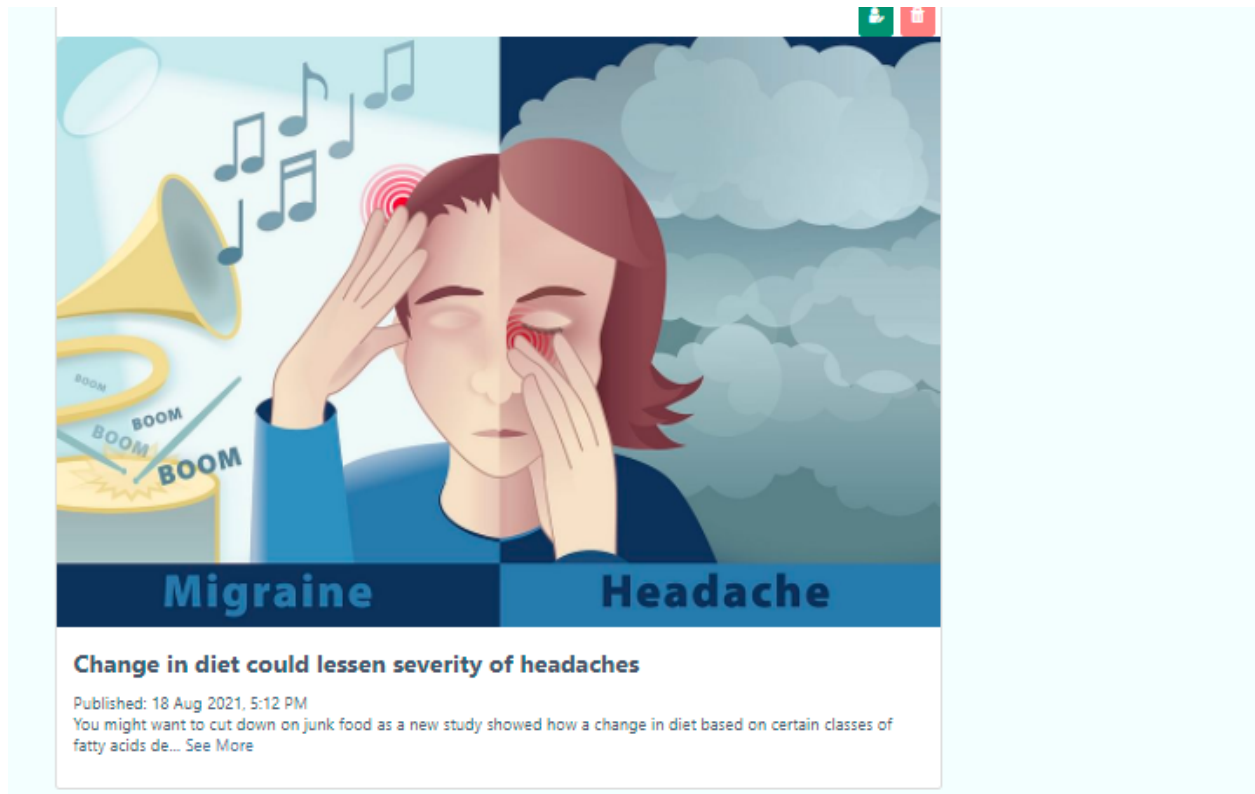


Calendar Settings

The user can customize the calendar with-

- **User Experience**
- **Add holiday/Block Room**
- **Calendar Cell Color Setup**
- **Calendar Event Color Setup**

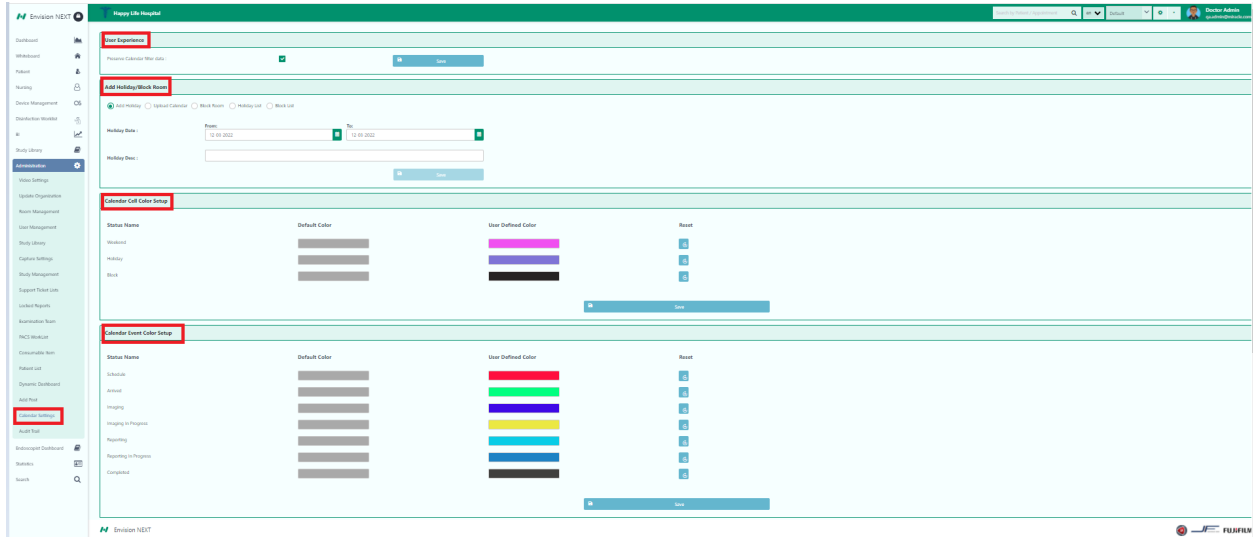

Example-

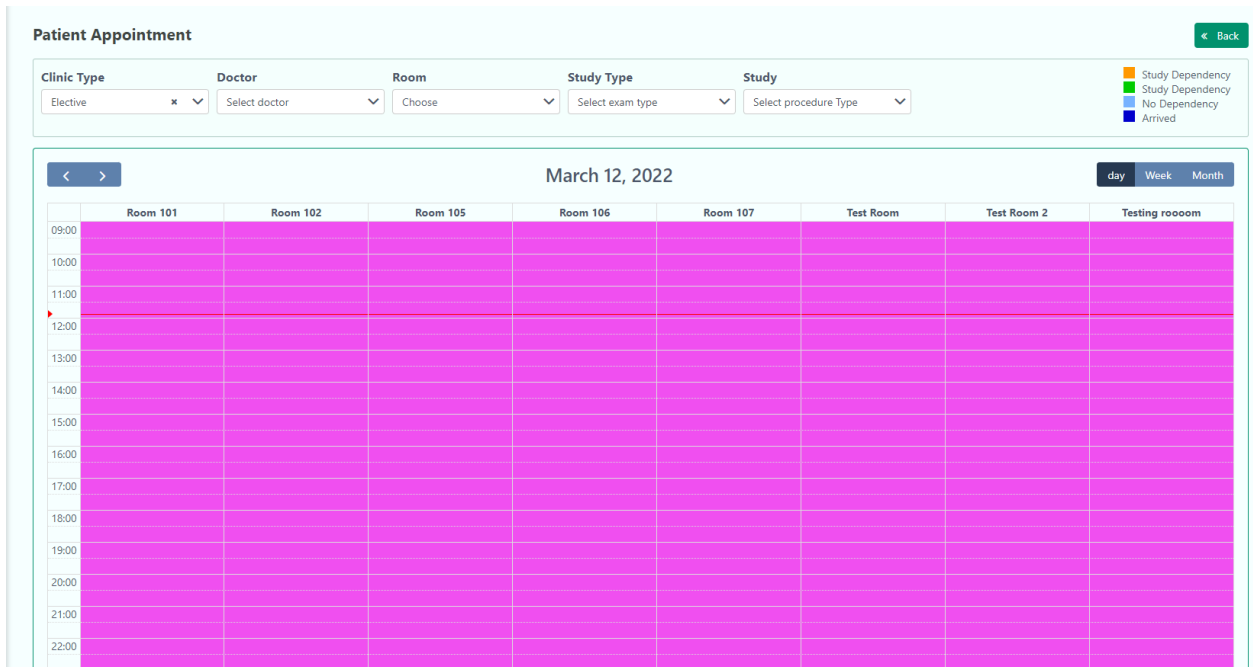

Weekend

Audit Trail

In this section, the user can know the details of usage of this application.

Audit Trail

Select A User (default all): Tasin Tazwar

From Date

01-03-2022

To Date

12-03-2022

Search

Today

Yesterday

Last X Day

0

IP Address	IP Port	Page Visited	Session On	Log-out Date & Time	Duration (d: h:m:s)
 Tasin Tazwar					
192.168.3.250	54989	Dashboard	01/03/2022 at 4:58 AM	01/03/2022 at 8:55 AM	0: 3:57:18
192.168.3.250	54989	login	01/03/2022 at 4:58 AM	01/03/2022 at 8:55 AM	0: 3:57:18
192.168.3.250	54989	whiteboard	01/03/2022 at 4:58 AM	01/03/2022 at 8:55 AM	0: 3:57:18
192.168.3.250	56799	Dashboard	01/03/2022 at 5:08 AM	01/03/2022 at 8:55 AM	0: 3:46:42
192.168.3.250	56799	endo-reporting-five	01/03/2022 at 5:08 AM	01/03/2022 at 8:55 AM	0: 3:46:42
192.168.3.250	56799	Logout	01/03/2022 at 5:08 AM	01/03/2022 at 8:55 AM	0: 3:46:42
192.168.3.250	56799	Patient Dashboard	01/03/2022 at 5:08 AM	01/03/2022 at 8:55 AM	0: 3:46:42

The above right corner, the user can view the data visually by clicking this  button.

Most Visited Page

Most Time Spent Pages

Most Time Spent Page (min)

Pie

Dashboard login whiteboard endo-reporting-five Logout Patient Dashboard Patient Register Post Scope
Digital Signature global-report Image Capture main-page patient-appointment-by-room search trending Update Organization
user-experience User Profile User Profile Update



Endoscopist Dashboard

In this section, the user will find the details of the endoscopy of a patient.


The screenshot displays the 'Endoscopist Dashboard' interface. The top header includes the 'Envision NEXT' logo, 'Happy Life Hospital' name, a search bar, and user information for 'Bipon Biswas'. The left sidebar lists various modules, with 'Endoscopist Dashboard' highlighted. The main content area features a date range filter (01-01-2022 to 01-03-2022) and a table of endoscopy studies. The table has columns for Study Date, Accession No., HN, Patient Name, Date of Birth, Gender, Study Id, Study Name, Status, and Action. The 'Action' column contains icons for Study Image, View Report, View Schedule, View Preparation Data, View Recovery Data, PACS, Resources, and Vital Sign. The 'Endoscopist Dashboard' option is highlighted in the left sidebar.

Study Date	Accession No.	HN	Patient Name	Date of Birth	Gender	Study Id	Study Name	Status	Action
28-02-2022	20220228146	P2022022812	Jack Beer	28-02-2002	Male	21172	Colonoscopy	Reporting	[Icons]
28-02-2022	20220228146	P2022022812	Jack Beer	28-02-2002	Male	20171	EGD	Reporting	[Icons]

The user can operate the dashboard by clicking these buttons.

- **Study Image**
- **View Report**
- **View Schedule**
- **View Preparation Data**
- **View Recovery Data**
- **PACS**
- **Resources**
- **Vital Sign**

Study Image- The user can view the images of the study. They can also download the selected files or can download all files from here.

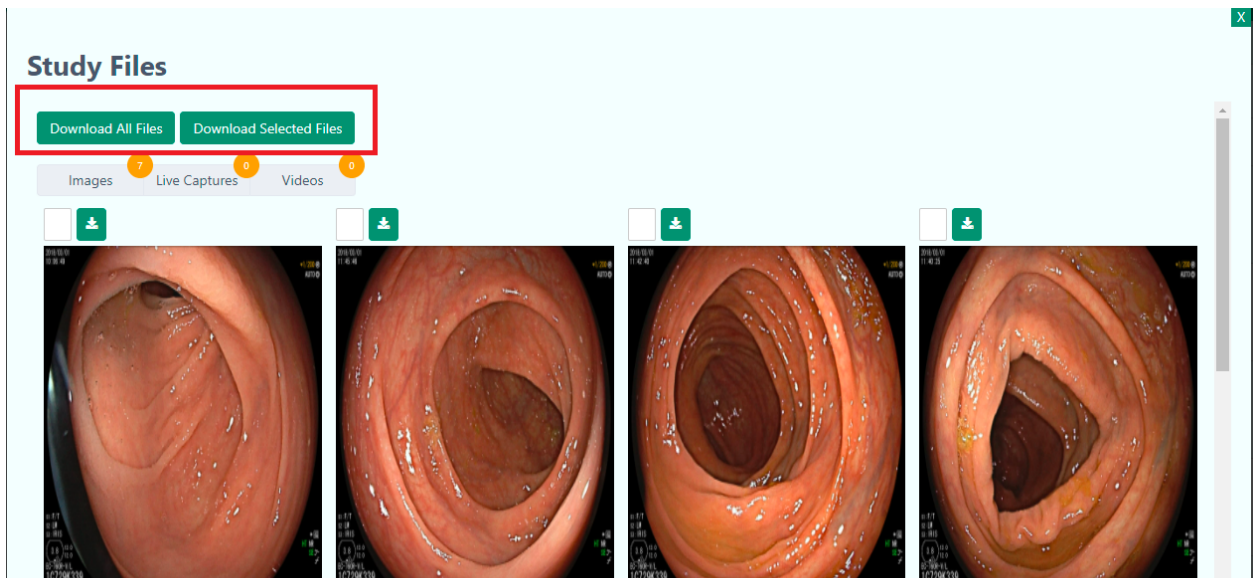

View Report

The user can view the report of the patient by clicking  this button.

View Schedule

The user can view the patient appointment and study details from here.

Patient Appointment

Patient Information

Name : Jack Beer
Phone :
Address :



Doctor Information

Name : Dr. Doctor Admin
Phone : +9999999999
Department : Internal Medicine



Study

Date* 28-02-2022
Duration (In Minutes)*
Doctor* Dr. Doctor Admin
Study Type* Endoscopy
Clinic Type* Elective
Chief Complaints*

Start Time* 22:00
End Time* 23:00
Room* Room 105
Study* Colonoscopy
Request Type* Routine

Study Type	Study Name	Start Time	End Time
Endoscopy	EGD	10:00 PM	10:30 PM
Endoscopy	Colonoscopy	10:30 PM	11:00 PM

View Preparation Data

The user can view preparation data from here.

View Recovery Data

The user can view the recovery data by clicking  this button. It will redirect to the **Post Scope** page.

Post Scope  [Back](#)

Patient :	Jack Beer	Age :	20 years	
HN :	P2022022812287	SSN :		
Allergy  :	Yes			

Study Information

Code :	20171.21172
Name :	EGD,Colonoscopy
Datetime :	28-02-2022, 28-02-2022

Doctor Information

Name :	Dr. Doctor Admin	
Phone :		
Department :		

The user can also add **consumable** items from the **Post Scope** page.

Post Scope  [Back](#)

Patient :	Jack Beer	Age :	20 years	
HN :	P2022022812287	SSN :		
Allergy  :	Yes			

Study Information

Code :	20171.21172
Name :	EGD,Colonoscopy
Datetime :	28-02-2022, 28-02-2022

Doctor Information

Name :	Dr. Doctor Admin	
Phone :		
Department :		

Time in Recovery Room	Recorded Nurse Bipon Biswas	Recorded Date 03/03/2022	Recorded Time 11:39
------------------------------	---------------------------------------	------------------------------------	-------------------------------

PACS(PICTURE ARCHIVING & COMMUNICATION SYSTEM)

PACS is a web-based archiving and communication system that improves an imaging centre's productivity and profits by offering correct patient care. ENVISION NEXT PACS is used to assist hospitals in the capture, management, storage, and viewing of diagnostic images.

Following the scan, diagnostic pictures are generated from modalities and transferred to PACS in electronic form. After receiving diagnostic pictures via PACS, technologists can verify the following imaging with the operations done. After reviewing and validating the images, radiologists can begin reading them and making diagnoses.

Resource Allocation

The user can view endoscopy preparation from here which includes- examination team, scope list, and endoscopy machine.

Vital Sign

The user can add vital signs from here.


Vital Sign

×

Patient Name : HN :

Gender :

DOB :

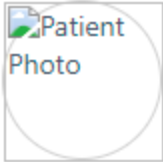



+ Add

Save

Height (cm)		Weight (kg)	
BMI		BP (mmHg)	Systolic / Diastolic
SPO ₂		Heart Rate (BPM)	
Temperature (F)		Temperature (C)	
Respiration Rate (RPM)		Pulse rate (BPM)	
Capillary Refill Time (Sec)		BCS	

Statistics

The user can view all reports from **Reports** dropdown.

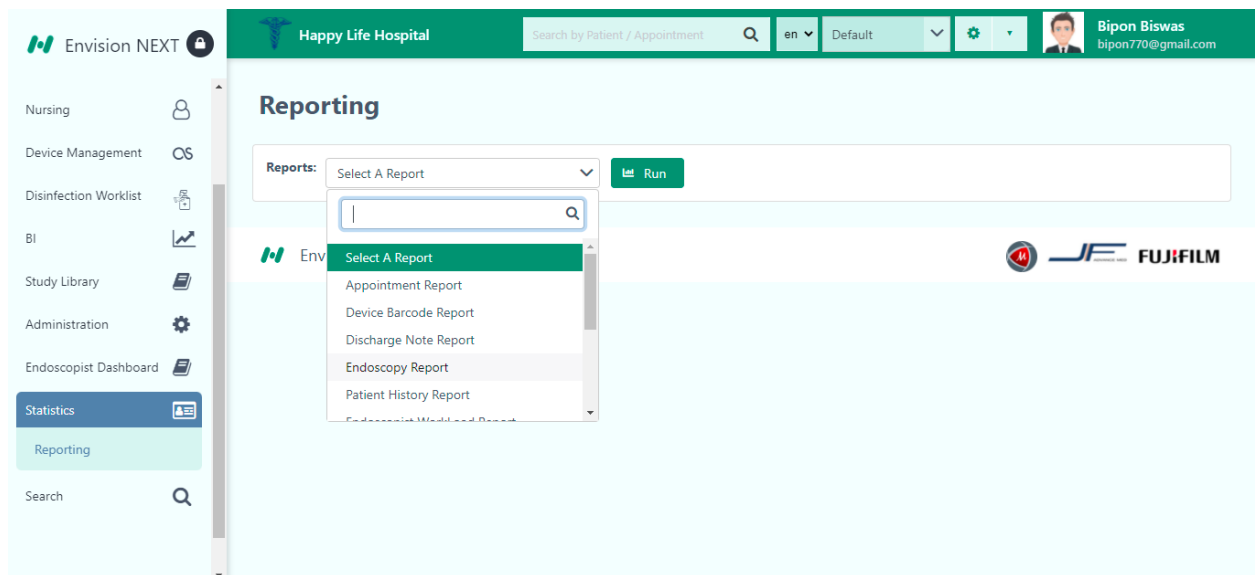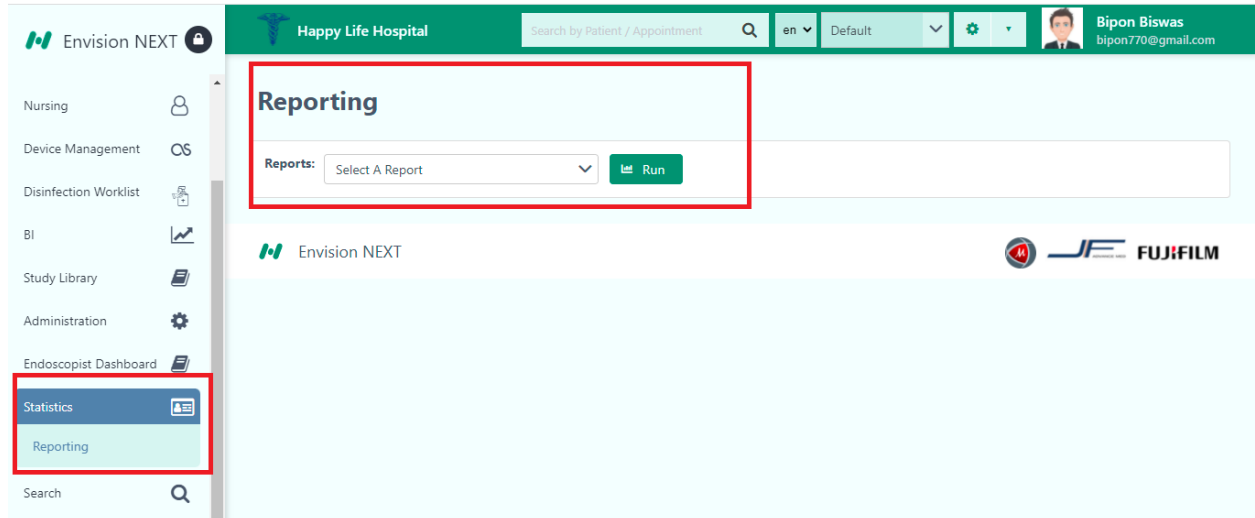


Search

The user can search anything regarding study from here.


ENVISIONNEXTBack

[All](#) [Images](#) [Videos](#)

EGD/March 02, 2022/Miss. NookNick Maha/ อภิสิทธิ์ อักษรระ
Findings : Oropharynx and Hypopharynx: Normal Esophagus: 1) One Barrett's esophagus at Diffused esophagus 3 cm from incisor. Intervention : EMR Technique: Banding, Device: **Cold** snare, Result: complete piece-meal resection, Specimen size: 2 mm, Defect closure: Yes.Clipping Type: TTS, Number: 1, Result: Success.2) Caustic injury at Middle esophagus, Lower esophagus, Diffused esophagus 3 cm from incisor. Zargar's classification -II B (Circumferential lesions), 3) Esophagitis at Middle esophagus, Lower esophagus, Diffused esophagus 3 cm from incisor. Characters **-non-erosive**. Stomach: 1) Three Cratered Ulcer size 0.2,0.9 cm at Cardia stomach, Body stomach, Anterior wall of body stomach. Forrest classification -Forrest II b (Adherent clot). Intervention : Banding: Number of band: 3, Result: Success.2) One Scarring size 0.5 cm at Cardia stomach. Body stomach, Anterior wall of body stomach.....

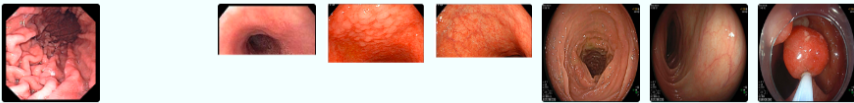


ERCP/May 27, 2022/ Error2 Test2/Dr. Doctor Admin
Findings : Anatomy: Post-surgical. Post surgical type: Billroth 1, Billroth 2. Complications: Yes: Complication Type: Cardiopulmonary: Event: Bradycardia, Cardiac ar.....

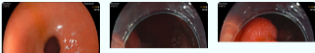

They can also search images from here only.

ENVISIONNEXTBack

[All](#) [Images](#) [Videos](#)

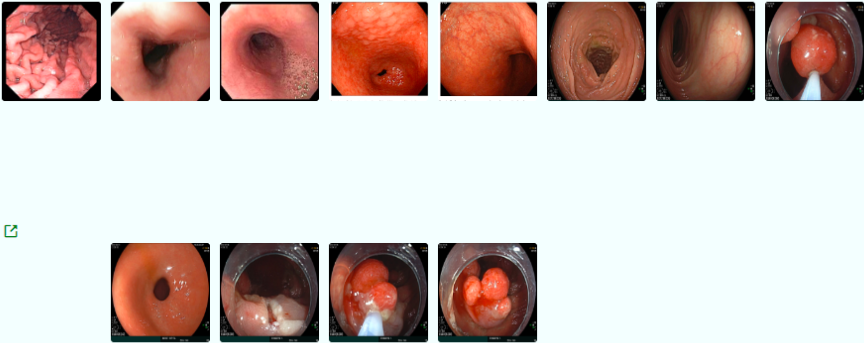

Finally, they can also find video contents from advance search.

