



Share and access medical imaging data securely.



System Administrator Manual



A Product of



**MIRACLE ADVANCE
TECHNOLOGIES CO. LTD.**

Integrate. Accelerate. Differentiate.



Version 1.0

Document Control

Release No.	Date	Purpose/Change	Author
V 1.0	Sep 8, 2025	Initial Release	Shahaneaz Shaikat





Table of Contents

1. Onboarding	3
1.1. Login to the Patient Portal	3
1.2 How to configure AD/LDAP for New Site	3
1.3 Application Versions and Release Date	5
2. Dashboard	5
2.1 Overview	5
2.2 Get Started	5
3. Patient Portal Whiteboard	8
3.1 Create New Order	8
3.2 View List	11
3.3 Filter	11
3.4 Search	12
3.5 View and Update Patient Information and	13
3.6 Shareable link or QR code generator buttons on the whiteboard	14
4. Patient Portal Settings	16
4.1 Setup and Update Organization	16
4.1.1 Update organization information	16
4.1.2 General information	17
4.1.3 API Settings page	19
4.1.4 Workflow settings(MRay settings)	20
4.1.5 General Setting	21
4.1.6 mRay Setting	22
4.1.7 Password Policy	23
4.1.8 Package Subscription	24
4.1.9 Email Settings	24
4.1.10 Set Up Email	25
4.2 User List of Patient Portal	26
4.2.1 Adding New Users	26
4.2.2 Update Existing Users	28
4.3 Exam Information	29
4.3.1 Exam List	29
4.3.2 Adding New Exam	30
4.3.3 Edit an exam	31
5. View Patient Information	32
6. Technical Documentation	33
6.1 API Security	33



1. Onboarding

1.1 Login to the Patient Portal

Log in using your login credentials (User Name and Password)

Patient Portal

Login to your **Account**

👤 | IdapUser

🔒 |

Login

Copyright 2025 © Miracle Advance Technologies Co. Ltd.

Now, you should see the Whiteboard after logging in to the application.



1.2 How to configure AD/LDAP for New Site


To configure a user with role-based access control, we have to set up AD/LDAP settings. You have to configure whether this site is on-premises or not. Now, you need to enter the server IP, Port number, and Domain Name for LDAP. Once you've adjusted these settings, click the **"Save"** button to finalize the configuration.

After saving the configuration, you can configure RBAC. To configure, click on the RBAC button. Here, you can create groups and decide who can access the Patient Portal.



1.3 Application Versions and Release Date

You can find the Patient Portal version and release date in the footer of the Patient Portal App. The most recent version is 1.13.0, and the latest release date is Aug 08, 2025.

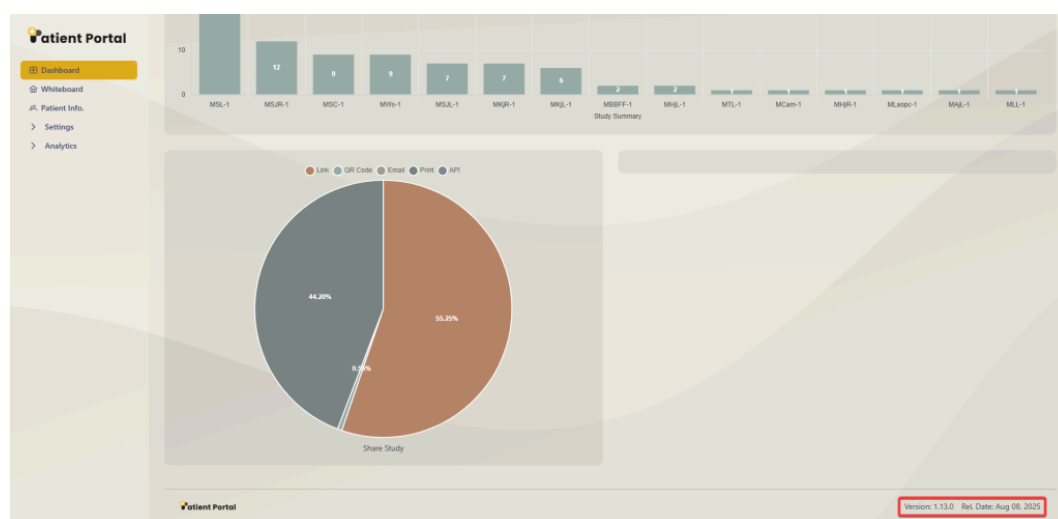


2. Dashboard

2.1 Overview

The Dashboard gives you an at-a-glance view of study volume and sharing activity across your sites. Use it to track the creation and sharing of studies (via Link, Email, QR Code, Print, API) and the timing of activity.

2.2 Get Started


Open the Dashboard after signing in to the Patient Portal.

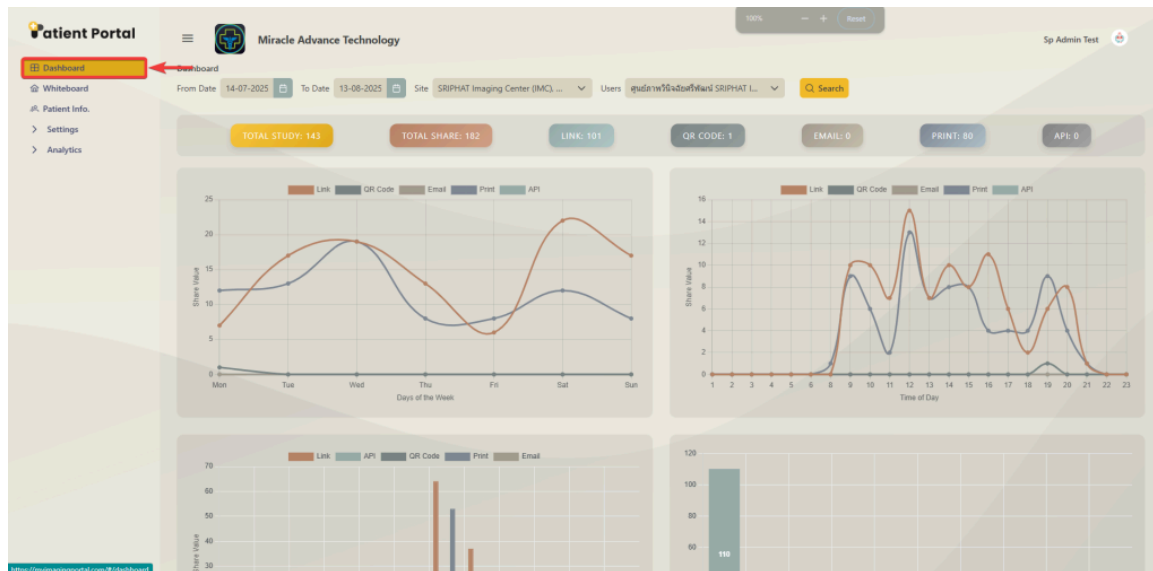

You can apply filters like (date range, site, users) to focus on the data you need.



In the next section, you can review the Total Study, Total Share Count, Shares by method (Share count): Link, Email, QR Code, Print, API.

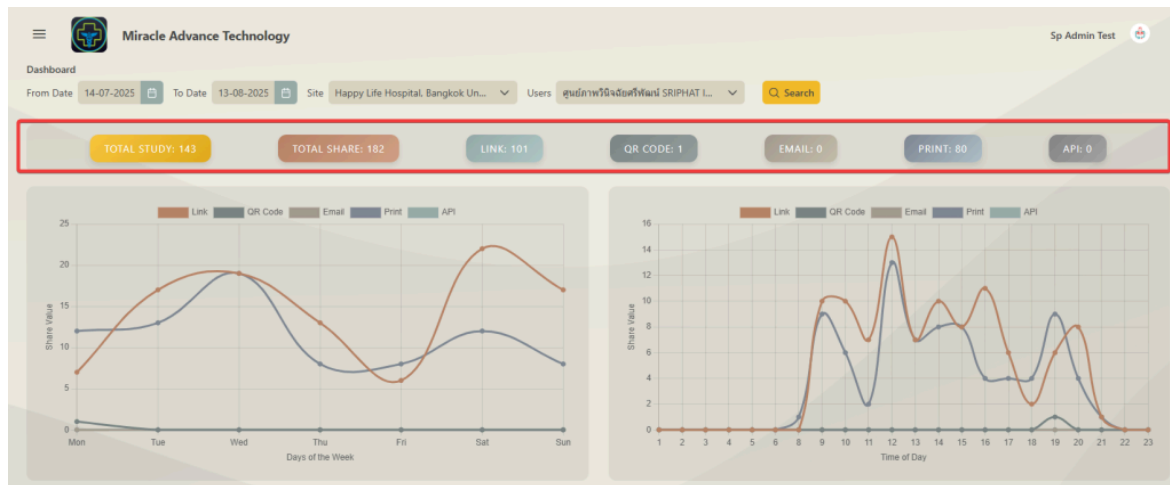

In the next section, you can review the KPI graph cards and charts to understand volume and trends.




3. Patient Portal Whiteboard

3.1 Create New Order

In this section, we will see the workflow for ordering a new case in the Patient Portal. The orders are created using the [Patient Portal Public API](#).

To create a new order, you have to use the Patient Portal API.

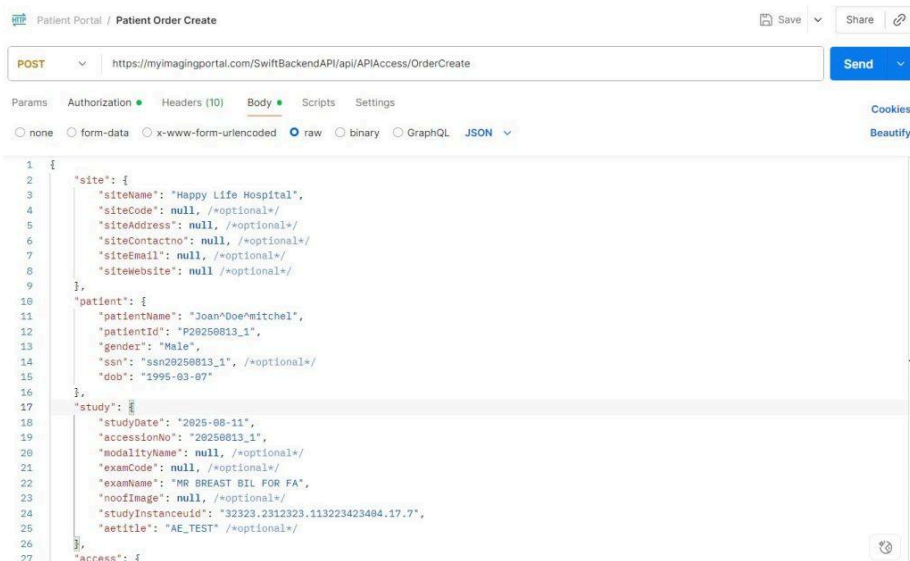


REST client interface showing a POST request to `https://myimagingportal.com/SwiftBackendAPI/api/APIAccess/OrderCreate`. The request body is a JSON object:

```

1 {
2   "site": {
3     "siteName": "Happy Life Hospital",
4     "siteCode": null, /*optional*/
5     "siteAddress": null, /*optional*/
6     "siteContactno": null, /*optional*/
7     "siteEmail": null, /*optional*/
8     "siteWebsite": null /*optional*/
9   },
10  "patient": {
11    "patientName": "JoanDoe@mitchel",
12    "patientId": "P20250813_1",
13    "gender": "Male",
14    "ssn": "ssn20250813_1", /*optional*/
15    "dob": "1995-03-07"
16  },
17  "study": {
18    "studyDate": "2025-08-11",
19    "accessionNo": "20250813_1",
20    "modalityName": null, /*optional*/
21    "examCode": null, /*optional*/
22    "examName": "MR BREAST BIL FOR FA",
23    "nooImage": null, /*optional*/
24    "studyInstanceid": "32323.2312323.113223423484.17.7",
25    "aetitle": "AE_TEST" /*optional*/
26  },
27  "access": {


```

The APIs are secured and require authentication. Use your username and password to submit an order. If you don't have credentials, contact the Sales team.

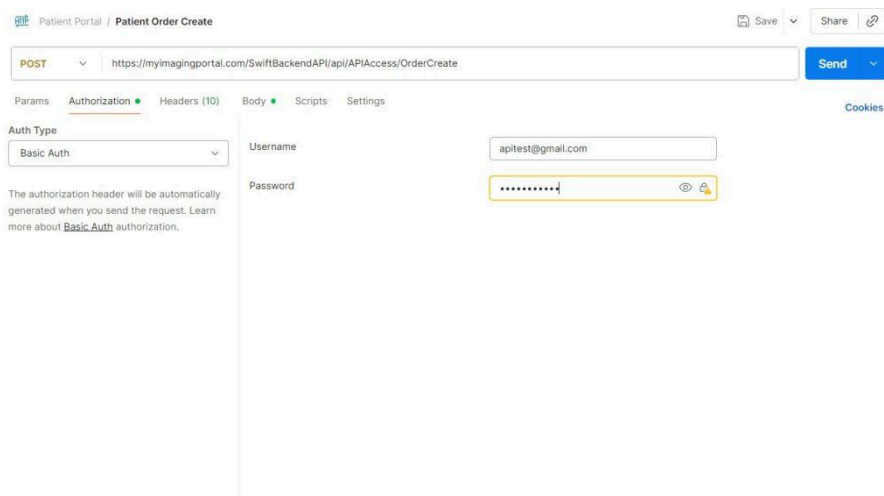


REST client interface showing the Authorization tab. The Auth Type is set to Basic Auth. The Username field is filled with `apitest@gmail.com` and the Password field is masked with dots. A note states: "The authorization header will be automatically generated when you send the request. Learn more about [Basic Auth](#) authorization."

After configuring your credentials, provide the required information and send the request.

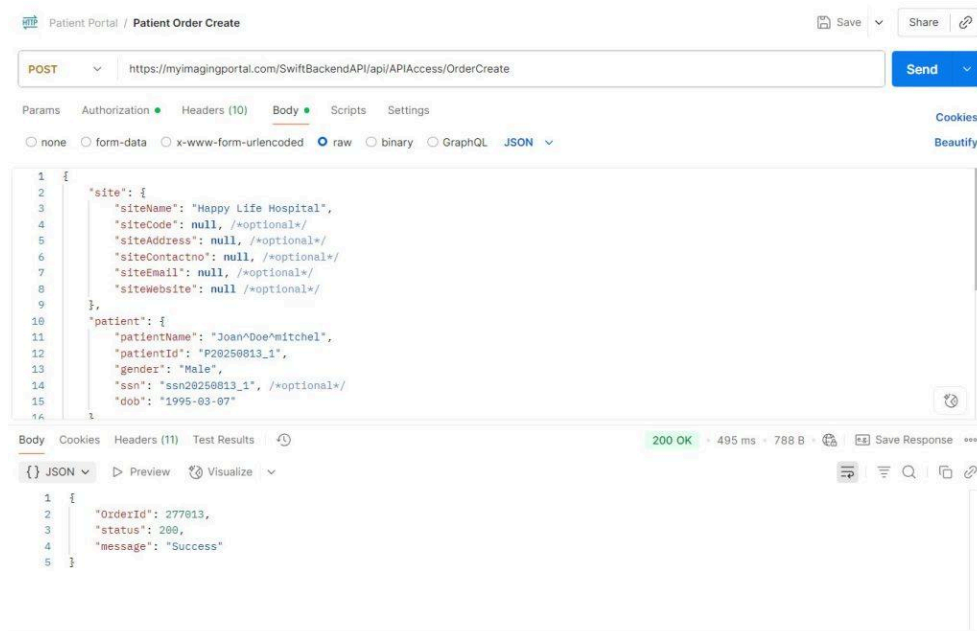

If the request is successful, the new order will appear in the Patient Portal, where you can track its status through the workflow.

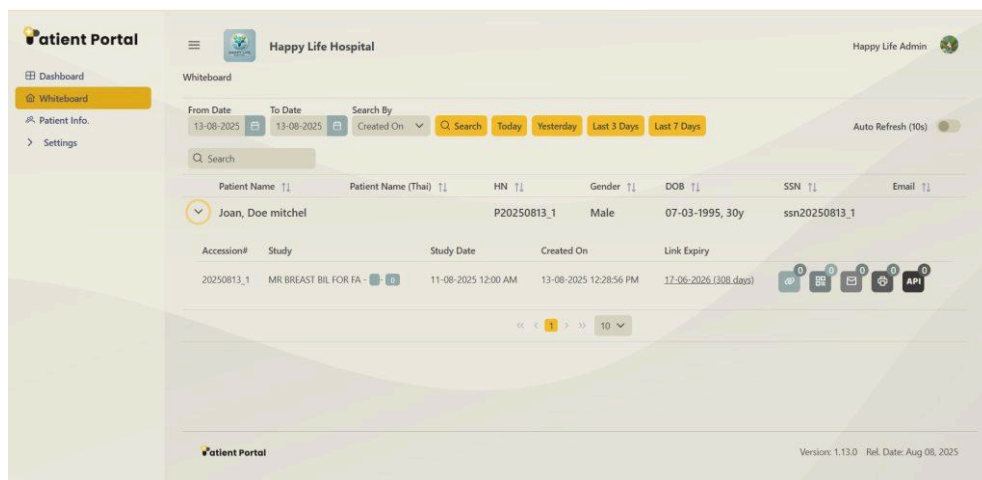



3.2 View List

The studies are listed on the whiteboard. Once a study is uploaded through DICOM or API, it will appear there.

Patient Name	Patient Name (Thai)	HN	Site	Gender	DOB	SSN	Email
>		4020175		Male	07-03-1951, 74y		
>		2115711		Female	02-10-1969, 55y		
>		3142526		Female	14-04-1957, 68y		
>		2862153		Female	18-11-1959, 65y		
>		3108376		Male	12-11-1990, 34y		


3.3 Filter

On the whiteboard, you can filter studies by:

- **Date:** Select a From Date and To Date, or use presets such as Today and Last 7 Days.
- **Site:** Choose a location from the Sites dropdown.
- **Keyword:** Type a patient's name or Hospital Number (HN) into the Search field.


Patient Name	Patient Name (Thai)	HN	Site	Gender	DOB	SSN	Email
>		4020175		Male	07-03-1951, 74y		
>		2115711		Female	02-10-1969, 55y		
>		3142526		Female	14-04-1957, 68y		
>		2862153		Female	18-11-1959, 65y		
>		3108376		Male	12-11-1990, 34y		

Auto-Refresh: Enable the Auto-Refresh toggle to view new studies as they become available.

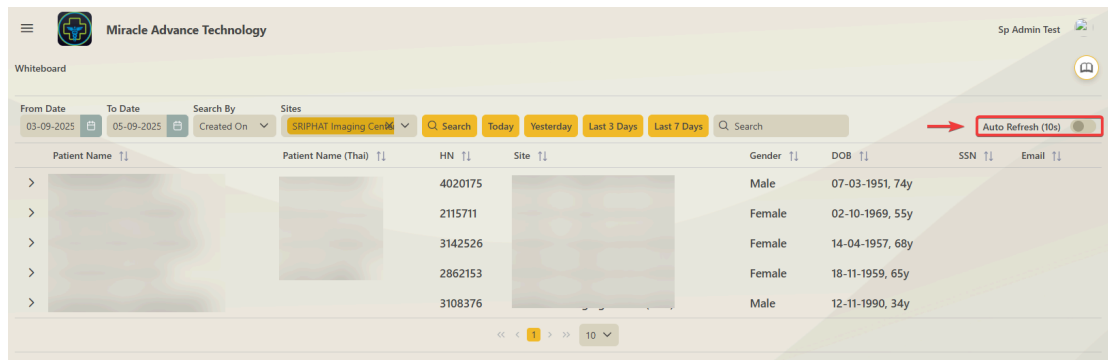


3.4 Search

The search function enables you to locate patient information using the following criteria: Patient Name, Patient Name (Thai), Site, Gender, Date of Birth (DOB), Social Security Number (SSN), and Email address.

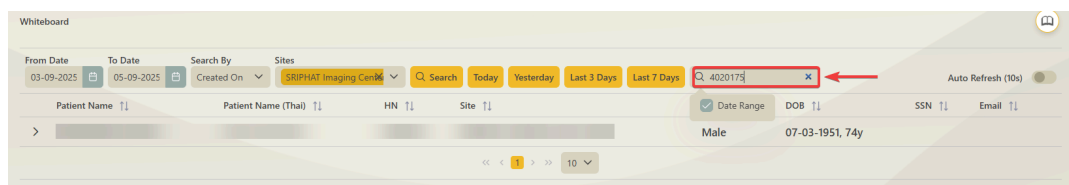

3.5 View and Update Patient Information and

Click on the patient's name to view their details.

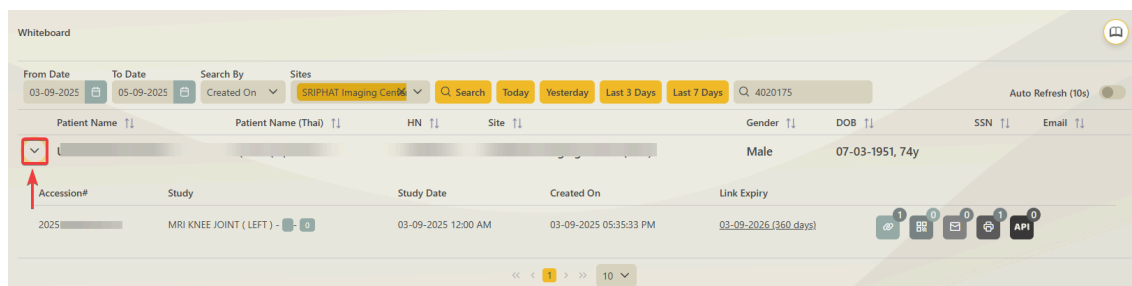

Once clicked, you can modify the patient's details.




To save changes, click on the **Update** button.


3.6 Shareable link or QR code generator buttons on the whiteboard

On the Patient Portal whiteboard, you have the option to share reports with the patient. You can share reports with patients by:


- Copying a secure link
- Copying, emailing, or printing a QR code
- Retrieving access link via the API

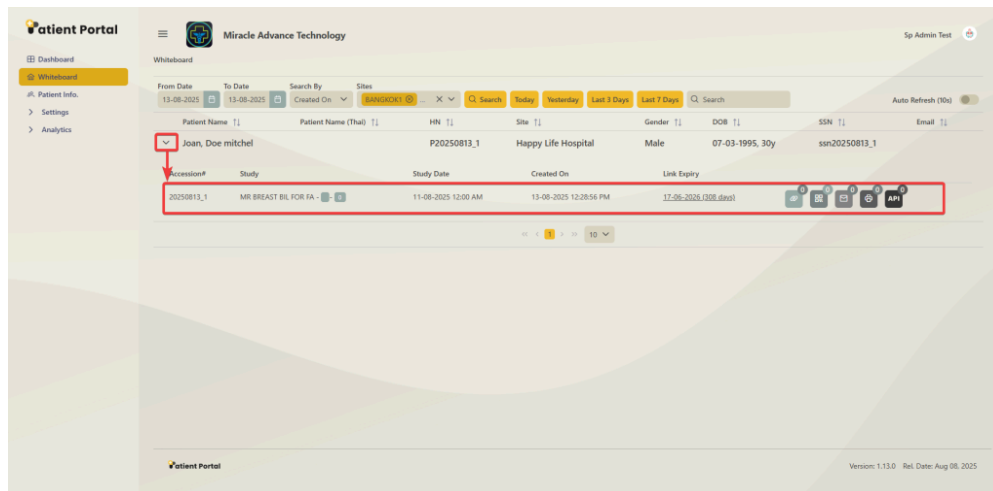

You can copy links directly. You can also copy, mail, or print the QR code. There is another way to track how many times the api has been used to retrieve links.

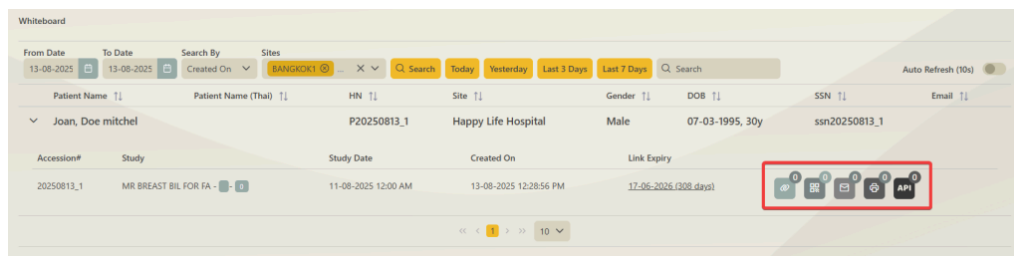

A counter next to each method shows how many times it's been used. After you click the Copy button for the link, the counter updates immediately (for example, from 0 to 1). Each action, copy link, email QR code, and print QR code increments its respective counter. API usage is tracked separately and increases each time the API successfully retrieves the link.

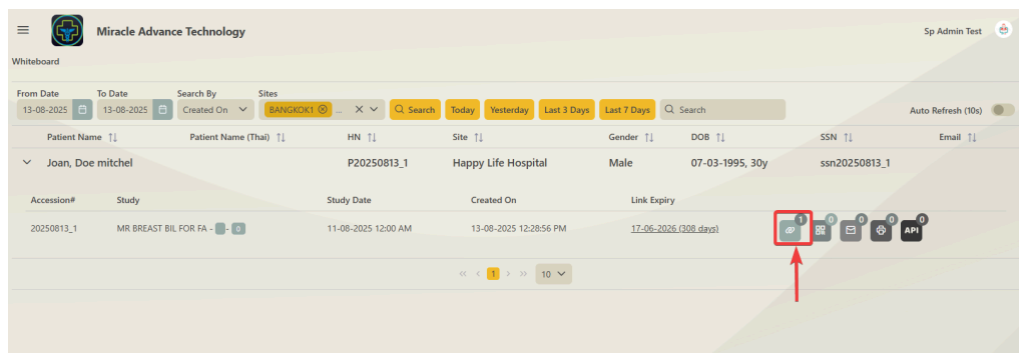




4. Patient Portal Settings


4.1 Setup and Update Organization

This page displays all the existing organizations. You can also update organization information on this page.

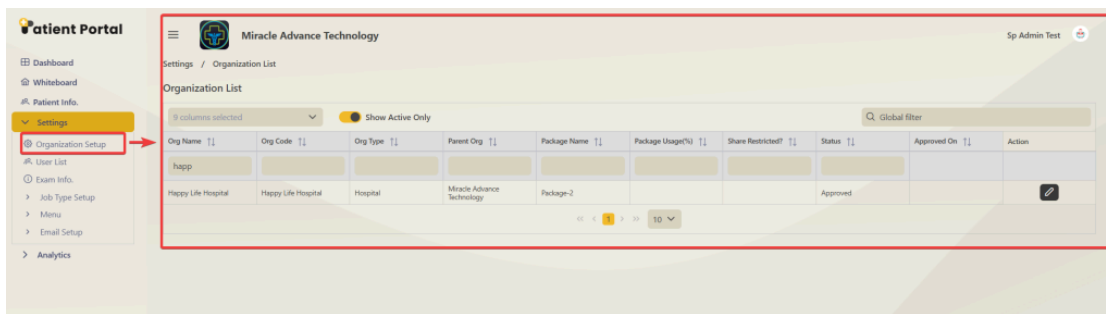

4.1.1 Update organization information

To update an organization, click on the pen icon.

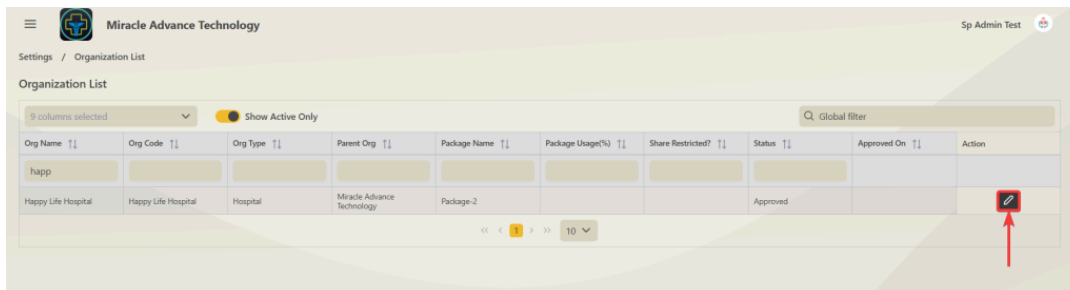

When the edit button is clicked, a pop-up page appears, allowing you to update the general information, API settings, workflow, Package information, and email settings for that specific organization.

Client Portal

Update Organization Information

General Information API Setting Workflow Package Subscription Email Setting

Org Type* Hospital

Org Name* Happy Life Hospital

Tel. No.

Line ID

Org Suffix

Logo Size: 50*50 + Choose

Logo (Bottom Left) Size: 130*220 + Choose

App Guide QR Code Size: 50*50 + Choose

Parent Org* Miracle Advance Technology

Web URL

Email

Org Prefix

Active? ☒

Logo (Top Left) Size: 32*140 + Choose

Logo (Bottom Right) Size: 30*160 + Choose

Organization Address

4.1.2 General information

The general information page allows you to update the organization's information, Parent Organization information, logo, and address.

Client Portal

Update Organization Information

Org Type* Hospital

Org Name* Happy Life Hospital

Tel. No.

Line ID

Org Suffix

Logo Size: 50*50 + Choose

Logo (Bottom Left) Size: 130*220 + Choose

App Guide QR Code Size: 50*50 + Choose

Parent Org* Miracle Advance Technology

Web URL

Email

Org Prefix

Active? ☒

Logo (Top Left) Size: 32*140 + Choose

Logo (Bottom Right) Size: 30*160 + Choose

Organization Address

Country Select Country Province Select Province

You can also set the organization as inactive by clicking on the **Active?** check mark.

Update Organization Information

General Information API Setting Workflow Package Subscription Email Setting

Org Type* Hospital

Org Name* Happy Life Hospital

Tel. No. Tel. No.

Line ID Line ID

Org Suffix Org Suffix

Logo Size: 50*50 + Choose

Logo (Bottom Left) Size: 130*220 + Choose

App Guide QR Code Size: 50*50 + Choose

Parent Org* Miracle Advance Technology

Web URL Web URL

Email

Org Prefix Org Prefix

Active? ☒

Logo (Top Left) Size: 32*140 + Choose

Logo (Bottom Right) Size: 30*160 + Choose

Organization Address

4.1.3 API Settings page

All incoming and outgoing APIs for a specific organization can be managed here.

Update Organization Information

General Information API Setting Workflow Package Subscription Email Setting

Submit API Settings

Auth API Auth API

Auth Type Select Auth Type

Submit API Submit API

Submit API Type Select Submit API Type

Base URL Base URL

User Name UserName

Password Password

Secret Key

Receive API Settings

Base URL Receive Base URL

User Name Receive User Name

Password Receive Password

Receive API Receive API

Receive API Type Select Receive API Type

Update Close

Here you will find the **Secret Keys** to be used as the bearer token.

Patient Portal

Update Organization Information

General Information | **API Setting** | Workflow | Package Subscription | Email Setting

Submit API Settings

Auth API	Auth API	Base URL	Base URL
Auth Type	Select Auth Type	User Name	UserName
Submit API	Submit API	Password	Password
Submit API Type	Select Submit API Type	Secret Key	Secret Key

Receive API Settings

Base URL	Receive Base URL	Receive API	Receive API
User Name	Receive User Name	Receive API Type	Select Receive API Type
Password	Receive Password		

Update Close

4.1.4 Workflow settings(MRay settings)

Use Workflow Settings to control how the Patient Portal behaves (refresh, timeouts, display formats), configure the mRay integration for launching studies, and enforce password rules. Changes here affect all users in your organization.

Patient Portal

Update Organization Information

General Setting | Workflow | Package Subscription | Email Setting

General Setting

Web Link Report Paper Size	A5	Web Link Report Site Info. in Top?	☐
Auto Refresh?	☐	Auto Refresh Second	10
Display Name Format	Last, First Middle	Idle Time Out (Second)	180 Default Time Out 180

mRay Setting

mRay Full URL	https://start.mray.app/?serverurl={serverurl}&login={login}&studyuid={studyuid}	Server URL	portal.jfav.co.th
Login User	napichead	Pre Shared Key	*****
Link shared to Patient Level	☐		

Password Policy

Enforce Policy?	☐	Upper Case?	☐
Lower Case?	☐	Number?	☐
Special Character?	☐	Length	Length

Update Close



4.1.5 General Setting

General Setting	
Web link Report Paper Size	A5
Auto Refresh?	☐
Display Name Format	Last, First Middle
Web link Report Site Info. in Top?	☐
Auto Refresh Second	10
Idle Time Out (Second)	180
	Default Time Out 180

- Web link Report Paper Size
- Web link Report Site Info. in Top?
 - When enabled, it adds your site information to the top header of the web-link report.
- Auto Refresh?
 - Enables automatic refresh for pages that support live updates (e.g., whiteboard/queues).
- Auto Refresh Second
 - Interval in seconds between automatic refreshes. Example: 10.
 - Only used when Auto Refresh is enabled.
- Display Name Format
 - Controls how patient names are shown across the portal and on printed reports.
 - Example options: “Last, First Middle” or “First Middle Last.”
- Idle Time Out (Second)
 - Logs users out after the specified inactivity period. Example: 180 (3 minutes).
 - Use a longer timeout for workstations with frequent but brief inactivity; keep security in mind.

4.1.6 mRay Setting

mRay Setting	
mRay Full URL	https://start.mray.app?serverurl={serverurl}&login={login}&studyuid={studyuid}
Login User	napichead
Link shared to Patient Level	☐
Server URL	portal.jfav.co.th
Pre Shared Key	*****

Use these fields to integrate with the mRay viewer/app.

- mRay Full URL



- Template URL used to launch mRay with runtime values.
- Tokens are replaced automatically when links are generated.
- Common tokens:
 - {serverurl} – your mRay server base URL
 - {login} – the service/login user
 - {studyuid} – the study identifier
- Example pattern:

<https://start.mray.app?serverurl={serverurl}&login={login}&studyuid={studyuid}>
- Server URL
 - Your mRay server host (e.g., portal.example.com). Must be reachable by end users.
- Login User
 - Service account used for the integration (e.g., mray_service). Ensure it has the minimal required permissions.
- Pre Shared Key
 - Secret used by the integration. Keep it confidential; rotate it periodically.
 - Use the eye icon to show/hide while typing (if available).
- Link shared to Patient Level
 - When enabled, shared links open at the patient level (may include multiple studies).
 - When disabled, links are scoped to a single study.
 - Enable only if this aligns with your privacy policy and user permissions.

4.1.7 Password Policy

Password Policy	
Enforce Policy?	☐
Upper Case?	☐
Lower Case?	☐
Number?	☐
Special Character?	☐
Length	Length

These settings apply to user passwords in the Patient Portal.


- Enforce Policy?



- Turns on password complexity enforcement for new passwords and resets.
- Lower Case? / Upper Case? / Number? / Special Character?
 - Require inclusion of each selected character type.
- Length
 - Minimum required password length (enter a number).
 - Recommendation: 10–14 characters with multiple character types.

4.1.8 Package Subscription

You can easily find all the details about your package subscriptions here.

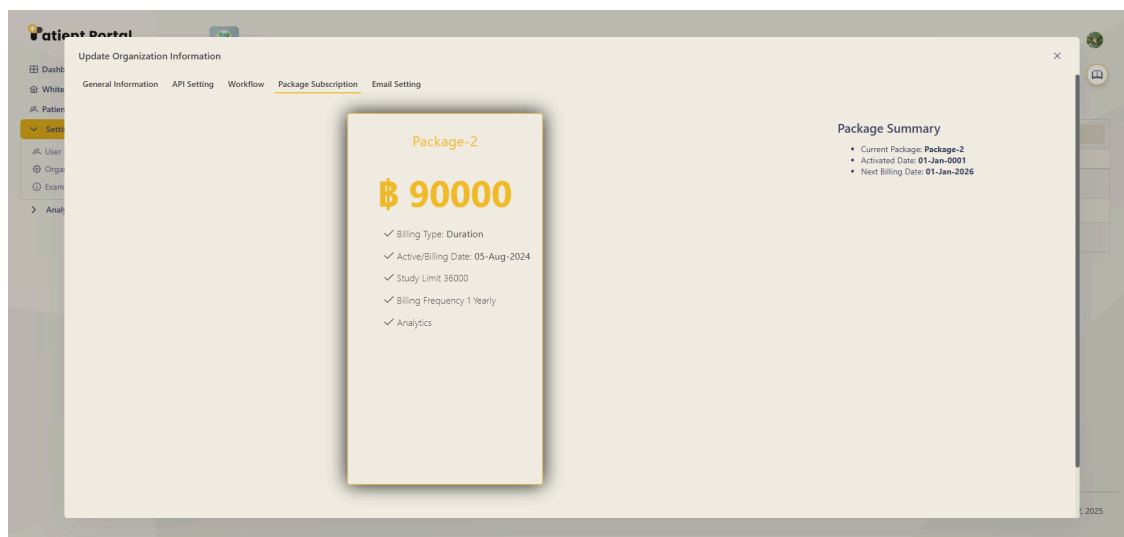

4.1.9 Email Settings

Use Email Settings to configure how the Patient Portal sends transactional emails (for example, emailing a QR link/report to a patient). Each setup ties your email provider, authentication, event mapping, and template.



4.1.10 Set Up Email

1. Click Add Email Setup.
2. Choose your provider.
3. Enter API credentials and base URL.
4. Map the Email Event to a Template.
5. Set the Queue Interval (how often queued emails are sent).
6. Click Configure Email Setup, then Validate Email Setup.
7. Save, then Update to apply changes to the organization.

Example (SendGrid)

- Sender Email: SendGrid
- Email Send From: noreply@yourdomain.com
- API Auth Type: Bearer Token
- API Base URL: <https://api.sendgrid.com/v3>
- API Bearer Token: SG.xxxxxx
- Email Event: Email QR Link Report to patient
- Email Template: Your Imaging Study Information
- Queue Interval: 30



4.2 User List of Patient Portal

All the users or roles can be reviewed, added, removed, or updated on this page.

Full Name	Email	Org Name	Org Type	Parent Org	Mobile No	Created On	Status	Action
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology			Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology	1234567890		Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology			Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology		Jun 27, 2025	Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology		Jun 27, 2025	Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology		May 1, 2025	Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology	9876543210	Apr 30, 2025	Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology			Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology	1234	Mar 19, 2025	Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology		Mar 19, 2025	Approved	

4.2.1 Adding New Users

A new user is added by clicking on the Add New button.

Full Name	Email	Org Name	Org Type	Parent Org	Mobile No	Created On	Status	Action
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology			Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology	1234567890		Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology			Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology		Jun 27, 2025	Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology		Jun 27, 2025	Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology		May 1, 2025	Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology	9876543210	Apr 30, 2025	Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology			Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology	1234	Mar 19, 2025	Approved	
John Doe	john.doe@miracleadvancetech.com	Test Org	Health	Healthcare Technology		Mar 19, 2025	Approved	

After clicking, an "Add New User" page appears. Before submitting the form, you must fill in all the information marked with a Red Asterisk(*).



Add User

Org Type*	Hospital	Organization*	Happy Life Hospital
User Role*	Select Role	User Level	Select Level
First Name*	First Name	Last Name*	Last Name
Phone No	Phone No	Email*	Email
Password*	Password	Confirm Password*	Confirm Password
Line Id	Line Id	Medical License	Medical License
Country	Select a country	Province	Select a Province
District	Select a District	Sub-District	Select a Sub-District
Subspecialties	Select Subspecialties	Zip Code	Select a Zip Code

4.2.2 Update Existing Users

To update a user, click on the pen icon

User List

[+ Add New](#)

[Clear](#)

Full Name	Email	Org Name	Org Type	Parent Org	Mobile No	Created On	Status	Action
Dr. John Doe	john.doe@happy-life-hospital.com	Happy Life Hospital	Hospital	Happy Life Hospital	0123456789		Approved	
Dr. Jane Smith	jane.smith@happy-life-hospital.com	Happy Life Hospital	Hospital	Happy Life Hospital	0987654321		Approved	
Dr. Alex Lee	alex.lee@happy-life-hospital.com	Happy Life Hospital	Hospital	Happy Life Hospital		Jun 27, 2025	Approved	
Dr. Emily White	emily.white@happy-life-hospital.com	Happy Life Hospital	Hospital	Happy Life Hospital		Jun 27, 2025	Approved	
Dr. Michael Brown	michael.brown@happy-life-hospital.com	Happy Life Hospital	Hospital	Happy Life Hospital		May 1, 2025	Approved	
Dr. Sarah Green	sarah.green@happy-life-hospital.com	Happy Life Hospital	Hospital	Happy Life Hospital		Apr 30, 2025	Approved	
Dr. David Black	david.black@happy-life-hospital.com	Happy Life Hospital	Hospital	Happy Life Hospital			Approved	
Dr. Lisa Gray	lisa.gray@happy-life-hospital.com	Happy Life Hospital	Hospital	Happy Life Hospital		Mar 19, 2025	Approved	
Dr. James Blue	james.blue@happy-life-hospital.com	Happy Life Hospital	Hospital	Happy Life Hospital		Mar 19, 2025	Approved	

« < 1 2 3 4 5 > »

Once done, click on the **Update** button to save the changes.

Update User

User Information | Documents

Org Type* Hospital X

User Role* Admin X

First Name* TestLN

Phone No. Phone No.

Password

Line Id Line Id

District Select a District

Sub-District Select a Sub-District

Organization* Test-01 X

Auth Type* Database Credentials X

Last Name* TestLN

Email

Status Approve X

Country Select a country

Province Select a Province

Zip Code Select a Zip Code

Update

4.3 Exam Information

4.3.1 Exam List

Use the Exam Information page to view, search, and edit exam definitions for a selected site. Each exam includes a unique ID/UID, display name, and optional subspecialty.

Patient Portal | Miracle Advance Technology | Sp Admin Test

Settings / Exam Information

Exam Information

4 columns selected | Sites: SRUPHAT Imaging Center (L... | Global filter


Exam ID	Exam UID	Exam Name	Subspecialty	Action
981	MAJ-1	MRI Ankle joint (Left)		
980	MCen-1	MRI Chest and/or mediastinum		
979	MTL-1	MRI Thigh (Left)		
978	MLL-1	MRI Leg (Left)		
977	MLasp-1	MRI Lower abdomen (or pelvic cavity)		
976	MR-1	MRI Neck		
973	MLA-3	MRI LOWER ABDOMEN		
972	MLA-2	MRI Lower Abdomen		
970	MSpec-1	MRI Spine : Screening whole spine		
969	MSTj-1	MRI Spine : Thoracolumbar junction		

Version: 1.13.0 | Rel. Date: Aug 08, 2025




4.3.2 Adding New Exam

A new exam can be added. After clicking the "Add Exam" button, you should see a pop-up window where you can enter the name, UID, and Subspecialty for the new exam.

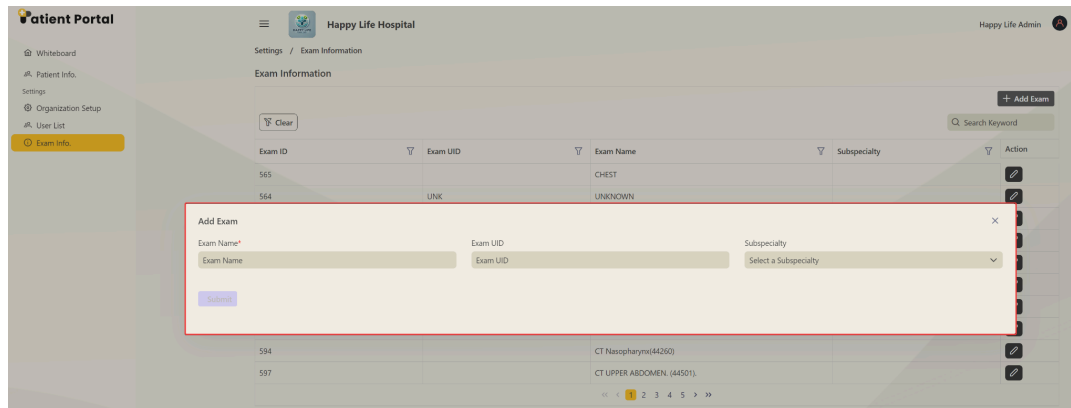

4.3.3 Edit an exam

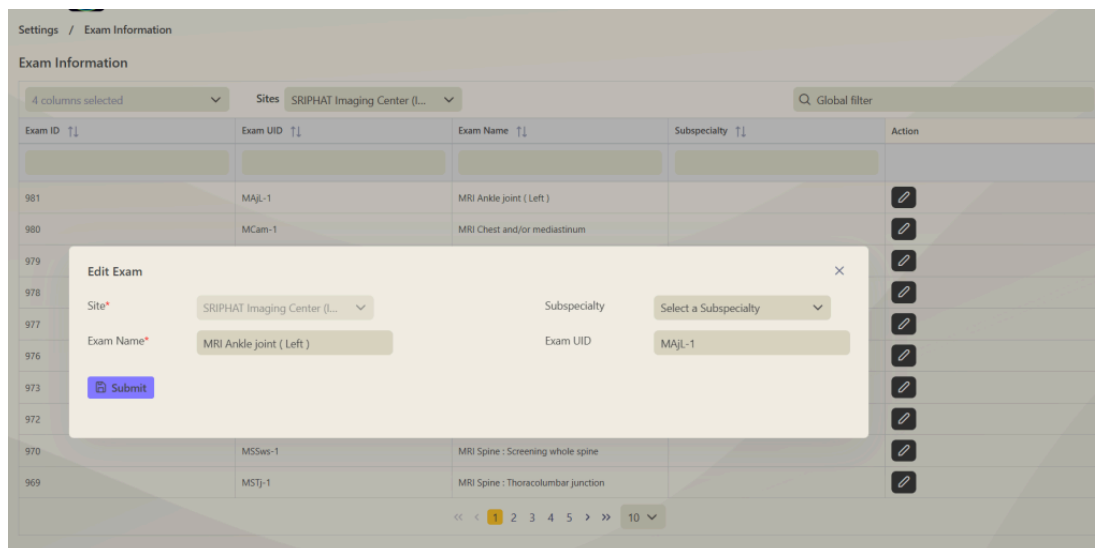


1. Find the exam (use Sites and Global filter).
2. Click the Edit (pencil) icon.
3. In the Edit Exam modal:
 - Site: Shows the site this exam belongs to (often read-only).

- Exam Name (required): Enter the display name (e.g., "MRI Ankle joint (Left)").
- Subspecialty: Select from the dropdown, if applicable.
- Exam UID: Unique identifier; may be read-only depending on your permissions. Avoid changing it—integrations may rely on this value.


4. Click Submit to save.

5. View Patient Information

Use the Patient Info page to find and manage patient demographics and identifiers.

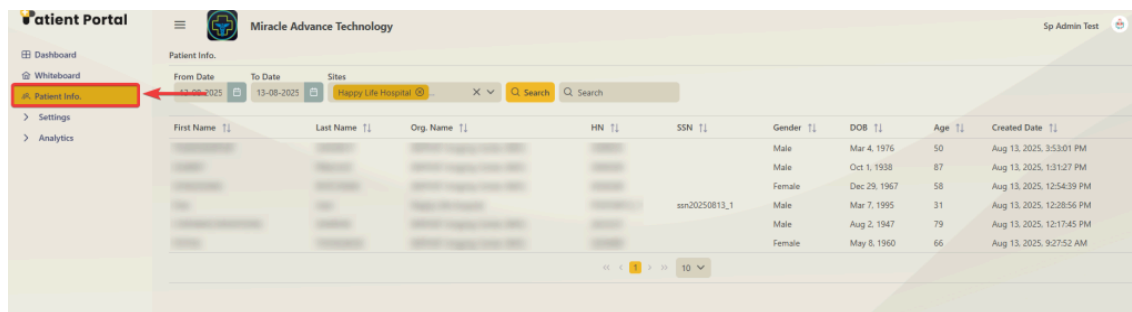

You can filter by date range, site, and organization, and quickly search by HN or SSN to locate a specific patient.

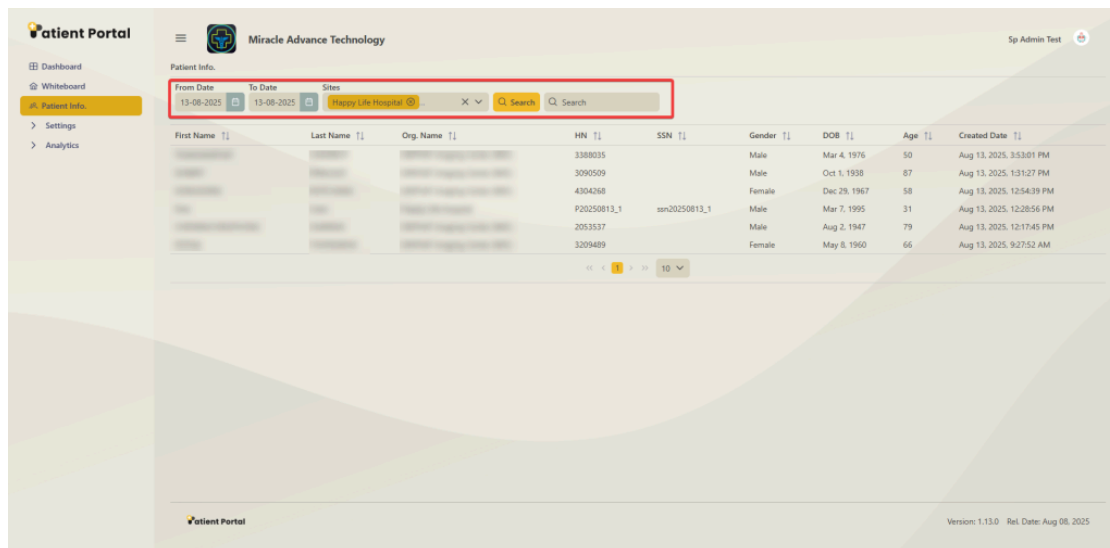



6. Technical Documentation

This section provides detailed API documentation, outlining how to create imaging orders and share links and QR codes through the Patient Portal API. It also covers the security measures implemented to protect these interactions.

Please refer to the document for API documentation on creating imaging orders and sharing links and QR codes - [☰ Patient Portal API Documentation](#)

6.1 API Security

The Patient Portal API is secured using Basic Authentication (Basic Auth), which requires API users to provide a valid username and password combination to access API functions. This security measure ensures that only authorized users can execute API calls, providing a simple yet effective layer of security.

- **Basic Authentication:** Each API request includes a Basic Auth header. This header consists of a **Base64**-encoded string representing the username and password. This method helps prevent unauthorized access to the API and ensures that interactions are secure.