



Technical Manual

RadPortal

Developed By -

Miracle Advance Technology

Last Update Date -

12/22/2022

This Document Is for the RadPortal Application

Version – 1.1.1

Written By -

S. M. Al Shahaneaz Shaikat

Revision Sheet

Release No.	Date	Revision Description	Author Name
Rev. 1.1.1	22/12/22	Add Medical licence while signing up, Radiologist Dashboard	S M Al Shahaneaz Shaikat
Rev. 2	DD/MM/YY	Short description of Major changes	Name
Rev. 3	DD/MM/YY	Short description of Major changes	Name
Rev. 4	DD/MM/YY	Short description of Major changes	Name

USER'S MANUAL OF RADPORTAL

Table of Contents

1. General Information about RadPortal.....	5
1.1 System Overview	5

1.2	Project References.....	6
1.3	Organization of the Manual	6
1.4	Acronyms and Abbreviations	6
1.5	Additional Help	7
2.	System Summary.....	7
2.1	Minimum Hardware Requirements, Supported Browsers List, And Operating System List for RadPortal	7
2.1.1	Hardware Requirements Personal Computers and Laptops:	7
2.1.2	Browsers:	8
2.1.3	Operating System:	8
2.2	User Roles and Permission Levels.....	8
2.2.1	Service Provider(Admin)	8
2.2.2	Customer	9
2.2.3	Hospitals(Healthcare Service Provider)	9
2.2.4	Radiologist(Subject Matter Experts)	9
2.2.5	Viewer	9
3.	Getting Started.....	9
3.1	Signing Up	10
3.2	Logging In	14
3.3	System Menus	16
3.4	How to Change Password, Change Profile Picture, and Add Your Digital Signature.....	17
3.4.1	Changing Password	17
3.4.2	Changing Profile Photo	18
3.4.3	Adding Digital Signature	19
3.5	Sign Out from The System.....	21
4.	Exploring Dashboard.....	23
5.	Exploring Worklist Menu	25
5.1	How to Search Studies.....	26
5.2	How to Sort Studies	27
5.3	How to view Images, Decline Cases, write Reports, and View Documents Additional Documents of Studies	28
5.3.1	Review Images of a Specific Study.....	28
5.3.2	How to Decline Assigned Studies	30
5.3.3	How to View Additional Documents of a Study.....	31
5.3.4	How to Write Reports and Finalize Reports	32
5.3.5	View Study Reports	37
5.3.6	Adding an Addendum	38



1. General Information about RadPortal

Radportal is a Secured cloud-based radiology platform. The main problem it solves is that there are many hospitals in remote places worldwide with all the modern medical equipment for medical tests such as (MRI, CT scan, x-rays, nuclear medicine, ultrasound, and more). Still, hospitals may not have expert radiologists/doctors readily available to write the reports for these studies when patients come to these remote hospitals and complete their tests.

Sometimes the Patient has to physically go to hospitals/Other more developed cities where radiologists/Doctors are available, but this process is difficult and time-consuming.

Using the RadPortal System, the hospital can easily save patient test results. Then the hospitals can forward the test results to an expert doctor/radiologist who may not live nearby.

1.1 System Overview

RadPortal is a secured cloud-based radiology platform. It helps hospitals to send their Patient's radiology test images to doctors/radiologists worldwide and receive written reports from remote doctors/radiologists.

RadPortals Core functionalities are

- All hospitals/sites that have already configured modality devices can automatically send DICOM files to their routers, which are then forwarded to the RadPortal cloud database.
- RadPortal system reads the DICOM headers and collects Information from the header, and saves it as an entity called "study".
- Hospitals also can upload single or multiple DICOM files manually.
- Hospitals can then access DICOMs from their RadPortal profiles and assign specific radiologists to write reports on that particular study or keep the study open for other radiologists.
- Suppose any open study is available on RadPortal. In that case, the independent radiologists can log in to their RadPortal profiles and independently select any open study. Then the individual Radiologist can write reports for that particular study.
- The Radiologist has a roaming profile, so they do not need specific devices to interact with the RadPortal system. The Radiologist can interact with the RadPortal system with laptops, desktops, or tabs and a reliable Internet connection.

- RadPortal also helps Government Institutes such as the National Health Security office to track the number of studies and reports the hospitals and the radiologists completed and disburse fees and payments accordingly.

1.2 Project References

References used in preparing this document are listed here. Project References may help you understand the complex topics and terms used in the RadPortal system.

<http://dicomlookup.com/modalities.asp>

1.3 Organization of the Manual

This document is patterned to familiarize you with the core functionalities of the RadPortal application. This document explains how the application is intended to use by properly explaining its functionalities and limitations.

You will also know how the RadPortal application saves Information and how the Information is shared across multiple user roles within the RadPortal system. We will discuss the most common elements so readers can later relate to the terms as they look deeper into the application's features and functionalities.

The functionalities will vary depending on user roles, but the role **Admin** will have access to all the functionalities of RadPortal. User roles are extensively explained in section [2.2, User Roles and Permission Levels](#).

1.4 Acronyms and Abbreviations

Common acronyms and abbreviations are introduced in this section so the reader can get comfortable with the acronyms and abbreviations used in this documentation and the RadPortal application.

DICOM	Digital Imaging and Communications in Medicine.
PACS	Picture archiving and communication system.

1.5 Additional Help

Common terms are introduced throughout the user manual to familiarize the reader with the common terms used in the application. Sections are indexed (numbered, bullet points, and roman numbers). Sections are also linked to common keywords so the reader can go directly to any relevant section to find more information. For easier navigation, keywords are marked in [blue](#).

Here are the styles used for lists

Bulleted lists:

- **Text 1**
- **Text 2**
- **Text 3**

Numerical numbered lists:

- 1. Text 1**
- 2. Text 2**
- 3. Text 3**

Roman numbered lists:

- I. Text 1**
- II. Text 2**
- III. Text 3**

Here are the styles for Notes and warnings.

| **Note:** Notes are marked in the color blue.

| **Warning:** Warning messages are marked in the color red.

2. System Summary

2.1 Minimum Hardware Requirements, Supported Browsers List, And Operating System List for RadPortal

2.1.1 Hardware Requirements Personal Computers and Laptops:

- Processor (CPU) with 2 gigahertz (GHz) frequency or above

- A minimum of 2 GB of RAM
- Monitor Resolution 1024 X 768 or higher
- A minimum of 20 GB of available space on the hard disk
- Stable internet connection with a speed of 4 megabits per second or higher
- Keyboard and mouse or compatible pointing device

2.1.2 Browsers:

- Chrome 36+
- Edge 20+
- Mozilla Firefox 31+
- Internet Explorer 11+ (Windows only)
- Safari 6+ (macOS only)

2.1.3 Operating System:

- Windows 7, Windows 8, Windows 10, or Windows 11
- Mac OSX 10.8, 10.9, 10.10, or 10.11

2.2 User Roles and Permission Levels

Before we look deeper into the RadPortal system, you must understand the user roles properly. The user roles are listed below.

- I. Service Provider (Admin)
- II. Customer
- III. Hospitals (Health Service Providers)
- IV. Radiologists (Subject Matter Experts)
- V. Viewer

2.2.1 Service Provider(Admin)

Miracle Advance Technology is the service provider, and Miracle Advance Technology plays the Administration role inside the RadPortal System. The main functionality of an administrator role is to manage the system so that the system runs smoothly and without errors. For the RadPortal system, the Service Provider will also install, configure and manage DICOM routers for hospitals.

2.2.2 Customer

This role is very important in the RadPortal system. The only customer role can Configure Turnaround Time and manage and verify Radiologists. They also have full control over the RadPortal system. The Customer is also in charge of disbursing Fees and payments to the Hospitals and Radiologists.

2.2.3 Hospitals(Healthcare Service Provider)

Hospitals are also called "**sites**" in the RadPortal system. Hospitals play an significant role in the RadPortal System. They have to perform medical tests such as (MRI, CT scan, x-rays, nuclear medicine, ultrasound, and more) and manually upload files and images from browsers or the modality devices can automatically send **DICOM** files to RadPortal servers to create **cases** which are also called **studies**. The hospital can access the files, images, and DICOMs uploaded to the RadPortal server.

2.2.4 Radiologist(Subject Matter Experts)

The Radiologist also plays a significant role in the RadPortal System. There are two types of radiologist roles in the RadPortal system.

- I. Normal Radiologists
- II. Supervisor Radiologist

The main difference between the Supervisor and a normal Radiologist is that a Supervisor Radiologist can view all cases or studies, and a normal radiologist can only see the studies/cases they are assigned to. Only the verified normal Radiologist can write reports of submitted cases.


2.2.5 Viewer

The Viewer role is a bit unique. The viewer does not interact with any aspects of the RadPortal system. The viewers only view how many cases/studies the hospital/sites completed and how many reports the radiologists submitted and finalized in the RadPortal system.

This role is given to people/organizations with the authority to monitor and audit the RadPortal system. This role is also given to the government institutes that might want to monitor the whole RadPortal system.

3. Getting Started


Signup and Login Workflow

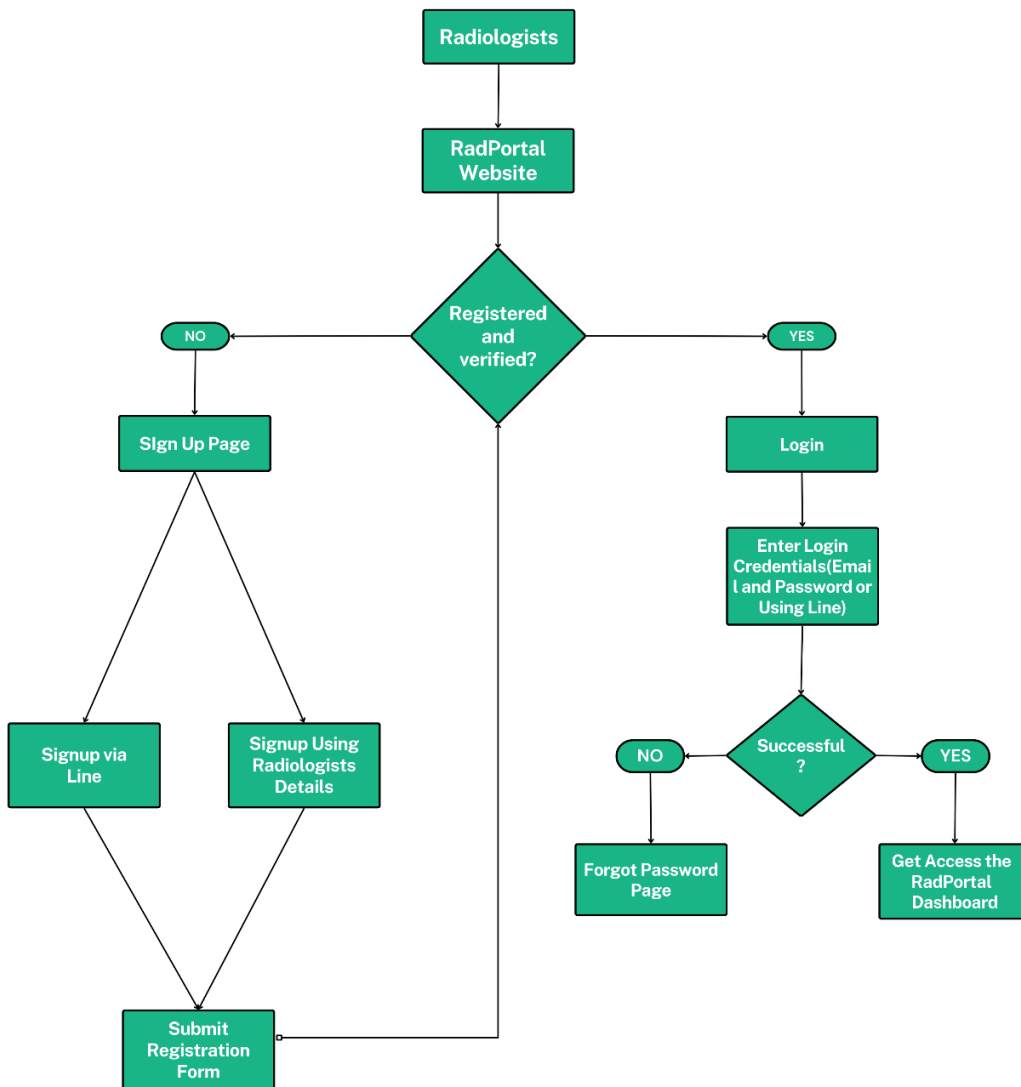


3.1 Signing Up

Before using the RadPortal system, you have to sign up as a Radiologist. To sign up as a radiologist, you have to go to this URL - <https://radportal.net>.


After the website loads, you should see a page like this.

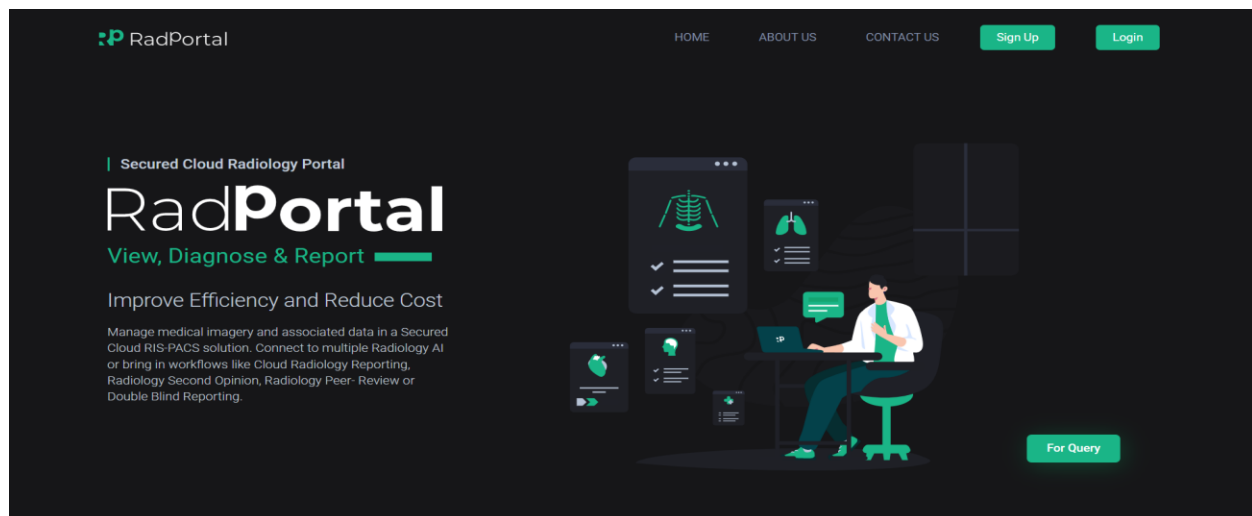

To sign up as a radiologist, click on the **Sign-Up** button.

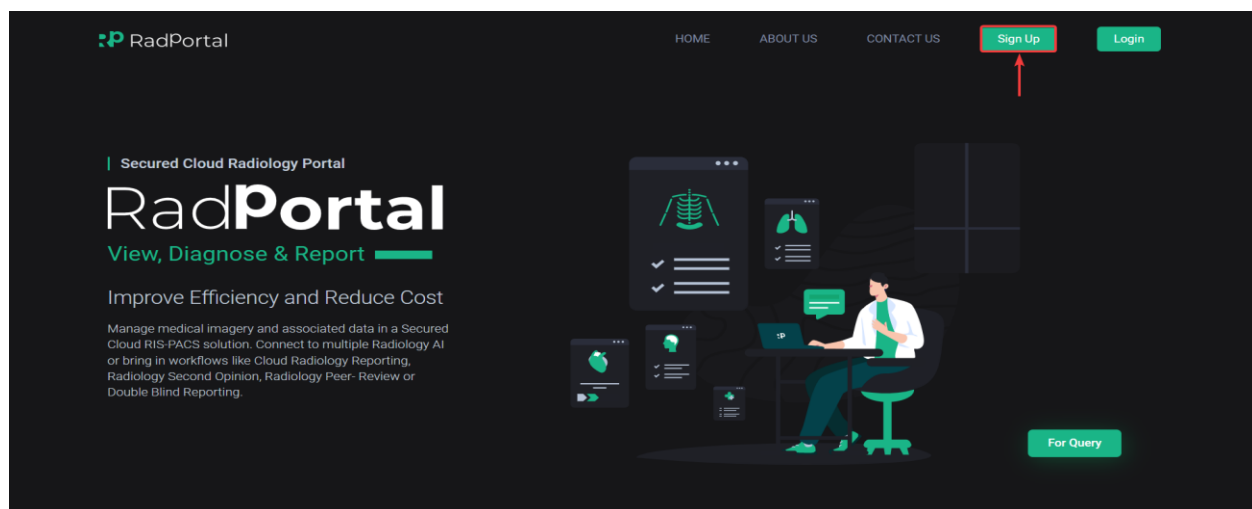

After clicking, you should see a page where you have to click on the **Sign up as a Radiologist** button.

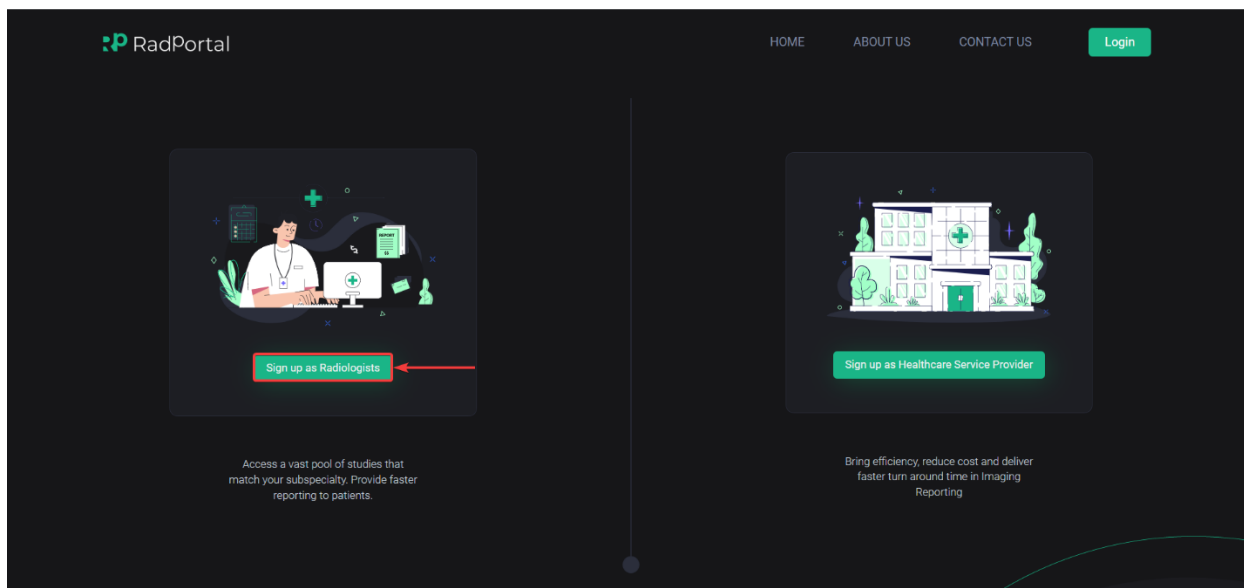

Now you will be forwarded to the Radiologist sign-up page, where you have to sign up using your **Email [1]** and **Phone number [2]** and agree to our Account Terms of Use by clicking on the **Checkbox [3]**. Finally, click on the **Send Verification code [4]** button.

The image shows the "Sign up as Radiologist" form. The title "Sign up as Radiologist" is at the top center. Below it is a form box titled "Radiologist Sign up". Inside the form box, there are four numbered callouts: 1 points to the "Email" input field, 2 points to the "Phone No" input field, 3 points to the checkbox for "I have read and agree to the RadPortal Account Terms of Use", and 4 points to the "Send Verification" button. The "Send Verification" button is green and located at the bottom right of the form box.



Sign up as Radiologist

Radiologist Sign up

Email

Phone No

PhoneNo

Verify

☒ I have read and agree to the [RadPortal Account Terms of Use](#)


After you complete your verification process, you can move forward to the next section, where you have to fill in your **Password [1]**, **Confirm Password [2]**, **First Name [3]**, **Last Name [4]**, **Line ID [5]**, **Country [6]**, **District [7]**, **Province [8]**, **Sub-District [9]**, **Zip Code [10]**, **Subspecialties [11]**, upload your profile picture by clicking on **Choose [12]**, Upload any relevant documents by clicking on **Upload Document [13]**. To start the signing-up process from the start, click on the **Start Over [14]** button. Finally, click on **Finish [15]** to complete your sign-up process.

The screenshot shows the 'Sign up as Radiologist' form on the RadPortal website. The form is dark-themed with white text. It includes fields for Email, Phone No., Password, Confirm Password, First Name, Last Name, Line Id, Country, District, Province, Sub-District, Zip Code, and Subspecialties. There are also buttons for '+ Choose', 'Upload Document', 'Start Over', and 'Finish'. Red numbered arrows point to specific elements: 1 points to the Password field, 2 to the Confirm Password field, 3 to the First Name field, 4 to the Last Name field, 5 to the Line Id field, 6 to the Country dropdown, 7 to the District dropdown, 8 to the Province dropdown, 9 to the Sub-District dropdown, 10 to the Zip Code field, 11 to the Subspecialties dropdown, 12 to the '+ Choose' button, 13 to the 'Upload Document' button, 14 to the 'Start Over' button, and 15 to the 'Finish' button. At the bottom left, there is a checkbox for 'I have read and agreed to the RadPortal Account Terms of Use'.

After you successfully signed up to the RadPortal system. You have to wait a certain period to get verified by our system. All of your documents will be manually verified. After your details are verified, you should receive a welcome email from our end.

3.2 Logging In

You will be able to log in. To log in, you have to go to this URL – <https://app.radportal.net/#/login>. After the login page loads, you should see a page like this

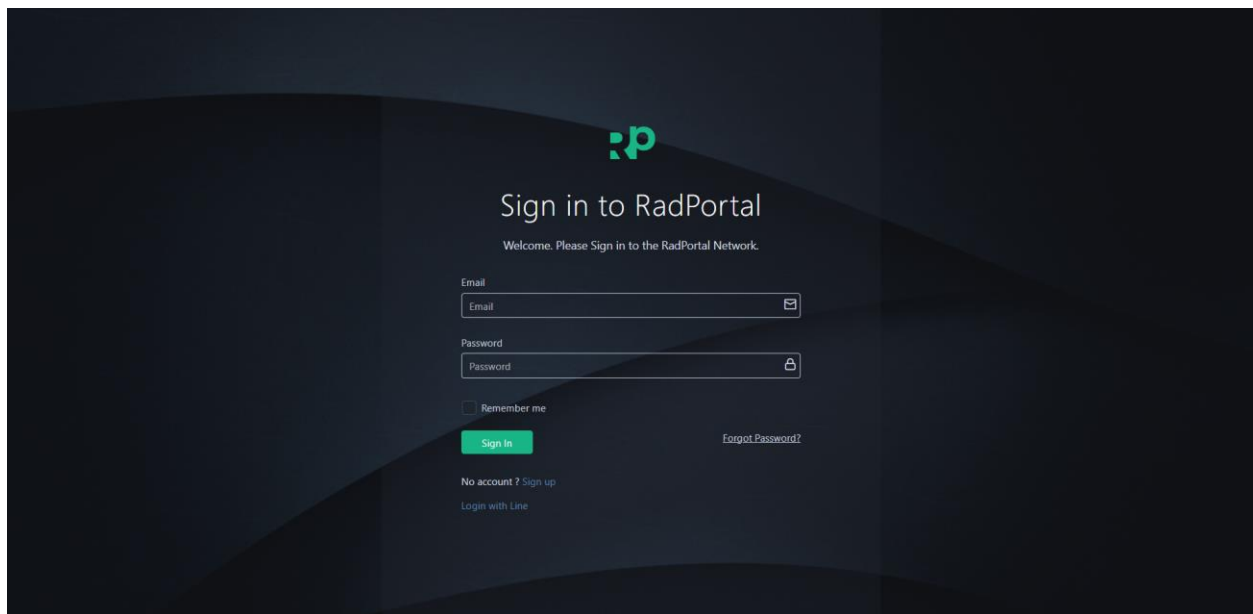

After you access the login page, you can use your **Email [1]** and **Password [2]** to log in to the RadPortal. You have to use the login credentials you used to create your account. Finally, click on the **Sign In [3]** button.

The image shows the RadPortal sign-in interface. At the top is the RadPortal logo (a green 'P' with a colon). Below it is the heading "Sign in to RadPortal" and a welcome message: "Welcome. Please Sign in to the RadPortal Network." The form contains two input fields: "Email" and "Password", each with a red box around it and a red arrow pointing to it labeled with a number (1 and 2 respectively). Below the password field is a "Remember me" checkbox. A red box around the "Sing In" button is labeled with a number 3. To the right of the button is a link "Forgot Password?". At the bottom, there are links for "No account? Sign up" and "Login with Line".

Sign in to RadPortal

Welcome. Please Sign in to the RadPortal Network.

Email

Password

☐ Remember me

Sing In [Forgot Password?](#)

No account? [Sign up](#)

[Login with Line](#)

After successfully logging in, you should see the dashboard page of RadPortal.

The image shows the RadPortal dashboard. At the top is the RadPortal logo and the text "Miracle Advance Technology". On the right is a user profile icon and the name "Sp Admin Test". The main section is titled "Study Worklist". It contains filters for "From Date" (21-11-2022), "To Date" (21-11-2022), "Customer" (Customer 2, Customer 1), and "Site" (HOSPITAL 3, Hospital 1, Hospital 2). There are buttons for "Search", "Today", "Yesterday", and "Last 3 Days". An "Upload DICOM" button is also present. Below the filters are radio buttons for "Completed", "Submitted", "Finalized", "Submitted / Finalized", and "Show All". A table with 11 columns is shown, including "Accession# ↑↓", "HN ↑↓", "Patient Name ↑↓", "Status ↑↓", "Date of Birth ↑↓", "Gender ↑↓", "Exam Name ↑↓", "Image ↑↓", "Created On ↑↓", "Updated On ↑↓", "Assigned To ↑↓", and "Action". The table is empty. At the bottom is the RadPortal logo and the text "RadPortal".

Miracle Advance Technology

Sp Admin Test

Study Worklist

From Date: 21-11-2022 To Date: 21-11-2022 Customer: Customer 2, Customer 1 Site: HOSPITAL 3, Hospital 1, Hospital 2

Search Today Yesterday Last 3 Days

Upload DICOM

☐ Completed ☐ Submitted ☐ Finalized ☐ Submitted / Finalized ☒ Show All


11 columns selected

Global filter

Accession# ↑↓	HN ↑↓	Patient Name ↑↓	Status ↑↓	Date of Birth ↑↓	Gender ↑↓	Exam Name ↑↓	Image ↑↓	Created On ↑↓	Updated On ↑↓	Assigned To ↑↓	Action
---------------	-------	-----------------	-----------	------------------	-----------	--------------	----------	---------------	---------------	----------------	--------


RadPortal

3.3 System Menus

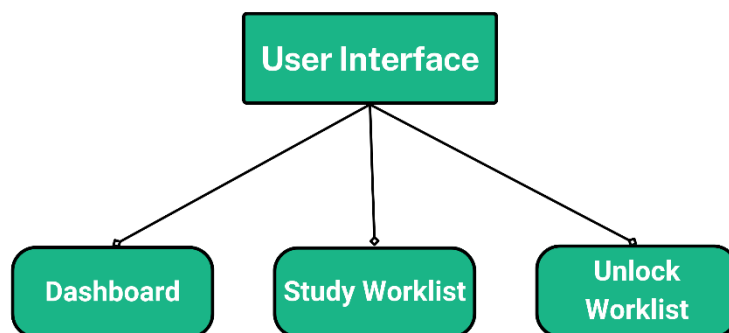

The RadPortal system has these menus currently available for Radiologists


- 1 Dashboard
- 2 Worklist
- 3 Unlock Worklist

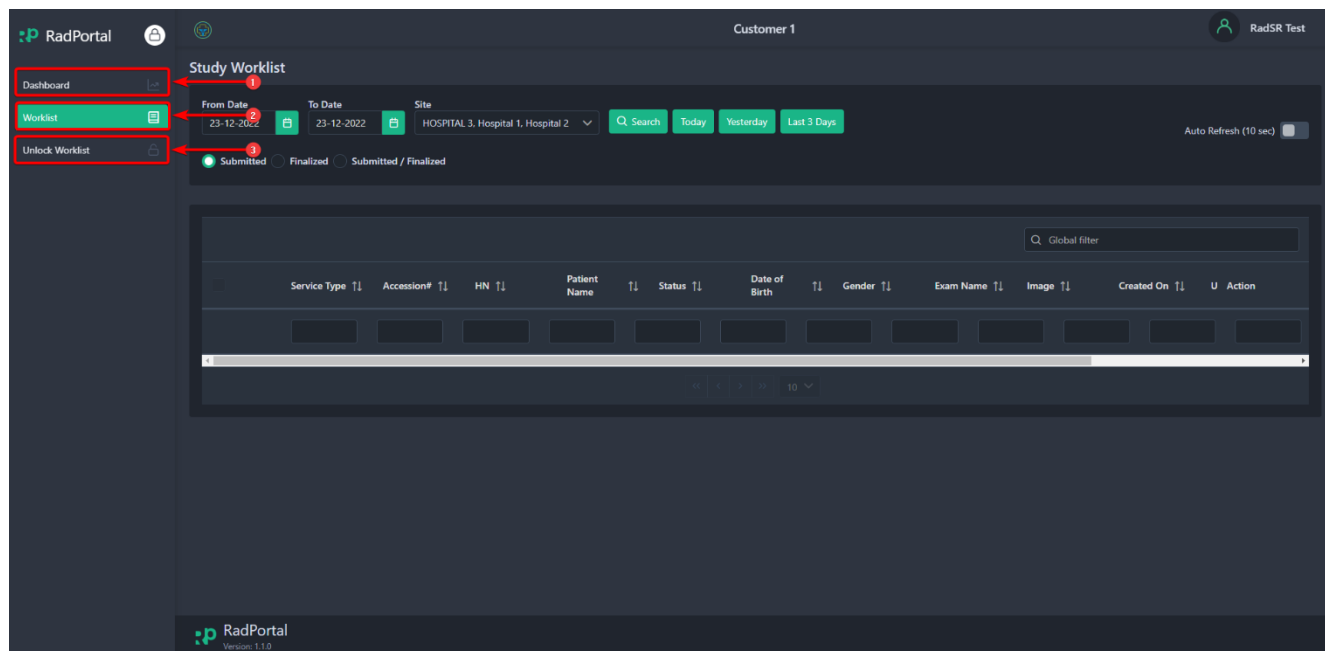

3.4 How to Change Password, Change Profile Picture, and Add Your Digital Signature

3.4.1 Changing Password

You can easily change your RadPortal login password from your profile. To access your profile, click on your **username [1]** in the top right corner and click on **Profile [2]**.

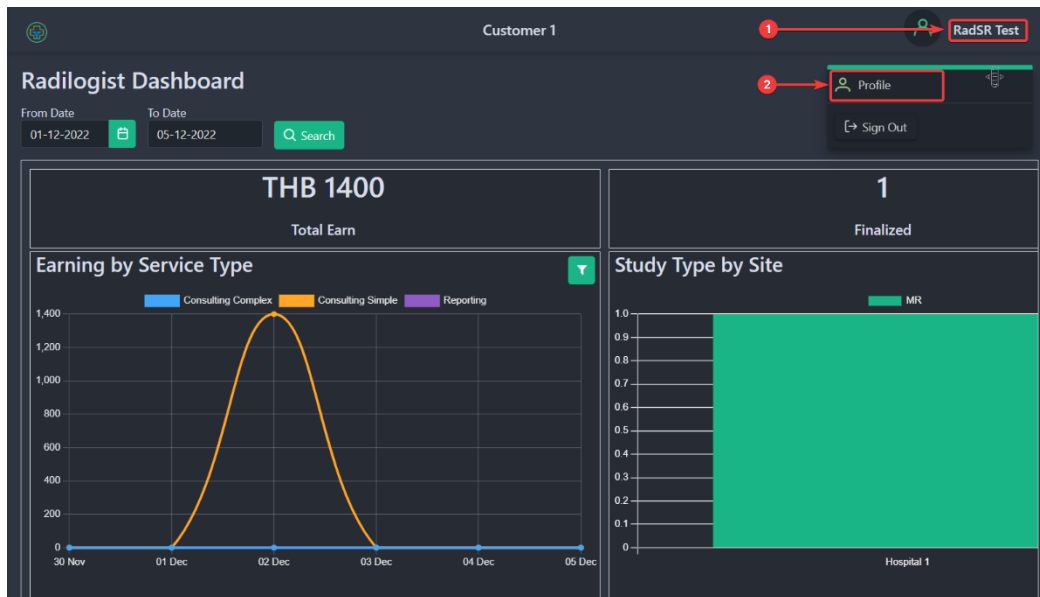

After clicking on The Profile, the system will redirect you to your profile page. On this page, you can change your Password. To change your user Password, click on the **Change User Password** button.

The screenshot shows the 'User Profile' page. On the left, there is a profile picture placeholder with a '+ Choose' button. On the right, there is a table of user information:

Name	: RadSR Test
Email	: RadSR@gmail.com
Phone	: 1234
Address	:

Below the table, there are two buttons: 'Change User Password' (highlighted with a red box and an arrow) and 'Digital Signature'.

After clicking, you should see a new page where you have to confirm your **Old Password [1]**, enter your **New Password [2]**, and **confirm your new Password [3]**. Click on the **Update button [4]** to finalize your new password setup.

The image shows a 'Update User Password' form. It has three input fields: 'Old Password', 'New Password', and 'Confirm New Password', each containing six asterisks. Red boxes and arrows with numbers 1, 2, and 3 point to these fields respectively. Below the fields are two buttons: 'Update' (with a checkmark icon) and 'Cancel' (with an 'X' icon). A red box and arrow with the number 4 point to the 'Update' button.

3.4.2 Changing Profile Photo


To change your profile photo, go back to your profile again and click on the Choose button.

The image shows a 'User Profile' page. On the left, there is a profile picture placeholder with a '+ Choose' button. A red arrow points to this button. To the right, there is a list of user details: Name (RadSR Test), Email (RadSR@gmail.com), Phone (1234), and Address. Below these details are two buttons: 'Change User Password' and 'Digital Signature'.

After clicking on the button, you can choose a photo for your profile. You can select your photos directly from your computer. After you upload your profile photo, you will see the preview in the section below.


The image shows the 'User Profile' page after a photo upload. The profile picture placeholder now shows a preview of a file named 'dummy-pro file-pic.png' with a size of 11.928 KB. A red box highlights this preview area. A green cross icon is visible next to the preview. The '+ Choose' button is still present. The user details and other buttons remain the same.

You can remove the photo by clicking on the green cross icon.

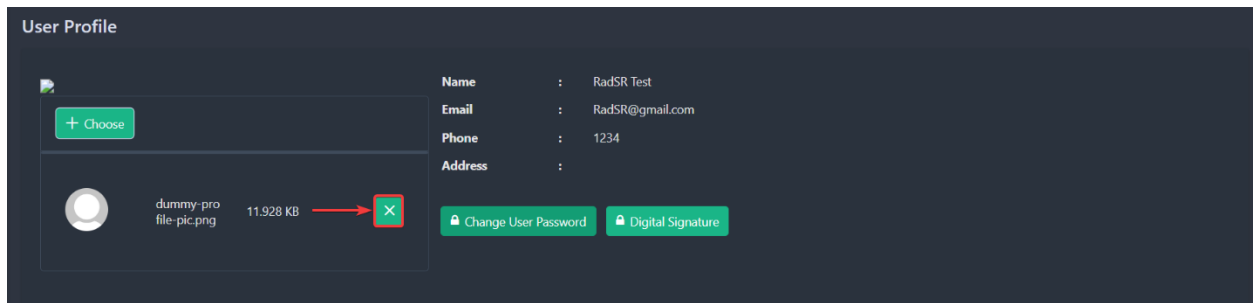


3.4.3 Adding Digital Signature


To add your Digital Signature, go back to your profile and click the **Digital Signature** button.

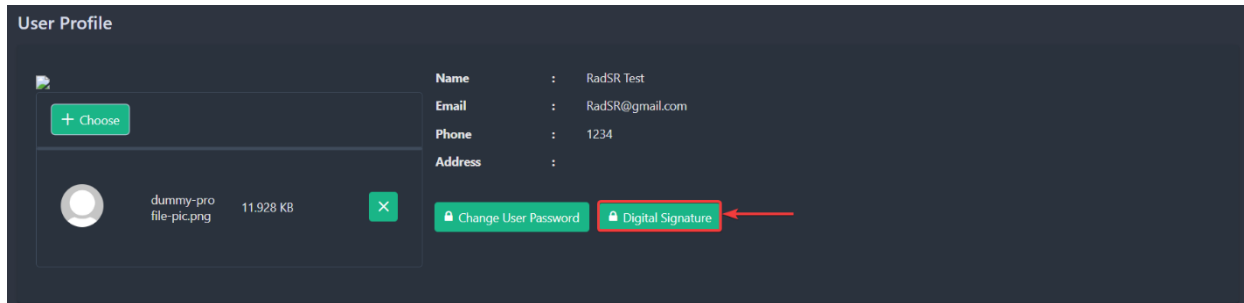

After clicking, you should see a new page where you will be able to upload your digital signature

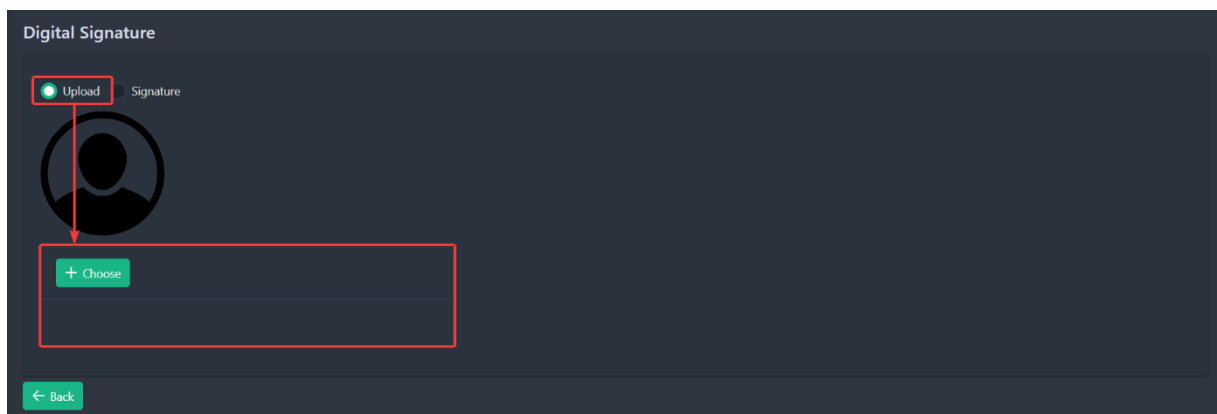

Or you can directly sign in the white blank section.

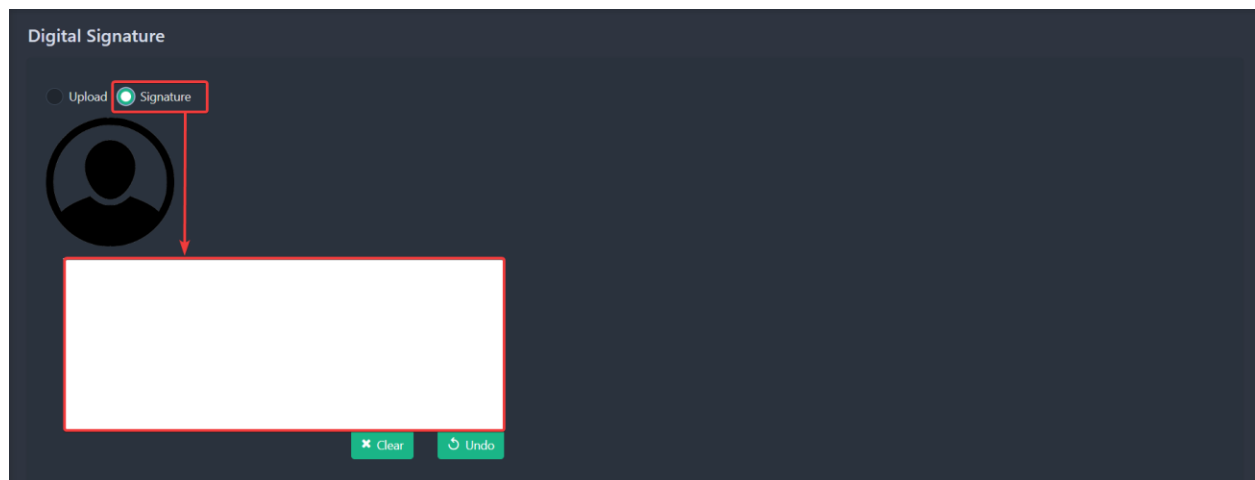


You can Clear your signature or undo a specific part of the signature by clicking on the **Clear** or **Undo** Button.




3.4.4 Exploring Finalized Workload

To add your Digital Signature, go back to your profile and in the bottom you should see a section called **My Finalized Workload**.

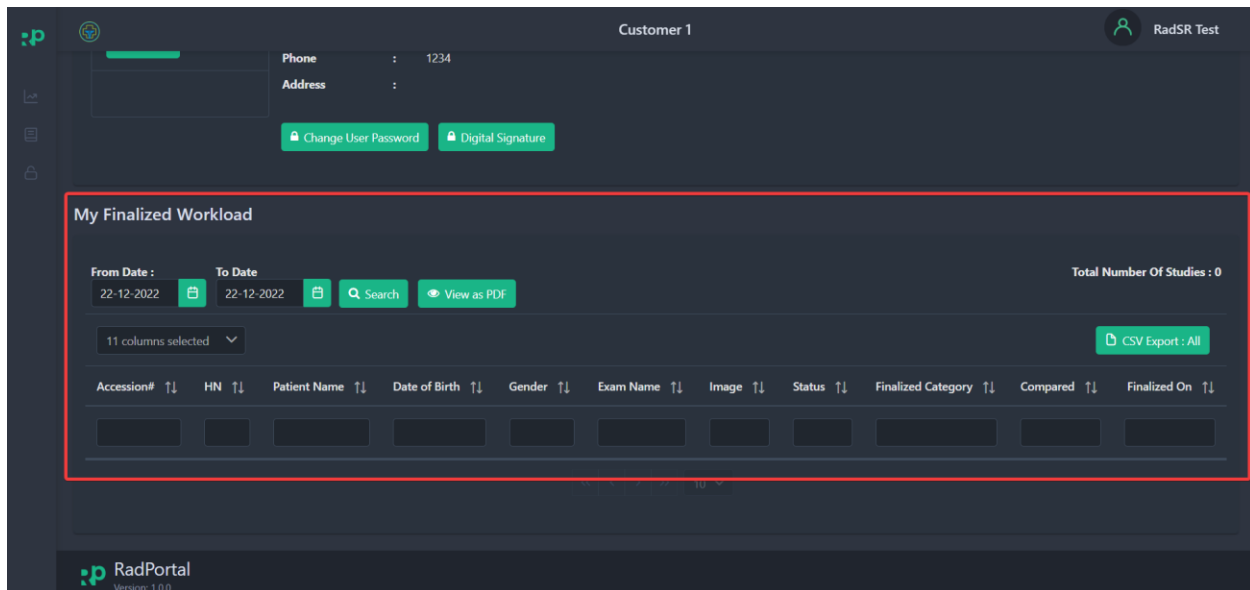

Now you can review a custom workload report by filtering the **From Date [1]** and **To Date [2]**. After selecting custom date range click on the **Search [3]** button.

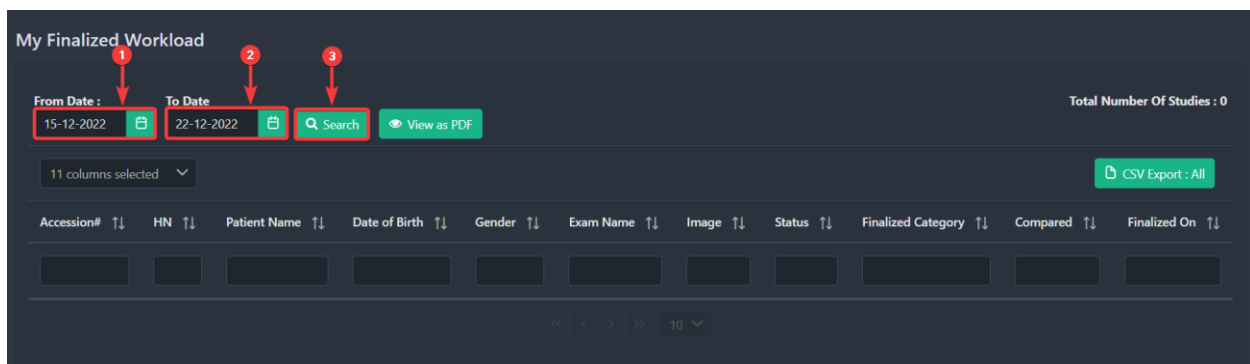

After clicking you should see your custom finalized workload list.

3.5 Sign Out from The System

To Sign Out from the RadPortal system, click on your **username [1]** in the top right corner and click on **Sign Out [2]**.

Sp Admin Test

Profile

Sign Out

Auto Refresh (10 sec)

Miracle Advance Technology

Study Worklist

From Date

To Date

Customer

Site

Search

Today

Yesterday

Last 3 Days

Upload DICOM

Completed

Submitted

Finalized

Submitted / Finalized

Show All

11 columns selected

Global filter

Accession#	HN	Patient Name	Status	Date of Birth	Gender	Exam Name	Image	Created On	Updated On	Assigned To	Action

RadPortal






In the Radportal Dashboard, you will see a quick summary

1. Total Earnings
2. Total reports you finalized
3. Pie chart of service type
4. Total earnings by service type
5. Pie chart of study type
6. Chart of service type by hospitals
7. Chart of Service type by a study
8. Pie chart of Study distribution by modality

5. Exploring Worklist Menu

Study Workflow

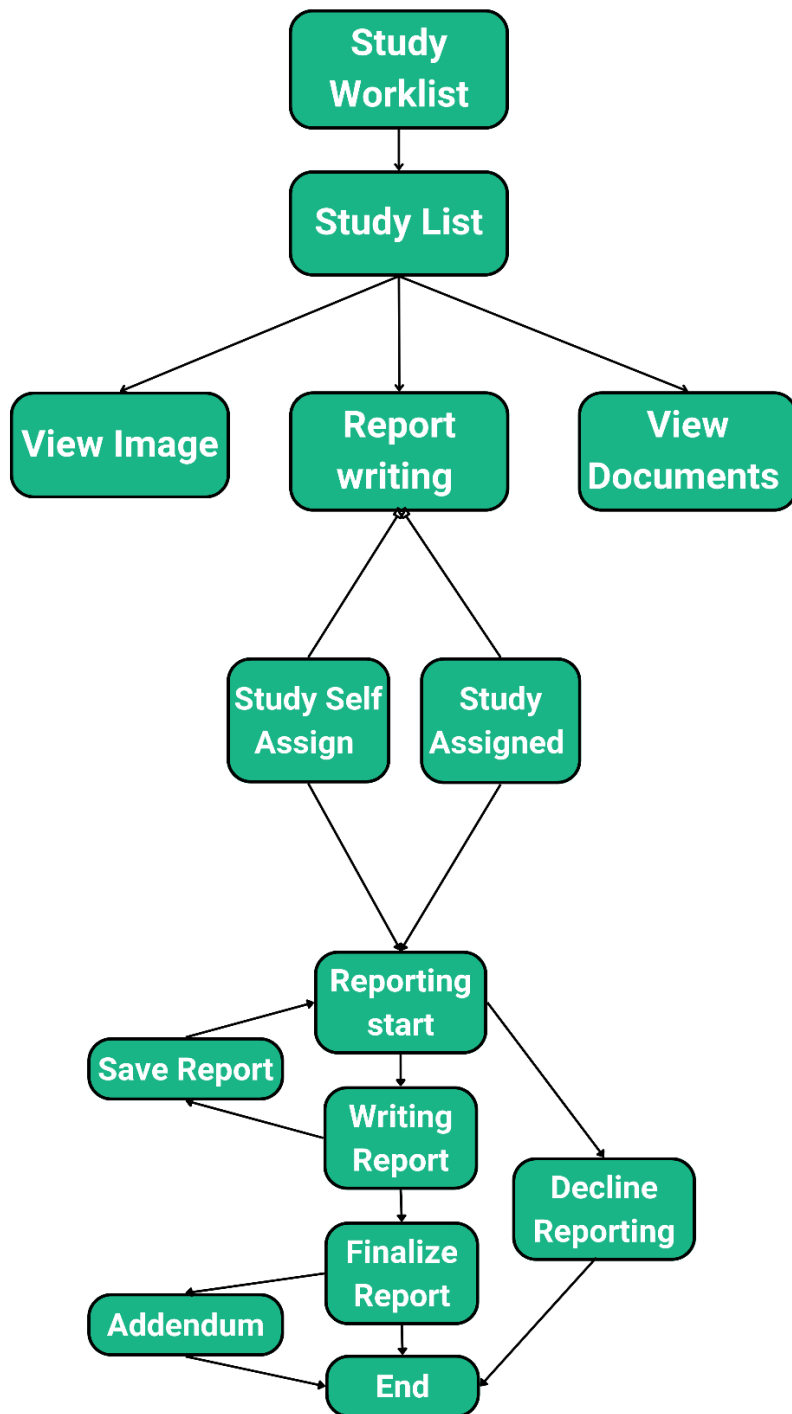

This menu in the RadPortal system allows you to keep track of studies that are assigned to you and all the studies that are kept open for any radiologist on the system. You can also view all modality test images in the Study worklist option. You can easily write reports for the studies that are assigned to you.

The worklist Menu page should look like the screenshot below.

The screenshot shows the RadPortal Study Worklist interface. At the top, there's a header with the RadPortal logo, a user profile icon labeled 'RadSR Test', and the text 'Customer 1'. Below the header, the title 'Study Worklist' is displayed. The main area contains search filters: 'From Date' (22-12-2022), 'To Date' (22-12-2022), and 'Site' (HOSPITAL 3, Hospital 1, Hospital 2). There are buttons for 'Search', 'Today', 'Yesterday', and 'Last 3 Days'. Below these are radio buttons for 'Submitted', 'Finalized', and 'Submitted / Finalized'. An 'Auto Refresh (10 sec)' toggle is on the right. A table with columns: Service Type, Accession#, HN, Patient Name, Status, Date of Birth, Gender, Exam Name, Image, and Action is shown. A 'Global filter' search bar is above the table. The footer shows 'RadPortal Version: 1.0.0'.


5.1 How to Search Studies

You can custom search based on criteria such as **From Date [1]**, **To Date [2]**, or you can also search by selecting a specific **Site or Hospital [3]**. Finally, click on the **Search [4]** button.

This screenshot is identical to the one above but includes red arrows and boxes highlighting the search criteria. Arrow 1 points to the 'From Date' field (05-12-2022). Arrow 2 points to the 'To Date' field (05-12-2022). Arrow 3 points to the 'Site' dropdown menu (HOSPITAL 3, Hospital 1, Hospital 2). Arrow 4 points to the 'Search' button. The table columns are: Service Type, Accession#, HN, Patient Name, Status, Date of Birth, Gender, Exam Name, Image, Created On, Updated On, Assigned To, Modality, and Action. The footer shows 'RadPortal Version: 1.0.0'.


After clicking, you should see all studies/cases that match your custom criteria.

You can even see all studies/cases created **Today [1], Yesterday [2], or in the last 3 days [3]**.

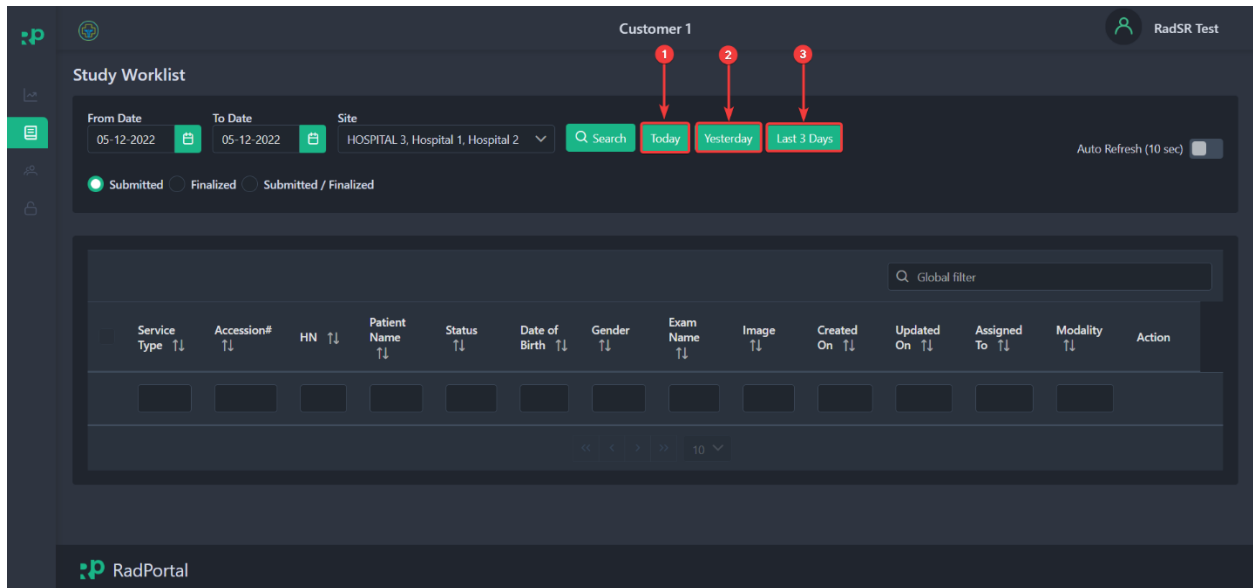


5.2 How to Sort Studies

You can easily sort study-list based on status, such as **Submitted [1], Finalized [2], or both Submitted and Finalized [3]**. To sort studies, click on the radio button that you see before each status section.

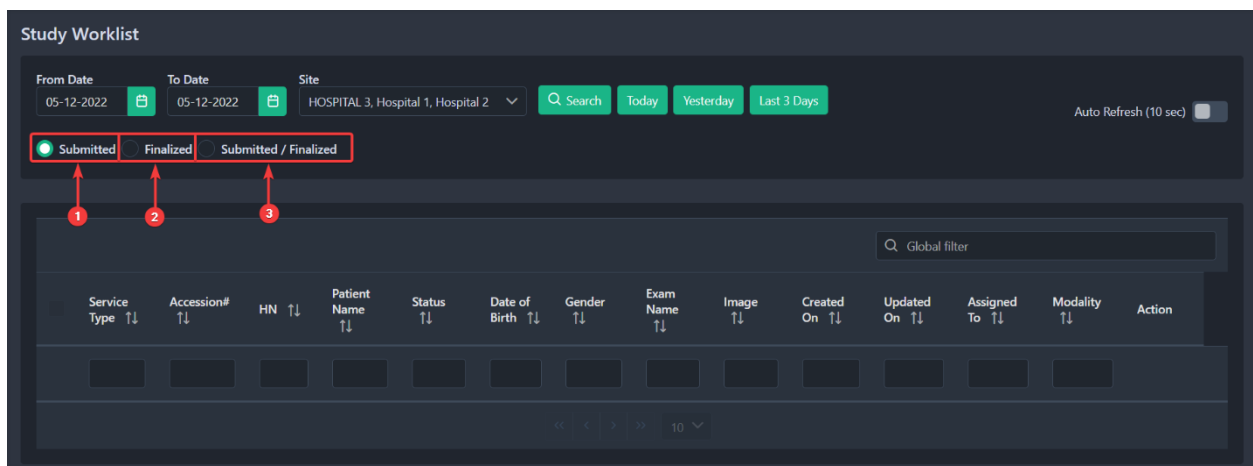

You can enable auto-refresh by clicking on the **auto-refresh slide button**. After enabling auto-refresh, the page will refresh continuously after every 10 seconds.

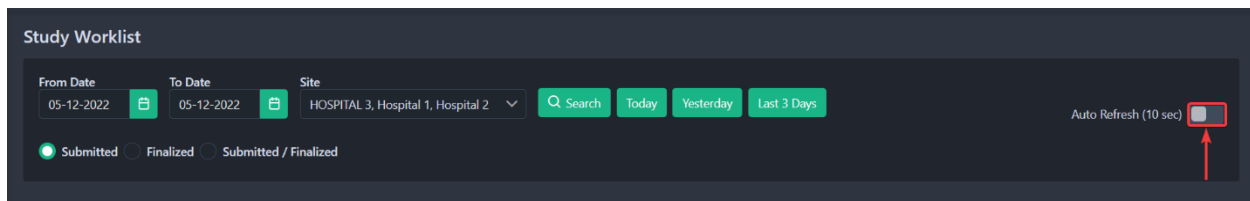

5.3 How to view Images, Decline Cases, write Reports, and View Additional Documents of Studies


On the Study Worklist page, you can only access the studies that are left open for all radiologists or specifically assigned to you. You can review Images of these studies, add reports, Finalize your report, and add Addendums after finalizing studies or you can decline assigned cases before you finalize your report. You can view additional documents that are attached to the studies.

5.3.1 Review Images of a Specific Study

You can easily access all images of a specific study by clicking on the **Image Icon** in the screenshot below.

Global filter												
☐	Service Type ↑↓	Accession# ↑↓	HN ↑↓	Patient Name ↑↓	Status ↑↓	Date of Birth ↑↓	Gender ↑↓	Exam Name ↑↓	Image ↑↓	Created On ↑↓	Updated On ↑↓	Action
☐												
☐	Consulting Complex Case	993	2000051	Michelle Mosconia	Submitted	27-09-1959 Age: 63	Female	MRI BREAST BILATERAL	665	06-12-2022 01:52 PM	06-12-2022 04:31 PM	 
☐	 Reporting	01509091655054	20150909162215	Kirstie Allie	Submitted	31-03-1971 Age: 51	Female	BREAST^Routine	12	02-12-2022 06:29 PM	06-12-2022 04:31 PM	R   
☐	 Reporting	1015	2000024	Melissa Dayton	Submitted	12-12-1955 Age: 67	Female	MRI CARDIAC FUNCTION	12	02-12-2022 06:10 PM	06-12-2022 04:31 PM	R   
☐	Reporting	AA20019740	1097693	RAMNAH BINTI MOHD RAMLI	Finalized	27-02-1969 Age: 53	Female	PET CT WHOLE BODY WITH RPT FFF	10	29-11-2022 06:35 PM	06-12-2022 02:03 PM	R   

After clicking, you should see a popup window where you must sign in with a specific **Username [1]** and **Password [2]**. This step ensures that only authorized personnel can view the medical test images. Finally, click on the **Sign In[3]** button to sign in.

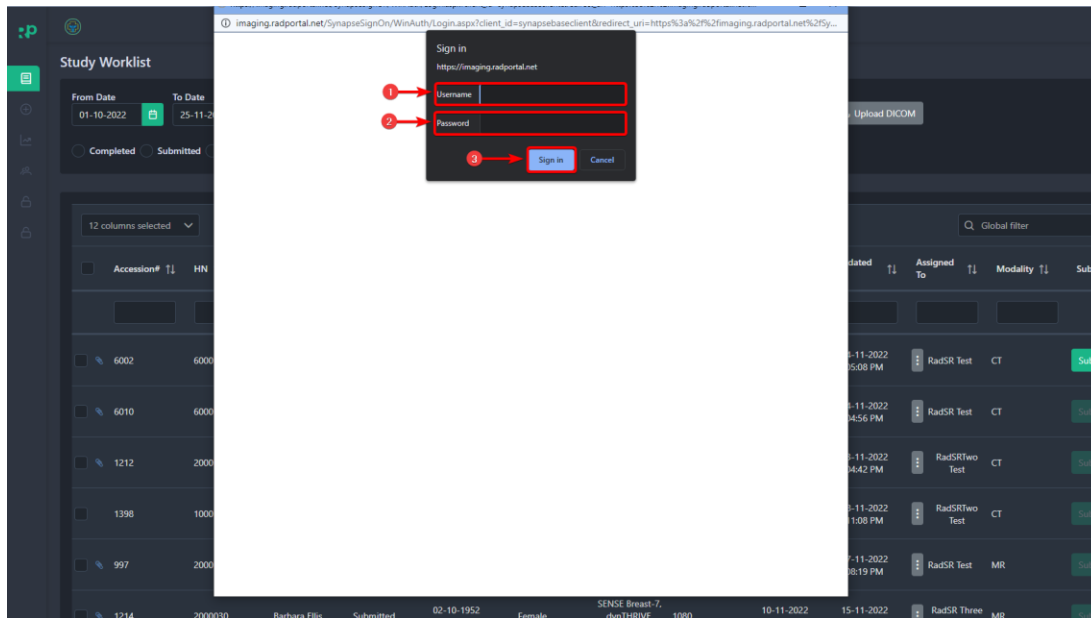

After successfully logging in to the system, you should be able to see all the test images of that specific Patient's medical test.



5.3.2 How to Decline Assigned Studies

If a hospital has assigned a study/case to you, then you should see the study/case on the study worklist page. To decline an assigned case, click on the **Cross Button**.


Global filter												
☐	Service Type	Accession#	HN	Patient Name	Status	Date of Birth	Gender	Exam Name	Image	Created On	Updated On	Action
☐												
☐	Consulting Complex Case	993	2000051	Michelle Mosconia	Submitted	27-09-1959 Age: 63	Female	MRI BREAST BILATERAL	665	06-12-2022 01:52 PM	06-12-2022 04:31 PM	  
☐	 Reporting	01509091655054	20150909162215	Kirstie Allie	Submitted	31-03-1971 Age: 51	Female	BREAST^Routine	12	02-12-2022 06:29 PM	06-12-2022 04:31 PM	R   
☐	 Reporting	1015	2000024	Melissa Dayton	Submitted	12-12-1955 Age: 67	Female	MRI CARDIAC FUNCTION	12	02-12-2022 06:10 PM	06-12-2022 04:31 PM	R   
☐	Reporting	AA20019740	1097693	RAMNAH BINTI MOHD RAMLI	Finalized	27-02-1969 Age: 53	Female	PET CT WHOLE BODY WITH RPT FEE	10	29-11-2022 06:35 PM	06-12-2022 02:03 PM	R   
☐	 Consulting	6002	6000002	Michael Cadens	Finalized	18-06-1961 Age: 61	Male	MR BREAST BIL FOR CA	11	25-11-2022 09:45 PM	05-12-2022 05:26 PM	R   

After clicking, you should see a confirmation popup window. You can click the **No [1]** button to stop the decline process and take you back to the study worklist page. Click on **Yes[2]** button to decline the case.

Global filter													
☐	Service Type	Accession#	HN	Patient Name	Status	Date of Birth	Gender	Exam Name	Image	Created On	Updated On	As To	Action
☐	Consulting Complex Case	993	2000051	Michelle Mosconia	Submitted	27-09-1959 Age: 63	Female	MRI BREAST BILATERAL	665	06-12-2022 01:52 PM	06-12-2022 04:31 PM		 
☐	Reporting	01509091655054	20150909162215	Kirstie Allie	Submitted	31-03-1971 Age: 51	Female	BREAST^Routine	12	02-12-2022 06:29 PM	06-12-2022 04:31 PM	Rad	 
☐	Reporting	1015	2000024	Melissa Dayton	Submitted	12-12-1955 Age: 67	Female	MRI CARDIAC FUNCTION	12	02-12-2022 06:10 PM	06-12-2022 04:31 PM	Rad	 
☐	Reporting	AA20019740	1097693	RAMNAH BINTI MOHD RAMLI	Finalized	27-02-1969 Age: 53	Female	PET CT WHOLE BODY WITH RPT FEE	10	29-11-2022 06:35 PM	06-12-2022 02:03 PM	Rad	 
☐	Consulting	6002	6000002	Michael Cadens	Finalized	18-06-1961 Age: 61	Male	MR BREAST BIL FOR CA	11	25-11-2022 09:45 PM	05-12-2022 05:26 PM	Rad	 
☐	Reporting	1212	2000098	Kevin Jeffers	Finalized	29-05-1956 Age: 66	Male	UNKNOWN	3	02-12-2022 10:07 PM	05-12-2022 05:25 PM	Rad	 
☐	Consulting	1400	1000156	Donna Bella	Finalized	01-02-1943 Age: 79	Female	BREAST^Routine	26	02-12-2022 04:56 PM	02-12-2022 05:53 PM	Rad	 
☐	Consulting	206	206206	Carotids	Finalized	06-11-1966 Age: 56	Female	UNKNOWN	10	25-11-2022 09:45 PM	28-11-2022 09:45 PM	Rad	 

Another way you can decline a study is by double-clicking on your preferred study/case. After double-clicking, you should see two popup windows.

In one of the windows, you will see a **Decline** button. Click on the **Decline** button to decline that specific study.

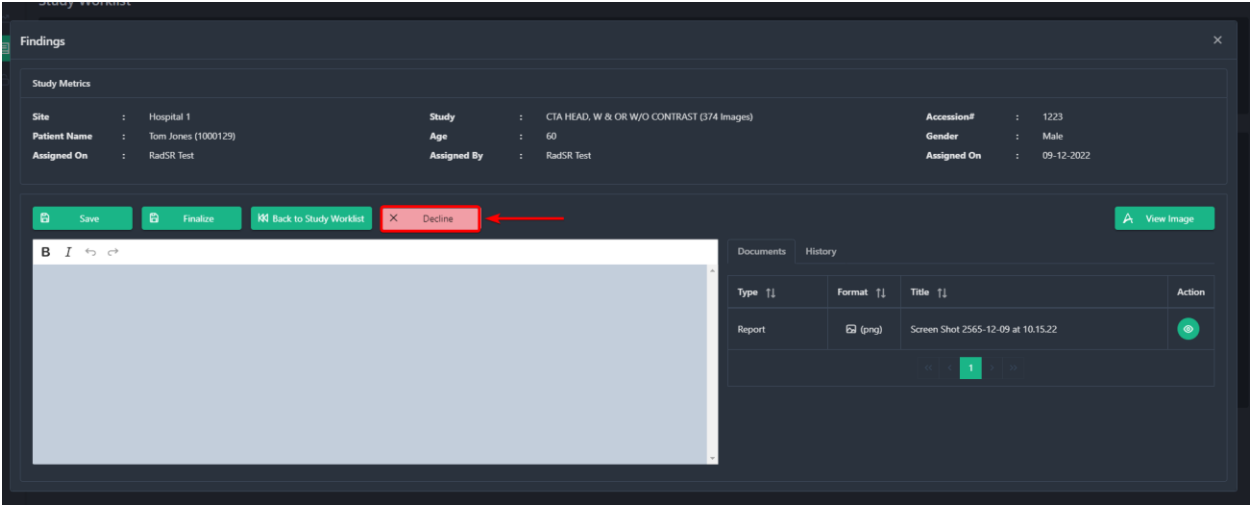

5.3.3 How to View Additional Documents of a Study

You can view any additional documents attached to open studies or specifically assigned to you. To review documents, you have to click on the **Paperclip Icon**.

☐	Service Type	Accession#	HN	Patient Name	Status	Date of Birth	Gender	Exam Name	Image	Create On	Action
☐	Reporting	1364	1000140	Rutledge Wood	Submitted	03-12-1964 Age: 58	Male	CT - CT CHEST ANGIO-PE	660	06-12-02:2	
☐	Reporting	330-CT-11-031670	2011092107	Bangkok Liver Tumor Us	Submitted	01-01-1970 Age: 52	Male	CT Abdomen and Pelvis w + w/o Contrast	720	06-12-02:1	
☐	Consulting Complex Case	1389	1000145	Sean Penn	Finalized	27-06-1960 Age: 62	Male	CT - CT ABD PELVIS W	730	06-12-02:1	
☒	Consulting Complex Case	180820SDOCR21	5399378	Photichanprasert Phumiphat	Submitted	18-04-1989 Age: 33	Male	Chest	1	06-12-01:5	
☐	Reporting	1024	2000104	Knox Pitt	Finalized	01-10-1982 Age: 40	Male	MRA HEAD	672	06-12-01:4	
☐	Consulting Complex Case	993	2000051	Michelle Mosconia	Submitted	27-09-1959 Age: 63	Female	MRI BREAST BILATERAL	665	06-12-01:5	
☒	Reporting	01509091655054	20150909162215	Kirstie Allie	Submitted	31-03-1971 Age: 51	Female	BREAST^Routine	12	02-12-06:2	
☒	Reporting	1015	2000024	Melissa Dayton	Submitted	12-12-1955 Age: 67	Female	MRI CARDIAC FUNCTION	12	02-12-06:1	

This paper **Paperclip** icon indicates that this specific study has additional documents attached to it.

RadPortal

Customer 1

RadSR Test

Dashboard

Worklist

User list


Unlock Worklist

Service Type	Accession#	HN	Patient Name	Status	Date of Birth	Gender	Exam Name	Image	Create On	Action
Reporting	1364	1000140	Rutledge Wood	Submitted	03-12-1964 Age: 58	Male	CT - CT CHEST ANGIO-PE	660	06-11-2018 02:2	  
Reporting	330-CT-11-031670	2011092107	Bangkok Liver Tumor Us	Submitted	01-01-1970 Age: 52	Male	CT Abdomen and Pelvis w + w/o Contrast	720	06-11-2018 02:1	 
Consulting Complex Case	1389	1000145	Sean Penn	Finalized	27-06-1960 Age: 62	Male	CT - CT ABD PELVIS W	730	06-11-2018 02:1	  
Consulting Complex Case	180820SDOCR21 5399378		Photichanprasert Phumiphat	Submitted	18-04-1989 Age: 33	Male	Chest	1	06-11-2018 01:5	  
Reporting	1024	2000104	Knox Pitt	Finalized	01-10-1982 Age: 40	Male	MRA HEAD	672	06-11-2018 01:4	  
Consulting Complex Case	993	2000051	Michelle Mosconia	Submitted	27-09-1959 Age: 63	Female	MRI BREAST BILATERAL	665	06-11-2018 01:5	  
Reporting	01509091655054 20150909162215		Kirstie Allie	Submitted	31-03-1971 Age: 51	Female	BREAST^Routine 12	12	02-11-2018 06:2	  
Reporting	1015	2000024	Melissa Dayton	Submitted	12-12-1955 Age: 67	Female	MRI CARDIAC FUNCTION	12	02-11-2018 06:1	  

5.3.4 How to Write Reports and Finalize Reports

After a hospital has assigned a study/case to you or you see an open case/study, you can write a report for these cases. To write a report, you have to double-click on your preferred study/case. After double-clicking on the case, you should see two popup windows.

The first window lets you review all the images of that particular study/case.

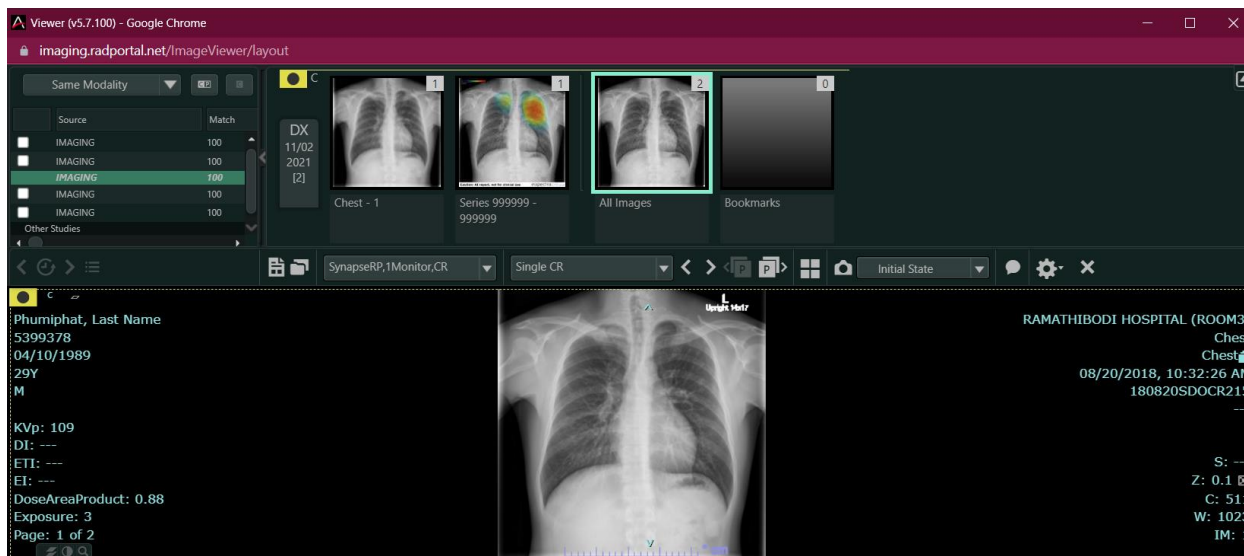

And on the second window, you can write a report in the blank section.

Findings

Customer 1 RadSR Test

Study Metrics

Site	: Hospital 1	Study	: Chest (1 Images)	Accession#	: 1808205DOGR215
Patient Name	: Photichanprasert Phumiphat (5399378)	Age	: 33	Gender	: Male
Assigned To	: RadSR Test	Assigned By	: RadSR Test	Assigned On	: 29-11-2022

Scoring

Complexity Rating: 0 1 2 3 4 5

Documents History

Type	Format	Title	Action
Report	(png)	Screenshot (2)	
Report	(png)	School Map On Nut 3F	
Report	(pdf)	Eric HSK2 Admission Ticket	


1 2 3 4 5 10

After writing your report, you can just save the report and come back to it later by clicking on the **Save** [1] button, or you can finalize the report by clicking on the **Finalize** [2] button. After finalizing your report, you can not change it afterwards. To go back to the Study worklist page, click on the **Back to Study Worklist** [3] button.

If you do not want to take on the study/case, then you can decline it by clicking on the **Decline** [4] button. On this page, you are also allowed to view all images.


To view images, click on **View Image** [7] button.

On this same window, you are allowed to provide Scores and Ratings based on the type of the study/case. You will not be able to see these two options if the case/study is

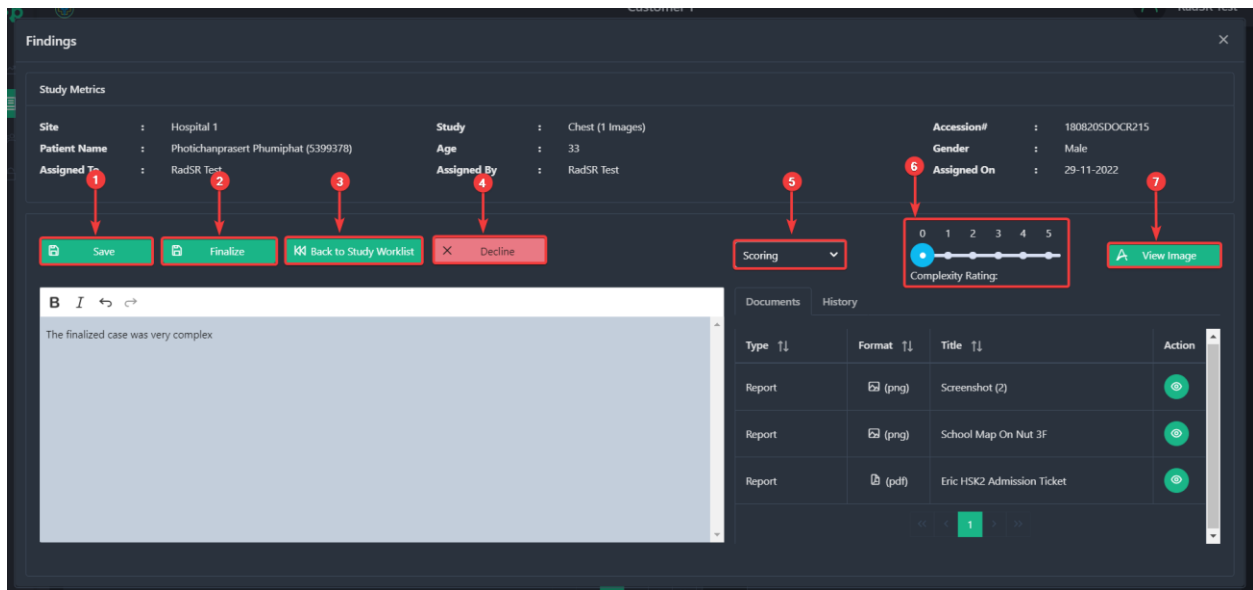

NOTE: You only can provide scores for Consulting Simple studies. You can provide both scores and complexity ratings for Consulting Complex studies.

On this window, you can access all additional documents attached by hospitals.

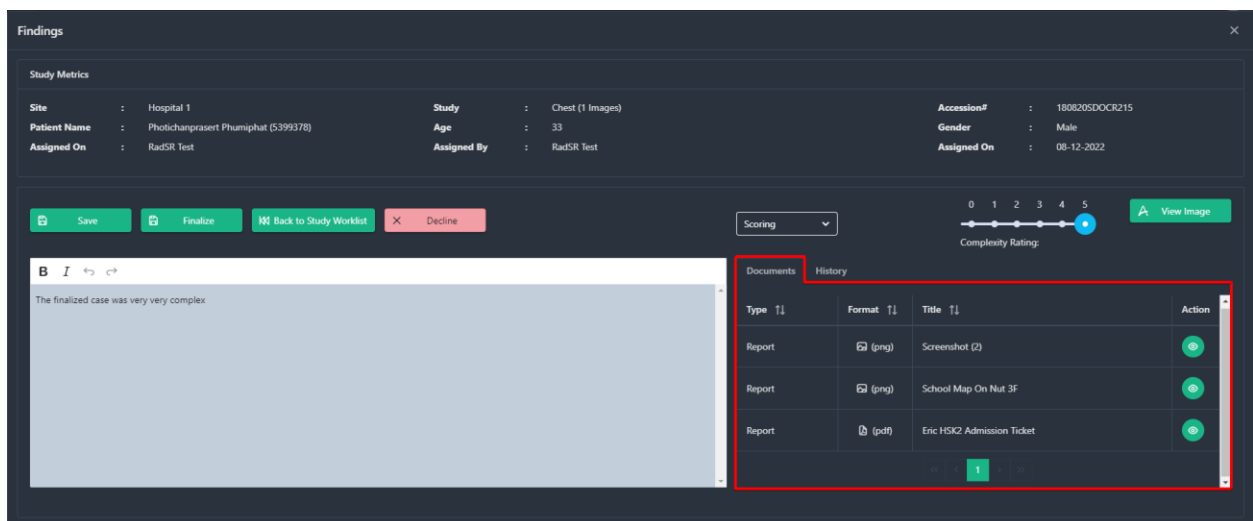

To review the document, click on the **Eye** Icon.

Findings

Study Metrics

Site	: Hospital 1	Study	: Chest (1 Images)	Accession#	: 180820SDOCR215
Patient Name	: Photichanprasert Phumiphat (5399378)	Age	: 33	Gender	: Male
Assigned On	: RadSR Test	Assigned By	: RadSR Test	Assigned On	: 08-12-2022

Save Finalize Back to Study Worklist Decline

Scoring 0 1 2 3 4 5 Complexity Rating: View Image

B I ↔ ⇄

The finalized case was very very complex

Documents History

Type	Format	Title	Action
Report	(png)	Screenshot (2)	View
Report	(png)	School Map On Nut 3F	View
Report	(pdf)	Eric HSK2 Admission Ticket	View

You can also see patient test history by clicking on the **History** tab.

Findings

Study Metrics

Site	: Hospital 1	Study	: Chest (1 Images)	Accession#	: 180820SDOCR215
Patient Name	: Photichanprasert Phumiphat (5399378)	Age	: 33	Gender	: Male
Assigned On	: RadSR Test	Assigned By	: RadSR Test	Assigned On	: 08-12-2022

Save Finalize Back to Study Worklist Decline

Scoring 0 1 2 3 4 5 Complexity Rating: View Image

B I ↔ ⇄

The finalized case was very very complex

Documents History

Patient History

Accession#	Exam Name	Exam Date	Actions
180917SDOCR211	Chest	05-12-2022 01:14 PM	View Report
180430SDOCR431	Chest	29-11-2022 06:35 PM	View Report
190513SDOCR360	Chest	07-12-2022 06:10 PM	View Report

In the history option, you can review Information such as Accession, Exam Name, and Exam Date, and under Actions, you can check history study images and previous Finalized reports.

To review Images of previous studies, click on the File icon.

The screenshot shows the 'Findings' interface with the following details:

- Study Metrics:**
 - Site: Hospital 1
 - Patient Name: Photichanprasert Phumiphat (5399378)
 - Assigned On: RadSR Test
 - Study: Chest (1 Images)
 - Age: 33
 - Assigned By: RadSR Test
 - Accession#: 1808205DOCR215
 - Gender: Male
 - Assigned On: 08-12-2022
- Actions:** Save, Finalize, Back to Study Worksheet, Decline.
- Scoring:** A dropdown menu and a complexity rating slider (0-5) with a blue dot at 4.
- Text Input:** A field with 'B I' and a text area containing 'The finalized case was very very complex'.
- Patient History Table:**

Accession#	Exam Name	Exam Date	Actions
1809175DOCR211	Chest	05-12-2022 01:14 PM	 
1804305DOCR431	Chest	29-11-2022 06:35 PM	 
1905135DOCR369	Chest	02-12-2022 06:10 PM	 

To review Finalized Reports of previous studies, click on the File icon.

The screenshot shows the 'Findings' interface with the following details:

- Study Metrics:**
 - Site: Hospital 1
 - Patient Name: Photichanprasert Phumiphat (5399378)
 - Assigned On: RadSR Test
 - Study: Chest (1 Images)
 - Age: 33
 - Assigned By: RadSR Test
 - Accession#: 1808205DOCR215
 - Gender: Male
 - Assigned On: 08-12-2022
- Actions:** Save, Finalize, Back to Study Worksheet, Decline.
- Scoring:** A dropdown menu and a complexity rating slider (0-5) with a blue dot at 4.
- Text Input:** A field with 'B I' and a text area containing 'The finalized case was very very complex'.
- Patient History Table:**

Accession#	Exam Name	Exam Date	Actions
1809175DOCR211	Chest	05-12-2022 01:14 PM	 
1804305DOCR431	Chest	29-11-2022 06:35 PM	 
1905135DOCR369	Chest	02-12-2022 06:10 PM	 

If you save your unfinished study you will see a status **Draft** in the upper section of the status. The studies that have the draft status can be edited again

Customer 1 RadSR Test

Global filter

	Service Type	Accession#	HN	Patient Name	Status	Date of Birth	Gender	Exam Name	Image	Action
	Consulting Complex Case	1623	1000184	Recist Liver	DRAFT SUBMITTED	10-12-1958 Age: 64	Male	Chest	578	
	Consulting Complex Case	0221214104332292	JFTEST04	TEST JFTEST04	DRAFT SUBMITTED	17-06-1952 Age: 70	Male	Chest	1	
	Consulting Complex Case	0221214105832620	JFTEST05	Unknown1 TEST JFTEST05	FINALIZED	16-05-1951 Age: 71	Male	CHEST	2	
	Reporting	6010	6000010	Mary Cadens	SUBMITTED	14-01-1962 Age: 60	Female	CT Chest W/ IV Contrast	10	
	Consulting Complex Case	0221209122823859	JFTEST01	TEST01 JFTEST01	FINALIZED	05-11-1961 Age: 61	Female	CHEST-AP	1	
	Reporting	1212	2000098	Kevin Jeffers	FINALIZED	28-05-1956 Age: 66	Male	Chest	10	
	Reporting	180302DIACD001	5399378	Photichanprasert Phumiphat	SUBMITTED	19-04-1989 Age: 33	Male	UNKNOWN	1	
	Reporting	180430SDOCR431	5399378	Photichanprasert Phumiphat	SUBMITTED	19-04-1989 Age: 33	Male	Chest	1	

5.3.5 View Study Reports

After you have finalized a report of a study, you can view it by clicking on the **File** Icon

Customer 1 RadSR Test

01-09-2022 08-12-2022 HOSPITAL 3, Hospital 1, Hospital 2 Search Today Yesterday Last 3 Days Auto Refresh (10 sec)

Submitted Finalized Submitted / Finalized


Global filter

	Service Type	Accession#	HN	Patient Name	Status	Date of Birth	Gender	Exam Name	Image	Created On	Updated On	Assigne	Action
	Reporting	330-CT-11-031670	2011092107	Bangkok Liver Tumor Us	FINALIZED	01-01-1970 Age: 52	Male	CT Abdomen and Pelvis w + w/o Contrast	720	06-12-2022 02:18 PM	06-12-2022 08:10 PM	RadSR	
	Reporting	1364	1000140	Rutledge Wood	FINALIZED	03-12-1964 Age: 58	Male	CT - CT CHEST ANGIO-PE	660	06-12-2022 02:26 PM	06-12-2022 06:08 PM	RadSR	
	Consulting Complex Case	1389	1000145	Sean Penn	FINALIZED	27-06-1960 Age: 62	Male	CT - CT ABD PELVIS W	730	06-12-2022 02:11 PM	06-12-2022 05:37 PM	RadSR	
	Reporting	1024	2000104	Knox Pitt	FINALIZED	01-10-1982 Age: 40	Male	MRA HEAD	672	06-12-2022 01:48 PM	06-12-2022 05:37 PM	RadSR	
	Reporting	AAJ20019740	1097693	RAMNAH BINTI MOHD RAMLI	FINALIZED	27-02-1969 Age: 53	Female	PET CT WHOLE BODY WITH RPT FILE	10	29-11-2022 06:35 PM	06-12-2022 02:03 PM	RadSR	
	Consulting	6002	6000002	Michael Cadens	FINALIZED	18-06-1961 Age: 61	Male	Cardiac^26_CoroCTA (Adult)	11	25-11-2022 09:45 PM	05-12-2022 05:26 PM	RadSR	
	Reporting	1212	2000098	Kevin Jeffers	FINALIZED	29-05-1956 Age: 66	Male	UNKNOWN	3	02-12-2022 10:07 PM	05-12-2022 05:25 PM	RadSR	
	Consulting	1400	1000156	Donna Bella	FINALIZED	01-02-1943 Age: 79	Female	BREAST^Routine	26	02-12-2022 04:56 PM	02-12-2022 05:53 PM	RadSR	
	Consulting	206	206206	Carotids Vincent	FINALIZED	06-11-1966 Age: 56	Female	UNKNOWN	10	25-11-2022 09:45 PM	28-11-2022 07:55 PM	RadSR	


5.3.6 Adding an Addendum

You can add information to your finalized report by writing an addendum. To write an addendum double, click on any finalized study/case.


After the study loads, you will see the **Addendum** button.

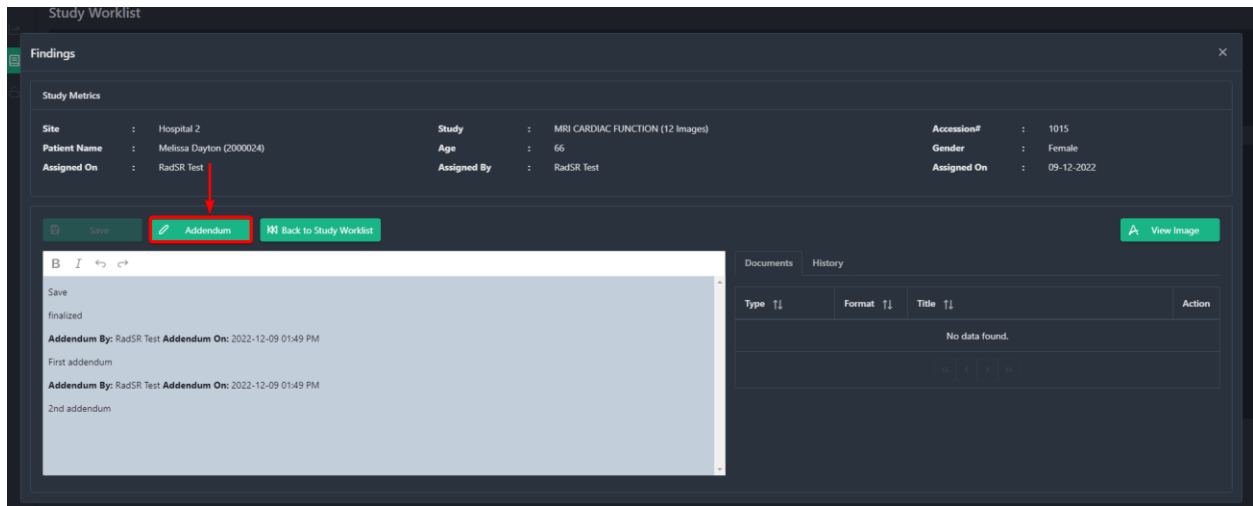


After clicking on the Addendum button, you should see a popup when the system shows you the previous **Finalized Reports and Addendums [1]** and a **Blank [2]** section where you can add the new Addendum.

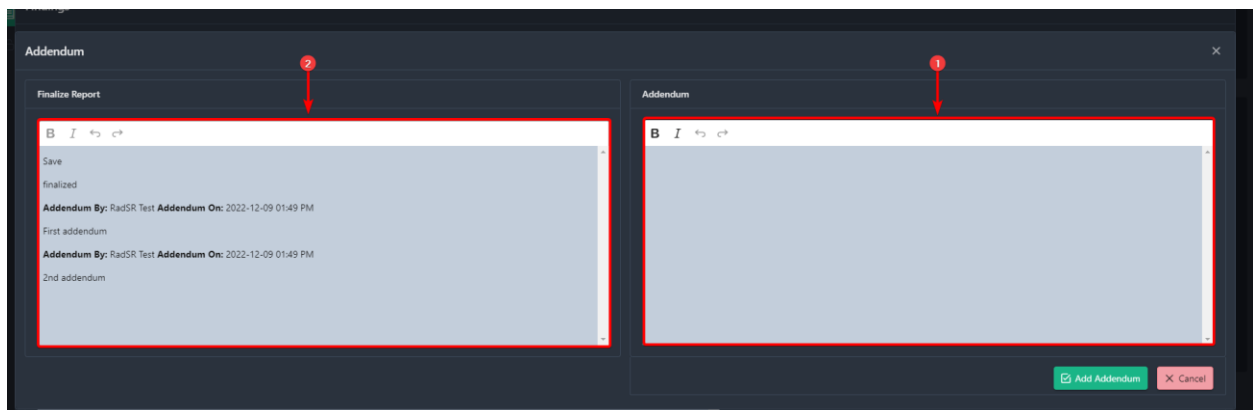


You can **bold** and *italicize* your addendum texts by clicking on the **B [1]** and *I [2]* icons.

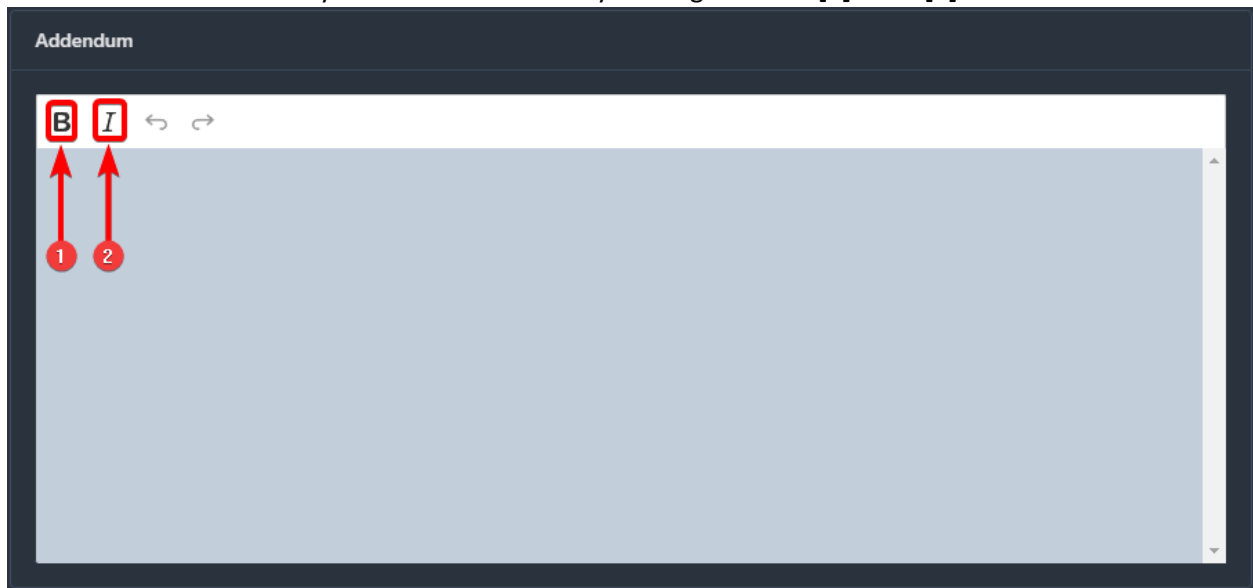

You can undo your edits by clicking on the **Undo [1]** icon and redo edits by clicking on the **Redo[2]** icon.

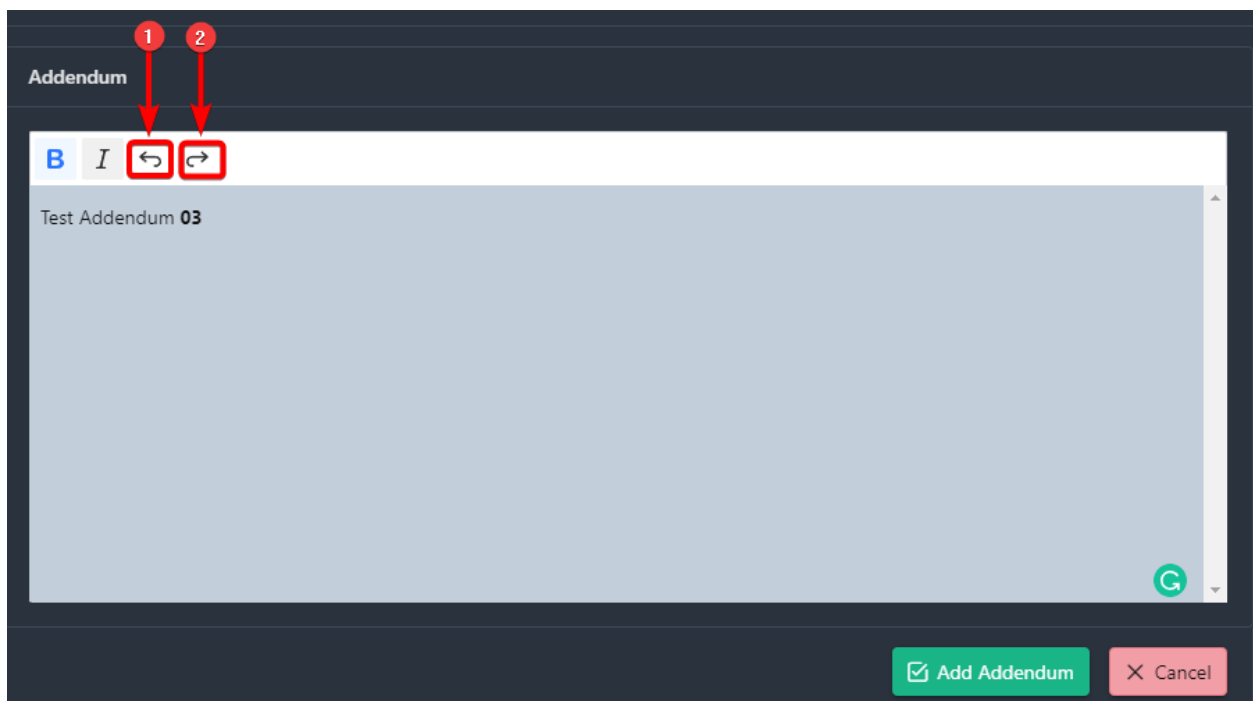

Finally, when you are done with the Addendum, click the **Add Addendum** button.

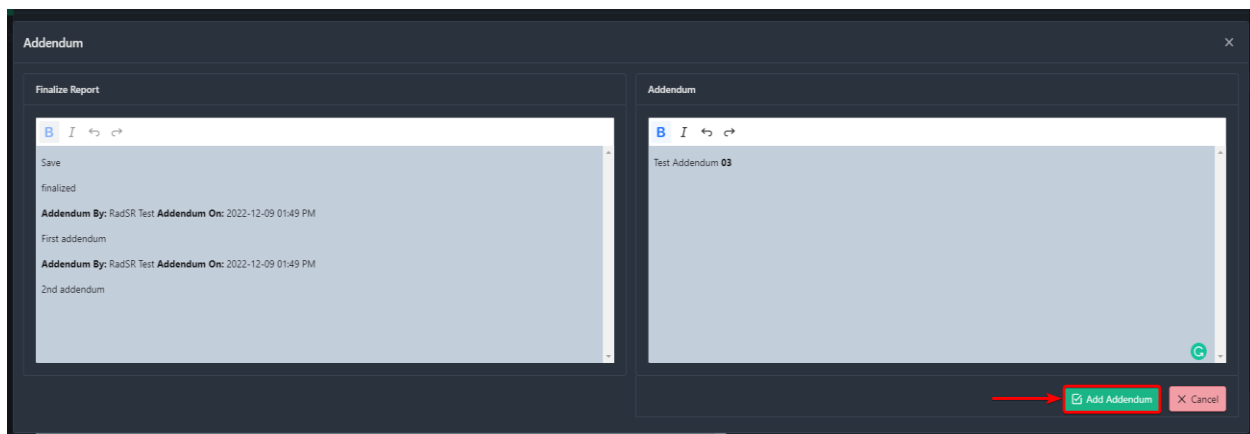

NOTE: You can not edit the Addendum again after clicking the **Addendum** button.

Now the system will take you back to the previous page, where you can review all addendums along with the Information about who added them and the date they were created.

